iluminage[®]

YOUTH ACTIVATOR

USER MANUAL
MODE D'EMPLOI

CONTENTS

CHAPTER 1: INTRODUCTION

CHAPTER 2: GETTING TO KNOW YOUR YOUTH ACTIVATOR

CHAPTER 3: WARNINGS

- 8 Contraindications
- 10 Precautions

CHAPTER 4: USING YOUR YOUTH ACTIVATOR FOR THE FIRST TIME

- 13 Setting Up Your **YOUTH ACTIVATOR**
- 13 Device Controls
- 14 Skin Sensitivity Test
- 14 How to Use Your **YOUTH ACTIVATOR**
- 16 Treatment Schedule

CHAPTER 5: TAKING CARE OF YOUR YOUTH ACTIVATOR

- 18 How to Clean and Store
- 18 Disposal

CHAPTER 6: TROUBLESHOOTING

CHAPTER 7: WARRANTY

CHAPTER 8: MANUFACTURER INFORMATION

INTRODUCTION

Congratulations!

1

NEED HELP? Please visit us at
iluminagebeauty.co.uk
or call +44 (0) 2033 020 606 (UK)

INTRODUCTION

Powered by breakthrough technology, the iluminage **YOUTH ACTIVATOR** is proven to effectively contour and tone the skin while reducing the appearance of fine lines and wrinkles. The unique combination of Radio-Frequency, Red LED and Infrared lights is proven to safely deliver energy deep within the dermal layer where collagen and elastin form, tightening and toning the face. When used with the **YOUTH ACTIVATOR** Contouring Serum, this technology results in smoother, younger, fresher-looking skin.

GETTING TO KNOW YOUR YOUTH ACTIVATOR



2

NEED HELP? Please visit us at
<http://www.iluminagebeauty.co.uk>
or call +44 (0) 2033 020 606 (UK)

YOUTH ACTIVATOR PACKAGE INCLUDES

- **YOUTH ACTIVATOR** device
- **YOUTH ACTIVATOR** Contouring Serum
- USB Cord
- Quick Reference Guide
- User Manual

DEVICE COMPONENTS

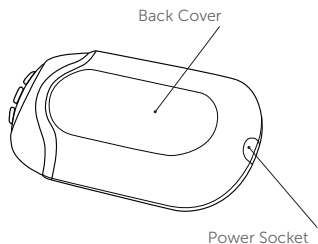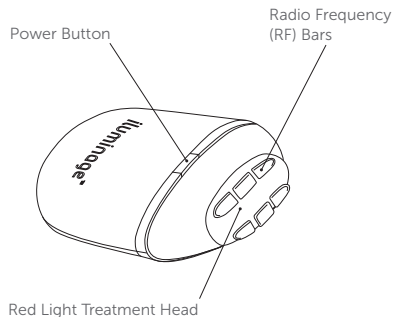WARNINGS

Read the following important safety information before using your **YOUTH ACTIVATOR**.

3

NEED HELP? Please visit us at
<http://www.iluminagebeauty.co.uk>
or call +44 (0) 2033 020 606 (UK)

WHO SHOULD NOT USE THE YOUTH ACTIVATOR (CONTRAINDICATIONS)

DO NOT use if you have an implanted pacemaker, arrhythmias or any other severe known heart disorder.

DO NOT use if you have any implantable metal device or body piercing in the treatment area.

DO NOT use if you have any blood clotting disorders or are taking blood thinning medications.

DO NOT use if you have a history of skin cancer or areas of potential skin malignancies. Use of the device has not been tested in areas with skin cancer and therefore the effects are unknown.

DO NOT use over thyroid gland (neck area).

DISCONTINUE USE if you encounter any discomfort or skin redness that lasts for more than 24 hours. You may have sensitivity to light or another condition. Consult with your healthcare provider.

WARNINGS

DO NOT use the **YOUTH ACTIVATOR** over any areas that are not recommended in this user manual.

DO NOT use the **YOUTH ACTIVATOR** on eyelids, eyebrows, lips or ears.

DO NOT use the **YOUTH ACTIVATOR** without the **YOUTH ACTIVATOR** Contouring Serum. The device will not operate properly if used without the serum or with another product.

DO NOT use if you are pregnant or breastfeeding.

DO NOT use if you have experienced a migraine in the last two years.

DO NOT use if you are sensitive to light:

- Some diseases can be related to photosensitivity including porphyria, polymorphic light eruption, solar urticarial and lupus.

- Some medications can be related to photosensitivity including: NSAIDS, steroids, antibiotics, phenothiazine, thiazide, diuretics, sulfonylurates, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha-Hydroxy Acids (AHAs), Beta-Hydroxy Acids (BHAs), Retin-A, Accutane, topical retinoids and/or Azelaic acid.

- Undergoing radiation or chemotherapy treatments.

These conditions can make the skin more sensitive to light and more susceptible to injury. Consult with your doctor.

DO NOT use if you have an active skin condition in the treatment area including sores, psoriasis, eczema, sunburns, herpes simples, open wounds or a rash. These conditions can make the skin more sensitive to light and more susceptible to injury. Consult with your doctor before use.

DO NOT use over dark brown or black spots such as large freckles, birthmarks, moles or warts on the area to be treated.

DO NOT use over or near a steroid injection received within the last six months. Doing so may result in a photosensitivity reaction.

Keep away from children. The device may contain small parts that may present as a choking hazard for young children.

In order to reduce the risk of electric shock that may cause injury or death:

- **KEEP AWAY** from water (including bath or shower). If the unit accidentally falls into water, unplug immediately.
- **DO NOT** place or store the device where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.
- **NEVER OPERATE** this device if it has a damaged cord, if it is not working properly, if it has been dropped or damaged, or submerged in water. Return the device to Illuminage Beauty Inc. for examination.

- DO NOT DISASSEMBLE or modify the device as this may cause unit damage, malfunction, electrical shock, fire or personal injury. There are no user-serviceable parts inside the unit.
- NEVER drop or insert any object into any openings on the device.
- KEEP CORD AWAY from heated surfaces.

To prevent fire or serious burns:

- TURN OFF immediately if device begins to overheat.
- KEEP AWAY from flammable aerosol products or where oxygen is being administered.

DO NOT operate this device in an environment in which other devices that radiate electromagnetic energy in an unshielded manner are being used. Portable or mobile Radio-Frequency communication equipment can affect medical electrical equipment.

PRECAUTIONS

STOP IMMEDIATELY if you experience any adverse effect or have a concern regarding use and consult your doctor.

READ ALL INSTRUCTIONS prior to use.

USE this device only for its intended use as described in this manual.

FAILURE TO USE and maintain the **YOUTH ACTIVATOR** in accordance with the instructions outlined in this manual will invalidate the warranty.

DO NOT look directly at the light source of the treatment head while it is in operation.

DO NOT shine the light source into someone else's eyes.

DO NOT use if you have increased sensitivity to light or you are unsure whether you have increased sensitivity

to light. We recommend that you perform the Light Sensitivity Test described on page 14.

Prolonged use beyond the recommended treatment time may cause the treatment surface of the device to overheat and may cause adverse effects such as discomfort or prolonged redness in the area.

DO NOT operate this device in an environment in which other devices radiate electromagnetic energy in an unshielded manner are being used. Portable or mobile Radio-Frequency communication equipment can affect medical electrical equipment.

DO NO scratch, mark or otherwise damage the light-emitting surface.

INSPECT this device prior to each use. Check for damage.

OPERATE, transport and store this device per the manner specified in this user manual.

DISPOSE of all products in accordance with local and national regulations and codes.

USING YOUR YOUTH ACTIVATOR FOR THE FIRST TIME

4

NEED HELP? Please visit us at
<http://www.iluminagebeauty.co.uk>
or call +44 (0) 2033 020 606 (UK)

SETTING UP YOUR YOUTH ACTIVATOR

1. Unpack the device and USB cord.
2. Remove the back cover of the device and remove the battery tab underneath it.
3. Connect the USB cable to the charging port at the bottom of the device and plug the other end into any certified USB charging source.
4. A green battery indicator will blink to indicate the device battery is charging.
5. Keep it plugged in and charge for up to 4 hours or until the battery light indicator is steady green.
6. The device is now ready for use. Unplug the USB cable from the device and remove the treatment cap to prepare the device for treatment.

YOUTH ACTIVATOR DEVICE CONTROLS

The **YOUTH ACTIVATOR** has a treatment indicator that appears on the front of the device, above the on/off button. It also has a battery indicator near the USB charging port.

The device indicators provide the following information:

TREATMENT INDICATOR LIGHT	
Steady Blue Light	The device has been activated and ready for treatment 1 Blue Light – Device is on LOW treatment setting 2 Blue Lights – Device is on HIGH treatment setting
Flashing Blue Light	Device error. Refer to Troubleshooting Table for resolution.
No light	The device is off

BATTERY LIGHT INDICATOR

Steady Green Light	The device battery is fully charged
Blinking (during treatment)	The device battery is low and requires charging
Blinking (during charging)	The device battery is charging
No light	There is enough battery for at least one full treatment

VIBRATION INDICATOR

Quick double vibration	Treatment time for a zone has elapsed (4 minutes). Move device to new treatment zone to continue treatment.
Constant vibration	Device error. Refer to Troubleshooting Table for resolution.

SKIN SENSITIVITY TEST


If you are using the **YOUTH ACTIVATOR** for the first time, the skin sensitivity test should be performed on the LOW energy level prior to your first treatment.

- Prepare your skin for the test and follow the treatment instructions.
- Complete a test treatment on your inner forearm for 4 minutes.
- Wait 30 minutes and check your skin for any reaction. If the skin appears normal, then you are ready to perform a full treatment.
- If the treated spot appears red or irritated, please wait 24 hours and repeat the test. If it persists, consult with your doctor.

HOW TO USE YOUR YOUTH ACTIVATOR

STEP 1: PREP YOUR SKIN

Wash and dry your face. Clean, dry skin is critical for optimal results. Ensure that all creams, cleansers and moisturisers etc. are removed - they will stop the device working.

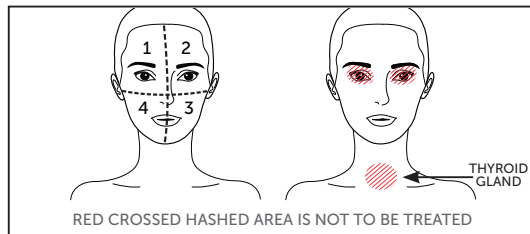
STEP 2: PLANNING YOUR TREATMENT

Divide your face into four treatment areas as shown in the above illustration. Press the power button briefly until the device turns on. One blue light will illuminate, indicating that the device is in LOW energy setting. To increase the setting to HIGH, press the power button once more. A second blue light will illuminate. It is recommended to start your first treatment on the low level and move up to high if you do not feel any discomfort.

STEP 3: APPLY CONTOURING SERUM

Apply a generous layer of the **YOUTH ACTIVATOR** Contouring Serum to the area of skin that you plan to treat. Do not use any other product with the **YOUTH ACTIVATOR** device during treatment as it will not work. Reapply the serum as needed throughout the treatment process, as the device will not function properly without a sufficient amount of serum.

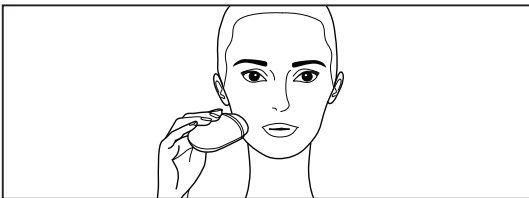STEP 4: TREAT YOUR SKIN

Place the tip of the applicator in full contact with your skin. The treatment head will illuminate red, indicating that the energy is being delivered. You will feel a gentle warming sensation. Move the device in a constant circular motion to treat the selected zone. The device will vibrate after four minutes, signaling you to move to the next treatment zone. After all four zones have been treated, the device will automatically power off.

STEP 5: POST TREATMENT CARE

Once your treatment has ended, you may either remove the serum with a warm, damp cloth or massage it in to your skin for added hydration. If your skin looks slightly red after treatment, this should disappear within 24 hours.

WHAT TO EXPECT DURING TREATMENT

During treatment with the **YOUTH ACTIVATOR**, you may feel a gentle warming sensation. You should always start your treatment on the LOW energy level and if you are comfortable, increase to HIGH for greater efficacy. You may experience some small heat sensations. These are completely normal and are a sign the product is working.


TREATMENT SCHEDULE

Skin rejuvenation is a gradual process. Actual results will vary from person to person. For best results, follow the full protocol of treating the targeted areas 5 times a week for 4 weeks and twice weekly for an additional 4 weeks. After this 8 week protocol, you can use your **YOUTH ACTIVATOR** to maintain results as you wish (up to maximum of 5 times a week on going).

TAKING CARE OF YOUR YOUTH ACTIVATOR

5

NEED HELP? Please visit us at
<http://www.iluminagebeauty.co.uk>
 or call +44 (0) 2033 020 606 (UK)

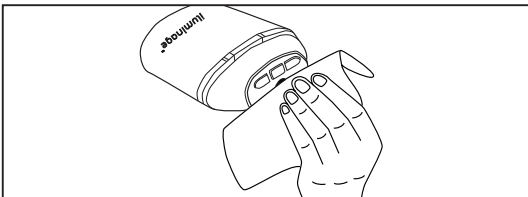HOW TO CLEAN AND STORE YOUR YOUTH ACTIVATOR

Clean the treatment window and outside surface of the **YOUTH ACTIVATOR** using ONLY non-alcoholic based wet wipes or a soft tissue after every treatment. Ensure all serum has been wiped off the treatment surface and device. Keep the cap on the treatment window when device is not in use.

Store the **YOUTH ACTIVATOR** in a cool, dry place; keep out of reach of children.

TRAVELLING WITH YOUR YOUTH ACTIVATOR

The **YOUTH ACTIVATOR** is designed to be used globally and includes a USB cord for charging. Make sure that the **YOUTH ACTIVATOR** does not get wet while traveling, and do not use it in a wet environment.

DISPOSAL/RECYCLING

Please follow your local regulations for disposal of electric systems when recycling or disposing of your **YOUTH ACTIVATOR**. If you have questions, please contact us at www.iluminagebeauty.co.uk.

TROUBLESHOOTING

6

NEED HELP? Please visit us at
<http://www.iluminagebeauty.co.uk>
 or call +44 (0) 2033 020 606 (UK)

In the unlikely event that the device fails to perform as intended, the following notes will help to identify potential problems with the device and its setup.

IF	DO THIS
Battery light indicator is off while the device is plugged into a USB charging source.	Ensure the DC input is connected to the charging port of the device and the other end is connected into a USB charging source. It is recommended to use a 1A charger for best charging performance.
Device turns off after 16 minutes.	This indicates the end of your treatment.
Device will not start when the power button has been pressed.	Ensure you press the on/off button OR ensure the device is fully charged. Charge the device for at least 4 hours by connecting it to a charger with the provided USB cable.
Treatment indicator light is glowing steady but treatment head is not illuminated.	The device treatment surface must be in full contact with the skin with generous amount of serum. Re-apply the serum and ensure you make full skin contact between the device and skin.
Device shuts off in mid-treatment.	Low Battery: check that the device battery is not low.
Device vibrates and the blue lights begin flashing during treatment.	The device is over heated. Shut down the device for 5-10 minutes to allow the device to cool off before using again. If problem persists, please contact our Customer Service for further assistance.
Battery indicator light is blinking.	Device battery is low and cannot complete a full treatment. Charge the device before attempting to treat.

WARRANTY

7

NEED HELP? Please visit us at <http://www.iluminagebeauty.co.uk> or call +44 (0) 2033 020 606 (UK)

ILUMINAGE YOUTH ACTIVATOR WARRANTY

NOTE: If your **YOUTH ACTIVATOR** is under warranty with your local distributor it will prevail over the below warranty terms.

Your **YOUTH ACTIVATOR** is warranted to be free from defects in materials and workmanship for a period of 2 years (or otherwise stated by your local distributor/retailer) from the original purchase date under normal use. This warranty extends only to the original purchaser with original proof of purchase and only when purchased from an authorized **YOUTH ACTIVATOR** retailer or reseller. If the product should become defective within the warranty period, please contact your local retailer/distributor or go to www.iluminagebeauty.com for further assistance.

iluminage Beauty, Inc. ("iluminage") reserves the right to replace a defective product with the most comparable product currently available. Your exclusive remedy for any defective product covered by this warranty is limited to the repair or replacement of the defective product. iluminage's entire liability for any defective product shall in no event exceed the original purchase price of the defective product.

TO THE EXTENT PERMITTED BY LAW, iluminage SHALL NOT BE RESPONSIBLE FOR ANY DIRECT OR INDIRECT DAMAGES RESULTING FROM THE USE OR MISUSE OF THIS PRODUCT INCLUDING ANY GENERAL, INCIDENTAL, SPECIAL, OR CONSEQUENTIAL DAMAGES—TO THE MAXIMUM EXTENT PERMITTED BY THE LAW. THE WARRANTY PROVIDED IN THIS DOCUMENT IS THE SOLE AND EXCLUSIVE WARRANTY PROVIDED BY iluminage TO THE PRODUCT AND IT REPLACES ANY REPRESENTATION AND/OR WARRANTY INCLUDING FITNESS FOR ANY PARTICULAR USE AND MERCHANTABILITY.

This warranty does not cover products determined by iluminage to be damaged by any of the following: accident, misuse, abuse or alteration, servicing by unauthorized personnel, use with unauthorized accessories or use other than as instructed, connecting to an incorrect current or voltage and/or any other conditions beyond iluminage's control.

MANUFACTURER INFORMATION

8

NEED HELP? Please visit us at
<http://www.iluminagebeauty.co.uk>
or call +44 (0) 2033 020 606 (UK)

LABELS AND SYMBOLS

The **YOUTH ACTIVATOR** device shows information pertaining to the device's part number, serial number, manufacturer information and the electrical and regulatory standards with which the **YOUTH ACTIVATOR** complies.

This device is not suitable for use in the presence of a flammable anesthetic mixture with air or with Oxygen or Nitrous Oxide.

WEEE - Waste Electrical and Electronic Equipment

Follow Operating Instructions

Class II Equipment

ETL Mark for USA and Canada

Do Not Use in Wet Environments

Degree of protection against electric type BF applied part

CE Mark

Illuminage Ltd. Yokneam Illit 2069200, IL

mdi Europa GmbH
Langenhagener Str. 71
D-30855 Hannover-Langenhagen

YOUTH ACTIVATOR SPECIFICATIONS

Radio Frequency Output: 1 MHz

RF Power: Up to 10 W

Light Spectrum: 660 nm; 850 nm

Light Output Power Density: Up to 50 mw/cm²

Maximal Optical Radiation: 100mW

Rating: 5VDC, 0.5-4A

Dimensions: 115 x 66.5 x 33mm (without cap)

Weight: 150g

ENVIRONMENT CONDITIONS

Operating/between-use-storage Specification:

Ambient Temperature: 5°C To 30°C; Relative Humidity: 15% To 93%, Non-Condensing; 700-1060 hPa

Transport/Store Environmental Specifications:

Ambient Temperature: 0°C To 60°C; Relative Humidity: 10% To 90%, Non-Condensing ;700-1060 hPa

IP22 Protected against access to hazardous parts with a finger, protected against vertically falling water drops when enclosure tilted up to 15 degrees.

APPLICABLE STANDARDS

SAFETY STANDARDS		EMC STANDARDS
EN60335:1 2002 (A1:2004; A2:2009)	IEC60601-2- 57:2011	EN61000-6- 1:2007
EN60335-2- 27:2009	ETL LISTED	EN61000-6- 3:2007
IEC62471:2006	CONFORMS TO UL STD	IEC60601-1- 2:2007
IEC60601-1, 3rd edition	UL60601-1 CERT TO CSA STD.	
IEC60601-1-11	C22.2 No. 601.1	

YOUTH ACTIVATOR CONTOURING SERUM

INGREDIENTS

Aqua (Water), Glycerin, Carbomer, Phenoxyethanol, Sodium Chloride, Sodium Hydroxide, Ethylhexylglycerin, Chondrus Crispus Extract, Nicotiana Sylvestris Leaf Cell Culture FOR EXTERNAL USE ONLY.

EMC INFORMATION

Note: Portable and mobile RF communications equipment can affect MEDICAL ELECTRICAL EQUIPMENT.

Table 1

GUIDANCE AND MANUFACTURER'S DECLARATION – ELECTROMAGNETIC EMISSIONS		
The Youth Activator device is intended for use in the electromagnetic environment specified below. The customer or the user of the Youth Activator device should assure that it is used in such an environment		
EMISSION TEST	COMPLIANCE	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
RF emissions CISPR 11	Group 1	The Youth Activator device must emit electro-magnetic energy in order to perform its intended function. Nearby electronics equipment may be affected
RF emissions CISPR 11	Class B	The Youth Activator device is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

Table 2

GUIDANCE AND MANUFACTURER'S DECLARATION – ELECTROMAGNETIC EMISSIONS			
The Youth Activator device is intended for use in the electromagnetic environment specified below. The customer or the user of the Youth Activator device should assure that it is used in such an environment			
IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 2, 4, 6 kV contact ± 2, 4, 8 kV air	N/A

Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2 kV for power supply lines ± 1 kV for input/output lines	Mains power quality should be that of a typical household environment
Surge IEC 61000-4-5	± 1 kV line(s) to line(s) ± 2 kV line(s) to earth	± 0.5 and 1 kV line(s) to line(s) ± 0.5, 1 and 2 kV line(s) to earth	Mains power quality should be that of a typical household environment
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5 % 100V (>95 % dip in 100/240V for 0.5 cycle 40 % 220V (60 % dip in 100/240V) for 5 cycles 70 % 220V (30 % dip in 100V/240) for 25 cycles <5 % 100V (>95 % dip in 100/240V) for 5 s	<5 % 100V (>95 % dip in 100/240V for 0.5 cycle 40 % 220V (60 % dip in 100/240V) for 5 cycles 70 % 220V (30 % dip in 100V/240) for 25 cycles <5 % 100V (>95 % dip in 100/240V) for 5 s	N/A
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m 50Hz 3A/m 60Hz	

Table 4

GUIDANCE AND MANUFACTURER'S DECLARATION – ELECTROMAGNETIC EMISSIONS			
The Youth Activator device is intended for use in the electromagnetic environment specified below. The customer or the user of the Youth Activator device should assure that it is used in such an environment			
IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE

Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 V	Portable and mobile RF communications equipment should be used no closer to any part of the Model 006, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2,5 GHz	3 V/m	Recommended separation distance $d=1.2\sqrt{P}$ $d=1.2\sqrt{P}$ 80 MHz to 800 MHz $d=2.3\sqrt{P}$ 800 MHz to 2,5 GHz where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, a should be less than the compliance level in each frequency range. b Interference may occur in the vicinity of equipment marked with the following symbol: 
Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.			

Table 6

RECOMMENDED SEPARATION DISTANCES BETWEEN PORTABLE AND MOBILE RF COMMUNICATIONS EQUIPMENT AND THE YOUTH ACTIVATOR DEVICE			
The Youth Activator device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Youth Activator device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Youth Activator device as recommended below, according to the maximum output power of the communications equipment.			
RATED MAXIMUM OUTPUT POWER OF TRANSMITTER W	SEPARATION DISTANCE ACCORDING TO FREQUENCY OF TRANSMITTER M		
	150 KHZ TO 80 MHZ D=1.2√P	80 MHZ TO 800 MHZ D=1.2√P	800 MHZ TO 2,5 GHZ D=2.3√P
.01	12	12	23
1	38	38	73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

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© 2017 ALL RIGHTS RESERVED. THE **YOUTH ACTIVATOR** IS
MANUFACTURED BY ILUMINAGE LTD. (YOKNEAM ILLIT 2069200,