Copyright 2011 © Home Skinovations Ltd. All rights reserved.
Print date: March 2013

Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd. is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd.

Data is subject to change without notification.

Home Skinovations Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd.

Specifications are subject to change without notice.

Home Skinovations Ltd. Tavor Building, Shaar Yokneam POB 533, Yokneam 20692, ISRAEL www.silk'n.com, info@silk'n.com

Table of Contents

1. Before You Start	2
1.1. Silk'n FaceTite™ device description	2
1.2. Package contents	2
2. Silk'n FaceTite™ Intended Use	3
3. Safety with Silk'n FaceTite™	3
4. Counter indications	3
5. Do not use Silk'n FaceTite™	4
6. Get Started!	5
7. Treatment regimen	6
7.1. What to Expect when treating with Silk'n FaceTite™?	7
7.2. After treating with Silk'n Glide™	7
8. Maintenance of Silk'n FaceTite™	8
8.1. Cleaning Silk'n FaceTite™	8
9. Troubleshooting	8
10. Customer Service	9
11. Specifications	9
12. Labeling	118

1. Before You Start

Before using Silk'n FaceTite™ for the first time, please read this Instructions for Use booklet in its entirety, and watch the Silk'n FaceTite™ training video.

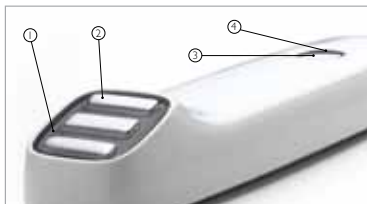
Pay particular attention to sections on device use procedures, device operation, and after-use procedures. We recommend you re-familiarize yourself with this Instructions for Use booklet before each use of Silk'n FaceTite™.

Silk'n FaceTite™ is a powerful electrical device. As such, it should be used with special attention to safety. Please read all warnings and safety precautions before use, and strictly follow them when using Silk'n FaceTite™

1.1 Silk'n FaceTite™ device description

Silk'n FaceTite™ is an anti-aging face treatment device which utilizes HT™ technology. This technology combines three kinds of energy: Infrared Radiation (IR) that affects the upper layer of skin for improved skin texture, Low Level Laser Therapy (LLLT) that reaches into the dermis for rejuvenation and collagen restoration, and focused Bi-Polar Radio Frequency energy (Bi-Polar RF) that penetrates deep into the sub-dermal fat for deep tissue tightening and wrinkle reduction

Silk'n® FaceTite has a TREATMENT SURFACE (1) with three ELECTRODES (2) and a CENTRAL BUTTON (3) with five surrounding INDICATOR LIGHTS (4)

- ① TREATMENT SURFACE
- ② ELECTRODES
- ③ CENTRAL BUTTON
- ④ INDICATOR LIGHTS

1.2. Package contents

Upon opening the Silk'n FaceTite™ package, you will find the following parts:

- Silk'n® FaceTite™ Device
- A bottle of Slider gel
- An AD/DC wall mount adaptor
- This Instruction for Use bookle
- A Quick Guide

2. Silk'n FaceTite™ Intended Use

Silk'n® FaceTite™ is intended for skin texture improvement and skin tightening, collagen restoration and wrinkle reduction.

Please note: Silk'n® FaceTite™ is intended both for self-treatment, and for the treatment of others.

3. Safety with Silk'n FaceTite™

Silk'n FaceTite™ protects your skin.

Silk'nFaceTite™ comes with a THERMAL SENSOR designed to measure skin temperature and stabilize the emitted heat energy.

To prevent overheating, a unique built-in heat stabilizer inside the Silk'n FaceTite™ measures the treated skin temperature during session. As a rise in the treated skin temperature is detected (42°C-43°C), the stabilizer will automatically stop emitting heat, until skin temperature lowers enough for the stabilizer to resume heating.

Note: Make sure the device is tightly pressed to your skin before starting treatment!

4. Counter Indications

Important Safety Information – Read Before Use!

Silk'n FaceTite™ is not designed for everyone. Please read and consider the information in the following section before use. For further information and personalized advice you may also visit www.silkn.com or your local Silk'n™ domain.

- Do not use Silk'n FaceTite™ if you have a Pacemaker or internal defibrillator, or any other active electrical implant anywhere in the body.
- Do not use Silk'n FaceTite™ if you have a permanent implant in the treated area such as metal plates and screws, silicone implants or an injected chemical substance.
- Do not use Silk'n FaceTite™ if you have a current history or condition of cancer, or pre-malignant moles.
- Do not use Silk'n FaceTite™ if you have received radiation therapy or chemotherapy treatments within the past 3 months.
- Do not use Silk'n FaceTite™ if you are pregnant or nursing.
- Do not use Silk'n FaceTite™ if you have a history of diseases stimulated by heat, such as recurrent Herpes Simplex in the treatment area. If so you may treat only following a prophylactic regime.
- Do not use Silk'n FaceTite™ if you suffer from poorly controlled endocrine disorders, such as diabetes and thyroid conditions.
- Do not use Silk'n FaceTite™ if you suffer from epilepsy
- Do not use Silk'n FaceTite™ if you have a disease related to photosensitivity, such as Porphyria, Polymorphic Light Eruption, Solar Urticaria, Lupus, etc.
- If you have any other condition which in your physician's opinion would make it unsafe for you to be treated.
- If your skin was treated with a physician-prescribed drug within the past 6 months, please consult with your physician.
- Do not use Silk'n FaceTite™ if you have any active condition in the treatment area, such as sores, psoriasis.
- Do not use Silk'n FaceTite™ if you had any surgery in the treated area within 3 months prior to treatment.

5. Do not use Silk'n FaceTite™:

With attachments or accessories are not recommended by Home Skinovations Ltd.

- If it has a damaged cord or plug.
- If you see or smell smoke when it is in use.
- If it is not working properly or if it appears damaged.
- If its THERMAL SENSOR is cracked, or broken.
- If the outer shell is cracked or is coming apart.
- While bathing or if it becomes damp or wet. Unplug immediately if it has fallen into water.

Warning: Keep Silk'n FaceTite™ away from water! Do not place or store Silk'n FaceTite™ where it can fall, be pushed or placed into a tub, sink or any other vessel containing water.

6. Get Started!

1. Remove Silk'n FaceTite™ DEVICE and other components from box.
2. Plug the AC/DC ADAPTOR cord into the Silk'n FaceTite™ socket.
3. Plug the AC/DC ADAPTOR into an electrical outlet. All five indicator lights will cycle repeatedly one after the other.
4. The device is now in STBY mode and is ready to start.
5. Press the CENTRAL BUTTON. The device will turn on to energy level 1 automatically. One INDICATOR LIGHT will turn on.

Setting Energy Levels:

- Silk'n® FaceTite™ offers 5 levels of energy represented by 5 ENERGY LEVEL INDICATOR LIGHTS.
- The selected energy level determines the intensity of Silk'n® FaceTite™ energy delivered to your skin.
- To set an energy level, press the CENTRAL BUTTON repeatedly to move between energies, from 1 to 5 and back to 1 again.
- When the desired energy is set, the ENERGY LEVEL INDICATOR LIGHTS will indicate the level chosen.
- To put FaceTite™ back in standby mode, give the CENTRAL BUTTON a long press. All five lights will cycle continuously.
- If you wish to change energy levels in mid-treatment, press the button again. This will change energy level mode on the spot and allow you to treat the same or different areas with different energy levels according to your needs.

6. After setting the energy level for the treatment, spread a thin layer of Silk'n Slider Gel over the area you intend to treat.

Please Note: Slider Gel (or any other water-based clear gel) facilitates smooth movement of the treatment surface over the skin and enables conductivity. To ensure safe and efficient treatment, always start by applying Slider Gel!

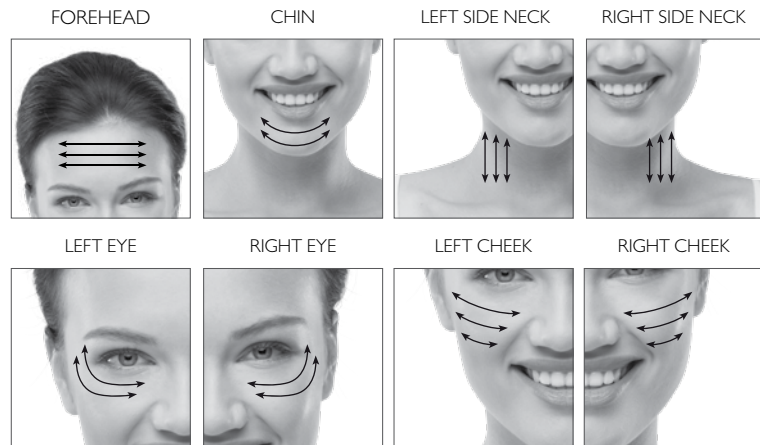
7. Press the device to the designated skin area. You will feel a mild sensation of warmth. Move the device in gentle massage- like linear motions on the area you wish to treat. Progress slowly and carefully.

Please Note: In order for the device to work properly, it must be firmly attached to the skin!

8. Immediately after treatment apply a facial lotion or mask.

7. Treatment regimen

Treating with Silk'n FaceTite™ is easy! Just follow the illustrations below, and you are on your way to a younger look! For ultimate results with FaceTite™, treat each of the 8 facial areas separately:

Please note: For your first treatment, you'll want to determine which FaceTite™ energy level is right for you. Start at energy level 1 for a few minutes and raise the energy level gradually, in order to find the highest energy level that you are still comfortable with.

Test each treatment area separately this way.

In following sessions, simply set FaceTite™ to the last comfortable energy level setting.

7.1. What to Expect when treating with Silk'n FaceTite™?

For many people, using Silk'n FaceTite™ may be their first experience with a heat-based device designed for home use. Silk'n FaceTite™ is simple to use! During a Silk'n FaceTite™ session it is normal to experience and feel:

- **A Sensation of Warmth** – During each session it is normal to feel a mild sensation of warmth from the RF and heat energy.
- **Some Mild Red or Pink Color** – During and just after your Silk'n FaceTite™ session it is not uncommon to see some very mild, pink-like coloration of the skin.
- However if you see full redness of the skin, or blistering stop using Silk'n FaceTite™ immediately.

For best results keep the following schedule

Start with sessions twice a week for 10 weeks to give your body an intensive treatment. By then you will be able to see a change in your skin texture. Maintain results by treating once every two weeks.

Please note: optimal results depend on accurately following the treatment regimen!

7.2. After treating with Silk'n FaceTite™

1. Turn Silk'n® FaceTite™ off by unplugging the AC/DC ADAPTOR cord from the electrical outlet.
2. Clean your Silk'n® FaceTite device, especially the TREATMENT SURFACE (See 8.1: "Cleaning Silk'n® FaceTite™"). Clean the device after each treatment.
3. Return your Silk'n® FaceTite™ device to its original box.

8. Maintenance of Silk'n FaceTite™

8.1. Cleaning Silk'n FaceTite™

- Unplug Silk'n® FaceTite™ before cleaning. After each treatment session, it is necessary to clean your Silk'n® FaceTite™ device, and especially the TREATMENT SURFACE.
- Use absorbent paper (such as toilet paper or hand drying napkins) to gently wipe the gel residue off the TREATMENT SURFACE.
- Never immerse Silk'n FaceTite™ or any of its parts in water!

9. Troubleshooting

“My Silk'n FaceTite™ does not start.”

- Make sure the POWER ADAPTOR is properly connected to the FaceTite™ device.
- Make sure the POWER ADAPTOR is plugged into an electrical outlet on the wall.

“Heat energy is not emitted and I don't feel any heat emitted from the device even though I pressed the CENTRAL BUTTON”

- Heat will not be emitted if the TREATMENT SURFACE is not attached to your skin. This is a safety feature. Make sure the device is tightly pressed to your skin before starting treatment. Do not attempt to open or repair your Silk'n FaceTite™ device. Only authorized Silk'n™ personnel are permitted to perform repairs. Trying to open Silk'n FaceTite™ may also damage the device and will void your warranty. Please contact Silk'n™ Customer Service if you have a broken or damaged device in need or repair.
- When Silk'n FaceTite™ reaches its End of Life stage, the device is to be orderly disposed of. Deliver it to the nearest electronic waste collection point for safe and efficient disposal.

10. Customer Service

For more information about Silk'n™ products please enter your local Silk'n™ website, www.silkkn.eu or www.silkkn.com.

If your device is broken, damaged, in need of repair, or for any other user assistance, please contact your local Silk'n™ Customer Service:

For USA and Canada customers: 1-877-DO-SILKN / 1-877-367-4556, contact@silkkn.com

For customers in other countries: info@silkkn.com

11. Specifications

Technology		HT™
Wavelength Red light		620nm±20nm;
Wavelength IR		850nm±20nm;
Radio Frequency		1 MHz ±30%; 12W max
Operation and Safety		Heat stabilizer preventing skin overheating
Package Size		16.6cm X 45.5cm X 30.2cm
System weight		125gr
Temperature	Operating	10°C to 35°C
	Storage	-40 to +70°C
Relative Humidity	Operating	30 to 75%rH
	Storage	10 to 90%rH
Atmospheric pressure	Operating	700 to 1060hPa
	Storage	500 to 1060hPa
Adaptor		Model: KSA24A 200150HE
Rated		100-240V; 50-60Hz; 0.5A
Output		12Vdc; 1.5A

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Home Skinovations Ltd. Tavor Building, Shaar Yokneam POB 533, Yokneam 20692, ISRAEL www.silkn.com, info@silkn.com

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Tous droits réservés. Copyright 2011 © Home Skinovations Ltd.
Copyright 2011 © Home Skinovations Ltd. Alle rechten voorbehouden.
Copyright © 2011 Home Skinovations Ltd. Με επιφύλαξη παντός δικαιώματος.
Copyright 2011 © Home Skinovations Ltd. Alle okeudet pidatetän.
Copyright 2011 © Home Skinovations Ltd. Alle rettigheder forbeholdes.