
waterpik®

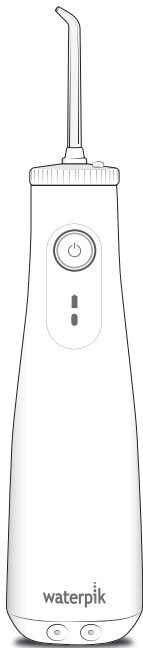Waterpik® Water Flosser Model WF-10 Series

Munddusche Waterpik® Modell WF-10
Hydropulseur Waterpik® Modèle WF-10
Idropulsore Waterpik® Modello WF-10
Irrigador bucal Waterpik® Modelo WF-10
Waterpik® Water Flosser Model WF-10
Irrigador oral Waterpik® Modelo WF-10
Ирригатор Waterpik® Модель WF-10

www.waterpik.com/intl

samcooper**design**

DO NOT PRINT

DATE: 04.07.2023
CLIENT: WATERPIK

FILENAME: 20032226-FAA_WF-10Series_EU_IM_R4.indd
DESCRIPTION: WF-10 EU Instruction Manual
LANGUAGES: 8 (Eng, Ger, Fr, Ital, Span, Dut, Port, Russ)
FILE TYPE: Adobe InDesign CC
SCD JOB# / PREPARED BY: 22WPOH7566/ JS
CLIENT ROUND: 4

SCALE: 100%
UPC: N/A
SIZE: 14.0 x 4.5 flat; 7.0 x 4.5 folded
FORM NUMBER: 20032226-F AA
INK: Process Black
NOTES:

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER:

To reduce the risk of electrocution:

- Do not handle the user supplied USB charger with wet hands.
- Do not place or drop the user supplied USB charger into water or other liquid.
- Do not use this product while bathing.
- Do not place or store this product where it can fall or be pulled into a tub or sink.
- Do not reach for this product if it has fallen into water. Unplug from electrical outlet immediately.
- Check the charging cord and user supplied USB charger for damage before the first use and during the life of the product.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- Use this product only as indicated in these instructions or as recommended by your dental professional.
- A user supplied USB charger with the following specifications is required to charge this product: USB Type-A, Rated Output: 5.0 VDC, 1.0 A.
- Do not plug this device into a voltage system that is different from the voltage system specified on the device or charger.
- Keep the charging cord and user supplied USB charger away from heated surfaces.
- Do not use the Water Flosser for more than five minutes in each two hour period.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into the nose or ear.
- If your physician advised you to receive antibiotic premedication before dental procedures, you should consult your dental professional before using this product or any other oral hygiene aid.

- Fill reservoir with warm water or other dental professional recommended solutions only.
- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce performance and will shorten the life of the product.
- Only use tips and accessories that are recommended by Water Pik, Inc.
- If the Pik Pocket™ tip (not included in all models) separates from the shaft for any reason, discard tip and shaft and replace with a new Pik Pocket™ tip.
- Do not drop or insert any foreign object into any opening on the product.
- Do not use this product outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Remove any oral jewelry prior to use of this product.
- Magnetic Field: Pacemakers and similar implanted medical devices may be affected by magnet.
- Do not use if you have an open wound on your tongue or in your mouth.
- There are no consumer serviceable parts in this product and it does not require maintenance.

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- If the supply cord is damaged, it must be replaced by a special cord or assembly available from the manufacturer or its service agent.
- Do not use this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into any liquid. Contact Water Pik, Inc. at 1-800-525-2774 or visit www.waterpik.com/intl.

SAVE THESE INSTRUCTIONS

Product Description	5
Getting Started	6
Cleaning and Troubleshooting Your Waterpik® Water Flosser	9
Battery Removal	10
Warranty Information	11

SYMBOLS


The following symbols may appear on your product and/or packaging. Some of these symbols may not be relevant in your region and are listed for information purposes only.

Symbol	Definition	Symbol	Definition	Symbol	Definition	Symbol	Definition
	Refer to instruction manual		Detachable Power Supply		Do not use this device in a shower, bathtub and do not immerse in any liquid.		Regulatory Compliance Mark
	Manufacturer		Direct Current (DC) Voltage		Authorized Representative in the European Community		Separate collection for waste electrical and electronic equipment.
IPX7 or IPX4	Ingress Protection Code		Alternating Current (AC) Voltage		CE marking of conformity		
	IEC Class II Equipment		Do not immerse this device in any liquid.		UK Conformity Assessed		

1. Tip knob
2. Power button
3. Mode/pressure button
4. Battery charge indicator
5. Charging contacts
6. Tip
7. Tip eject button
8. Reservoir fill door
9. Reservoir

USB charging cable
20026734-1 (White)
20026734-2 (Black)

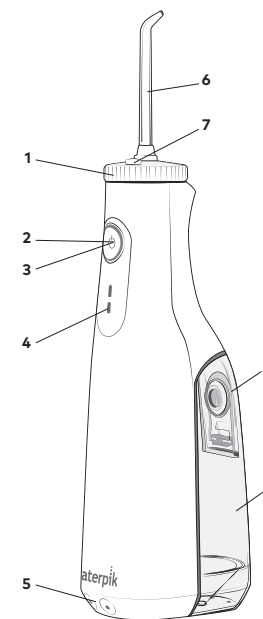TIPS

A. Orthodontic Tip*

- Braces/General use

B. Pik Pocket™ Tip*

- Periodontal pockets/Furcations

C. Plaque Seeker™ Tip*

- Implants/Crowns/Bridges/
Retainers/General use

D. Toothbrush Tip*

- General use


E. Classic Jet Tip**

- General use

F. Tongue Cleaner**


- Fresher breath

NOT ALL TIPS INCLUDED IN ALL MODELS.

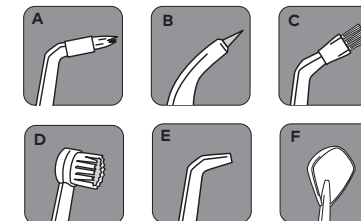*Replace Every 3 months **Replace Every 6 months

Replacement tips can be purchased from your local distributor, which can be found at www.waterpik.com/intl.

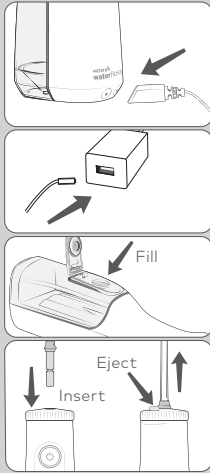Charging Unit Prior to First Use

Plug the USB charger into the wall outlet and connect the magnetic plug to the unit for approximately 4 hours (see image). A user supplied USB charger with the following specifications is required to charge this product: USB Type-A, Rated Output: 5.0 VDC, 1.0 A.

Ensure contacts on BOTH the cable and unit are clean and dry before connecting the charger. The LED light will display a single, rapid blink every 3 seconds when unit is fully charged and the charger is connected. A full charge will last for approximately 1 week of typical use. When the batteries are in need of charging (approximately 30% charge left) the charge indicator light will blink 8 times quickly after use to indicate that you should recharge the unit. A complete recharge will take 3-4 hours. If you use your Water Flosser once a day or less, it is recommended you charge it once a week. If you use it twice a day or more, charge the unit more frequently.

Filling the Reservoir

Remove the magnetic charger from the unit prior to filling reservoir and product use. Lift flip top on reservoir and fill the reservoir with lukewarm water.

Inserting and Removing Tips

Insert tip into the center of the knob at the top of the Water Flosser handle. Press firmly until the tip clicks into place. To remove the tip from the handle, press the tip eject button and pull tip from the handle.

CAUTION: Do not press tip eject button while unit is running.

Adjusting the Pressure Setting

- First press (from OFF) starts unit on "LOW".
- Second press turns unit to "HIGH".
- Third press turns unit "OFF".

The single LED exhibits the following behavior:

- Solid - If motor is running and battery is greater than half full.
- Single blink every 3 seconds - When charger is connected and battery is fully charged.

Recommended Technique

Lean low over sink and place tip in mouth. Aim the tip toward gumline, and turn the unit ON (I).

For best results, start with the back teeth and work toward the front teeth. Glide tip along gumline and pause briefly between teeth. Continue until you have cleaned the inside and outside of both the upper and lower teeth.


Direct the jet stream at a 90-degree angle to your gumline. **Slightly close lips to avoid splashing** but allow water to flow freely from mouth into the sink. Keep unit upright during use for best results.

When Finished

Turn the unit OFF. Empty any liquid left in the reservoir.

If desired, the reservoir may be removed from the power handle by simply pulling down on the reservoir to remove.

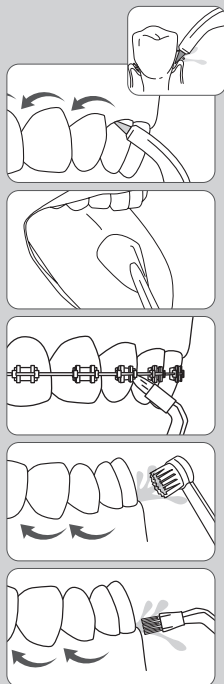NOTE: NOT ALL TIPS INCLUDED WITH ALL MODELS.

TIP USAGE

Pik Pocket™ Tip

The Pik Pocket™ tip is specifically designed to deliver water or anti-bacterial solutions deep into periodontal pockets. To use the Pik Pocket™ tip, **set the unit to the lowest pressure setting**. Place the soft tip against a tooth at a 45-degree angle and gently place the tip under the gumline, into the pocket.

Tongue Cleaner Tip

Set to lowest pressure setting, place tip in the center/middle of your tongue about half way back. Pull forward with light pressure. Increase pressure as you prefer.

Orthodontic Tip

To use the orthodontic tip, gently glide tip along gumline, pausing briefly to lightly brush area between teeth and all around orthodontic bracket, before proceeding to the next tooth.

Toothbrush Tip

Place the toothbrush tip in mouth with the brush head on the tooth at the gumline. The toothbrush tip can be used with or without toothpaste. Using a **light pressure** (bristles should not bend) massage the brush back and forth with very short strokes – much like you would with a manual toothbrush.

Plaque Seeker™ Tip

To use the Plaque Seeker™ tip, place the tip close to the teeth so that the bristles are gently touching the teeth. Gently glide tip along gumline, pausing briefly between teeth to gently brush and allow the water to flow between the teeth.

Using Mouthwash and Other Solutions

Your Waterpik® Water Flosser can be used to deliver mouthwash and antibacterial solutions. For ease of adding other solutions remove the reservoir (see "When Finished" section) and add from top opening on reservoir. If you are adding additional water you may replace the reservoir and fill reservoir fully with the fill door. After using any special solution, rinse unit to prevent clogging by partially filling reservoir with warm water and running unit with tip pointed into sink until unit is empty.

Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior.

It is recommended that the reservoir is removed and cleaned weekly.




Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance. Cleaning internal parts: add 1 tablespoon of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

This process should be done every 1 to 3 months to ensure optimal performance.

Battery Care and Disposal

If the Water Flosser is to be stored for an extended period of time (six months or more), be sure to fully charge it prior to storage. The unit contains a rechargeable Nickel Metal Hydride battery. At the end of the battery life, please recycle the unit at your nearest recycle center per local or state requirements. The battery is not replaceable.

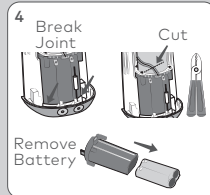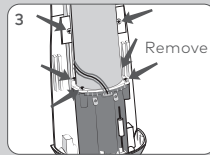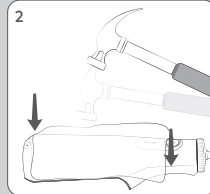Battery Removal

The product contains battery and/or recyclable electronic waste. In order to protect the environment, please do not dispose along with household waste. When battery life ends, please dispose or recycle according to local administrative regulations. Please contact local waste administrative department for any question.

If you intend to dispose the battery on your own, for example, for WF-10 model, please follow below steps to remove the battery:

1. Power off the device, remove reservoir.
2. Wrap in a towel and hit with a hammer to break the joint between the front, rear, and lower housings.
3. Remove 2 screws holding pump assembly to rear cover and 4 screws holding battery housing to pump assembly.
4. Carefully break the joint between the battery housing and rear cover.
Cut battery wires (careful not to short positive and negative wires together) and remove battery.

The removed battery should be disposed safely by special recycling center.

Warning: The battery of this product is not replaceable. Please do not disassemble the product unless you want to dispose it. Otherwise the warranty service will be void.

Service Maintenance

Waterpik® Water Flossers have no consumer - serviceable electrical items and do not require routine service maintenance. For all your service needs go to www.waterpik.com/intl to find your local distributor. Some replacement parts and accessories for the Water Flosser are also available from your local distributor which can be found at www.waterpik.com/intl.

Refer to the serial and model numbers in all correspondence. These numbers are located on the back of the handle and bottom of the reservoir respectively.

Still have questions? Go to www.waterpik.com/intl to find your local distributor.

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 2 years from the date of purchase (Spain Only: 3 years). Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion any part of the product, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. Repair or replacement services provided for by the guarantee will be carried out by the local distributor, whose details may be found at www.waterpik.com/intl. Upon request of the consumer, a copy of the guarantee complete with the local distributor's information will be provided by email or mail. This limited warranty excludes accessories or consumable parts such as tips, etc.

This warranty gives you specific legal rights. You may also have other rights that vary depending on the law in your place of residence.

To locate your local distributor, go to www.waterpik.com/intl.

WATER PIK, INC. DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

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©2023 Water Pik, Inc.
**Water Pik, Inc., a subsidiary
of Church & Dwight Co., Inc.**
SOFIBEL SAS