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Intended for the treatment of mild to moderate inflammatory acne.

YOUR SHARK CRYOGLOW JOURNEY

Choose the LED routine that best suits your skin goals and use it once daily.
Use Under-Eye Routine anytime to quickly depuff under eyes.

Powered by InstaChill Cold Technology. This cooling sensation helps reduce under-eye puffiness. Cold level can be adjusted or turned off using the remote.

LED Routines

with InstaChill Cold Technology

Select either the **Red + Infrared Light Therapy** or **Blended Blue Light Therapy** routine below and use your selection exclusively once daily for a minimum of **8 weeks**. For both the treatment of acne and to reduce redness and improve roughness, we recommend using the Blended Blue Light Therapy mode exclusively for 8 weeks. Do not stack. After 8 weeks, use Skin Sustain to maintain results.

FIRST 8 WEEKS

Red + Infrared Light Therapy (6 mins.)
Helps reduce redness and improve roughness.

OR

Blended Blue Light Therapy (8 mins.)
Helps reduce acne for a clearer, smoother complexion.

AFTER 8+ WEEKS

Skin Sustain (4 mins.)
Helps to boost brightening and maintains complexion with daily use.

Under-Eye Routine

with InstaChill Cold Technology

Use InstaChill without LEDs. Select length of time and use as often as you like to cool down under eyes.

Under-Eye Revive (5–15 mins.)

Firm, tighten, and soothe under eyes.

InstaChill will automatically start with every routine. You can adjust the level of coldness or turn it off with the remote.

Scan for how to videos and more, or visit us at sharkclean.eu

Please read the **Warnings, Instructions, and Contraindications**—and charge the mask fully*—before first use.

HOW TO USE YOUR SHARK CRYOGLOW MASK

1. Start with a clean, dry face. Put on mask and adjust straps for comfort.
2. Turn on the mask by pressing and holding the remote's center dial.
3. Select a routine and press "Start." Adjust InstaChill with the remote's chill button.
4. After the session, remove the mask. It will power down automatically.
5. Clean the mask with a cloth slightly dampened with water.

GET THE PERFECT FIT


Everyone's face is unique, so use these adjustments to make the mask more comfortable.

Top and back straps

Adjust the fit to your eye level so it's most comfortable.

Self-adjusting top strap

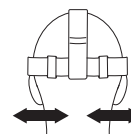Slide back strap

Chill pads (inside the mask)

with InstaChill Cold Technology

Correct positioning

Optional clip-on chill pads
(available in the box)

Position the mask on your face, then stretch the back strap over your head so it sits on the crown. The top strap can be adjusted so the mask sits at your desired eye level.

The mask's built-in chill pads should be positioned under the eyes. If the built-in chill pads do not touch your face, place the optional clip-on chill pads on the mask to achieve your best fit.

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SHARK and SHARK CRYOGLOW are registered trademarks in the European Union
© 2025 SharkShark Operating LLC. All rights reserved.
of SharkShark Operating LLC.
SharkShark Operating LLC 89 A St #100, Needham, MA 02494, USA
Emergo Europe B.V. Westervoortsedijk 60, 6827 AT Arnhem, The Netherlands