clear

Instructions for use
Mode d'emploi

Thank You...

For purchasing the **clear**.

No need to head to the dermatologist when you now have the same blemish fighting technology at home, in the palm of your hand. The **clear** delivers a unique combination of blue light, sonic vibration, and gentle warming that has been clinically proven to clear or fade a blemish within 24 hours without dryness or irritation. Rechargeable and portable, **clear** is safe and easy to use anywhere, anytime to clear blemishes fast. For important safety and how to use information, read this User Manual.

Contact the experts at mē to help you get started and set up your treatment routine. Our goal is to help you get the best results possible.

Customer care contact:

Toll free (NA): 1.877.651.1177

Email: info@mepower.com

Hours of operation: Mon-Fri, 9am-5pm EST

Manufacturer:

Tanda Beauty Canada Inc.

11-380 Jamieson Parkway, Cambridge ON, N3C 4N4

Contents

About the clear device	2	Battery life indicator	14
Device features	3	Light sensitivity test	14
Acne and blue light technology	4	Treatment with clear device	15
Contraindications	5	Device maintenance	16
Warnings	6	Troubleshooting	17
Precautions	8	Specifications	18
Risks	9	Device servicing	20
Benefits	11	Customer care	20
How to know that the clear device is right for you	12	Warranty	20
Charging the battery	13	Satisfaction Guarantee	22

EN About the clear device

What does it do?

The **clear** anti-blemish device is for the treatment of individual or small clusters of acne blemishes. It is rechargeable so it is always ready to treat. The device is indicated for the treatment of mild to moderate inflammatory acne.

Package Contents


- **clear** device
- Protective Cap
- USB Cable
- Quick Reference Guide
- Instructions for Use

Device components

- Protective cap
- Treatment head surface
- Power button
- Battery life indication lights

- **clear** device cordless rechargeable unit

- Plastic battery tab
- USB cable

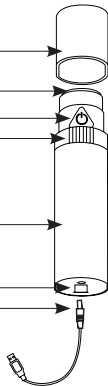
Device features

- **Single on/off power button** will activate **clear** device. Battery life indicator lights will turn white and a short audible beep will let you know that the device is ready for use.
- **Skin contact sensor** on the treatment surface will activate the device only when it is touching your skin.
- **Continuous mode** allows you to treat multiple blemishes at one time. Simply move the device to the next blemish every time you hear a short audible beep (2 minutes) without having to press the gray on/off button.
- **One short audible beep** will indicate that 2 minutes have passed and you have completed 1 treatment. You can move and treat another blemish or simply put your device down to terminate the treatment. For best results, treat your blemish 3 times a day with the **clear** device.
- **Built in temperature sensor** monitors the temperature of the treatment surface. The device can be used for 10 minutes continuously before two short audible beeps sound and the device shuts off automatically.
- **Smart device shut off** by simply removing the device from your treatment area. The LED indicator lights will glow for 20 seconds before complete shut off occurs.

EN Acne and Blue Light Technology

Blue light technology offers an effective solution for the treatment and prevention of mild to moderate inflammatory acne. Once only available at the dermatologist's office, this advanced technology is proven to work better than Benzoyl Peroxide* without any of the negative side effects like dryness, irritation or flakiness. The specific wavelength of professional strength blue light delivered by the **clear** device has a phototoxic effect on the acne-causing bacteria, causing it to self destruct, breaking the acne cycle simply and cleanly. New sonic vibration and gentle warming features improve results by opening pores, exposing even more of the *P. acnes* bacteria to the powerful blue light.

* Papageorgiou P., Katsambas, A. Chu, A. (2000). Phototherapy with blue (415nm) and red (660 nm) light in the treatment of acne vulgaris. British Journal of Dermatology.

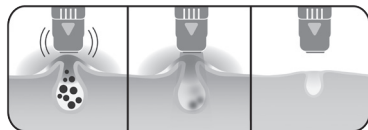Sonic vibration
and gentle warming
opens pores

Blue light destroys
acne-causing
bacteria

Blemish disappears in
just 24 hours

Contraindications

Never use this device over any suspicious or cancerous lesion. Failure to do so can result in delaying the best medical treatment. Consult your healthcare provider before using this device on any skin growth or lesion that appears abnormal.

Do not use this device on severe acne or on non-inflammatory acne. Failure to do so can result in delaying the best medical treatment.

Do not allow children to use the device unsupervised.

Do not use if you are pregnant or breastfeeding. **clear** has not been tested on pregnant or breast-feeding women. Hormonal changes occur during pregnancy or breast-feeding could increase the risk of skin sensitivity.

- Do not use if you are exposed to other non-drug “sensitizers” known to cause photosensitive reactions. Doing so can result in skin discoloration, discomfort or pain.

Common examples of non-drugs that may have made you sensitive to light are: deodorants, antibacterial soaps, artificial sweeteners, fluorescent brightening agents for cellulose, nylon and wool fibers such as in laundry detergent, naphthalene (mothballs), petroleum products, and in cadmium sulfide, a chemical injected into the skin during tattooing.

Do not use over or near a steroid injection received within the last six months. May result in photosensitive reaction which can result in skin discoloration, discomfort or pain.

Do not use if you have a history of light triggered seizures.

Do not use if you are sensitive to light. Doing so can result in skin discoloration, discomfort or pain. It is recommended that you conduct the light sensitivity test (page 14).

Do not use if you have experienced chronic migraines in the last two years.

Do not use over the thyroid gland (neck area). This may result in stimulation of the gland.

Warnings

If you experience any adverse effects, or have a concern regarding use, stop immediately. Contact your health care provider. You may be having a photosensitive reaction.

To reduce the risk of electric shock that may cause injury or death:

Always attach the AC/DC power adapter to device first, then into the outlet.

Always use a USB certified port when charging **clear**.

Do not place or store the device where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

Keep **clear** away from water (including bath and shower). If the unit accidentally falls in water unplug immediately.

Keep the cord away from heated surfaces.

Never drop or insert any object into any openings on the device.

Never operate this device if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or submerged in water. Return the device to Tanda Beauty Canada Inc. for replacement.

Do not use the device when the USB cable is plugged into the device when its battery is charging.

Do not disassemble or modify the device as this may cause unit damage, malfunction, electrical shock, fire or personal injury. There are no user-serviceable parts inside the unit.

To prevent fire or serious burns:

Turn **clear** off immediately if it begins to overheat.

Keep away from where flammable aerosol products are being used or where oxygen is being administered.

EN Precautions

Discontinue use if you experience any discomfort. You may have sensitivity to light or other condition which can result in skin discoloration, discomfort or pain. Consult with your health care provider.

Use this device only for its intended use as described in this manual. Failure to do so may result in injury to yourself or others.

Dispose of the device in accordance with local and national regulations and codes. Failure to do so may result in violation of the law.

Do not look directly at the light source of the treatment head of this device while it is in operation. Failure to do so may result in headache.

Do not operate this device in an environment where other devices, which intentionally radiate electromagnetic energy in an unshielded manner, are being used. Portable or mobile Radio Frequency communication equipment can affect the **clear** device. Doing so may result in the device not functioning properly.

Do not scratch, mark or otherwise damage the light emitting surface. To do so may diminish the effectiveness of the light.

Do not shine the light source into someone else's eyes. Doing so may result in giving the person a headache.

Failure to use and maintain the device in accordance with the instructions outlined in this manual will invalidate the warranty.

Read all instructions before using your **clear** device. Failure to do so may result in failure to understand safety information or how the device is used.

Prolonged or excessive use beyond the recommended treatment time may cause the treatment surface of the device to become warm. Failure to follow the recommended treatment time may cause you to experience adverse effects such as discomfort, or prolonged redness in the area or an ineffective treatment due to the automatic shut down of the device

Inspect this device prior to each use. Check for possible damage. Make sure that the outer shell of the device isn't cracked. Failure to do so may result in improper delivery of the light.

Operate, transport and store clear as specified in this manual. Failure to do so may result in damage to the unit.

Adjustments or performance of procedures other than those specified herein may result in injury.

Risks

Hazard	Likelihood of hazardous event	Effects to User
Possible electrical problems	Not likely to happen	Shock
Overheating	May occur from time to time	Discomfort, interrupted treatment
Inadequate power to the unit	Not likely to happen	Unit will not function
Operation outside prescribed environmental conditions	Not likely to happen	Damage to unit
User fails to follow maintenance instructions	Not likely to happen	Lack of product effectiveness

Hazard	Likelihood of hazardous event	Effects to User
Battery overheating	Not likely to happen	Handle can overheat resulting in discomfort
Battery explodes	Not likely to happen	Can cause serious burns or injuries
Battery leaking	Not likely to happen	Can cause serious burns
Short circuit wiring	Not likely to happen	Can cause overheating of handle, which may result in user discomfort
User may inadvertently shine light in eyes	May occur from time to time	Discomfort and possibly headache
User may attempt to shine light in eyes	Not likely to happen	Discomfort and possibly headache
User may have a photosensitive skin condition	Not likely to happen	May cause skin to temporarily turn red/darken.

Hazard	Likelihood of hazardous event	Effects to User
User may be taking a photosensitizing medication or herb	Not likely to happen	May cause skin to temporarily turn red/darken.
User attempts to use product past recommended maximum exposure per site per day	May occur from time to time	Ineffective treatment
Person with thyroid disease may attempt to use around neck area	Not likely to happen	Possible stimulation of disease growth
Person may attempt to use over recent steroid injection site	Not likely to happen	Unanticipated reactions

Hazard	Likelihood of hazardous event	Effects to User
Incorrect identification of "mild to moderate inflammatory acne"	Likely to happen	User may attempt to use on more severe forms of acne or on other forms that normally would require other forms of treatment or even physician care. The attempt would delay getting best treatment

Only the known risks of the device are listed. Users of the device can control the risks of the device only if they read and understand all parts of the user manual and follow all the instructions, contraindications, warnings and precautions in it exactly.

Benefits

EN

Benefit	Likelihood of Benefit	Effects to User
Reduction in mild to moderate inflammatory acne	Frequent	Improved skin condition because of the reduction of mild to moderate inflammatory acne
Users of this device can obtain benefit from using it only if they read and understand all information in the user manual and follow exactly all the instructions in the user manual.		

EN How to know that the clear device is right for you

There are different forms of acne. Some forms are best treated by a doctor. Others may be suited for **clear**. You can use **clear** for red pimples and pimples with white or yellow matter in them known as pustules.

The following chart provides descriptions of the different types of acne and indicates whether the **clear** device can be used to treat the condition.

Acne type	Description of condition	Is clear device right for you?
Mild Acne	Several inflamed pimples - red in color Less than 20 whiteheads, or less than 15 inflammatory red pimples or less than 30 total pimples, not all inflamed	Yes
Moderate Acne	Any inflamed pimples (red in color) and pustules (visible accumulation of pus in skin) 20 to 100 whiteheads, or 15 to 50 inflammatory pimples, or 30 to 125 total pimples, not all inflamed	Yes

Acne type	Description of condition	Is clear device right for you?
Severe Nodular Acne	Inflamed pimples and pustules (visible accumulation of pus in skin) with a few deep nodular lesions (solid mass can be felt under skin – can sometimes be raised) Greater than 5 cysts (solid mass of skin like a knot, can be raised or felt under the skin), or total whitehead count greater than 100, or total inflammatory count greater than 50, or greater than 125 total lesions	No, consult your dermatologist for the best prescribed treatment
Severe Cystic Acne	Many nodular cystic lesions (with signs of scarring)	No, consult your dermatologist for the best prescribed treatment
Comedonal Acne	Several whiteheads and blackheads	No

Charging the battery

EN

1. Plug the DC input directly into the bottom of the device and plug the other end into any certified USB charging source. Keep plugged in and charge for 4 hours.
2. The indicator lights will flash one by one from left to right while the device is charging (see illustration); they will turn off when charging is complete.
3. The device is now ready for use. Unplug and begin treating.

EN Battery life indicator

Battery capacity:

Indicator lights illuminate to show 100%-25% battery capacity when device is turned on.

Battery charging:

Each indicator light flashes one by one from left to right when charging.

Light sensitivity test

Follow these instructions to see if you are sensitive to light:


1. Press the gray power button.

2. Place the device near the area that you intend to treat that is free from a blemish. Hold the device in place for one 2-minute treatment cycle or until you hear a beep.
3. Wait 24 hours. You should NOT experience any adverse reaction during the next 24 hours.
4. An adverse reaction is pain or discomfort in the treatment area or any prolonged skin redness.

Device Set-up

If this is your first time using the device, be sure to conduct the Light Sensitivity Test described below.

1. Unpack the device and all the accessories.
2. Remove the plastic battery tab at the bottom of the unit near the DC input receptacle prior to charging the unit.

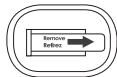Treatment with Clear device

Step 1: Prepare skin for treatment

Before use, cleanse your skin to remove makeup or creams that may contain reflective materials.

Step 2: Activate device

Press the gray button.

A short audible beep and illuminated battery indicator lights will signal the activation of **clear**.

Step 3: Treat your skin

Hold the device steady over your

blemish. Ensure you maintain full contact between the treatment surface area and your skin, while sonic vibration and blue light are delivered. Treat until you hear a single audible beep (2 minutes).

To treat another blemish, move the device to next area for treatment or simply put your device down to terminate the treatment.

For best results, treat your blemish 3 times a day with the **clear** device.

Step 4: After treatment

Apply your favorite non-comedogenic skincare products.

EN Device maintenance

The following maintenance instructions are important to ensure that your device continues to work as it was designed. Failure to follow these instructions may cause your device to stop delivering the required dose or to stop functioning altogether.

Environment for storage and use

It is recommended to store and use the device in a controlled environment with constant ambient temperature of 5°C to 35°C; relative humidity of 15% to 93% non-condensing; atmospheric pressure of 700 – 1060 hPa.

Battery care

It is recommended to keep your battery capacity at 50% or higher (2 Indicator lights illuminated). Refer to the Battery Life Indicator Section for details on how to read the battery capacity on your device.

Cleaning

The treatment head and outer surfaces of the device can be cleaned using a damp, soft, lint-free cloth and a mild alcohol-free detergent following each treatment or when required. Dry the surface using a clean dry cloth. Your device is not waterproof. Do not submerge the device in water. Do not put in dishwasher, as doing so will result in permanent damage. Should the unit accidentally become submerged, do not attempt to use it. Contact the retailer or Tanda Beauty Canada Inc.

Please visit mepower.com

Troubleshooting

In the event that the **clear** device fails to perform as intended, the following notes will help to identify potential problems with the device and its setup.

If	Do this
You see all 4 indicators lights are flashing.	Ensure you are in the ideal environment for treatment (ie temperature, humidity). If not, move to an appropriate environment (ie indoors). OR Wait 5-10 minutes and attempt to re-activate the device after moving to an appropriate environment (ie. indoors) OR Charge your device as the battery capacity may be too low.

If	Do this
During treatment, you hear 2 short audible beeps.	This means the device has delivered five 2 minute cycles of treatment and must be re-activated by pressing the on/off button to continue use.
During treatment, you hear 3 short beeps and see all 4 indicators lights are flashing	Ensure you are in the ideal environment for treatment (ie temperature, humidity). If not, move to an appropriate environment (ie indoors). OR Wait 5-10 minutes and attempt to re-activate the device after moving to an appropriate environment (ie. indoors)
Indicator lights are on, but the treatment head is not illuminated.	Ensure the treatment surface is in full contact with the skin. Otherwise, the smart shut off will be activated and the device will turn off in 20 seconds.

EN

EN

If	Do this
During treatment, the light continues to illuminate even without skin contact.	Turn off the device and re-activate by pressing the on/off button once.

If you still have questions please visit mepower.com for more tips and troubleshooting or contact customer care directly.

Specifications

clear device = Series 14

TYPE BF APPLIED PART

Attention, consult
ACCOMPANYING DOCUMENTS

TUV Certification Mark

Follow instructions for use

FC FCC - DECLARATION OF CONFORMITY:

Model: **clear** device Series 14 conforms to the following specification: FCC Part 15, Subpart B: 2006.

Supplemental Information: This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

INDUSTRY CANADA NOTICE:

This Class B digital apparatus complies with Canadian ICES-003.

DISPOSAL

Waste Electrical and Electronic Equipment (WEEE)

For information and instructions regarding the proper disposal of this product in your location please visit mepower.com

EN Device servicing

The device is not serviceable. If **clear** needs to be replaced, please contact the selling retailer, local distributor or Tanda Beauty Canada Inc. for more information.

THIS HOTLINE CAN BE USED IF YOU HAVE ANY QUESTIONS REGARDING USE OF THE DEVICE

Email: info@mepower.com

Hours of operation:

Mon-Fri, 9am-5pm Eastern Standard Time

Warranty

Register your **clear** device warranty at mepower.com

Your mē product is warranted to be free from defect in material and workmanship for a period of one year after its original purchase date under normal use. This warranty extends only to the original retail purchaser with original proof of purchase and only when purchased from an authorized mē retailer or reseller.

If the product should become defective within the warranty period, contact customer care team toll free at 1.877.651.1177 for repair or replacement. If you are outside of North America, contact your local distributor. Contact information can be found at mepower.com. Tanda Beauty Canada Inc. reserves the right to replace a defective product with the most comparable product currently available. Your exclusive remedy for any defective product covered by this limited warranty is limited to the repair or replacement of the defective

product. Tanda Beauty Canada Inc.'s entire liability for any defective product shall in no event exceed the purchase price for the defective product.

This warranty does not cover products determined by Tanda Beauty Canada Inc. to be damaged by any of the following:

- Accident, misuse, abuse or alteration
- Servicing by unauthorized persons
- Use with unauthorized accessories or use other than as instructed
- Connecting to incorrect current and voltage
- Any other conditions beyond our control

TO THE EXTENT PERMITTED BY LAW, TANDA BEAUTY CANADA INC. SHALL NOT BE RESPONSIBLE FOR ANY DIRECT OR INDIRECT DAMAGES RESULTING FROM THE USE OR MISUSE OF THIS PRODUCT INCLUDING ANY GENERAL, INCIDENTAL, SPECIAL, OR CONSEQUENTIAL

DAMAGES. ALL IMPLIED REPRESENTATIONS, WARRANTIES OR CONDITIONS INCLUDING BUT NOT LIMITED TO ANY IMPLIED REPRESENTATIONS, WARRANTIES OR CONDITIONS OF FITNESS AND MERCHANTABILITY ARE LIMITED IN DURATION TO ONE YEAR FROM DATE OF ORIGINAL PURCHASE.

This warranty gives you specific legal rights, and you may also have other rights that vary depending on your state, province or other jurisdiction. Some jurisdictions do not allow the exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so the above limitations and exclusions may not apply to you.

EN Satisfaction Guarantee

We want you to be 100% happy with your **clear** device. If you are not completely satisfied please return the product within 30 days for a full refund. Please note that this guarantee applies to all mē products purchased at mepower.com or directly from Customer care. All returns are subject to the return policy at point of purchase.

Copyright© Tanda Beauty Canada Inc. 2015, mē is a trademark of Tanda Beauty Canada Inc. All rights reserved.

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Copyright © Tanda Beauty Canada Inc.. 2015, mē et **clear** device sont des marques de commerce de Tanda Beauty Canada Inc. Tous droits réservés.
Copyright © Tanda Beauty Canada Inc. 2015,
All rights reserved. PUBLICATION: HU-PB00243 (November, 2015)