Therabody

Theragun Mini

3rd Generation

User Manual & Unit Warnings

Table of Contents

EN	1.	Product Overview	EN	3-15
	2.	What's in the Box	DE	16-23
	3.	Download the Therabody App	FR	24-31
	4.	Getting to Know the Device	IT	32-39
	5.	Using the Device	ES	40-48
	6.	Attachments	RU	49-56
	7.	Care, Cleaning, and Charging	FRCA	57-65
	8.	Troubleshooting Guide	ESMX	66-73
			SC	74-81
			TC	82-89
			KR	90-99

Theragun Mini

EN

1. Product Overview

The redesigned Theragun Mini is smaller, lighter and more travel-friendly than ever. Scientifically proven to relieve aches and pains, ease tension and reduce stress, Theragun Mini makes it easy to get the relief you need anytime, anywhere. Comfortably treat different areas of the body with a new, more secure grip made for all hand sizes and 3 attachments scientifically designed to target specific pain points. And unlike other massage guns, Theragun Mini is made with a quiet, anti-stall motor, so you can massage even the toughest knots without stopping your treatment. Perfect for long flights and daily commutes, lunch breaks or a round of golf – Theragun Mini tucks right into your travel bag for quick, effective relief wherever life takes you.

Intended Use

The Theragun Mini's percussive massage therapy is intended to relieve minor daily aches and pains, reduce muscle tension, tightness, and stress, and ease painful knots.

2. What's in the box

- A** Theragun Mini
- B** USB-C charging cable
- C** Three attachments:
 - I** Dampener
 - II** Standard Ball
 - III** Thumb
- D** Travel pouch
- E** Instruction manual

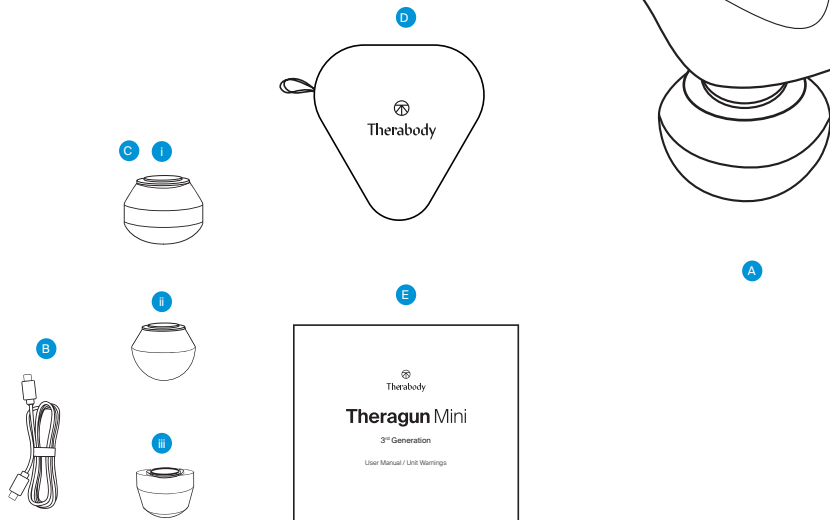
Figure 1. What's in the box

3. Download the Therabody App

If you currently have the Therabody app, proceed to step 4.

EN

1. Download the Therabody app from the App Store or Play Store.
2. Follow the in-app instructions to create a new app account and verify your email.
3. Once logged in, follow the onboarding steps to connect your device to the app. Be sure that your phone's Bluetooth settings are enabled.
4. Select "Theragun Mini" as your Therabody device when prompted.
5. Press and hold the power button on your device to turn it on. The device must be powered on to locate and pair to the app.
6. Select "Theragun Mini" to confirm the device and app connection.
7. Once connected, you can customize your treatments based on your preferences. Explore the library of guided pain and stress relief routines that show you exactly which muscle to treat and which attachment to use.
8. You can also select from an expansive catalog of TheraMind sound therapy, science-backed audio to listen to during your treatment for a fully immersed experience.

Figure 2. QR code to download the app

4. Getting to Know the Device

- A Ergonomic grip
- B Multi-Functional Power/Speed Button
- C LED light battery indicator light
- D LED light speed indicator light
- E USB-C charging port
- F Interchangeable percussive attachments

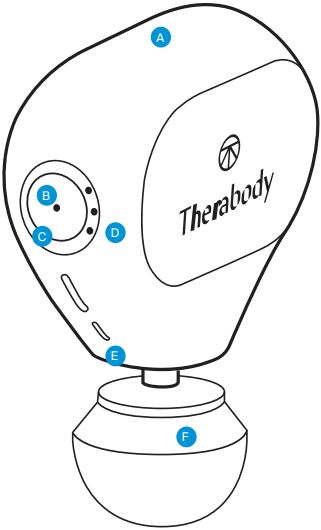
Figure 3. Corresponding pictures/diagrams of the device

5. Using the Device

- 1. Turn on the device.** To turn on the device, firmly press and hold the button on the left side of the device.
- 2. Choose your treatment speed.** Press the button to toggle between the three speeds (low, medium, and high). The white LED light speed indicator will light up depending on the selected speed setting (low: 1 light, medium: 2 lights, high: 3 lights).
- 3. Use the device.** Move the device back and forth across your intended treatment areas. Apply light to moderate pressure as needed. Never apply excessive pressure or any level of pressure that causes pain, discomfort, or impedes the normal function of the device. Do not start by putting pressure directly on an area of discomfort. Instead, work slowly around the area to reduce tension.
- 4. Turn off the device.** To turn off the device, firmly press and hold the button until the device turns off. Warning: The device does not automatically shut off. To turn off your device, firmly press and hold the button until the device turns off.
- 5. Lock the device.** When not in use or during transport, lock the device to prevent turning it on accidentally. Quickly press the button three times to initiate the lock function. Each LED light will flash sequentially to indicate that the device is locked. Quickly press the button three times again to unlock. Each LED light will flash sequentially (in the opposite direction) to indicate that the device is unlocked.

6. Attachments

The device comes with three attachments scientifically designed to target aches, pains, and tension all over the body.

Dampener

For tender or sensitive muscles.

The Dampener attachment is specifically designed for use on tender or sensitive muscles. With lower density foam and a wide base that maximizes contact with the body, this attachment delivers a soothing, low-impact treatment to effectively reduce aches and pains and keep you feeling your best.

Standard
Ball

For large & small muscle groups.

The Standard Ball attachment is specifically designed for all-over use on both large and small muscle groups. This attachment offers a more targeted percussive therapy treatment with a firm, smaller surface area and medium impact level. Perfect for quick re-activation and easing soreness, the Standard Ball is your go-to during a workout, delivering effective relief while keeping you moving.

Thumb

For lower back and trigger points.

The Thumb attachment is specifically designed to deliver the precision of a human thumb, making it ideal for trigger point therapy and targeting the lower back. Designed with a higher impact level and a shape that mimics the feel of an actual thumb, this attachment offers a dense, focused treatment for those stubborn areas. Perfect for anywhere you'd use a real thumb, the Thumb attachment provides the deep, targeted relief you need, all from the comfort of your home.

How to Connect and Remove Attachments

1. Ensure that the device is powered off before connecting or removing any attachments.
2. Align the attachment with the tip of the arm and push to connect.
3. To remove the attachment, hold the device with one hand, and the attachment with the other hand. Place one finger on each side of the plastic connector and firmly pull off the attachment.
4. **Warning:** Do not place your finger between the attachment and the end of the device, as this may cause pinching or injury.

**Attachments can be used across other Theragun Mini device models.
It is recommended that you swap out the attachments if you start to notice signs of wear.*

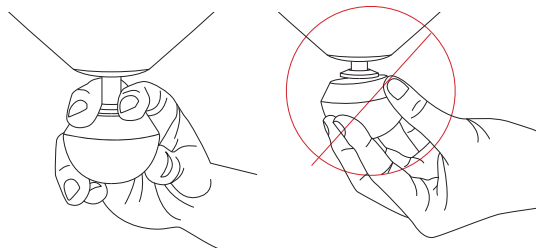
Figure 4. Corresponding diagram of connecting attachments

7. Care, Cleaning, and Charging

Device Maintenance

The following maintenance instructions are important to ensure that your device continues to work as it was designed. Failure to follow these instructions may cause your device to stop working.

Cleaning the Device

Read the following instructions carefully before cleaning the device to ensure longevity of the device and its components. For the most hygienic device experience, sanitize it after each use.

1. Wipe the device with a damp cloth or alcohol-free cleansing wipe.
2. Use a new wipe to disinfect the surface area by wiping all sides in a downward motion.

Note: This device is not waterproof. Only clean with a damp cloth or alcohol-free cleansing wipe by wiping the device. Do not submerge the device in water or clean it under running water. Do not allow the device to come into contact with corrosive solutions, which can damage the device's appearance and function.

Charging the Device

1. The device is USB-C enabled and comes with a charging cable. Before charging, ensure that the device is powered off.
2. Connect one end of the charging cable to the USB-C charging port on the device. Connect the other end of the charging cable to a power adapter.
3. The LED light will indicate the battery status when the device is on or charging. Turn the device on to check the battery status.
 - Flashing Orange: Low battery
 - Flashing Blue: Half battery
 - Flashing Green: High battery

The battery life of the device is up to 120 minutes, depending on the attachment and treatment intensity used.

****Note:** If using a charger from a third party, ensure that it is from a certified manufacturer and has not suffered any structural damage. Do not try to use the device while charging.

8. Troubleshooting Guide

SCENARIO	ERROR CODE DISPLAYED	POTENTIAL CAUSE	POTENTIAL SOLUTIONS
The motor stops working	The power button LED indicator blinks orange 10 times and then turns off	The device is not charged	Use a power adapter and charging cable to connect to the device to a power outlet to charge.
	The power button LED indicator blinks orange, blue, and green	The battery pack has overheated during charging or use	Ensure the device is powered off, remove from charger, and let the device rest for 60 minutes.
	The 3 speed LED lights blink white 5 times then turn off	The motor has overheated	Ensure the device is powered off, remove from charger, and let the device rest for 60 minutes.
			If the issue still persists after verifying all of the above, the device may be damaged. Contact Therabody for further assistance. Visit https://www.therabody.com/us/en-us/support/support.html for the available contact methods.
Cannot charge the device	N/A	The power adapter cable is not well connected to the adapter socket on the console	Verify that the power adapter cable is well connected to the adapter socket on the console, and the power adapter is connected to the 100-240 Volt power outlet. If this doesn't work, contact Therabody for further assistance. Visit https://www.therabody.com/us/en-us/support/support.html for the available contact methods.
		The power outlet connected is beyond the range of 100-240 Volt	
The device is charging very slowly	Slow charging or incompletely charging	There may be an adapter issue, or the adapter power is too low	Ensure you are using the USB-C to USB-C cable provided, or one from a certified manufacturer.
			If the issue still persists, Contact Therabody for further assistance. Visit https://www.therabody.com/us/en-us/support/support.html for the available contact methods.

SCENARIO	ERROR CODE DISPLAYED	POTENTIAL CAUSES	POTENTIAL SOLUTIONS
Mechanical noise	There is a rattling noise when the device is turned on	Internal malfunction	Contact Therabody for further assistance. Visit https://www.therabody.com/us/en-us/support/support.html for the available contact methods.
Bluetooth connectivity issue	Cannot connect to the app or the device isn't working properly once connected	The firmware is not updated to the latest version	Update the device firmware and try to connect to Bluetooth again.
		The phone version or system cannot support the connectivity	Update your phone to the latest version. If this doesn't work, try to install the app on another phone device.

Store and transport the device

Store and transport the device in the travel pouch. Store the device in a cool and dry place (Temperature: -20°C - 45°C. Relative humidity: 10~85% RH. Air Compression: 70 – 106 kPa). Do not store the device where temperatures may exceed 40°C/104°F, such as in direct sunlight or in a vehicle.

Warnings and Guidance (Precautions and Contraindications)

Background

Therabody products are designed to unlock the body's natural ability to achieve health and well-being. Through science and technology, the Therabody portfolio of products allows people to access the therapeutic benefits of different natural phenomena to meet both their needs and preferences. There will be times when it is advisable to modify how devices are used (precautions) or times when it is not appropriate to use certain devices (contraindications). Read the following safety information for the Theragun Mini in its entirety prior to use.

Important Safety Information

General Device Use

Read the full Warnings and Guidance prior to using the device. This device is intended for use by people in good health. This device is contraindicated against and should not be used by or on anyone with a history of epilepsy, seizures or cardiopathy.

The device is not recommended for anyone with an electronic implanted device (such as a pacemaker), cardiac arrhythmia, tumors, or acute episodes of inflammatory diseases. The device is not recommended for those who have arteriosclerosis, thromboses, or implants in the body region being treated.

Please consult your physician prior to using the device if you are pregnant and/or nursing. **Immediately stop using the device at the first sign of discomfort.** If you have any medical considerations, please consult your doctor before using the device.

Safety, Precautions, and Contraindications

These recommendations are derived from consultation with medical experts and published research regarding precautions and contraindications as of the printing date. For up-to-date information, please visit us online at <https://www.therabody.com/us/en-us/precautions-and-contraindications.html>.

Precautions:

Due care is required in these circumstances and device use may need to be modified. Consult with a medical professional if you currently have or suspect you may have any of the following conditions or if you have any questions.

- Hypertension (controlled)
- Osteopenia
- Osteoporosis
- Pregnancy
- Diabetes
- Varicose veins
- Bony prominences or regions
- Abnormal sensations (e.g., numbness)
- Sensitivity to pressure
- Recent injury or surgery
- Scoliosis or spinal deformity
- Medications that may alter client sensations

Contraindications:

The following are circumstances where the potential risks may outweigh the benefits.

Consult a medical professional before use.

- Skin rash, open wounds, blisters, local tissue inflammation, bruises, or tumors
- Deep vein thrombosis, osteomyelitis
- Bone fracture or myositis ossificans
- Hypertension (uncontrolled)
- Acute or severe cardiac, liver, or kidney disease
- Neurologic conditions resulting in loss or altered sensation
- Direct application to the Face, Throat, or Genitalia
- Bleeding disorders
- Recent surgery or injury
- Connective tissue disorders
- Peripheral vascular insufficiency or disease
- Medications that thin the blood or alter sensations
- Direct pressure over surgical site or hardware
- Direct pressure over eyes or throat
- Extreme discomfort or pain felt by client
- Severe scoliosis or spinal deformity
- Pacemaker, ICD, or history of embolism
- Allergy to device material: plastic, silicone, aluminum, stainless steel, copper alloy, paint

Limited Warranty and Customer Service

For full warranty information, please visit www.therabody.com/warranty. Customers who are in need of product support should visit <https://www.therabody.com/us/en-us/support/support.html> for the available contact methods.

FDA-Registered
Limited Warranty Only With Authorized Retailer Purchase

© 2025 Therabody, Inc. All Rights Reserved.
Patents at www.therabody.com/patents

Manufactured and Distributed by:
Therabody, Inc.
1640 S Sepulveda, Suite 300
Los Angeles, CA 90025
USA

In the Australian market, the Therabody Theragun Mini device warranty is 1 year. For Australian Market Only: The purchase of this device is accompanied by a limited warranty. For complete details about the rights afforded to you under this warranty, visit <https://therabody.com/anz/en-au/warranty-support/>.

Importers:
Therabody Australia Pty Ltd
Waterman Business Centre
175 Maroonah Highway
Suite 16, C307, Level 3
Ringwood Vic 3134 Australia

Therabody International Limited
5th Floor, 40 Mespil Road,
D04C2N4 Dublin, Ireland

Therabody UK Limited
Portsoken House 155-157
The Minorities, London,
United Kingdom, EC3N1LJ

Capacity: 1200mah
Chemistry: Li-Ion
Weight of the battery: 90.5 g/0.20lbs

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UNIT WARNINGS

READ ALL INSTRUCTIONS AND CAUTIONARY MARKINGS IN THIS MANUAL, ON THE CHARGER, AND ON THE DEVICE BEFORE USING OR CHARGING THE DEVICE AND ITS ACCESSORIES.

When using the device, the following precautions should always be adhered to:

1. **USE ONLY AS INSTRUCTED.** Use only Therabody recommended accessories and replacement parts. This device is not intended to diagnose, mitigate, or treat diseases. The device is an advanced mechanical tool with electric components. If the device and its accessories are not used or maintained properly, there is a risk of fire, electric shock, or injury. Failure to use and maintain the device and corresponding attachments in accordance with the instructions outlined in this manual will void the warranty and may result in product damage or physical injury. The device is intended for commercial and home use.
2. **NOT FOR CHILDREN.** The device is not intended for use by young children or persons with reduced physical, sensory, or reasoning capabilities or lack of experience and knowledge, unless given supervision or instruction concerning use of the device by a person responsible for their safety. The device is not a toy. Do not play with, bend, or pull the electrical components. Advise children not to play with the device or accessories. Adult supervision should be provided for those under 18 using the device.
3. **CHARGING.** If your device doesn't turn on or the battery indicator displays a low battery level, please charge before first use. Charge the battery using the USB-C cable (included). When you use the charging cable (included), make sure its USB-C connector is fully inserted into a compatible power adapter before you plug the adapter into a power outlet. You can also charge the battery with third party cables that are compliant with USB-C and with applicable country regulations and international and regional safety standards. Other adapters may not meet applicable safety standards, and charging with such adapters could pose a risk of death or injury. Do not force the cable into place. Fully charge the device at least once every six months to prevent device and battery damage. Unplug the power adapter and charging cable when not in use.
4. **DO NOT OVERCHARGE.** Do not leave the device connected to the power adapter for more than one hour after the battery has been fully charged. The battery includes a protection system to avoid the risk of overcharging. However, overcharging may reduce its life over time.
5. **CHARGING LOCATIONS.** Ensure that the battery status LED is OFF before charging. Insert the USB-C cable into the port (behind the power button) to charge the device. Do not use the device while charging. Avoid skin contact when the USB-C is plugged into a power source as it may cause discomfort or injury. Do not sleep or sit on the USB-C connector. The device should be charged indoors in a well-ventilated, dry location. Do not charge the device outdoors, in a bathroom, or within 10 feet (3.1 meters) of a bathtub or pool. Do not use the device or charger on wet surfaces, and do not expose the charger to moisture, rain, or snow. Do not use the device or its compatible charger in the presence

of explosive atmospheres (gaseous fumes, dust, grain, metal powders, or flammable materials). Sparks may be generated, possibly causing a fire.

6. **UNPLUGGING THE CHARGER.** Pull the plug, not the cord, to reduce the risk of damage to the electrical plug and cord. Never carry the charger by its cord. Keep the cord away from heat, oil, and sharp edges. Do not stretch the charger cord or place the cord under strain. Do not handle the charger, terminals, or the device with wet hands.

Always unplug this appliance from the electrical outlet immediately after using and before cleaning. Store cable to ensure it is not stepped on, tripped over, or otherwise subjected to damage or stress. Do not use a power adapter or charging cable that has received a sharp blow, been dropped, run over, or damaged in any way. For long-term storage, store with a fully charged battery.

7. **DO NOT CRUSH, DROP, OR DAMAGE THE DEVICE OR CHARGER.** Do not use a power adapter or charging cable that has received a sharp blow, been dropped, run over, or damaged in any way. Do not puncture or damage the device. Puncturing the device can lead to incorrect operation.

8. **POWER ADAPTER AND CHARGING CABLE CARE.** Unplug the power adapter and charging cable when not in use. Pull the plug, not the cable, to reduce the risk of damage to the electrical plug and cable. Store cable to ensure it is not stepped on, tripped over, or otherwise subjected to damage or stress. Keep the cable away from heated surfaces, oil, and sharp edges. Never operate this device if it has a damaged cable or plug if it is not working properly, if it has been dropped or damaged, or dropped into water. Do not stretch the charger cable or place the cable under strain. Do not handle the cable with wet hands. For long-term storage, store with a fully charged battery. Therabody is not responsible for damages that may occur due to the use of third-party chargers.

9. **USAGE.** Discontinue use of the device if you are feeling severe pain. Make sure the Battery status LED is OFF when the device is not in use. While in use, hold the device in your hands only.

10. **DO NOT OPERATE UNDER BLANKET AND PILLOW.** Excessive heating can occur and cause fire, electric shock, or injury.

11. **SERVICE.** If the device is not working properly, has received a sharp blow, or has been dropped, damaged, left outdoors, or dropped into water, then do not use it. Do not attempt to repair or disassemble the device which may result in an electric shock or fire.

12. **DEVICE CARE.** Handle the device with care. The device contains electronic components that can be damaged if dropped, burned, punctured, or crushed. Do not disassemble the device or attempt to repair it yourself. Disassembling the device may damage it, result in loss of water resistance, and may cause injury. If the device is damaged or if malfunctions occur, contact Therabody Customer service. Do not use a damaged device, for example, visible water intrusion, as these may cause injury. Avoid heavy exposure to dust or sand.

13. **CLEANING INSTRUCTIONS.** Visually inspect the device for any obvious signs of debris build-up. Wipe your device with a damp cloth or alcohol-free cleansing wipe. After cleaning, allow the device to dry thoroughly before storing or beginning another treatment. A properly cleaned device should have no visible signs of debris or moisture. Note: This product is not waterproof. Only clean with a damp cloth or alcohol-free cleansing wipe by wiping the device. Do not submerge the device in water or clean it under running water. Do not allow the device to come into contact with corrosive solutions, which can damage the device's appearance and function.

14. **DO NOT DISASSEMBLE.** Disassembly or incorrect reassembly may result in the risk of electric shock, fire, or exposure to battery chemicals. The warranty will be void if the device, batteries, or charger are disassembled or if any parts have been removed.

15. **STORING THE DEVICE AND CHARGING CABLE.** Store in a cool, dry place. Only charge the device when the ambient temperature is between 40°C/104°F and 0°C/32°F. Do not store the device or charging cable where temperatures may exceed 70°C/158°F, such as in direct sunlight, in a vehicle, or in a metal building during the summer.

Reporting adverse events to FDA

MedWatch is the Food and Drug Administration's (FDA) program for reporting serious reactions, product quality problems, therapeutic inequivalence/failure, and product use errors with human medical products, including drugs, biologic products, medical devices, dietary supplements, infant formula, and cosmetics.

If you think you or someone in your family has experienced a serious reaction to a medical product, you are encouraged to take the reporting form to your doctor. Your health care provider can provide clinical information based on your medical record that can help the FDA evaluate your report. However, we understand that for a variety of reasons, you may not wish to have the form filled out by your healthcare provider, or your health care provider may choose not to complete the form. Your health care provider is not required to report to the FDA. In these situations, you may complete the Online Reporting Form yourself. You will receive an acknowledgement from the FDA when your report is received. Reports are reviewed by FDA staff. You will be personally contacted only if we need additional information.

Submitting adverse event reports to FDA

Use one of the methods below to submit voluntary adverse event reports to the FDA:

1. Report Online at: www.accessdata.fda.gov/scripts/medwatch/index.cfm?action=reporting.home
2. Consumer Reporting Form FDA 3500B. Follow the instructions on the form to either fax or mail it in for submission. For help filling out the form, see MedWatchLearn. The form is available at: www.fda.gov/downloads/aboutFDA/reportsmanualsforms/forms/ucm349464.pdf
3. Call FDA at 1-800-FDA-1088 to report by telephone.
4. Reporting Form FDA 3500 commonly used by health professionals. The form is available at: www.fda.gov/downloads/aboutFDA/reportmanualsforms/forms/ucm163919.pdf

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Medical EMC + Safety
IEC 60601-1 and IEC60601-11 Requirement (Warning)—For Medical device apply for US, CA

This device is Class III equipment with type BF applied. It complies with Medical Electrical Safety Standards (IEC 60601-1).

This device also complies with Medical EMC Standard (IEC 60601-1-2).

The has been tested and found to comply with the electromagnetic compatibility (EMC) limits for medical devices to IEC 60601-1-2: 2014. These limits are designed to provide reasonable protection against harmful interference in a typical medical installation.

1* WARNING: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

2* WARNING: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

3* WARNING: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the ME EQUIPMENT, including cables specified by the manufacturer. Otherwise, degradation of the performance of this equipment could result.

Please scan QR code for a copy of the Declaration of Conformity under Regulation (EU) 2019/1020

EN

No.	Icon	Description	Location of Label
1		CE Mark	On rating label
2		Please read the user manual before use.	In packaging
3		Separate collection for waste electrical and electronic equipment. Note: For more information about disposal of equipment, its parts and accessories, please contact your local distributor.	On rating label
4	IP22	IP Classification: The first number 2: Protected against access to hazardous parts with a finger, and the jointed test finger of 12 mm Φ , 80 mm length; shall have adequate clearance from hazardous parts, and protected against solid foreign objects of 12.5 mm Φ and greater.	On rating label
5		RCM Mark	On rating label and packaging
6		Type BF apply part	On rating label
7		Manufacturer information	On rating label
8	UK CA	UKCA Mark for UK	On rating label
9		Bluetooth	On packaging
10	IC	IC Canada	On rating label
11	 UDI : (01) 00840295405810 (10) 0000 (21) 00000	Unique Device Identification (UDI)	On rating label and packaging
12	 UDI : (01) 00840295405810 (10) 0000 (21) 00000	Date of Manufacture	On rating label

**GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC EMISSION -
FOR ALL EQUIPMENT AND SYSTEMS**

Guidance and manufacturer's declaration – electromagnetic emission		
The Theragun Mini is intended for use in the electromagnetic environment specified below. The customer or the user of Theragun Mini should assure that it is used in such an environment.		
EMISSIONS TEST	COMPLIANCE	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
RF emissions CISPR 11	Group 1	The Theragun Mini uses RF energy only for its internal function. There for, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The Theragun Mini suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage function / flicker emissions IEC 61000-3-3	Not applicable	

EN

GUIDANCE AND MANUFACTURER'S DECLARATION - ELECTROMAGNETIC IMMUNITY - FOR ALL EQUIPMENT AND SYSTEMS

Guidance and manufacturer's declaration – electromagnetic immunity			
The Theragun Mini is intended for use in the electromagnetic environment specified below. The customer or the user of the Theragun Mini should assure that it is used in such an environment.			
IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	± 8 kV contact ± 2 kV, ± 4 kV, ± 8 kV, ± 15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrostatic transient / burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	N/A	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±1kV differential mode ±2kV common mode	N/A	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0 % UT; 0.5 cycle g) At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0 % UT; 1 cycle and 70 % UT; 25/30 cycles Single phase: at 0° 0 % UT; 250/300 cycle 0 % UT; 250/300 cycle	N/A	Mains power quality should be that of a typical commercial or hospital environment. If the user of the Theragun Mini requires continued operation during power mains interruptions, it is recommended that the Theragun Mini be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz magnetic field IEC 61000-4-8	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE U _i is the a.c. mains voltage prior to application of the test level.			

Guidance and manufacturer's declaration – electromagnetic immunity			
The Theragun Mini is intended for use in the electromagnetic environment specified below. The customer or the user of the Theragun Mini should assure that it is used in such an environment.			
IMMUNITY TEST	IEC 60601 TEST LEVEL	COMPLIANCE LEVEL	ELECTROMAGNETIC ENVIRONMENT - GUIDANCE
Conducted RF	3 Vrms	N/A	Portable and mobile RF communications equipment should be used no closer to any part of the Theragun Mini, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.Recommended separation distance: $d = \left[\frac{3.5}{V_1} \right] \sqrt{P}$ $d = \left[\frac{12}{V_2} \right] \sqrt{P}$ $d = \left[\frac{3.5}{E_1} \right] \sqrt{P}$ $d = \left[\frac{7}{E_1} \right] \sqrt{P}$ where p is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m).^a Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey,^b should be less than the compliance level in each frequency range.^c Interference may occur in the vicinity of equipment marked with the following symbol: 
IEC 61000-4-6	150kHz to 80MHz		
Radiated RF	10 V/m	10 V/m	
IEC 61000-4-3	80MHz to 2.7GHz	80MHz to 2.7GHz	
	385MHz-5785MHz Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communication equipment (Refer to table 9 of IEC 60601-1-2:2014)	385MHz-5785MHz Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communication equipment (Refer to table 9 of IEC 60601-1-2:2014)	
NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies. NOTE 2 These guidelines may not apply in all situations. Electromagnetic is affected by absorption and reflection from structures, objects and people.			

a The ISM (industrial, scientific and medical) bands between 150 kHz and 80 MHz are 6,765 MHz to 6,795 MHz; 13,553 MHz to 13,567 MHz; 26,957 MHz to 27,283 MHz; and 40,66 MHz to 40,70 MHz. The amateur radio bands between 0,15 MHz and 80 MHz are 1,8 MHz to 2,0 MHz, 3,5 MHz to 4,0 MHz, 5,3 MHz to 5,4 MHz, 7 MHz to 7,3 MHz, 10,1 MHz to 10,15 MHz, 14 MHz to 14,2 MHz, 18,07 MHz to 18,17 MHz, 21,0 MHz to 21,4 MHz, 24,89 MHz to 24,99 MHz, 28,0 MHz to 29,7 MHz and 50,0 MHz to 54,0 MHz.

b Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the **Theragun Mini** is used exceeds the applicable RF compliance level above, the **Theragun Mini** should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the **Theragun Mini**.

c Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3V/m.

Recommended separation distances between portable and mobile RF communications equipment and the EQUIPMENT or SYSTEM -
for EQUIPMENT and SYSTEMS

Recommended separation distances between portable and mobile RF communications equipment and the Theragun Mini

The **Theragun Mini** is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the **Theragun Mini** can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the **Theragun Mini** as recommended below, according to the maximum output power of the communications equipment

Rated maximum output of transmitter W	Separation distance according to frequency of transmitter m			
	150 kHz to 80 MHz outside ISM and amateur radio bands	150 kHz to 80 MHz in ISM and amateur radio bands	80 MHz to 800 MHz	80 MHz to 800 MHz
	$d = \left[\frac{3.5}{V_1} \right] \sqrt{P}$	$d = \left[\frac{12}{V_2} \right] \sqrt{P}$	$d = \left[\frac{12}{V_2} \right] \sqrt{P}$	$d = \left[\frac{7}{E_1} \right] \sqrt{P}$
0.01	0.12	0.20	0.035	0.07
0.1	0.38	0.63	0.11	0.22
1	1.2	2.00	0.35	0.70
10	3.8	6.32	1.10	2.21
100	12	20.00	35	70

For transmitters rated at a maximum output power not listed above the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

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