
MANUEL D'UTILISATION – INSTRUCTION MANUAL
THERMO SMART

CE0197

Terraillon®
www.terraillon.com

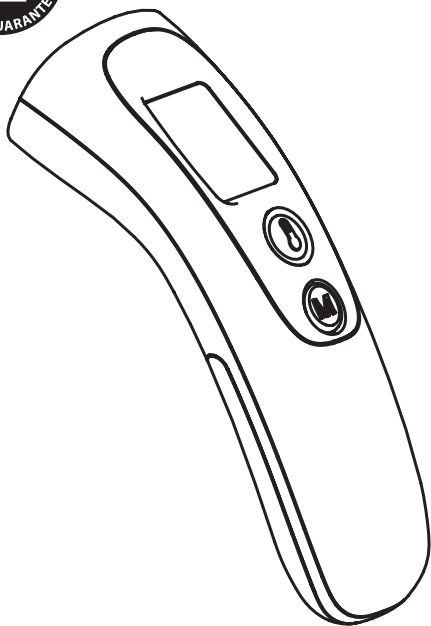THERMO SMART

FR	THERMOMÈTRE SANS CONTACT CONNECTÉ NOTICE D'UTILISATION	4
EN	CONNECTED NON-CONTACT THERMOMETER INSTRUCTION MANUAL	15
DE	KONTAKTLOSES VERBUNDENES THERMOMETER BEDIENUNGSANLEITUNG	26
NL	VERBONDEN CONTACTLOZE THERMOMETER GEBRUIKSAANWIJZING	37
IT	TERMOMETRO SENZA CONTATTO CONNESSO LIBRETTO DI ISTRUZIONI	48
SV	BERÖRINGSFRI INFRARÖD TERMOMETER BRUKSANVISNING	59
PL	TERMOMETR BEZDOTYKOWY INSTRUKCJA OBSŁUGI	70

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Terraillon®

CONNECTED NON-CONTACT THERMOMETER

Dear customer,

You have just bought this product and we thank you for your purchase. We hope you enjoy using the product. In order to be completely satisfied, we recommend that you read these directions carefully.

This thermometer is intended for home use by customer only.

The Non-contact Forehead Thermometer is a non-sterile, reusable clinical thermometer intended for the intermittent determination of human body temperature on the centre of the forehead as the measurement site.

If you have any difficulties or need more information, you can find our help centre by scanning the QR code opposite with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: serviceconsommateurs@terraillon.fr

RECOMMENDATION

- Store and use the thermometer in the same environment than the one in which the device is stored.
- Do not use in areas where it may be subject to room temperature variations.
- Determine your normal temperature range by taking your temperatures when healthy. These readings can be used as a form of reference to compare with future readings when ill.
- Before taking your temperature, make sure your forehead is free of sweat and hair, and make sure the infrared sensor is clean.
- Do not take temperature within 30 minutes of eating, exercising, or taking a bath.
- Wait one minute between each temperature reading. Two consecutive readings may give slightly different measurements. Should this occur, take the average of the two.
- Do not attempt to disassemble the thermometer.
- Keep thermometer away from direct sunlight and store well away from any source of heat in a clean, dry location.
- Do not drop or break the thermometer. Do not use if broken. Do not disassemble.
- If these recommendations are not respected, the readings may be inaccurate.
- Temperature will vary depending on the part of the body. This is why you should not compare rectal and forehead temperature readings, for example.
- Temperature will also vary depending on the person, age, and time of day.
- Don't store the thermometer to extremes temperature or extremes humidity. Failing to do so may cause inaccuracy. Don't operate this thermometer under extremes temperature or extremes humidity. Failing to do so may cause inaccurate measurement.
- When the thermometer is used to measure the temperature of a child, the temperature measurement must be carried out by an adult. Only use the thermometer as described in this user manual. Any other use is deemed improper and may result in damage to property or even personal injury. Do not use any accessories which have not specified in this instruction manual. The time required for the device to warm from the minimum storage temperature (-10°C) to ambient temperature (20 °C) until the device is ready for use: about 2 hours. The time required for the device to cool from the maximum storage temperature (55°C) to ambient temperature (20 °C) until the device is ready for use: about 2 hours.

CONNECTED NON-CONTACT THERMOMETER


You should contact the manufacturer or the manufacturer's representative:

- for assistance, if needed, in setting up, using or maintaining the thermometer; or
- to report unexpected operation or events.

To ensure an optimal temperature reading

- Take a reference temperature.
- Make sure the product is stored in the same environment as that in which the temperature will be taken.
- If skin has recently been exposed to hot or cold conditions, this will disrupt the measurement.
- Avoid bathing or showering just before taking the temperature.
- Do not take the temperature in a draught (fan, air conditioning).

PRODUCT DESCRIPTION

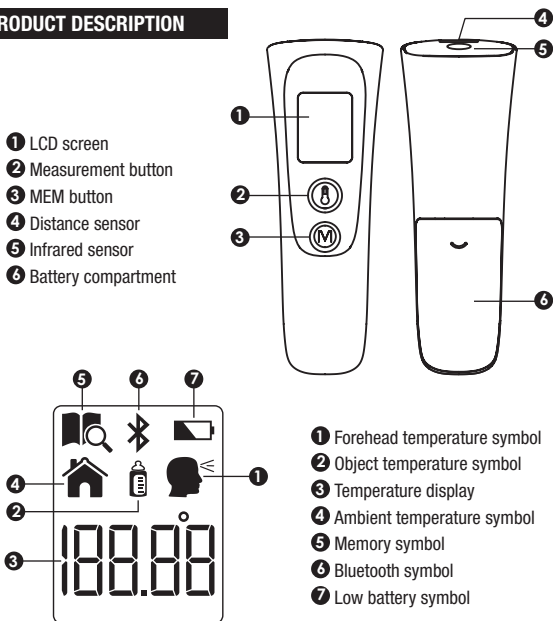
Terraillon®

CONNECTED NON-CONTACT THERMOMETER

INSERTING & CHANGING BATTERIES

- To insert or change the 2 LR3 batteries (supplied with the thermometer), open the battery compartment by sliding off its cover, insert the batteries ensuring that the negative and positive ends are in the right place, and replace the cover before use.
- When the low battery indicator  on the LCD, or the thermometer does not function at all, you should replace new battery as soon as possible.

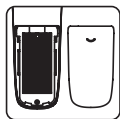CHOOSING THE TEMPERATURE MEASURING UNIT (C° OR F°)

- Press one of both buttons to turn on thermometer. Then, to switch from one measuring unit to the other, press around 5 seconds [M].


FIRST USE

- A** Download the MyHealth – Wellness Coach app.

- B** You must have an account to be able to use all the features of your Thermo Smart. To create your account, simply follow the instructions in the app.

- C** Activate the Bluetooth function on your smartphone (Settings > Bluetooth > ON).

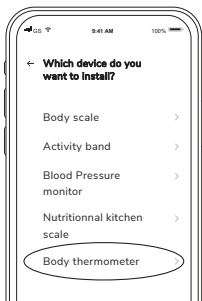If you need to create a MyHealth account:

When prompted by the app, select Non-contact smart thermometer.

If you already have a MyHealth account:

You can use your Thermo Smart straightaway – see chapter “HOW TO ADD A TEMPERATURE READING IN THE APP”, p.18


CONNECTED NON-CONTACT THERMOMETER

Follow the instructions on the mobile application.

Did you know?

You can create and assign a temperature to a family member (up to a maximum of 8 family members) in the MyHealth app.

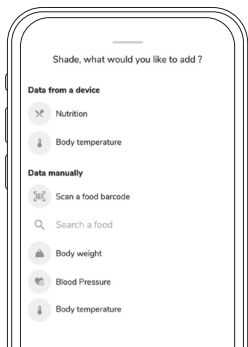HOW TO ADD A TEMPERATURE READING IN THE APP

There are two options:

1. Take your temperature with your Thermo Smart (the thermometer must not be too far away from your smartphone): Follow this path: **+** - Data from a device - Thermo Smart. Now you can take your temperature with your Thermo Smart and follow the steps.
2. Take your temperature without your Thermo Smart (example scenario: you are at the doctor's surgery, you take your temperature without using your Thermo Smart but you want to enter the reading in your MyHealth app): Follow this path: **+** - Manual data - Thermo Smart.

Did you know?

1. The Thermo Smart is not added to the list of "devices" (no need to pair the product with the app); the thermometer can be used immediately by adding a temperature reading using the **+** button at the bottom centre of the app.
2. You can use your Thermo Smart without the app being open.

CONNECTED NON-CONTACT THERMOMETER

USING THE PRODUCT

FOREHEAD TEMPERATURE

- Place the thermometer in your hand as shown with your thumb on the [measure] button.
- Press the [measure] button  to turn on the thermometer.
- Place the thermometer at a distance of 0 to 5 cm from the forehead (specifically at the temples) and press the "take temperature" button  again.
- A double beep sound indicates that the distance is correct.
- You can read the temperature on the display and the color indicator (red or green) around the [measure] button helps to analyze the temperature, if fever or not.

Green: 34.0°C (93.2°F) ~ 37.9°C (100.3°F)

Red: 38.0°C (100.4°F) ~ 42.9°C (109.3°F)

- Temperature is automatically memorised.

Did you know?

1. The last temperature reading can be transferred to the MyHealth app at any time. Just open your app, go to  and follow the chapter "HOW TO ADD A TEMPERATURE READING IN THE APP" on p.18.
2. To comprehensively track your temperature in your application, you can add a note (e.g. symptoms) or medication taken at a specific time.

AMBIENT TEMPERATURE TAKEN

- Place the thermometer in your hand as shown with your thumb on the [measure] button  the thermometer automatically switches to room temperature mode .

CONNECTED NON-CONTACT THERMOMETER

CAUTION:

After press the measure button, the distance sensor will detect the distance between the forehead and thermometer. If the distance is between 1 - 5 cm, it will directly test the forehead temperature.

If the distance is more than 5cm or less than 1cm, short successive beep sounds will be issued indicating that the thermometer should be moved to the forehead in 1 - 5 cm. When the correct distance is reached a double beep signals it.

If it is not the proper distance for forehead temperature, the thermometer will measure ambient temperature.

CAUTION:

External factors may influence forehead temperature particularly when an individual has:

- been exposed to very hot or very cold temperatures.
- been recently swimming or bathing.
- had their forehead covered.

In these cases, let the individual stay away from the above scenario and station in a room for 30 minutes prior to talking a temperature.


OBJECTS OR LIQUID TEMPERATURE

To measure the temperature of the surface of an object or liquid, put the thermometer in object mode.

- Place the thermometer in your hand as shown with your thumb on the [measure] button  the thermometer automatically switches to room temperature mode .
- Press the button  and [M] at the same time to enter object mode and display the screen opposite:

- Hold the thermometer and ensure the probe is around 0-5 cm from the liquid or surface whose temperature you want to measure.
- Press the [measure] button  to start the measurement. A beep sounds signals the end of measurement. Remove the thermometer from the object. The LCD displays the measured temperature.

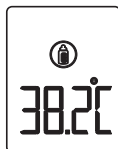- To return to room temperature mode, press [M].
- The thermometer will automatically return back to standby mode and to body temperature mode after 1 minute.

CONNECTED NON-CONTACT THERMOMETER


CAUTION:

- Make sure you do not use this mode to take your body temperature.
- No Bluetooth transmission to the APP for this data.
- With object measurement, the current surface temperature of the object is displayed. It can be different from its internal temperature, especially if the surface is exposed to direct sunlight or a draught.

CHECK THE MEMORY

- There are a total of 30 memory data for recording forehead temperature measurements.
- The current measurement is always stored in the last storage space. When all storage spaces have been occupied, the oldest measurement is always deleted from the memory.
- Press [M] button to turn on thermometer and [M] button again to enter in memory data .

The LCD displays memory number followed by measured temperature like figures. Press the [MEM] button repeatedly to go to the memory number you need (30 is the youngest, 1 the oldest).

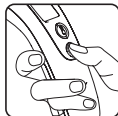PLEASE NOTE:

If you remove or replace the batteries, the temperatures taken with your Thermo Smart will not be saved.

MAINTENANCE

- Clean the thermometer with a soft, dry cloth.
- Clean both lenses with a cotton swab moistened with alcohol and immediately wipe dry with a clean cotton swab. After cleaning, allow at least 20 minutes drying time before taking temperatures.

CONNECTED NON-CONTACT THERMOMETER

- Do not touch either lens with your fingers.
- Do not submerge the thermometer in water or any other liquid.
- Do not rinse the thermometer in water.
- Do not use water directly on the LCD screen.
- Store the thermometer in the dry, clean place.
- If the thermometer is not used for a long time, remove the batteries.

SAFETY

- Keep out of reach of children and to avoid inhalation or swallowing of small parts. Do not allow children to take their temperatures unattended.
- Do not use this appliance for any other purposes than those for which it was designed (human body temperature).
- Use of this thermometer is for measuring purposes only and is not intended as a diagnostic tool. Make sure that any unusual readings are communicated to your physician. Under no circumstances should you change your prescribed medicine doses.
- Use of this thermometer is not intended as a substitution for consultation with your physician. Measurement results are for reference only. Contact your physician if you have or suspect any health.

PROTECT THE ENVIRONMENT



Dispose of used batteries in a designed container so that they can be collected and recycled. Do not mix different types of batteries. Do not mix new and used batteries. After life, dispose of this product at a designated waste recycling point.

GUARANTEE

This product is guaranteed against defects in materials or manufacturing. During the guarantee period, any such defects will be replaced free of charge (proof of purchase must be shown in the case of a complaint under guarantee). The guarantee does not cover damage resulting from improper use or negligence. In the event of a claim, first contact the shop where you purchased the product.

CONNECTED NON-CONTACT THERMOMETER

TROUBLESHOOTINGS

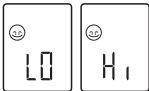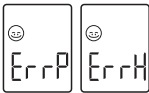Troubles	Checklists	Counter measures
	The measured temperature is lower than 32.0°C/89.6°F or higher than 42.9°C/109.3°F Please check the operation method.	Follow User's Manual for proper measurement.
	Hardware problem.	Contact your distributor.
	Operating temperature is out of the range.	Use the thermometer in the range of operating conditions.
The measurement is not accurate or if there is any doubt on the measured result.	Please check if the infrared sensor is clean or not.	Clean the infrared sensor with cotton swab according to user manual.
	Please check the measuring way is correct or not.	Ensure you have read the manual and know how to use the thermometer properly. Please keep the thermometer and patient in the measuring room least 30 minutes before using. Please take the measurement indoor. Do not held the thermometer in your hand too long. It could affect the accuracy.
	Please check if there is low battery symbol on LCD.	Change new battery.

CONNECTED NON-CONTACT THERMOMETER

SPECIFICATIONS

Power supply:	DC 3V (2 x AAA battery)
Mode of operation:	Adjusted mode
Reference body site:	Core
Body measuring range:	32.0°C ~ 42.9°C (89.6°F ~ 109.3°F)
Laboratory accuracy:	±0.2°C(±0.4°F) for 35.0°C (95.0°F) ~ 42.0°C (107.6°F) ±0.3°C (±0.5°F) for other range
Object measuring range:	-22°C ~ 80.0°C (-7.6°F ~ 176.0°F)
Object measuring accuracy:	±2.0°C(±4.0°F)
Displayed room temperature range:	5.0°C ~ 59.9°C (41.0°F~139.8°F)
Room temperature accuracy:	±1°C (±2°F)
Measure time:	1 second
Measure interval:	1 second
Dimension:	130.1(L)x45.0(W)x50.3(D)mm
Weight:	About 62.0g (without battery)
Forehead and object operating condition:	15.0°C~40.0°C (59.0°F~104.0°F) with a relative humidity of 15%~85% Atmosphere pressure: 700hPa-1060hPa
Storage and transport condition:	-20°C~55°C (-4°F~131°F) with a relative humidity of 15%~90% Atmosphere pressure: 700hPa-1060hPa

Description of the Wireless Functions and Technology:

RF Frequency	2402 MHz - 2480 MHz
Modulation type	GFSK
Number of channel	40
Chanel spacing	2 MHz
Bluetooth version	V4.0
Maximum Output Power	8,31 dBm

Service life: 2 years

Service life of the battery: With a new battery (Carbon) approx. 6 months based on the use frequency at 5 times/day

CONNECTED NON-CONTACT THERMOMETER

Item		Group A1	Group A2	Group B	Group C
Clinical bias	Δ_{cb}	-0,185	-0,175	-0,173	-0,182
Uncertainty	L_A	0,357	0,349	0,408	0,391
Clinical repeatability	σ_r	0,102	0,081	0,118	0,212

**Caution:**

This thermometer has been thoroughly tested and inspected to assure proper performance and operation!

NORMALIZED SYMBOLS

	Applied part of type BF (The applied part is the probe tip)
	Disposal in accordance with Directive 2012/19/EU (WEEE)
	Complies with the European Medical Device Directive (93/42/EEC) and amended by directive 2007/47/EC requirements. Notified body TÜV Rheinland (CE0197)
	The name and the address of the manufacturer
	The name and the address of the Authorized EC-representative in Europe
	The name and the address of the authorised representative in Swiss
	Caution
	Refer to Instruction Manual
	The first number 2: Protected against solid foreign objects of 12,5 mm Φ and greater. The second number: Protected against vertically falling water drops when enclosure tilted up to 15°. Vertically falling drops shall have no harmful effects when the enclosure is tilted at any angle up to 15°, on either side of the vertical.
	Date of Manufacture

DECLARATION OF CONFORMITY

Hereby, Terraillon SAS declares that Thermo Smart complies with the essential health and safety conditions. This product meets the essential requirements and other relevant provisions of the directive RED 2014/53/ UE. The complete declaration of conformity is available on [https:// terraillonhelp. zendesk.com/hc/ en-us](https://terraillonhelp.zendesk.com/hc/en-us) (Declaration of Conformity section). E.I.R.P : 8,31 dBm.

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Shenzhen Dongdixin Technology Co., Ltd.
Shanghai International Holding Corp. GmbH (Europe)
MedNet SWISS GmbH
Terraillon SAS France & Headquarters
Terraillon UK Limited
Terraillon Asia Pacific Ltd