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
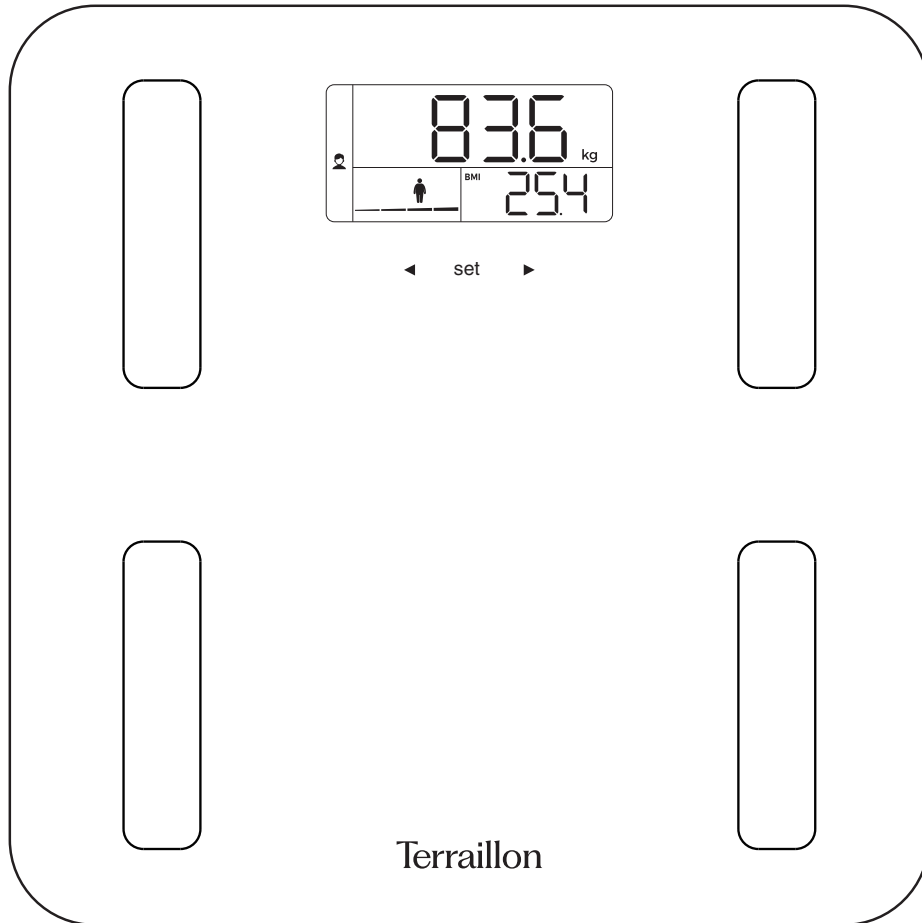
E.I.R.P: <=0 dBm.

Thank you very much for selecting the Terrailon Glass Body Fat Analyzer BODY CONNECT.


Don't use this product if you wear a pacemaker or implantable electronic device • Any information provided by this device is in no way meant to treat, cure or prevent any disease or illness from happening.

PRODUCT OVERVIEW

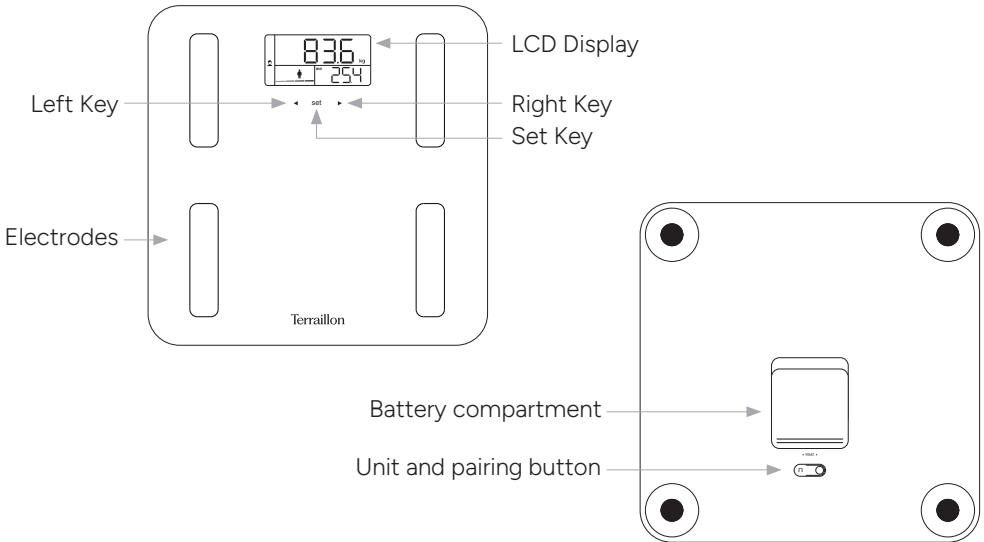LCD DISPLAY

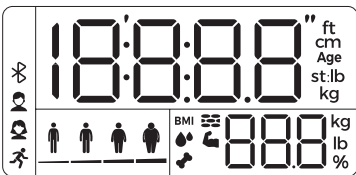Age	Age	cm	Height unit
	Gender	ft	
	Body muscle %	kg	Weight unit
	Body bone kg	lb	
	Body hydration %	st:lb	
	Body fat %		Underfat
	Athlete mode		Healthy
BMI	Body mass index		Overfat
	Bluetooth		Obese

If you have any difficulties or need more information, you can find our help centre by scanning the QR code opposite with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com

FIRST USE

- Insert the batteries (4xAAA).
- You will need your smartphone to install your bathroom scales.
- Download the Wellness Coach - MyHealth app, then open it.

- Follow the instructions that appear on your smartphone. Ensure you enter all the information requested.
- Your smartphone must remain near to the bathroom scales throughout the installation process.
- If you encounter any problems, please consult the help centre:

- It's important to bind your account with the scale for automatic recognition of user during weighing. You can pair the scale with 8 smartphones maximum. Before your first weighing, please select your profile number on the scale.
- If the user is not recognised, only the weight will be displayed.

START MEASURING


Measurement unit

The default measurement unit is "Kilogram". To select the desired unit (kg, lb or st), briefly press the button until the desired unit appears. The unit is then selected.

Initialization

Press the platform centre and remove your foot. "0.0" will be displayed.


The scale will switch off and is now ready for use. This initialisation process must be repeated if the scale is moved. At all other times step straight on the scale.

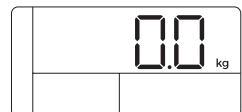To weigh and measure

1. You must have bare feet to be able to analyse your body mass.
2. Open the MyHealth application. If the application is not opened, the data will be saved in the scale.
3. Get on your scale. Stand still. Maintain contact with the electrodes until the "👤" signal stops moving from left to right.
4. Display the measurement


The LCD screen then displays BMI and body masses (fat percentage, water percentage, muscle mass and bone mass) for three times.

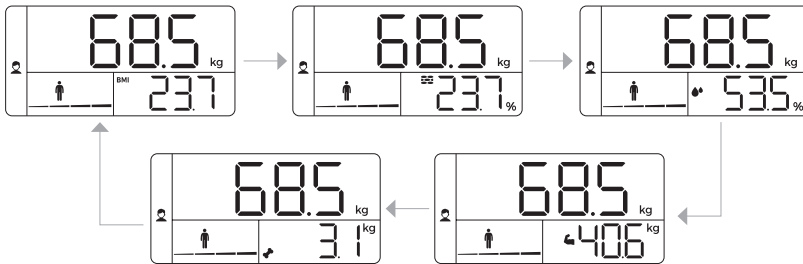When your scale is successfully paired with your smartphone and the Bluetooth is ON, the scale will process data transmission.

Automatically identify the user measurement.


1. When the user stand on the scale directly, the scale will turn on and lock the weight.
2. When the weight is stable, display moving "👤" means it is measuring. If it can identify the user successfully, it will display the measurement. If it detect several similar weights, the profile number flashes. Use the ◀ and ▶ buttons to select the user number, then press SET. Once you have selected the user, the measurements will be displayed and then the unit will turn off. If you do not select a user within 15 seconds, the weight will be displayed for a few seconds and then the unit will shut off (the weight will not be assigned to a user).
3. If it occurred error or it can't identify the user, it will display the present weight and turn off.

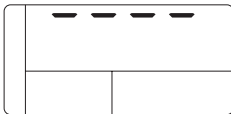WARNING INDICATORS

Overload:
remove weight immediately.

Battery low,
replace the batteries.

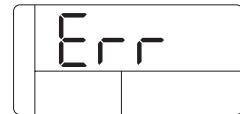Failure of pairing up
your scale
with your smartphone.

MAINTENANCE

Use water only. DON'T use propellant or other chemicals.

WARRANTY & ENVIRONMENT PROTECTION

This product is guaranteed for 2 years against material and manufacturing defects. During this period, such defects will be repaired free of charge (proof of purchase must be presented if the scale is under warranty). This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint please contact the store where you purchased the scale first.

Waste electrical products should not be disposed of with household waste. Please recycle where facilities exist. Check with your local authority or retailer for recycling advice.

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iPhone is registered trademark of Apple Inc.
Android is a trademark of Google Inc.
The Bluetooth word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.
Other trademarks and trade names are those of their respective owners.
Terraillon SAS France & Headquarters
Terraillon Asia Pacific Ltd