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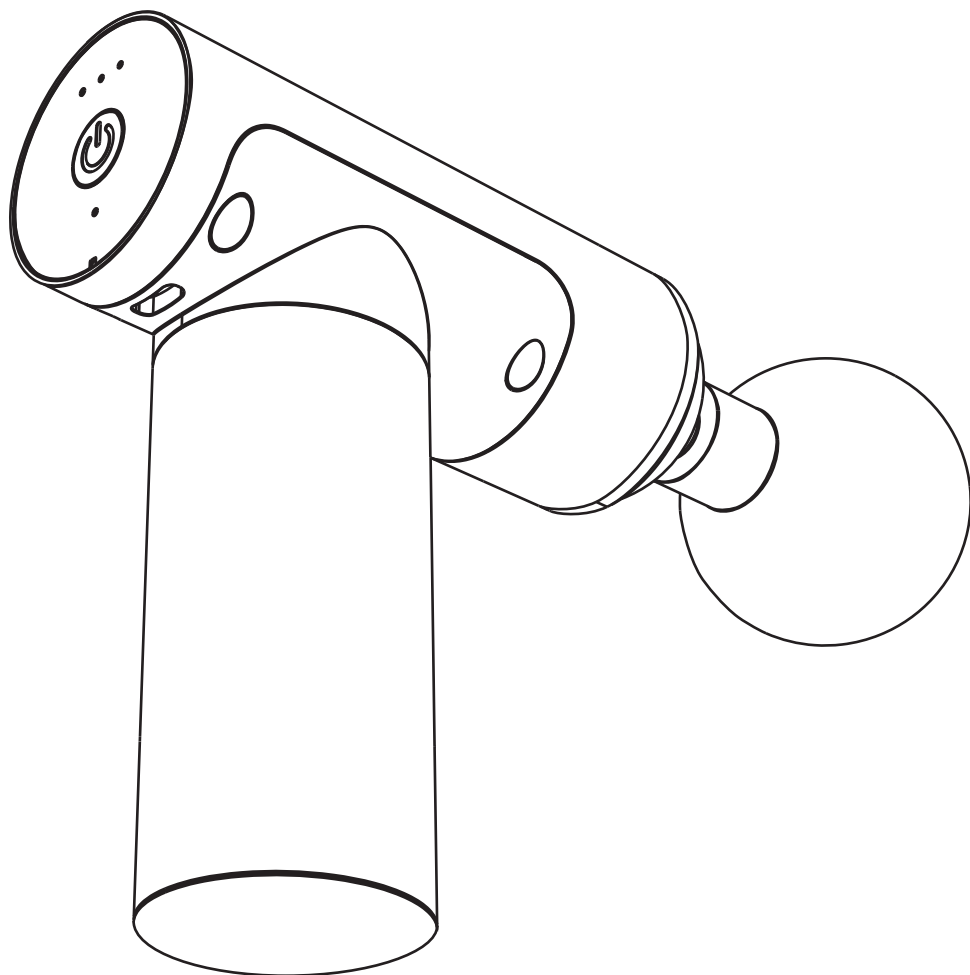
Dear Customer,

Thank you for purchasing the Muscle Massager Please read the instructions before using the massager as this manual provides additional information including safety precautions and warnings.

If you have any difficulties or need more information, you can find our help centre by scanning the QR code opposite with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com

DISCLAIMER

CAUTION

If you have the following conditions or any health issues, please consult with your doctor before use.

1. Pregnancy, diabetes including complications such as neuropathy and retinal damage.
2. Use of cardiac pacemaker, recent surgery, epilepsy or migraines, protrusion of lumbar intervertebral discs and anterior displacement of bone vertebrae.
3. Recent joint replacement surgery, intrauterine device placement, metal implants and other concerns of the body.

WARNING: consult doctor before using this product to reduce the risk of injury.

PRECAUTIONS

1. Minors are not allowed to use this product. This product is not a toy, please store it out of reach of children.
2. Do not immerse the product in water. If the product falls into the water, please submit it for inspection in time.
3. Do not place the product near the fire source or into the fire.
4. Do not use the gun massager on the head, near the heart or on the genitals.
5. Do not use the gun massager on the bones and joints. The massager is designed for use on soft tissues.
6. Do not drop the product intentionally.
7. Do not remove screws or attempt to disassemble the product.
8. It is forbidden to replace or supplement medical treatment with the device, if you feel pain or are receiving treatment, please consult before using the device your doctor.
9. Do not put fingers, hair or other parts of the body near the axis and back of the massage head to avoid crushing injuries.
10. If the equipment is damaged, stop operation immediately.
11. Do not place fingers, hair or other body parts near the axis or back of the massager to prevent injury.

PACKAGE CONTENTS

Gun Massager x 1

Massage Head x 4

Heated Massage Head x 1

Charging Cable x2

Carrying case x 1

Instruction Manual x 1


PRODUCT OVERVIEW

Speed LED
Indicator

On/Off
Speed Button


Charging LED
Indicator

USB-C Charging Port

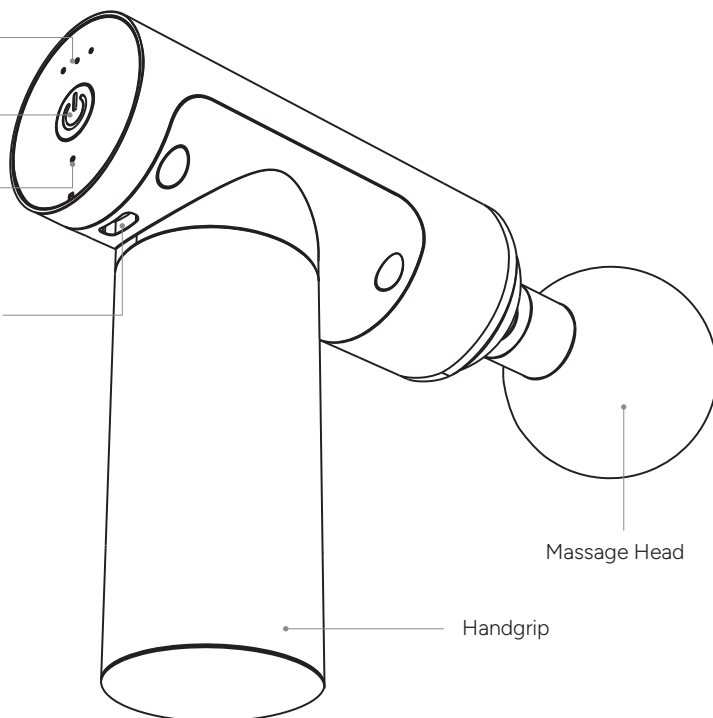Speed LED Indicator Diagram

○	Slow Speed
○ ○	Normal Speed
○ ○ ○	Fast Speed

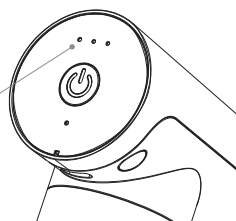CHARGING THE PRODUCT


 For your safety, always use an appropriate 5V input charger/adaptor to charge the massager. DO NOT USE a quick-charging USB adapter as this may damage the product.

Important: The massage gun must not be turned on.

1. Use the USB-C cable included with this massager.
2. Connect the USB-C plug to the USB-C charging port on the underside of the massager.
3. Connect the USB plug to a USB port on a computer or standard wall adapter to begin charging.

NOTE:

- You must select a qualified charger from a reputable manufacturer.
- You must charge the massager if it has not been used for more than 12 months.

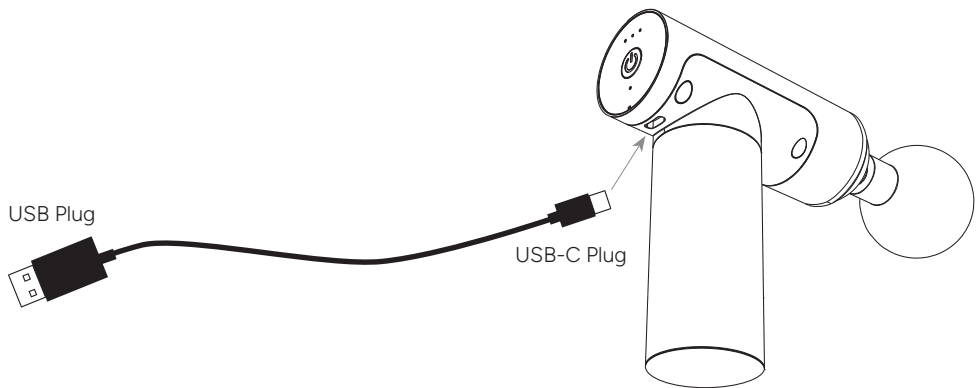 **Safety instructions for charging the product's lithium battery:**

- The input voltage is DC 5V 2A.
- In order to charge the product safely and reliably, please use a wall adapter that is certified and in good condition.
- Only charge the product in ambient temperatures between 10°C /40°F and 40°C /104°F.
- Do not heat above 70°C /158°F. Do not dispose of in a fire.
- Keep out of direct sunlight.
- Turn off the product after use.
- Please note: we will not take any responsibility for the misuse of the product, such as using voltage higher than indicated. Misuse may result in serious injury or loss of property.
- In extreme cases, abuse or misuse of a lithium battery can lead to an explosion, heat generation, fire or smoke development.

MESSAGE HEADS

- **Spherical Head:** suitable for Relaxing Large Muscle Groups, including your arms, waist, back, buttocks, and thighs.

- **Flat Head:** suitable for Relaxation and Shaping of Various Muscle Groups and especially the abdominal belt, the pectoral muscles of the chest or even the biceps and triceps femoris in the leg.

- **Cone Head:** suitable for deep tissue impact recovery sessions. Ideal for targeting trigger points and small specific areas: such as the soles of the feet, palms, Achilles heel and rhomboid muscles in the back.

- **U-Shaped Head:** suitable for trapezius muscles, areas around the achilles tendon muscle, ankle muscle and calf muscles.

- **Heating Massage Head:** for muscles that contract and helps to speed up recovery after physical exertion and relaxes the muscles more quickly.

CONNECTING THE ATTACHMENT

Place the massage head you want on the end of the gun and snap it on. Then turn on the gun to start the massage.

If you wish to use the heated tip:

- Charge the Heating Massage Head beforehand by connecting its included USB-C cable to an electrical outlet. A red light will flash to indicate the tip is charging. When the tip's battery is charged, the red light will appear continuously.

- Turn on the Heating Massage Head using the on/off button on the back of the head.
- Once the Heating Massage Head is turned on, press the on/off button once to access the first temperature level, 45 degrees. Press the on/off button a second time to access the second temperature level i.e. 55 degrees.
- After selecting the temperature, attach the Heating Massage Head to the end of the message gun and turn on the message gun with a long press of the on/off button on the gun.
- Turn off the heating tip using the on/off button on the back of the tip before removing it from the message gun. If you do not turn off the heating tip, it will discharge its battery: it will not turn off automatically.

To change the Massage Head, turn off the gun by holding down the on/off button: you can remove the used head and insert the new one.

USING THE PRODUCT

1. Long press the on/off speed button of the gun massager to start the massager. The massage will start at the slowest speed and will display the slow indicator light.
2. Short press the on/off speed button to change the intensity. One short press will start the normal speed and a second short press will activate the fast speed.
3. The Gun Massager switches off automatically after each 10-minute session, corresponding to the recommended massage time. To restart, simply long press the on/off speed button.
4. The gun massager features three speeds: slow (automatically selected when powered on), normal (first short press), fast (second short press), off (third short press).
5. If you wish to stop the session before the 15-minute timer, simply long press the on/off button until the massager powers off.

FAQ

• Why is my Muscle Massager not turning on?

If the Muscle Massager will not power on, please ensure the battery is fully charged.

• Why is my Muscle Massager is not charging?

If your Muscle Massager will not charge, please ensure each end of the charging cable is securely attached to the Gun Massager and the power source. Check to ensure the wall outlet has power. If using a computer USB port, please ensure that the computer has power.

TECHNICAL SPECIFICATION

Product Name	Muscle Massager
Battery capacity	2000 mAh
Vibration mode	3 levels
Speed	1800 RPM 2500 RPM 3200 RPM
Charging time	3.5 Hours
Charging interface	USB-C
Usage time	2-3 hours
Rated Power	20W ± 10%
Product size	127x153x42 mm

Product Name	Heating Massage Head
Battery capacity	900 mAh
Heating temperature	45°C / 55°C
Charging time	3 hours
Charging interface	USB-C
Usage time	3.5-4 hours (45°C setting) 2.5-3 hours (55°C setting)
Product Size	50x50x46 mm

CARE AND MAINTENANCE

To reduce the risk of failure, electric shock, injury, fire and equipment damage, please read all the instructions and warnings in this manual.

1. Do not dispose of, store or use the massager near or in fire or open flames.
2. Do not remove screws or attempt to disassemble the massager.
3. Please consult with your doctor before use. This massager is not intended to replace or supplement prescribed medical treatments.

WARRANTY

This product is guaranteed against defects in materials and workmanship. This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint, first contact the store where you made your purchase.

ENVIRONMENTAL PROTECTION

When the product is no longer usable, please take it to a collection point for the recycling of waste electrical and electronic equipment.

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Terraillon SAS France & Headquarters
Terraillon Asia Pacific Ltd
Terraillon Corp USA