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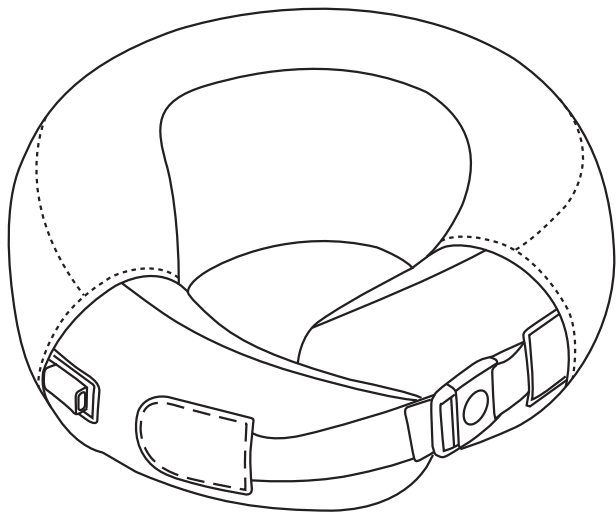
Dear Customer,

Thank you for purchasing the Neck Massager. Please read the instructions before using the massager as this manual provides additional information including safety precautions and warnings.

If you have any difficulties or need more information, you can find our help centre by scanning the QR code opposite with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: [support.uk@terraillon.com](mailto:support.uk@terraillon.com)

## DISCLAIMER

### CAUTION

**If you have the following conditions or any health issues, please consult with your doctor before use.**

1. Diabetes including complications such as neuropathy and retinal damage.
2. Use of cardiac pacemaker, recent surgery, epilepsy or migraines, protrusion of lumbar intervertebral discs and displacement of bone vertebrae.
3. Recent surgery, metal implants and other concerns of the body.

**WARNING: consult doctor before using this product to reduce the risk of injury.**

## PRECAUTIONS


1. Minors are not allowed to use this product. This product is not a toy, please store it out of reach of children.
2. Do not immerse the product in water. If the product falls into the water, please submit it for inspection in time.
3. Do not place the product near the fire source or into the fire.

4. Do not use the massager on the head.
5. Do not use the massager on the bones and joints. The massager is designed for use on neck.
6. Do not drop the product intentionally.
7. Do not remove screws or attempt to disassemble the product.
8. It is forbidden to replace or supplement medical treatment with the device, if you feel pain or are receiving treatment, please consult before using the device your doctor.
9. If the equipment is damaged, stop operation immediately.
10. Do not place fingers, hairs between the neck and the massager.

## PACKAGE CONTENTS

Neck Massager x 1  
Charging Cable x 1  
Instruction Manual x 1

## PRODUCT OVERVIEW

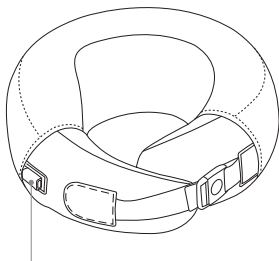USB-C Charging Port

On / Off

Direction of rotation

Intensity levels

Heating function

## CHARGING THE PRODUCT


 For your safety, always use an appropriate 5V input charger/ adapter to charge the massager. **DO NOT USE** a quick-charging USB adapter as this may damage the product.

**Important: The Neck Massager must not be switched on.**

1. Use the USB-C cable included with this massager.
2. Connect the USB-C plug to a USB-C charging port on the side of the massager.
3. Connect the USB plug to a USB port on a computer or standard wall adapter to begin charging.
4. When the massager is charging, the indicator light will be red. When fully charged, the indicator light will be green.

### NOTE:

- You must select a qualified charger from a reputable manufacturer.
- When the battery is not used for a long time, the service life of the battery may be shortened. It is recommended to charge the battery 1 or 2 times per month.

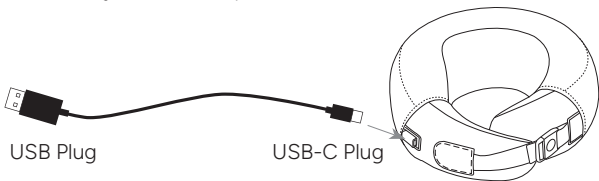
 **Safety instructions for charging the product's lithium battery:**

- The input voltage is DC 5V 2A.
- In order to charge the product safely and reliably, please use a wall adapter that is certified and in good condition.

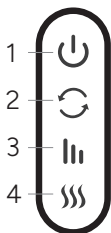- Only charge the product in ambient temperatures between 10°C /40°F and 40°C /104°F.
- Do not heat above 70°C /158°F. Do not dispose of in a fire.
- Keep out of direct sunlight.
- Turn off the product after use.
- Please note: we will not take any responsibility for the misuse of the product, such as using voltage higher than indicated. Misuse may result in serious injury or loss of property.
- In extreme cases, abuse or misuse of a lithium battery can lead to an explosion, heat generation, fire or smoke development.

## USING THE PRODUCT

Place the Neck Massager around your neck with the fastener facing forward, close the strap by snapping the 2 parts of the strap together and sit back to start using the product.

1. Turn on the button , press and hold for 3 seconds to start the machine, and press and hold for 3 seconds to shut down the Neck Massager.
2. If you wish to change the direction of rotation of your massage: short press the button  to change the massage direction between forward and reverse.
3. To adjust the intensity of the massage: short press the button  to change the massage speed between 3 different levels.
4. To benefit from the heating function: short press the button  to start and stop the heating function.

The fabric covering the cushion can be washed by hand or in the washing machine.

## TECHNICAL SPECIFICATION

Product Name	<b>Neck Massager</b>
Battery type and capacity	Li-Ion, 2 000 mAh
Charging interface	USB-C, DC 5V, 2A max
Charging time (empty battery)	About 4 hours
Fully charged battery usage	About 1h30
Dimensions	25 x 27 x 11.7 cm
Weight	~ 600 g

## CARE AND MAINTENANCE

**To reduce the risk of failure, electric shock, injury, fire and equipment damage, please read all the instructions and warnings in this manual.**

1. Do not use the product if the fabric is damaged and/or the rotating heads are visible. This could lead to severe injuries.
2. Do not remove fabric around the foam and/or do not try to disassemble the massager. This could lead to severe injuries.
3. Please consult with your doctor before use. This massager is not intended to replace or supplement prescribed medical treatments.

## WARRANTY

This product is guaranteed against defects in materials and workmanship. During the warranty period, such defects will be serviced free of charge (proof of purchase must be presented for warranty claims). This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint, first contact the store where you made your purchase.

## ENVIRONMENTAL PROTECTION

When the product is no longer usable, please take it to a collection point for the recycling of waste electrical and electronic equipment.

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**Terrailon SAS France & Headquarters**
**Terrailon Asia Pacific Ltd**
**Terrailon Corp USA**