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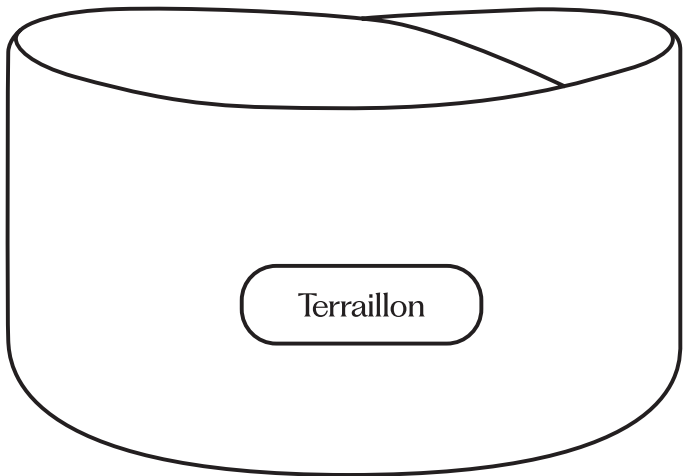
Dear Customer,

*Thank you for purchasing the Waist Massager Please read the instructions before using the massager as this manual provides additional information including safety precautions and warnings.*

If you have any difficulties or need more information, you can find our help centre by scanning the QR code opposite with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: [support.uk@terraillon.com](mailto:support.uk@terraillon.com)

## DISCLAIMER

### CAUTION

**If you have the following conditions or any health issues, please consult with your doctor before use.**

1. Pregnancy, diabetes including complications such as neuropathy and retinal damage.
2. Use of cardiac pacemaker, recent surgery, epilepsy or migraines, protrusion of lumbar intervertebral discs and anterior displacement of bone vertebrae.
3. Metal implants and other concerns of the body.
4. In vivo transplanted medical electronic devices such as pulse regulators.
5. Artificial heart, lungs and other life-sustaining medical electronic devices.
6. Wearable medical electronic devices such as electrocardiographs.
7. Patients with acute disease.
8. Patients who feel abnormalities in the body.

## PRECAUTIONS

1. Minors are not allowed to use this product. This product is not a toy, please store it out of reach of children.
2. The device has a heating function, and people who are not sensitive to heat must pay attention to it. Patients with infectious diseases should use this product alone. If you notice discomfort, discontinue use immediately.
3. This product is a personal care appliance and cannot replace a medical device.
4. Do not immerse the product in water. If the product falls into the water, please submit it for inspection in time.
5. Do not place the product near a heat source or into the fire.
6. Do not drop the product intentionally.
7. Do not attempt to disassemble the product.
8. It is forbidden to replace or supplement medical treatment with the device, if you feel pain or are receiving treatment, please consult before using the device your doctor.
9. If the equipment is damaged, stop operation immediately.
10. Do not use in places such as bathroom.
11. Do not use this product while intoxicated or unwell.
12. Do not fall asleep while the product is being handled.
13. It's forbidden to use product accessories from other suppliers

## PACKAGE CONTENTS


Waist Massager x 1

Instruction Manual x 1

Charging Cable x1

Fabric Extension x1

## PRODUCT OVERVIEW
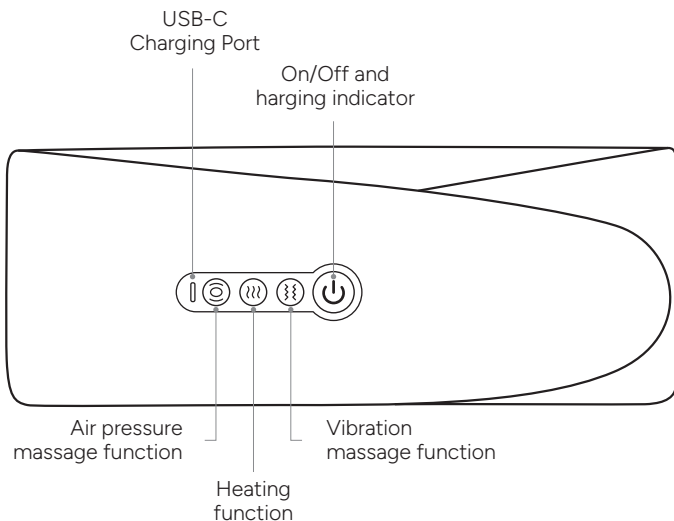
## CHARGING THE PRODUCT

 For your safety, always use an appropriate 5V input charger/ adapter to charge the massager.  
DO NOT USE a quick-charging USB adapter as this may damage the product.


**Important: the Waist Massager must not be turned on while charging.**

1. Use the USB-C cable included with this massager.
2. Plug the USB-C plug into the USB-C charging port located on the control panel.

3. Connect the USB plug to a USB port on a standard wall adapter to begin charging.
4. When the Waist Massager is charging, the indicator light is red. When fully charged, the indicator light will turn green.

**NOTE:**

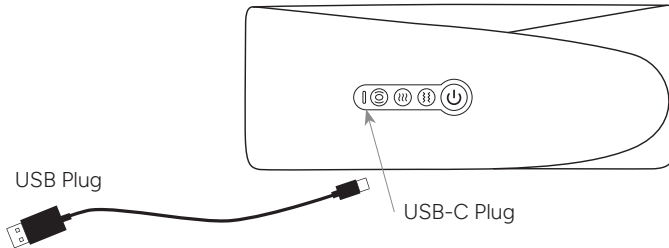- You must select a qualified charger from a reputable manufacturer.
- You must charge the massager if it has not been used for more than 12 months.

**⚠ Safety instructions for charging the product's lithium battery:**

- The input voltage is DC 5V 1A.
- In order to charge the product safely and reliably, please use a wall adapter that is certified and in good condition.
- Only charge the product in ambient temperatures between 10°C /40°F and 40°C /104°F.
- Do not heat above 70°C /158°F. Do not dispose of in a fire.
- Keep out of direct sunlight.
- Turn off the product after use.

- Please note: we will not take any responsibility for the misuse of the product, such as using voltage higher than indicated. Misuse may result in serious injury or loss of property.
- In extreme cases, abuse or misuse of a lithium battery can lead to an explosion, heat generation, fire or smoke development.

## USING THE PRODUCT

1. Press and hold down the button  on the Waist Massager to switch it on.
2. Briefly press the button  to start compression massage mode at the minimum intensity level. A short press starts the medium intensity and a second short press activates the high intensity. A third press stops the massage.
3. Press the button briefly  to start vibration massage at the minimum intensity level. A short press starts the medium intensity and a second short press activates the high intensity. A third press stops the massage.
4. Press the button briefly  to start the heating function at the minimum temperature. You will hear a beep to confirm that the heating function has been activated. A short press activates the medium temperature level (you will hear a beep) and a second short press activates the high temperature level (you will hear a beep). A third press stops the heating function.
5. A programme lasts approximately 15 minutes. The appliance will stop automatically at the end of the programme.

If you wish to stop the session before the end of the programme, press and hold the button  until the massager switches off. You can activate the 2 massage modes at the same time if you wish.

## FAQ

- **Why is my Waist Massager not turning on?**

If the Waist Massager will not power on, please ensure the battery is fully charged.

- **Why is my Waist Massager is not charging?**

If your Waist Massager will not charge, please ensure each end of the charging cable is securely attached to the Waist Massager and the power source. Check to ensure the wall outlet has power. If using a computer USB port, please ensure that the computer has power.

## TECHNICAL SPECIFICATION

Product Name	<b>Waist Massager</b>
Battery capacity	3 000 mAh
Air pressure massage function	3 levels
Vibration massage function	3 levels
Heating function	3 levels: 1: low temperature 2: medium temperature 3: high temperature
Charging time	4h30
Charging interface	USB-C
Usage time	Minimum 7 cycles of 15 minutes with the 3 modes activated (vibration + heat + compression), i.e. 1 hour 45 minutes.
Rated Power	6 W
Product size	118 cm + extension 45 cm

## CARE AND MAINTENANCE

**To reduce the risk of failure, electric shock, injury, fire and equipment damage, please read all the instructions and warnings in this manual.**

1. Do not dispose of, store or use the massager near or in fire or open flames.
2. Do not attempt to disassemble the massager.
3. Please consult with your doctor before use. This massager is not intended to replace or supplement prescribed medical treatments.
4. Only wipe clean with a damp towel, don't use corrosive detergents.

## WARRANTY

This product is guaranteed against defects in materials and workmanship. This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint, first contact the store where you made your purchase.

## ENVIRONMENTAL PROTECTION

When the product is no longer usable, please take it to a collection point for the recycling of waste electrical and electronic equipment.

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**Terrailon SAS France & Headquarters**
**Terrailon Asia Pacific Ltd**
**Terrailon Corp USA**