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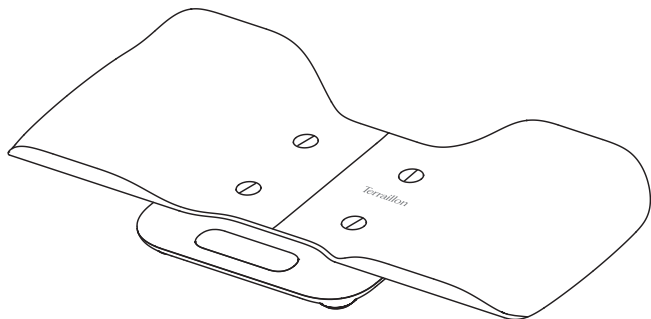
Dear Customer,

Thank you for purchasing *BABY FIT*. Please read the instructions before using the device as this manual provides additional information including safety precautions and warnings.

If you have any difficulties or need more information, you can find our help centre by scanning the QR code opposite with your smartphone. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com


READ BEFORE USE

1. This product is not a toy. Never use the appliance for any purpose other than its intended use. Any other use must be considered dangerous and inappropriate. Do not let children play with it. Use under parental supervision. For indoor use only. **WARNING:** to avoid the risk of suffocation, keep the packaging, screws and batteries out of reach of babies and young children.
2. Clean the scale with a slightly damp cloth. **DO NOT** immerse the scale in water or use chemical/abrasive cleaning agents.
3. All plastics parts should be cleaned immediately after contact with fats, spices, vinegar and strongly flavored / colored foods. Avoid contact with acids citrus juices.
4. Always use the scale on a solid, flat surface.
5. Do not knock, impact and drop products.
6. This product is a precision instrument, please be careful.

7. Keep away from high temperature, store in a cool dry place.
8. Check the batteries if the scale cannot be turned on.

BATTERY INSTALLATION

- Powered by 3 x 1.5V (AAA) batteries (excluded).
- Use a screwdriver to remove the screw from the battery compartment cover, open the battery cover and install the batteries, observing the correct polarity, as shown in the picture.
- Remark: Please use corresponding battery.
- Replace the cover and the screw.
- Please remove battery if it's not used for a long time (more than three months)

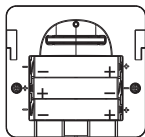
PRODUCT OVERVIEW

Infant weight
locked

Lock screen
mode

Comparison
mark

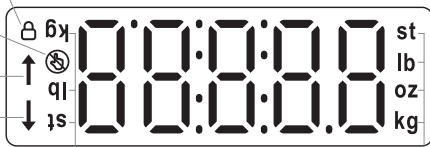Weight unit

Infant weight
unit

Terraillon.


WARNING INDICATORS

Err	OVERLOAD: REMOVE WEIGHT IMMEDIATELY
Lo	BATTERY FAILURE: REPLACE IT
	<p>LOCK SCREEN MODE</p> <p>In this mode, the touch-sensitive keys are disabled during weighing. Press and hold "padlock logo" to activate/deactivate this mode.</p>

BUTTON FUNCTIONS


	<ul style="list-style-type: none"> • Short press to switch ON and Tare function. • Long press to switch OFF.
	<ul style="list-style-type: none"> • Short press to convert units in on status. • In infant weight locked status, short press to clear and lock reading again. • Long press : activates/deactivates screen lock mode. Touchscreen keys are disabled while the child is being weighed.

TRAY ASSEMBLING

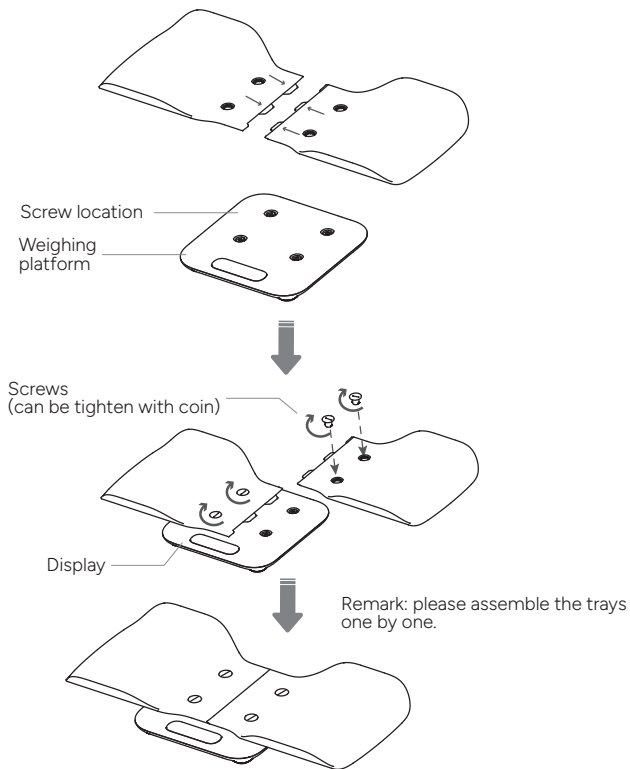REMOVING THE TRAY


Loosen the screws with a coin to remove the tray, then replace the screws on the scale.

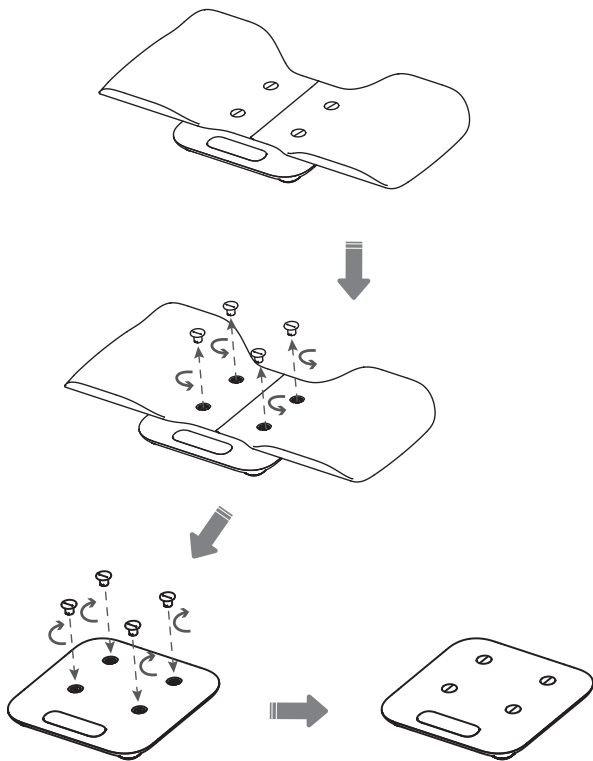Terraillon.

INFANT WEIGHING

0 - 15 kg

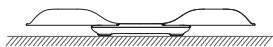
Place scale on flat surface and switch it on.
Infants are weighed using the tray.


Place infant
on scale

Remark : to avoid infant touching the buttons,
please select to lock screen mode before weighing infant.

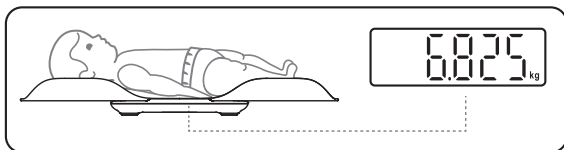Weigh reading is stable and locked
automatically.

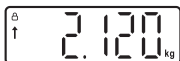If measurement result <25kg, the current weight flashes
three times and the difference in weight compared with the last
record is displayed.

Remark : only weight reading is recorded, no weight difference
for first usage after battery installed

Weighing is finished - Auto off

Terraillon.


CHILD/ADULT WEIGHING

15 - 180 kg

Weighing over 15 kg should be performed without the cradle. In this case, operation is as for a conventional scale. Note that the display rotates 180° in relation to the baby scale mode.

Place scale on flat surface.

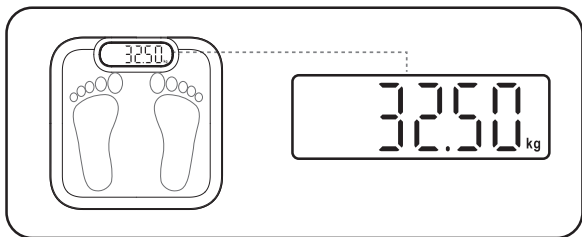Remark: when scale is placed in a new place, press firmly and briefly on the platform to switch on the scale, wait for it to display 0.0 and then switch off. You can now weigh yourself.

Weighing finished.

If no operation or 10 seconds auto off.

GUARANTEE

This product is guaranteed against defects in materials or manufacturing. During the guarantee period, any such defects will be replaced free of charge (proof of purchase must be shown in the case of a complaint under guarantee). The guarantee does not cover damage resulting from improper use or negligence. In the event of a claim, first contact the shop where you purchased the product.

PROTECT THE ENVIRONMENT

Dispose of used batteries in a designed container so that they can be collected and recycled. Do not mix different types of batteries. Do not mix new and used batteries. After life, dispose of this product at a designated waste recycling point.

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Terrailon SAS France & Headquarters
Terrailon Asia Pacific Ltd
Terrailon Corp USA