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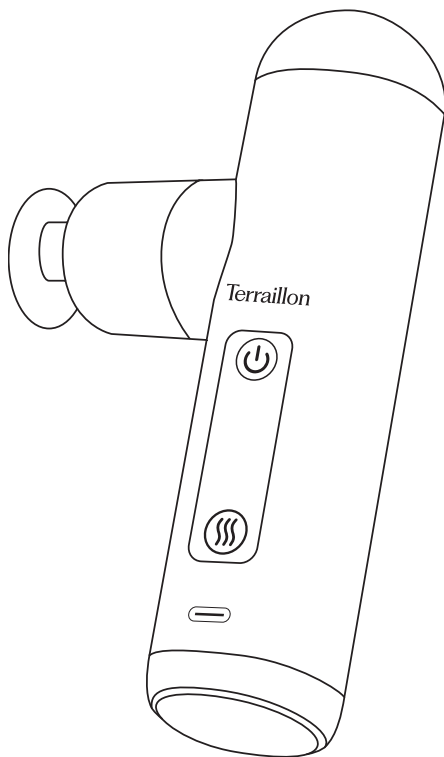
Dear Customer,

Thank you for purchasing the Mini Massager. Please read the instructions before using the massager as this manual provides additional information including safety precautions and warnings.

If you have any difficulties or need more information, you can find our help centre by scanning the QR code. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com

DISCLAIMER

CAUTION

If you have the following conditions or any health issues, please consult with your doctor before use.

1. Pregnancy, diabetes including complications such as neuropathy and retinal damage.
2. Use of cardiac pacemaker, recent surgery, epilepsy or migraines, protrusion of lumbar intervertebral discs and anterior displacement of bone vertebrae.
3. Recent joint replacement surgery, intrauterine device placement, metal implants and other concerns of the body.


WARNING: consult doctor before using this product to reduce the risk of injury.


PRECAUTIONS

1. Minors are not allowed to use this product. This product is not a toy, please store it out of reach of children.
2. Do not immerse the product in water. If the product falls into the water, please submit it for inspection in time.
3. Do not place the product near the fire source or into the fire.
4. Do not use the gun massager on the head, near the heart or on the genitals.
5. Do not use the gun massager on the bones and joints. The massager is designed for use on soft tissues.
6. Do not drop the product intentionally.
7. Do not remove screws or attempt to disassemble the product.
8. It is forbidden to replace or supplement medical treatment with the device, if you feel pain or are receiving treatment, please consult before using the device your doctor.
9. Do not put fingers, hair or other parts of the body near the axis and back of the massage head to avoid crushing injuries.
10. If the equipment is damaged, stop operation immediately.
11. Do not place fingers, hair or other body parts near the axis or back of the massager to prevent injury.

PACKAGE CONTENTS

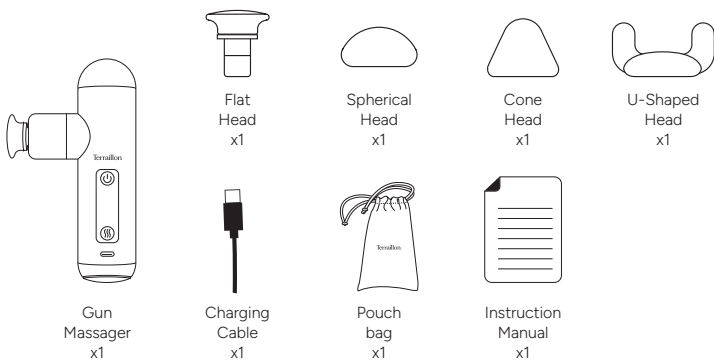PRODUCT OVERVIEW

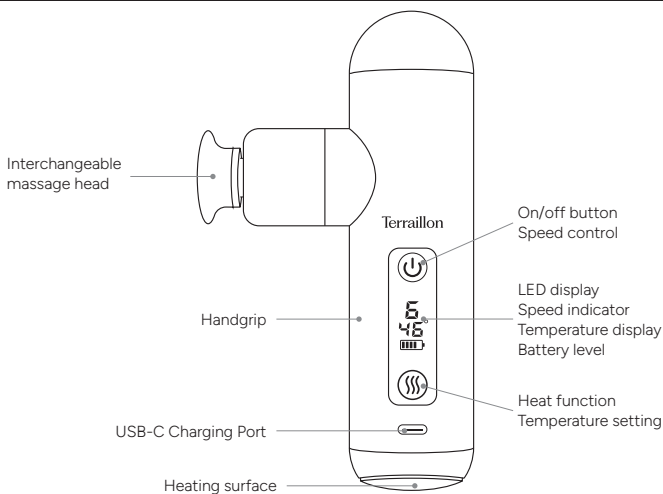Terraillon.

MESSAGE HEADS

- **Flat Head:** suitable for Relaxation and Shaping of Various Muscle Groups and especially the abdominal belt, the pectoral muscles of the chest or even the biceps and triceps femoris in the leg.

- **Spherical Head:** suitable for Relaxing Large Muscle Groups, including your arms, waist, back, buttocks, and thighs.

- **Cone Head:** suitable for deep tissue impact recovery sessions. Ideal for targeting trigger points and small specific areas: such as the soles of the feet, palms, Achilles heel and rhomboid muscles in the back.

- **U-Shaped Head:** suitable for trapezius muscles, areas around the achilles tendon muscle, ankle muscle and calf muscles.

For more specific programme tips, you can find an article detailing them in our help centre by scanning the QR code.


CHARGING THE PRODUCT

IMPORTANT: you must recharge the massager before first use and if it has not been used for more than 12 months. To recharge, the massage gun must be switched off.

1. Use the USB-C cable included with this massager.
2. Plug the USB-C plug into the USB-C charging port.
3. Connect the USB plug to a mains adapter.
4. During charging, the indicator flashes to show the current charge level.
When the battery is fully charged, the battery indicator is full and steady.
Charging time is approximately 3 hours.

Battery life is a maximum of 5 hours. This time will vary depending on the programmes and usage

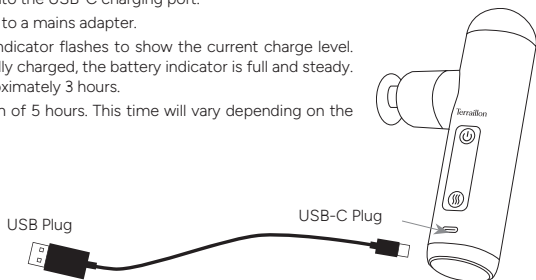CHECKING THE BATTERY LEVEL:

- The battery charge level is visible as soon as the product is switched on.

- The battery is discharged: a red light indicates that the battery needs recharging

- The battery is charging: the battery level flashes on the LED screen as charging progresses.

- The battery is fully charged: the battery symbol is complete.

**Safety instructions for charging and backing up the product's lithium battery:**

- For your safety, always use a 5 V input charger/adaptor (the Mini Massager's input voltage is 5 V DC 1 A).
- Select a qualified charger from a reputable manufacturer. It must be certified and in good condition.
- DO NOT use a fast charging USB adaptor as this may damage the product.
- Only charge the product at ambient temperatures between 10°C/40°F and 40°C/104°F.
- Do not heat the device above 70°C/158°F. Do not throw into a fire.
- Keep out of direct sunlight.
- Switch off the product after use.

PLEASE NOTE:

We will not take any responsibility for the misuse of the product, such as using voltage higher than indicated.

Misuse may result in serious injury or loss of property.


In extreme cases, abuse or misuse of a lithium battery can lead to an explosion, heat generation, fire or smoke development.

USING THE PRODUCT

SETTING-UP THE MESSAGE HEADS:

1. Place the flat message head at the end of the gun and clip it on.
2. The silicone tips: spherical head, conical head, U-shaped head clip directly onto the flat head.

To change the message head, the gun must be switched off.

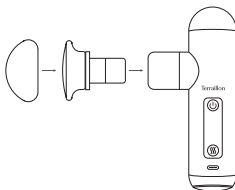**ON/OFF:**

Then switch on the gun to start the massage.

1. Press the massage gun's on/off button (⏻) for **3 seconds** to start the massage. Massage starts at the slowest speed.

The massage gun switches off automatically after every 15-minute session. If you wish to continue your massage, you can restart the unit.

2. If you wish to stop the session before the 15-minute timer expires, press the On/Off button (⏻) again for **3 seconds**.

SPEED SELECTION:

The Mini Massager has 6 massage speeds.

Briefly press the On/Off speed button  to change the intensity.

Each press increases the speed by 1.

The intensity selected will then be displayed on the unit's LED screen.

USING THE HEATING SURFACE:

To use the heating surface, the Mini Massager must be ON.

1. Press the heating surface button  for 3 seconds to start the programme. The heating surface takes approximately 1 minute to reach the temperature of the selected programme.

The Mini Massager has 3 heating programmes:

- 1: Temperature at 38°C - LED display + green light
- 2: Temperature at 42°C - LED display + blue light
- 3: Temperature at 46°C - LED display + red light

Briefly press the heating surface button  to change programme.

The heating function stops after 15 minutes of operation (the timer starts when the appliance is switched on).

2. If you wish to stop the heating function during use, press the heating surface button  for 3 seconds.

The appliance will also stop.

FAQ

- Why is my Mini Massager not turning on?

If the Mini Massager will not power on, please ensure the battery is fully charged.

- Why is my Mini Massager is not charging?

If your Mini Massager will not charge, please ensure each end of the charging cable is securely attached to the gun massager and the power source. Check to ensure the wall outlet has power. If using a computer USB port, please ensure that the computer has power.

TECHNICAL SPECIFICATION

Product Name	Mini Massager	Charging time	3 hours
Battery capacity	1800 mAh	Charging interface	USB-C
Vibration mode	6 levels	Usage time	3-5 hours
Speed	800 RPM	Rated Power	13,3 W
	1 300 RPM	Product size	90 x 40 x 155 mm
	1 800 RPM	Product weight	400 g
	2 300 RPM	Decibels	45 - 50 dB
	2 800 RPM	Excursion	6,5 mm
	3 200 RPM		
Temperature of heating programmes	38°C, 42°C, 46°C		

CARE AND MAINTENANCE

To reduce the risk of failure, electric shock, injury, fire and equipment damage, please read all the instructions and warnings in this manual.

1. Do not dispose of, store or use the massager near or in fire or open flames.
2. Do not remove screws or attempt to disassemble the massager.
3. Please consult with your doctor before use. This massager is not intended to replace or supplement prescribed medical treatments.

WARRANTY

This product is guaranteed against defects in materials and workmanship. This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint, first contact the store where you made your purchase.

ENVIRONMENTAL PROTECTION

When the product is no longer usable, please take it to a collection point for the recycling of waste electrical and electronic equipment.

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Terraillon SAS France & Headquarters
Terraillon Asia Pacific Ltd
Terraillon Corp USA