Notice d'utilisation
User guide
Bedienungsanleitung
Gebruiksaanwijzing

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Your product

Box contents

- 1 water flosser
- 1 charging station
- 2 tips
- 1 instruction leaflet

Technical characteristics

- 3 modes: 1) normal 2) soft 3) pulse
- On/Off button
- Mode selection button
- Charge & low battery indicator
- 1400mAh Li-ion battery
- Charge time: about 24 hours
- Battery life: about 45 minutes
- Removable and transparent 170 ml tank
- Power cord length: 1.2 m

Product Information for power consumption and maximum time to reach applicable low power mode	
Power consumption in off mode	NA
Power consumption in standby	NA
The maximum time needed to automatically reach the applicable low power mode or condition	NA
Power consumption in network standby if all wired network ports are connected and all wireless network ports are activated	NA

Battery maintenance mode: 0,4 W



The period after which the equipment reaches automatically battery maintenance mode: 0 minute

Before first use

- The water flosser must be charged before first use. Refer to the next section entitled "Charging".

Charging


1. Make sure that the charging base is placed on a flat, stable, and dry surface.
2. Place the water flosser on its charging base.
3. Plug the power cord into a mains socket. The charging indicator lamp will flash.
4. When charging is complete (about 24 hours), the indicator lamp will remain steadily lit. Unplug the power cord.

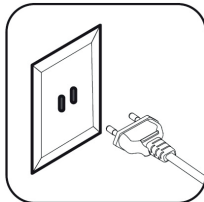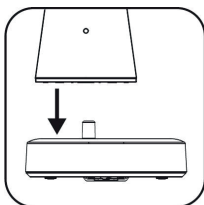Note: if, during or after use, the indicator lamp flashes, the water flosser needs to be charged.


Use

The water flosser is suitable for cleaning dental appliances and prostheses.

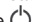1. Make sure the head of the tip is pointed towards the side with the  on/off button, and push it onto the handle until it clicks into place.


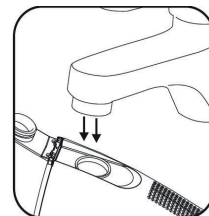
2. Open the cover of the tank and fill it with cold or warm water (max. 40°C).


EN

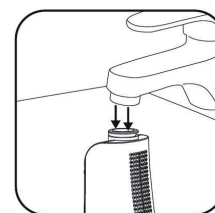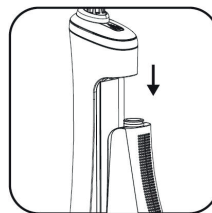You can also take the tank out of the handle to make it easier to fill, and then reassemble it on the water flosser.

3. Press the mode button  repeatedly to select the desired spray type:
 - Normal : for normal cleaning to eliminate bacteria that cause cavities and plaque.
 - Soft : to massage the gums and stimulate blood circulation. Ideal for sensitive gums.
 - Pulse : to effectively remove all food residue and clean interdental spaces.
4. Turn on the device with the on/off button  and start cleaning your teeth. Keep the water flosser in an upright position with the tip head facing the teeth.


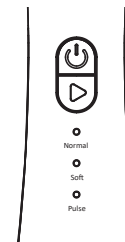
After use

1. The stream will automatically stop after **2 minutes**.
2. Switch off the device with the on/off button .
3. Open the cover of the water tank and empty the remaining water into a sink.
4. Close the cover, turn on the device again and drain the remaining water from the suction tube, then turn off the device.
5. Wipe the water flosser with a dry cloth and press the unlock button on the tip to remove it if you wish to replace it.

Note: for hygiene reasons and to avoid bacterial growth, make sure that the tank is completely empty after each use. If you do not intend to use the device for a long period of time, allow the inside of the tank to dry and leave the lid open.

Tips for use

- The water must be sprayed at a 90° angle to the teeth.
- Slowly move the water flosser along your teeth and gums.
- Pay particular attention to the area between the teeth and, if necessary, around dental appliances and prostheses.

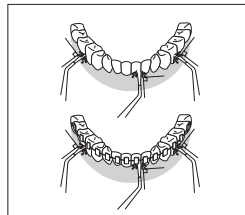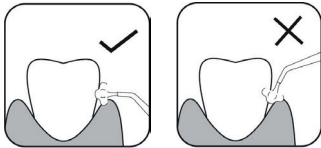Cleaning and maintenance

- Do not immerse the charging base or water flosser in water or any other liquid.
- Before cleaning, the device must be switched off.
- Never use alcohol, solvents or abrasives to clean your device. Only clean water or a neutral cleaning agent suitable for this device.

EN

Tip

- Press the unlock button on the tip to remove it, and clean it with clean water. Then wipe it with a soft, dry cloth.
- When used regularly, the tip should be replaced every 6 months.

Body of the water flosser

- Clean the water flosser's surfaces with a soft, damp cloth with liquid detergent and dry it thoroughly.

Suction tube

- Clean the suction tube with clean water and wipe it with a soft, dry cloth.
- Do not bend, pull, or twist the tube, as this may damage it.

Water tank

- Rinse the tank with clean water and let it dry completely and leave the lid open (especially if you do not plan to use the water flosser for a long time).

Charging base

- Clean the surface of the charging base with a soft, dry cloth.

Troubleshooting guide

The water flosser isn't working:

- You're using the water flosser for the first time or it hasn't been used for more than three months. Charge the water flosser for at least 24 hours.

The stream pressure is low:

- The tip is blocked or deformed. Replace it.

The water flosser doesn't work or stops working after a few minutes:

- The filter at the end of the suction hose is obstructed: clean it.
- The battery is flat or wasn't charged long enough. Charge the device for at least 24 hours before using it.
- If it still doesn't work after charging, contact your retail store's customer service department.

Water is not coming out of the tip:

- The tank is empty. Fill it.
- The water flosser is tilted too much during use. Place it in an upright position.

If the problems persist, contact your retailer's after-sales department.

Recycling the battery

This device has a built-in rechargeable battery that can be recycled but not replaceable. At the device's end of life, return it to your retail store or place it in a dedicated container for this type of waste.

EN

All information, designs, drawings and pictures in this document are the property of SOURCING & CREATION. SOURCING & CREATION reserves all rights to its brands, designs and information. Any copy and reproduction through any means shall be deemed and considered as counterfeiting.

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Protection of the environment

This symbol attached to the product means that it is an appliance whose disposal is subject to the directive on waste from electrical and electronic equipment (WEEE). This appliance may not in any way be treated as household waste and must be subject to a specific type of removal for this type of waste. Recycling and recovery systems are available in your area (waste removal) and by distributors. By taking your appliance at its end of life to a recycling facility, you will contribute to environmental conservation and prevent any harm to your health.

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* Tested in our laboratories

Warranty valid from the date of purchase (receipt as proof of purchase). This warranty does not cover defects or damage caused by improper set up, incorrect use, or normal wear and tear of this product.

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Art. 8009436

Ref. EJD2

SERVICE RELATION
CLIENTS / Customer
Relations Department /
Kundenservice /
Klantenafdeling

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