

essentiel**b**

Pistolet de massage / Massage gun /
Pistola de masaje / Massagepistole /
Massage gun
HyperB

Notice d'utilisation / User guide / Instrucciones de uso /
Bedienungsanleitung / Gebruiksaanwijzing

description de l'appareil

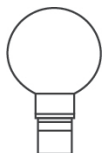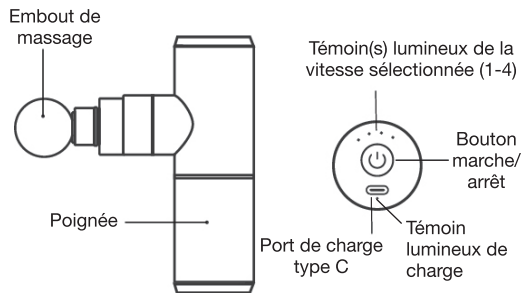Boule

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
appliance description

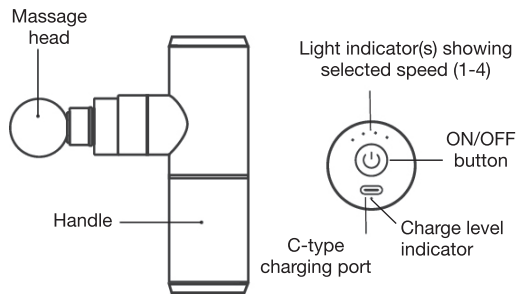Ball

For the abdomen, buttocks, legs and other muscle groups (thighs, calves, etc.)

Fork

For the neck, shoulders and Achilles tendon. **Avoid direct contact with the spine.**

Pointed

For less sensitive areas such as the soles of the feet or the palms of the hands.

Flat

Suitable for the relaxation of all muscles.

your product

Contents of the box

- 1 massage gun
- 1 AC adapter
- 1 USB-C cable
- 4 massage heads
- 1 carry bag
- 1 user guide

Technical features

- Massage by “percussion”
- 4 speeds
- 3.7V 2000mAh Li-ion battery
- Charging time: 300 minutes (5 hours)
- Battery life: 40 cycles of 10 minutes
- Low battery indicator
- Weight: 426g
- Dimensions (L x D x H): 14cm x 9cm x 4.5cm


User benefits

- Allows for the alleviation of muscle fatigue, muscle aches, sprains and strains.
- Relaxing or stimulating effect.
- Helps to reduce the pain related to sore muscles due to the build-up of lactic acid in the muscles after intense physical exertion.
- Relieves the discomfort and pain caused by repetitive injuries or a sedentary lifestyle.

Examples of use

1. I have a competition / marathon goal with a training plan in place. Do you have any products that could help me to prepare my body and my muscles?
2. I have just resumed sports after several years without any sports activity. The first sessions were painful, I had sore muscles, do you have a product to relieve the pain?
3. After a working day seated at my desk, in the evenings I have back pain and neck pain, how can I relieve this and relax for the evening? Watch programmes on applications with the keywords: “sore muscles, work out muscle knots, etc.”
4. If I don't do a lot of sport, is the product still suitable for me? Yes, absolutely, the product is as suitable for professional sportsmen and women as it is for sports enthusiasts doing running, cross fit, pilates, and other sports. The product is also suitable for daily use in order to relieve back pain or even as part of physiotherapy.

charging the battery

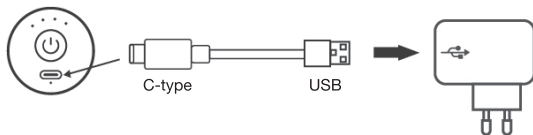1. Please charge the battery fully before first use. The **charging time** is approximately **300 minutes (5 hours)**. The gun massage can only be recharged to the voltage indicated on the charger CW0502000EU supplied.
2. Plug in the end of the power cord supplied to the appliance's charging port and the other end into the mains power socket. Then plug the AC adapter into a wall socket.
 - Red charge light: the appliance is busy charging.
 - Green charge light: the appliance is fully charged.

3. To prolong the battery life, please charge the appliance as soon as the red light starts flashing.
 - Full charge offers a **battery operating time** of approximately **40 cycles of 10 minutes**. The usage time varies according to the selected speed.
4. During use, the charge light flashes in various ways as follows depending on the battery charge level:
 - More than 70% battery: the green light flashes 3 times every 5 seconds.


- Less than 70% battery: the green light flashes 2 times every 5 seconds.
- Less than 30% battery: the green light flashes 1 time every 5 seconds.
- Less than 10% battery remaining: the red light flashes.

EN

Note: Do not use this appliance while it is charging.

insertion/removal of a head

1. Make sure the appliance is OFF.
2. To insert a massage head, slide it gently into its slot in the body of the appliance.
3. To remove a massage head, pull it gently towards you.

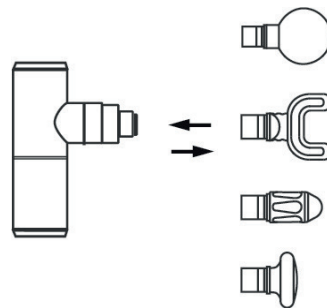use

1. **Before use, ensure a massage head is attached correctly and firmly to the appliance.**
2. Press the ON/OFF button  for 3 seconds to switch on the appliance.
 - **Always switch on the appliance before it enters into contact with the skin.**
3. The appliance has 4 massage speeds.
 - Always start with Speed 1.
 - You can then choose a massage speed according to your preference and the type of muscle.
 - Each short press on the ON/OFF button  increases the speed as follows: SPEED 1 > SPEED 2 > SPEED 3 > SPEED 4 > SPEED 1, etc.
 - The indicator lights switch on according to the selected speed:

Light indicator(s) showing selected speed

- Speed 1: 1 indicator light comes on
 - Speed 2: 2 indicator lights come on
 - Speed 3: 3 indicator lights come on
 - Speed 4: 4 indicator lights come on
4. Place the appliance in contact with the body and move it around gently. It is not necessary to apply additional pressure. Move the appliance over an area of approximately 3cm per second. Doing it more gently greatly enhances the massage effect.



- When you treat particularly sensitive areas or “knots”, work gently without pressing too hard.
 - For better results, treat an area for approximately 30 seconds, breathe in deeply, and then slowly move the appliance around the rest of the area to be treated.
5. The device automatically stops after **10 minutes**, which is the recommended usage time to treat a particular area.
 - It is recommended to wait for 30 to 40 minutes at least between two usage sessions.
 6. You can stop the appliance at any time by pressing on the ON/OFF button  for 3 seconds.

cleaning and storage

When necessary, you can clean the appliance and its accessories with a slightly damp, soft cloth.

You can store the appliance and its accessories in the carry bag supplied.

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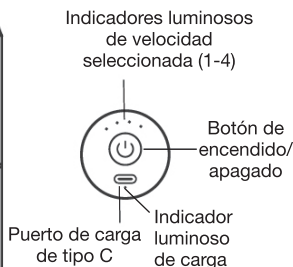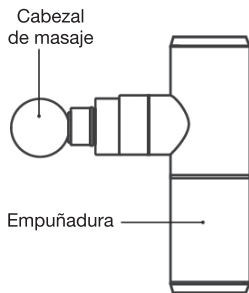Protection of the environment

This symbol attached to the product means that it is an appliance whose disposal is subject to the directive on waste from electrical and electronic equipment (WEEE). This appliance may not in any way be treated as household waste and must be subject to a specific type of removal for this type of waste. Recycling and recovery systems are available in your area (waste removal) and by distributors. By taking your appliance at its end of life to a recycling facility, you will contribute to environmental conservation and prevent any harm to your health.

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All information, designs, drawings and pictures in this document are the property of SOURCING & CREATION. SOURCING & CREATION reserves all rights to its brands, designs and information. Any copy and reproduction through any means shall be deemed and considered as counterfeiting.

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Tested in our laboratories

Warranty valid from the date of purchase (receipt as proof of purchase). This warranty does not cover defects or damage caused by improper set up, incorrect use, or normal wear and tear of this product.

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