
essentiel

Combiné dentaire : Brosse à dents électrique et
Hydropulseur / Dental set: Electric toothbrush and Water
flosser / Combinado dental: Cepillo de dientes eléctrico
e hidropulsador / Zahnpflegekombination Elektrische
Zahnbürste und Munddusche / Mondverzorgingsset:
Elektrische tandenborstel en Monddouche

B-HYDRO BRUSH

Notice d'utilisation / User guide / Instrucciones de uso /
Bedienungsanleitung / Gebruiksaanwijzing

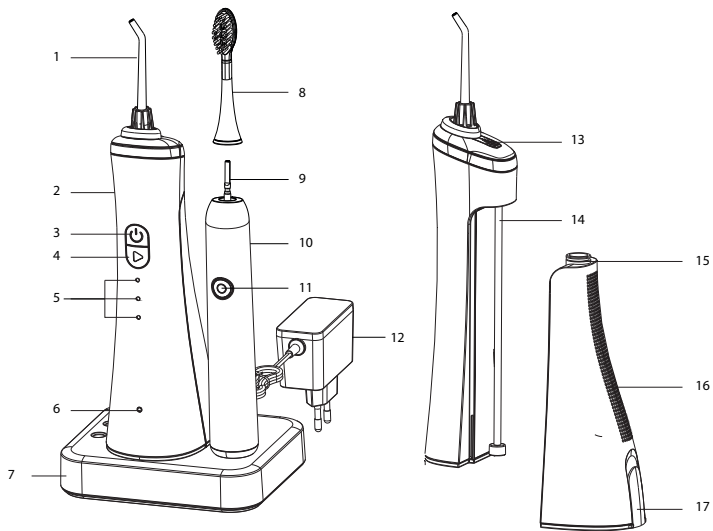Le produit ne doit être utilisé qu'avec le chargeur fourni portant la référence ZD006C050100EUE.
 The product should only be used with the charger supplied (reference ZD006C050100EUE).
 El producto solo se debe usar con el cargador incluido, que lleva la referencia ZD006C050100EUE.
 Das Gerät darf nur mit dem mitgelieferten Ladegerät mit der Bezeichnung ZD006C050100EUE verwendet werden.
 Het product mag alleen gebruikt worden met de meegeleverde lader met de referentie ZD006C050100EUE.

	FR	EN	NL	DE	ES
1	Canule	Tip	Canule	Kanüle	Cánula
2	Hydropulseur	Water flosser	Monddouche	Munddusche	Irrigador dental
3	Bouton marche/arrêt	On/Off button	Aan/uit-knop	Ein/Aus Knopf	Botón de encendido y apagado
4	Bouton de sélection du mode	Mode selection button	Knop voor het selecteren van de modus	Modus-Auswahlknopf	Botón de selección del modo
5	Indicateur lumineux du mode	Mode indicator lamp	Indicatorlampje van de modus	Leuchtanzeige für den Modus	Indicador luminoso del modo
6	Indicateur lumineux de charge	Charge indicator lamp	Oplaadindicatie-lampje	Leuchtanzeige für das Laden	Indicador luminoso de carga
7	Socle de charge	Charging base	Oplaadsokkel	Ladesockel	Base de carga
8	Brossette	Brush heads	Cepillo	Aufsteckbürste	Opzetborstel
9	Moteur	Motor	Motor	Motor	Motor
10	Brosse à dents	Electric Toothbrush	Cepillo de dientes	Zahnbürste	Tandenborstel
11	Bouton marche/arrêt et de sélection du mode	On/Off and Mode selection button	Botón de encendido/apagado y selector de modo	Ein/Aus-Knopf und Modusauswahl	Aan-/uitknop en standkeuzeknop
12	Adaptateur secteur	Mains adapter	Adaptador de red eléctrica	Netzadapter	
13	Bouton de déverrouillage de la canule	Tip unlock button	Ontgrendelingsknop voor de canule	Entriegelungsknopf der Kanüle	Botón de desbloqueo de la cánula
14	Tuyau d'aspiration	Suction tube	Aanzuigslang	Saugschlauch	Tubo de aspiración
	Filtre	Filter	Filter	Filter	Filtro
15	Joint étanche	Waterproof ring	Waterdichte pakking	Dichtungsring	Junta estanca
16	Réservoir d'eau	Water tank	Waterreservoir	Wassertank	Depósito de agua
17	Couvercle du réservoir d'eau	Water tank door	Deksel van het waterreservoir	Deckel des Wassertanks	Tapa del depósito de agua

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your product


Contents of the box

- 1 dental set: 1 electric toothbrush and 1 water flosser
- 1 charging station and its mains adapter
- 2 nozzles
- 2 brush heads
- 1 user guide

Technical features

- 2 brushing modes: Clean and Sensitive
- 3 water jet modes: Normal, Soft, Pulse
- Low battery indicator
- Weight of toothbrush with brush head: 110 g
- Weight of the filled water flosser with nozzle: 460 g
- Water jet timer: 2 minutes
- **Mains charger**
 - » Model: ZD006C050100EUE
 - » Input: 100-240 V~ 50/60 Hz, 0,3A
 - » Output: 5.0V $\equiv \equiv \equiv$ 1.0A; 5.0W
 - » IP44
- **Charging base**
 - » Input 5.0V $\equiv \equiv \equiv$ 1.0A; 5.0W
- **Water flosser**
 - » 3.7V $\equiv \equiv \equiv$ 1400mAh
 - » IPX7
 - » Clean water tank capacity: 170ml
- **Electric toothbrush**
 - » 3.7V $\equiv \equiv \equiv$ 800mAh
 - » Lithium Ion
 - » IPX7
 - » 31000 - 41000 vibrations per minute

installation

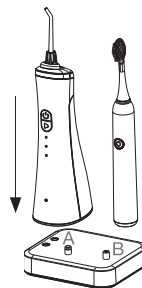1. Make sure the charging base is placed on a flat, stable and dry surface.
2. Place the water flosser and the toothbrush on the charging base, in the correct indicated position (A for the water flosser, B for the toothbrush).
3. Plug in one end of the power supply to the back of the charging station and the other end in to a mains socket.
4. The lights flash on both devices (red on the water flosser, white on the toothbrush).
5. When the devices are fully charged, the lights remain on without flashing. You can now unplug the power supply.

Note: if, during or after use the lights start flashing again, you need to recharge the water flosser or the toothbrush.


using the water Flosser

The water flosser is to be cleaned in the same way that is suitable for other dental devices and prostheses.

1. Ensure that the nozzle head is positioned on the same side as the On/Off button , and press it into the handle until it clicks locked into place.
2. Open the water tank lid and fill it with cold or warm water (40°C max.).


You can also remove the tank from the handle to make it easier to fill, and then fit it back onto the water flosser.

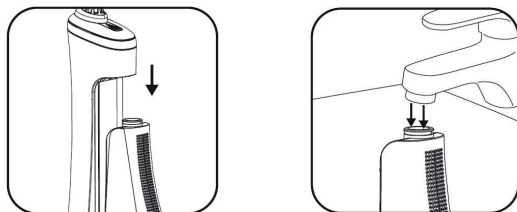3. Press the button mode  in succession to select the water jet type desired:

- Normal: allowing for a full brushing cycle to be completed.
- Soft: soft, gentle brushing. To massage the gums and stimulate blood circulation. Ideal for sensitive gums.
- Pulse: to eliminate any food residue between teeth more effectively.

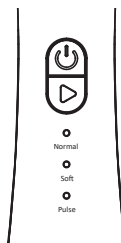
4. Turn the device on using the On/Off button  and start to clean your teeth. Hold the water flosser in an upright position, with the nozzle facing the teeth.

after use




1. The water jet will switch off automatically after **2 minutes**.
2. Turn the device off using the On/Off button .
3. Open the water tank lid and empty out the remaining water into the washbasin.
4. Close the lid again, switch on the device again and empty the remaining water into the suction hose, and then turn off the device.
5. Wipe the water flosser with a dry cloth and press the nozzle's unlock button to remove it if you want to replace it.

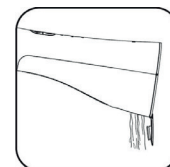Note: for hygiene reasons and to prevent the growth of bacteria, check that the tank is completely empty after each use. If you will not be using the device for a long time, leave the inside of the tank to dry out by leaving the lid off.

user tips

- The water needs to be sprayed at a 90° angle to the teeth.
- Slowly move the water flosser along all of your teeth and gums.
- Concentrate slightly longer on the spaces between the teeth and, where applicable, around dental devices and prostheses.

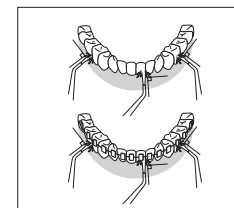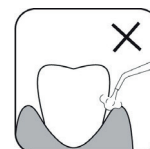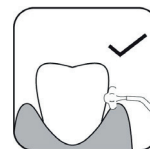using the electric toothbrush

1. Install the brush head on the body of the appliance.
2. Moisten the brush head by running it under tap water.
3. Apply a small amount of toothpaste to the brush head.
4. Place the toothbrush in your mouth and press the On/Off button . The toothbrush switches on in “Normal” mode. In “Normal” mode, the indicator light flashes light blue.

NOTE: In the case of using the electric toothbrush or water flosser after a long period of time, we advise you to gradually allow your gums to adjust to their use again in order to avoid any temporary irritation that may occur.

TIP: Do not press the toothbrush too hard on your teeth. Do not brush like a manual toothbrush: the toothbrush does it for you. Gently apply the electric toothbrush against your teeth.

EN

Your toothbrush has two modes of operation:

- “Clean” (Standard) mode: this is the standard mode for daily tooth cleaning.
- “Soft” (Sensitive) mode: mode suitable for sensitive teeth and gums.

1. To shift into “Sensitive” (Soft) mode, press the On/Off button  a second time. In “Soft” mode, the indicator light flashes dark blue.
2. To switch off the toothbrush when you have finished using it, press the On/Off button  a third time.
3. Remove the toothbrush from your mouth and then clean the brush head and toothbrush under tap water.

Smart reminder

Every 30 seconds, the toothbrush briefly stops to encourage you to change the brushing area (4 areas: top front/back and bottom front/back).

After two minutes, the toothbrush stops to indicate that you have reached the recommended brushing time.

cleaning and maintenance

- Do not immerse the charging base, the toothbrush or the water flosser in water or any other liquid.
- Always unplug the charging base from the mains before cleaning the devices.
- Do not use alcohol, solvent or abrasive products to clean the devices. Use only clean water or a neutral cleaning product specifically suitable for use on this device.

The water flosser's nozzle

- Press the release button of the nozzle to remove it, and clean it with clean water. Wipe it with a soft, dry cloth.
- During periods of regular use, the nozzle must be replaced every 6 months.

Cleaning the body of the water flosser

- Clean the surfaces of the water flosser using a soft, damp cloth with some liquid detergent, then dry completely.

The water flosser's suction hose


- Clean the suction hose with clean water and wipe it with a soft, dry cloth.
- Do not bend the hose, do not pull it or twist it, to prevent damaging it.

The water flosser's water tank

- Rinse the tank with clean water, then dry completely by leaving the lid open (especially when you don't plan to use the water flosser for a while).

Cleaning the brush head of the toothbrush

- Rinse the brush head under tap water after use. Squeeze the water out of the brush head using a towel and let it dry in the open air.
- To protect the brush head, place the protective cap on it.

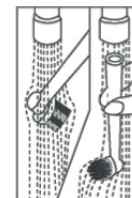Replacing the brush head of the toothbrush

- Replace the brush head when the blue bristles turn completely white and bend outwards. The discoloration of the brush bristles indicates that it is worn and you should replace it.
- Remove the brush head by pulling it upwards.
- When installing the new brush head, make sure that the coloured ring is positioned correctly.
- We recommend that you replace the brush head every 3 months. Replacement brush heads compatible with your electric toothbrush can be purchased from your retail store.

Cleaning the body of the electric toothbrush

- Clean the body of the appliance by running it under tap water first, and wipe it with a towel.
- To clean the body of the appliance, you can use a mild detergent if necessary.
- Make sure not to damage the coloured ring at the end of the brush head.

Cleaning the charging base

- Wipe the charging base with a soft cloth slightly moistened with water or possibly a mild detergent.
- Do not run tap water over the charging base.
- Make sure the power cord does not get twisted.

troubleshooting guide

The electric toothbrush or the water flosser is not working:

- You're using the electric toothbrush or the water flosser for the first time or they haven't been used for more than three months. Charge for at least 24 hours.

The water flosser water jet pressure is low:

- The nozzle is blocked or misshapen: replace it.

The water flosser or toothbrush only works slowly or stops working after a few minutes:

- The filter at the end of the water flosser suction hose is blocked: clean it.
- The battery is flat or wasn't charged long enough. Charge the device for at least 24 hours before using it.
- If it still doesn't work after charging, contact your retail store's customer service department.

Water does not spray through the water flosser's nozzle:

- The tank is empty: fill it.
- The water flosser is tilted at an excessive angle during use, use it in a more upright position.

If the problem persists, contact your retail store's customer service department.

recycling the battery

These two devices have a built-in rechargeable battery that cannot be removed. At the device's end of life, return it to your retail store or place it in a dedicated container for this type of waste.

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EN		
Information published	Value and precision	Unit
Manufacturer's name or trade mark, commercial registration number and address	GUANG DONG XINBAO ELECTRIC APPLIANCES HOLDINGS CO. , LTD. Zhenghe South Road, Leliu Town, Shunde District, Foshan City, Guangdong Province, P.R. China	
Model identifier	ZD006C050100EUE	
Input voltage	100-240	V
Input AC frequency	50-60	Hz
Output voltage	5.0	V $\equiv \equiv \equiv$
Output current	1.0	A
Output power	5.0	W
Average active efficiency	73.0	%
Efficiency at low load (10 %)	68.0	%
No-load power consumption	0,10	W

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All information, designs, drawings and pictures in this document are the property of SOURCING & CREATION. SOURCING & CREATION reserves all rights to its brands, designs and information. Any copy and reproduction through any means shall be deemed and considered as counterfeiting.

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Protection of the environment

This symbol attached to the product means that it is an appliance whose disposal is subject to the directive on waste from electrical and electronic equipment (WEEE). This appliance may not in any way be treated as household waste and must be subject to a specific type of removal for this type of waste. Recycling and recovery systems are available in your area (waste removal) and by distributors. By taking your appliance at its end of life to a recycling facility, you will contribute to environmental conservation and prevent any harm to your health.

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Tested in our laboratories

Warranty valid from the date of purchase (receipt as proof of purchase). This warranty does not cover defects or damage caused by improper set up, incorrect use, or normal wear and tear of this product.

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