essentiel

Notice d'utilisation
User guide
Bedienungsanleitung
Gebruiksaanwijzing

Epilateur électrique
Electric epilator
Elektrischer Epilierer
Elektrischepilineerapparaat

EEE2 EPILYS

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Your product

Box contents

- 1 electric epilator
- 1 epilation head
- 1 shaving head
- 1 cutting guide
- 1 USB-C charging cable
- 1 storage pouch
- 1 cleaning brush
- 1 user guide

Technical features

- Li-ion battery: 3.7 V $\equiv$ 800 mAh
- Power: 5 W
- Epilation and shaving heads washable with water
- 2 speeds
- Battery life: 60 minutes
- Full charging time: 90 minutes
- In-operation and low battery indicator lights
- On/Off button
- Release button

Charging the appliance

1. Connect the USB-C end of the charging cable supplied to the USB-C port on your epilator.
2. Connect the other end of the USB-A cable to the USB port of a 5 V 1,000 mA mains adaptor (not supplied).
3. Connect the mains adaptor to an electrical socket. Charging begins.
4. When charging, the indicator light glows solid red.
5. When charging is complete, the indicator light glows white. You can remove the adaptor from the power supply and remove the USB-C cable from the appliance.

- A **full charge** lasts for **approximately 90 min** and provides a **battery life** with operation at slow speed of **approximately 60 min**.

NOTES:

- **You cannot use the epilator while the battery is charging.**
- The epilator may heat up slightly during the charging process, but this does not affect its operation.
- When charging for the first time, or if the appliance has not been used for some time, the charging time will be slightly longer.


- If the battery level is low, the indicator light will flash red.

Optimising the lifespan of the rechargeable battery

- Unplug the AC adaptor from the wall socket when the battery is fully charged.
- Do not leave the adaptor plugged into the wall socket for more than 24 hours.
- Fully discharge the battery twice a year, leaving the shaver running until it stops. Then fully recharge the battery.
- Preferably place the appliance in a room where the ambient temperature is between 15°C and 35°C when charging.

Preparing the appliance

- Place the removable cutting head of your choice on the body of the epilator and press until you hear a “click”.
- To switch between the shaving head and the epilation head, press the **RELEASE button on the back of the appliance, then lift the cutting head into place.**

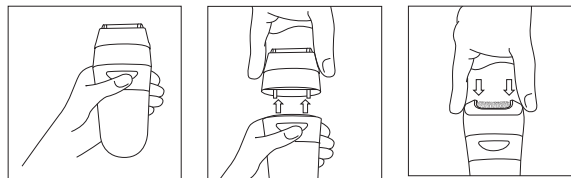- If you're going to be away (on a trip, for example), press and hold the On/Off button to lock the appliance. The indicator light flashes white three times and then goes out. Do the same to unlock your epilator.

Use

For hygiene reasons, the razor should only be used by one person.

Hair removal

1. Insert the epilation head.
2. Press the On/Off button to switch on the appliance.
The appliance starts at default speed 1 (fine hair). The indicator light glows white.
3. Press the On/Off button again to switch to speed 2 (thick hair).
4. Move the epilator by pressing lightly on the skin and in the opposite direction to hair growth. For the closest possible epilation, place the epilator at an angle of around 90° and stretch the skin slightly.
5. Press the On/Off button a third time to stop the epilator.
6. Clean the appliance after each use (see “Cleaning the appliance”) paragraph, using the cleaning brush supplied.

- Use your epilator on a maximum hair length of approximately 0.5 cm.
- Always use your epilator against the direction of hair growth and do not use circular movements when epilating.

EN

INFORMATION ON HAIR REMOVAL

- Hair removal at the root can cause skin irritation (itching, redness) depending on the condition of the skin and hair. This reaction is normal, but should disappear quickly. It may be more significant when an epilator is used for the first time, or if the skin is sensitive. If your skin remains irritated for more than 36 hours, we advise you to consult a doctor.
- Regular use of the epilator reduces the reaction of the skin and the sensation of pain. Inflammation can occur when bacteria penetrate under the skin. We recommend cleaning the epilation head with a disinfectant solution before each use.

- Before using the appliance, you should consult a doctor in the following cases: if you have eczema, cuts, inflammatory reactions of the skin such as folliculitis (infection of the hair follicles), varicose veins around a mole, reduced immunity of the skin, such as in cases of non-insulin-dependent diabetes, pregnancy, Raynaud's disease, haemophilia, immune deficiency or Candida.

AFTER HAIR REMOVAL

Once you've finished the hair removal, apply a soothing oil, a moisturising lotion or milk to limit hair regrowth.

To avoid ingrown hairs, remember to moisturise your skin regularly and to exfoliate, particularly when the hair grows back.

Exfoliation removes dead skin cells and refines your skin, making hair removal easier. We recommend that you exfoliate the day before using your electric epilator. However, avoid exfoliating just after waxing.

Shaving

BEFORE EACH USE

- The epilator must be free of hair and dirt. Use the cleaning brush supplied to clean the appliance. You can also remove the grille from the appliance and run it under a trickle of water. Then dry the grille before reattaching it to the appliance.
- Make sure the appliance is clean and in optimal working order.

PREPARING THE SKIN:

- Your skin should be **clean and dry**. Do not apply creams, moisturisers or alcohol-based lotions before shaving.
- Your appliance will be much more efficient and easier to use if you use it regularly.

WHILE SHAVING

1. Insert the shaving head, and the cutting guide if necessary.
2. Press the On/Off button to switch on the appliance. The appliance starts at default speed 1 (fine hair). The indicator light glows white.

3. Press the On/Off button again to switch to speed 2 (thick hair).
4. Shave by stretching the skin with your free hand. This lifts the hair and reduces the feeling of tightness.


For optimum performance, position the shaving head so that the shaving grille and trimmer are in direct contact with the skin. The On/Off button on the appliance must face the skin.

5. Move the appliance over the skin in a slow, continuous motion, against the direction of hair growth.


- To avoid any risk of irritation, do not press the appliance too hard against your skin while shaving.
 - In some areas of the body, hair grows in several directions. If this is the case, move the machine in several directions for best results.
 - If you sweat while shaving, pat your skin dry with a soft cloth.
6. When you have finished using it, press the On/Off button a third time to turn off the appliance.
 7. Clean the appliance after each use (refer to the “Cleaning and care” section) with the cleaning brush supplied.

- Your skin may need a few days to a few weeks to get used to the shaving system on your appliance.

The clipper:

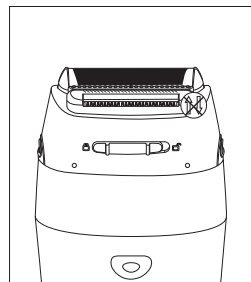To lock the clipper, turn the switch to the  position.

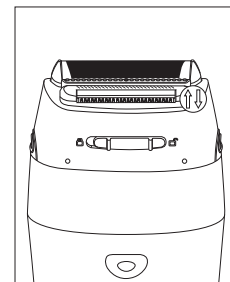To release it, turn the switch to the  position.

The cutting guide:

- The comb can be used as a cutting guide in the bikini area.
- Place the comb on the shaving head, with the teeth facing inwards, until it clicks into place. Position the shaving head at a 45° angle.
- To remove the comb, pull it back slightly.

AFTER SHAVING


- To relax the skin, we recommend you apply a moisturiser as soon as you've finished shaving.
- Use an exfoliating cream or sponge regularly (in the shower, for example) 24 hours after shaving to remove the top layer of dead skin.
- Clean the appliance immediately after use (see "Cleaning and care" section).

Cleaning the appliance

- Always switch off the appliance and unplug it from its power supply before cleaning it.
- Press the **RELEASE** button on the back of the epilator to remove the detachable cutting head.


THE EPILATION HEAD

- Run the cleaning brush over both the epilation head and inside the shaver head.
- Gently tap the shaver head against a flat surface (with the epilation head removed).

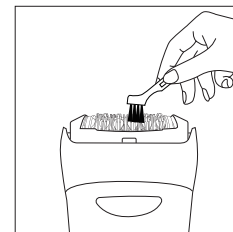SHAVING HEAD

- Run the cleaning brush over both the shaving head and the inside of the shaver head.
- Gently tap the shaver head against a flat surface (with the shaving head removed).
- Press the side eject buttons and remove the shaver head.
- Rinse the shaver head and cutter block separately under warm running water. Leave the cutter block and shaver head to air dry separately.
- **Lubrication:** for optimum shaving comfort, we recommend lubricating the shaving head about once a month: apply a drop of oil (supplied) to the shaving head, then run the appliance for a few seconds. Then switch off the shaver and dry off any excess oil using a soft, dry cloth.

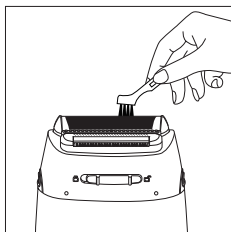BODY OF THE APPLIANCE

- To clean the body of the epilator, use a cloth moistened with a mild detergent. The body of the epilator should never be immersed in water or any other liquid. Never use a scouring sponge nor abrasive detergents to clean the appliance.
- Leave to air dry.
- Then replace the removable cutting head on the epilator.

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All information, designs, drawings and pictures in this document are the property of SOURCING & CREATION. SOURCING & CREATION reserves all rights to its brands, designs and information. Any copy and reproduction through any means shall be deemed and considered as counterfeiting.

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Protection of the environment

This symbol attached to the product means that it is an appliance whose disposal is subject to the directive on waste from electrical and electronic equipment (WEEE). This appliance may not in any way be treated as household waste and must be subject to a specific type of removal for this type of waste. Recycling and recovery systems are available in your area (waste removal) and by distributors. By taking your appliance at its end of life to a recycling facility, you will contribute to environmental conservation and prevent any harm to your health.

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* Tested in our laboratories

Warranty valid from the date of purchase (receipt as proof of purchase). This warranty does not cover defects or damage caused by improper set up, incorrect use, or normal wear and tear of this product.

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essentiel

Art. 8012107
Ref. EEE2 Epilys

SERVICE RELATION CLIENTS /
Customer Relations Department /
Kundenservice / Klientenafdeling

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