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FR APPAREIL ANTI-DEMANGEAISON
EN INSECT BITE HEALER
NL ANTI-PRIKAPPARAAT
ES DISPOSITIVO ANTIPICADURAS
DE ANTI-STICH-GERÄT

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Photo non contractuelle/ Foto is niet bindend / Imagen no contractual / Non contractual photo / Nicht-vertragliches Foto

Importer :

EURO-TECH DISTRIBUTION 37A RUE CESAR LORIDAN BONDUES France

info@eurotechdistribution.com

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Email: EU-Rep@share-info.com

AOJ HEALTH Technology Co., Ltd

Room 202, Building 5, Shuhe Industrial Park,

Sanwei Community, Hangcheng Street, Bao'an District

518126 Shenzhen, China

READ ALL INSTRUCTIONS IN THIS MANUAL CAREFULLY, FOLLOWING THE USAGE AND SAFETY GUIDELINES. KEEP THIS MANUAL AND ATTACH IT TO THE DEVICE IF ANOTHER PERSON IS TO USE IT.

SAFETY INSTRUCTIONS

The device provides relief from insect bites and reduces itching and swelling caused by insect bites or stings by producing heat (local hyperthermia).

Keep the device out of the reach of children.

Only use the insect bite relief device to treat insect bites or stings.

This device is intended for domestic use.

The device may be used by persons over 3 years of age under supervision and for self-treatment by persons over 12 years of age.

Children should be supervised to ensure that they do not play with the device.

Keep the device out of the reach of children and pets.

Keep packaging material out of the reach of children (risk of suffocation).

WARNING: There is a risk of injury if used incorrectly.

Caution! Never immerse the device in water or any other liquid.

Battery type: 2 x AAA (not included)

Used batteries must be disposed of in accordance with local regulations in your area.

If a battery leaks, stop using it. Place the batteries in a sealed container and dispose of them in accordance with local regulations in your area.

Do not expose the device or battery to excessive temperatures.

Be careful not to short-circuit the terminals of the battery-powered device or the battery with metal objects.

For information on cleaning, maintenance and disposal of the device, refer to the following sections of the manual.

For information on operating times, refer to the section below.

Before using the device, consult your doctor if you suffer from skin conditions, serious illness (especially if you have a predisposition to thromboembolic diseases or recurrent malignant tumours) or unexplained chronic pain in any part of your body.

This device is not a substitute for medical consultation. If symptoms persist despite using the insect bite treatment product, discontinue use. Consult a doctor immediately if symptoms worsen.

Consult a doctor immediately in the event of a tick bite. This bite can lead to the transmission of pathogens such as tick-borne encephalitis or borreliosis.

Consult a doctor immediately if the first signs of increased itching, skin redness, swelling, dizziness, nausea, shortness of breath, papule formation or a drop or rise in blood pressure appear.

Do not use near the eyes, on the eyelids, mouth or mucous membranes. In addition, the device should not be used on wet skin.

Do not use this device if you have diabetes. Pain sensitivity may be reduced in some circumstances and some people with diabetes may burn themselves.

Do not use the device on disabled persons, children under 3 years of age or persons with reduced sensitivity to heat. This includes persons with skin alterations due to disease or scarring, after taking painkillers, or after consuming alcohol or drugs.

Do not use the device on severely or chronically diseased skin (e.g. inflamed skin, whether painful or not, reddened skin, rashes, burns, frostbite, bruises, swelling, open and healing wounds, and surgical scars where the healing process could be affected).

Do not use the device if you have a fever.

Do not use the device on sensitive skin areas.

Do not use the device if you suffer from a sensory impairment that reduces your sensitivity to pain (e.g. metabolic disorders).

Do not use the device if you are simultaneously using lotions, creams or gels applied locally.

Do not use the device if you suffer from persistent skin irritation due to prolonged application of heat to the same area of skin.

If the device does not work properly, if you feel unwell or experience pain, stop using it immediately.

If users or patients have been involved in a serious incident related to the device, please inform the manufacturer and the competent authority of the Member State in which you are established.

The use of the device may be limited in the presence of electromagnetic interference. This could lead to problems such as error messages or device failure.


Avoid using this device in close proximity to other devices or stacked on top of other devices, as this may cause malfunction.

GENERAL PRECAUTIONS

- Do not touch the device with wet or damp hands. Do not allow water to splash onto the device. Only use the device when it is completely dry.
- Do not use detergents or solvents to clean the device. If liquid enters the device, it could damage the electronic components or other parts of the device and cause it to malfunction.
- The device should only be used by one person; its use by multiple people is not recommended.
- Protect the device from direct sunlight and store it in a dry place.
- Protect the device from impact, moisture, dust, chemicals, significant temperature fluctuations, and heat sources such as ovens.
- Ensure that no water enters the device. If this should happen, remove the batteries and only use the device again when it is completely dry.
- If you are storing the appliance, do not place any heavy objects on top of it.
- Before use, ensure that the device and its accessories show no visible damage. If in doubt, do not use the device and contact customer service.

- Store the device in a location protected from climatic influences. The device must be stored under the specified environmental conditions.
- When the ambient temperature is below -10°C or above 50°C , the time required between the minimum storage temperature and use is at least 4 hours, and the time required between the maximum storage temperature and use is also at least 4 hours.

DESCRIPTION OF THE DEVICE

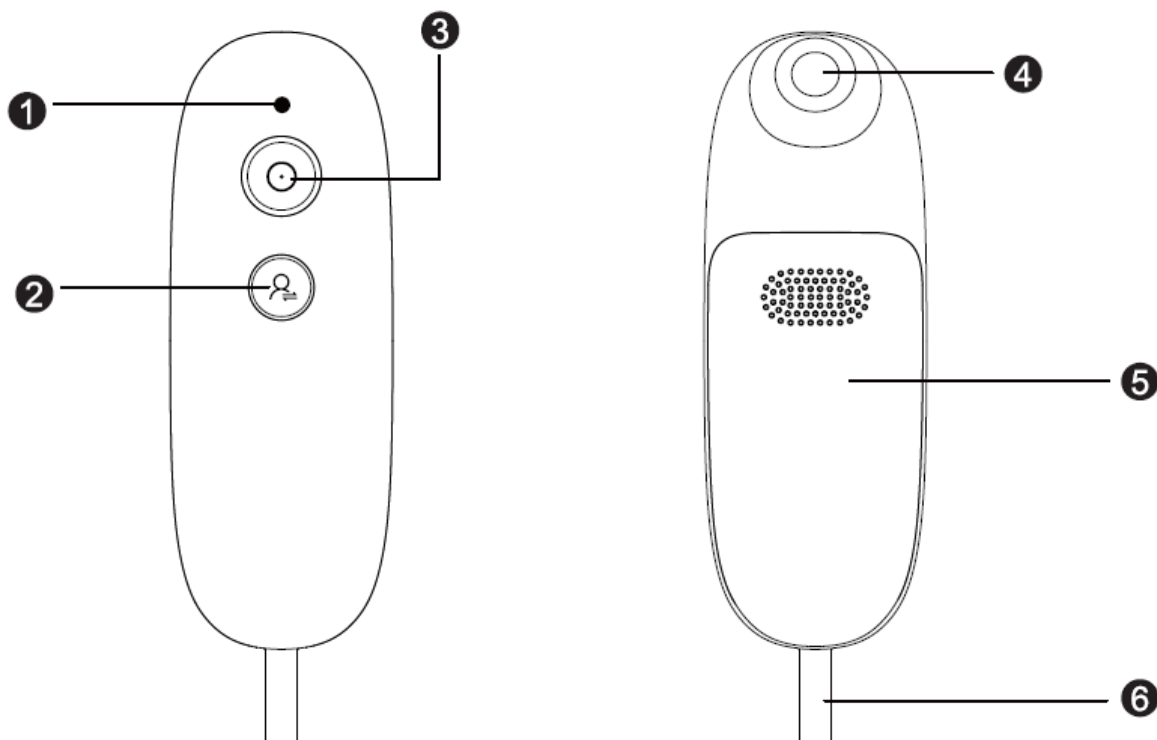1. LED indicator
2. Activation button
3. Power button
4. Heating disc
5. Battery compartment
6. Silicone hook

Before first use :

- Unpack the device and its various components and ensure that they are in good condition.
- Read all the instructions in this leaflet carefully, following the usage and safety guidelines. Keep this leaflet and advise potential users.

Inserting the batteries :

1. Slide the battery compartment cover (5) to open it.
2. Insert 2 AAA batteries into the battery compartment, ensuring that the polarity (+/-) is correct.
3. Slide the battery compartment cover until you feel and hear a click indicating that it is securely in place.

Important: remove the batteries if you do not plan to use the device for a long period of time (more than 3 months).

USE

Before using the device, ensure that there are no insect stingers left in the affected area.

Use the device immediately after an insect sting or bite for best results.

1. Press the power button (3) to start the device. The LED indicator (1) will light up green and a beep will sound to indicate that the device is ready for use.
2. **Mode 1 (for sensitive skin) :** After being stung or bitten, apply the heating disc (4) of the device to the affected area and press the activation button (2) once. The LED indicator light (1) will turn blue, indicating that the device is working on the bitten or stung area. The device maintains a treatment temperature of 47°C (give or take 2°C) for 3 seconds. This short application time is recommended for first-time use or on sensitive skin.

3. **Mode 2 (on normal skin)** : After being stung or bitten, apply the heating disc (4) of the device to the affected area and press the activation button (2) twice quickly. The LED indicator (1) will light up yellow, indicating that the device is working on the affected area.
The device maintains a treatment temperature of 50°C (give or take 2°C) for 6 seconds. This application time is recommended for regular use on normal skin.
4. Once the application is complete, an audible signal will sound, the heating disc (4) will automatically switch off and the LED indicator (1) will turn green again. Remove the heating disc (4) from the affected area of skin as soon as the audible signal sounds.

IMPORTANT : Depending on the sensitivity of the person using the device and the affected skin area, the treatment temperature may be uncomfortable and cause slight skin redness and, in rare cases, skin irritation.

To achieve this soothing effect on the skin nerves, a brief period of thermal discomfort is necessary. Only slight thermal discomfort is required for the treatment to be effective. If the heat function seems too hot, stop the treatment immediately.

5. If necessary, you can reuse the heat function of the device on the sting/bite after a 2-minute break. The maximum number of 5 applications per hour on the same treatment area must not be exceeded.
6. Once you have finished using the device, press the power button (3) to switch it off. The LED indicator (1) will turn off. If no further operation is performed within 5 minutes, the device will switch off automatically.
7. Before storing the appliance, ensure that it is switched off and that the heating disc (4) has cooled down.

Note : When the battery level is low, the LED indicator (1) will flash three times and an audible signal will sound three times. Replace the batteries before using the device.

Battery warnings

- The device operates with 2 AAA batteries, not included.
- Different types of batteries or new and used batteries should not be mixed.
- The batteries must be inserted with the correct polarity.
- Used batteries must be removed from the product.
- Batteries must be disposed of safely. Place them in designated collection bins (ask your retailer for details) to protect the environment.
- The terminals of a battery or accumulator must not be short-circuited.
- In the event of a battery or accumulator leak, avoid contact of the liquid with the skin or eyes.
- In case of contact, rinse immediately with plenty of clean water and seek medical advice.
- If your skin or eyes come into contact with the battery fluid, rinse the affected areas with water and consult a doctor.
- Choking hazard! Young children may swallow batteries and choke. Keep batteries out of reach of children.
- Risk of explosion! Do not dispose of batteries in fire.

- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.
- Do not dismantle, open or crush the batteries.
- Protect batteries from excessive heat.
- Do not recharge or short-circuit the batteries.
- If the device is not used for a long period of time, remove the batteries from the device.
- Only use identical or equivalent batteries.
- Always replace all batteries at the same time.
- Do not use rechargeable batteries!

SYMBOLS USED

The symbols shown below are indicated on the appliance, in this user manual, on the packaging or on the appliance's rating label.

WARNING

Consult the user manual

BF type device

Disposal in accordance with the European Directive on Waste Electrical and Electronic Equipment (WEEE)

Manufacturer

Date of manufacture

Batch number

IP 42 Device protected against solid foreign bodies of 1.0 mm or greater and against diagonally falling water droplets

Authorised representative in the European Community

⁰⁵⁹⁸ This product complies with the requirements of current European and national directives

RoHS This product complies with this European directive

Medical device

Indicates the entity importing the medical device

Unique device identifier

(01)16976313750028

(10)IBR1000001

(11)251021

CLEANING AND MAINTENANCE

- Ensure that the device is switched off, cooled down and that the batteries have been removed before each cleaning.
- It is recommended to clean the device before each use.
- Only use soft, dry cloths to clean the external parts of the appliance.
- Never use cleaning products that could be toxic if they come into contact with the skin or mucous membranes, or if they are swallowed or inhaled. Do not use abrasive cleaning products and never immerse the device in water or any other liquid.
- Do not clean the appliance or its accessories in the dishwasher.
- Avoid mechanical cleaning with brushes or similar tools, as this could cause irreparable damage.
- When using the device on multiple areas that have been stung or bitten, it is recommended that you clean the heating disc each time. This reduces the risk of transferring germs and bacteria from one area of skin to another. Use a commercially available alcohol-based disinfectant wipe to disinfect the device.

TROUBLESHOOTING

Problem	Possible causes	Solutions
The device does not switch on.	The batteries have been installed incorrectly / The batteries are flat.	Insert the batteries correctly, ensuring that the polarity is correct. Replace both batteries as indicated in the instructions.
The appliance no longer heats up sufficiently.	The battery level is too low.	Replace both batteries as indicated in the instructions.
The device emits 3 beeps/the LED indicator flashes, even after switching the device off and on again.	The battery level is too low, which means that the device may no longer function properly.	Replace both batteries as described in the instructions.
The device does not start up even though the LED indicator light turns green.	The device must be restarted.	Press the power button again.
The device emits 5 beeps/the LED indicator flashes even after switching the device off and on again.	The device may be faulty.	Please contact customer service.

TECHNICAL SPECIFICATIONS

Product reference	002956
Model	IBR10
Device class	Ila
Dimensions	109,5mm x 34mm x 24,5mm
Weight (without batteries)	25g
Power supply	2 x AAA batteries (LR03 1.5V) not included
Maximum temperature	50°C +/- 2°C
Functions	Application on sensitive skin: 3 seconds at 47°C +/- 2°C / Application on normal skin: 6 seconds at 50°C +/- 2°C
Lifespan	5 years
Expected operating time	Approximately 300 uses
Permitted storage and transport environment	-10 °C to 50 °C, relative humidity 30% to 85%, pressure 700 hPa to 1060 hPa
Permissible operating environment	10 °C to 40 °C, relative humidity 30% to 85%, pressure 700 hPa to 1060 hPa

This device complies with Regulation (EU) 2017/745 on medical devices and the following standards: EN 60601-1: 2006 + A1: 2013 + A2: 2021, EN 60601-1-1-2: 2015 + A1: 2021, EN 60601-1-11: 2015 + A1: 2021

and EN 60601-1-6: 2010 + A1: 2015 + A2: 2021 in accordance with :

Electromagnetic emission	
CISPR11	Group 1
CISPR11	Class B
Electromagnetic immunity	
IEC 61000-4-2	
IEC 61000-4-8	
IEC 61000-4-3	

The device meets the special safety requirements relating to electromagnetic compatibility.

WARRANTY

This product is guaranteed for a fixed period with the distributor from the date of purchase against any failure resulting from a manufacturing or material defect. This warranty does not cover defects or damage resulting from improper installation, incorrect use, or abnormal wear and tear.

The device must be returned in its original packaging and accompanied by the receipt, provided it has been handled correctly. Therefore, carefully follow the usage instructions and safety instructions stated above. Furthermore, we cannot guarantee your device if you or a third party have made any modifications or repairs to it.

SELECTIVE COLLECTION OF ELECTRICAL AND ELECTRONIC WASTE

Electrical products should not be disposed of with household products. According to European Directive 2012/19/EU on the waste of electrical and electronic equipment and its implementation in national law, used electrical products must be collected separately and disposed of at designated collection points. Contact your local authority or retailer for recycling advice.

Share Info GmbH

Am Schulzentrum 12, 41564 Kaarst, Germany

Tel: 0049 179 5666 508

Email: EU-Rep@share-info.com

AOJ HEALTH Technology Co., Ltd

Room 202, Building 5, Shuhe Industrial Park,

Sanwei Community, Hangcheng Street, Bao'an District

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