
Withings Body+

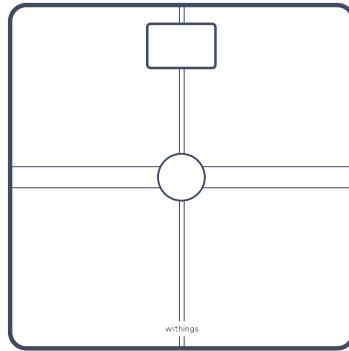QUICK INSTALLATION GUIDE

GUIDE D'INSTALLATION | INSTALLATIONSANLEITUNG | GUÍA DE INSTALACIÓN | GUIDA ALL'INSTALLAZIONE

EN | FR | DE | ES | IT

THANK YOU FOR CHOOSING THE WITHINGS BODY+


MERCI D'AVOIR CHOISI WITHINGS BODY+
WIR BEDANKEN UNS FÜR IHRE WAHL ZUR WITHINGS BODY+
TE AGRADECEMOS QUE HAYAS ELEGIDO WITHINGS BODY+
GRAZIE PER AVER SCELTO WITHINGS BODY+

This guide is available in additional languages at: www.withings.com/guides

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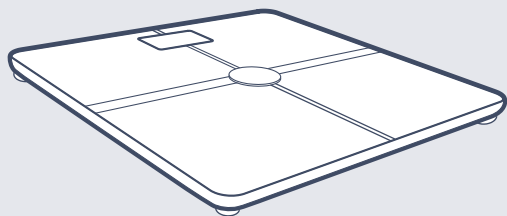
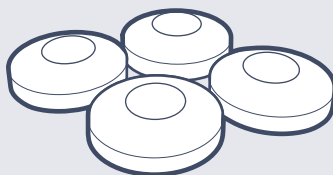
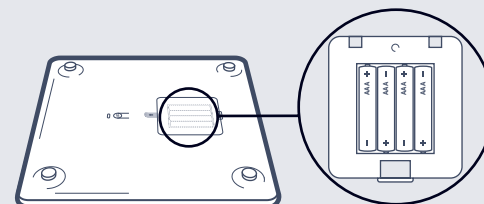
Installation of carpet feet (optional)

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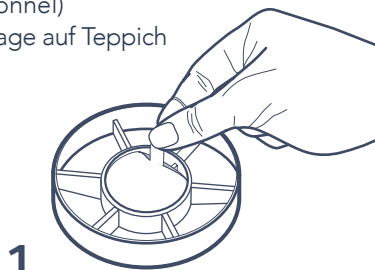
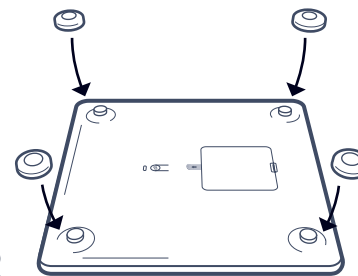
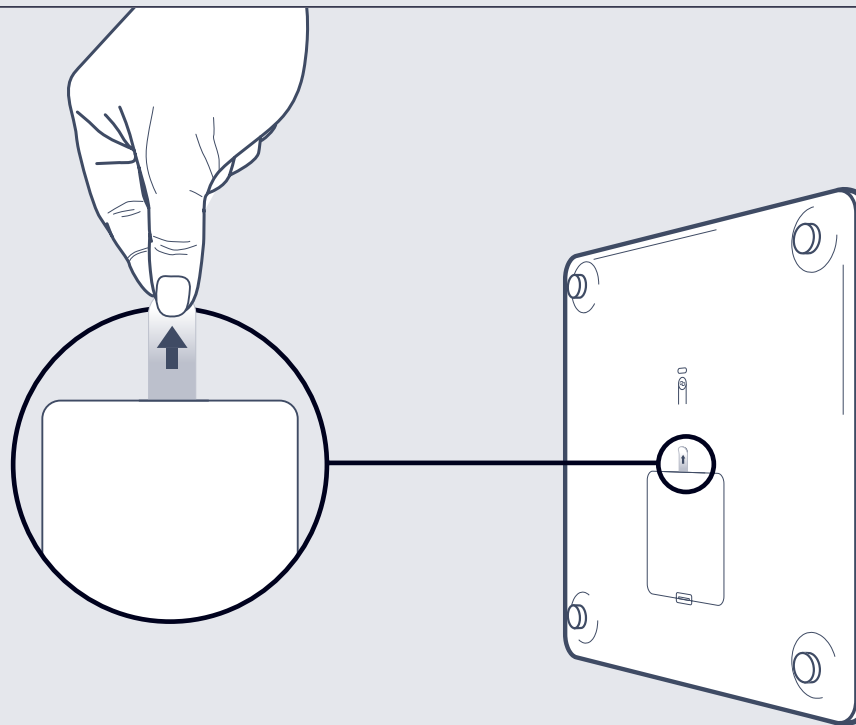
2

Visit go.withings.com to get the app and follow the instructions.

- FR Obtenez l'appli via go.withings.com et suivez les instructions.
- DE Holen Sie sich die App auf go.withings.com und folgen Sie den Anweisungen.
- ES Visita go.withings.com para conseguir la aplicación y sigue las instrucciones.
- IT Visita go.withings.com per ottenere l'app e segui le istruzioni.

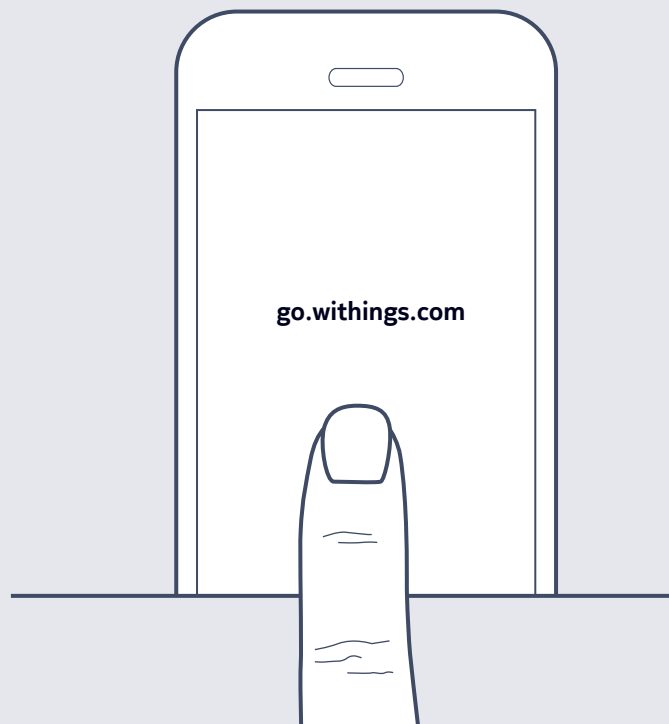Know more with Body+


En savoir plus sur Body+ | Mehr über Body+ erfahren | Descubre Body+ | Scopri Body+

Analyze your Body Composition, and be notified of Rain Forecast for the day, while weighing.

- FR** Analysez votre composition corporelle, et soyez avertis en cas d'alerte pluie ou météo, pendant la pesée.
- DE** Analysieren Sie Ihre Körperzusammensetzung und erhalten Sie beim wiegen den regenbericht für heute.
- ES** Analiza tu composición corporal y recibe el pronóstico de lluvia para ese día, todo ello mientras te pesas.
- IT** Analizza la tua composizione corporea e visualizzerai, mentre ti pesi, le previsioni meteo del giorno direttamente sullo schermo.

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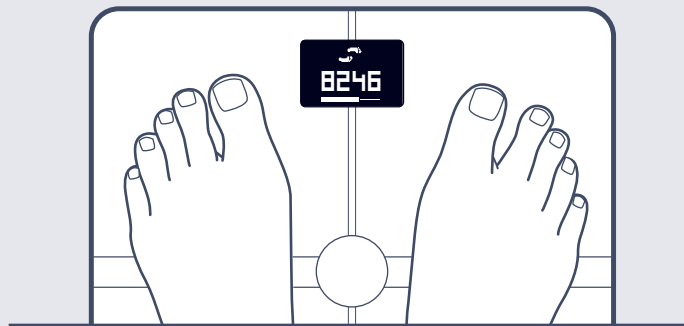Change the order of the screens and choose what matters most to you.

If you use a Withings watch, activity tracker or track your steps with the Health Mate™ app, activate the Steps screen to know how much you walked the previous day!

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

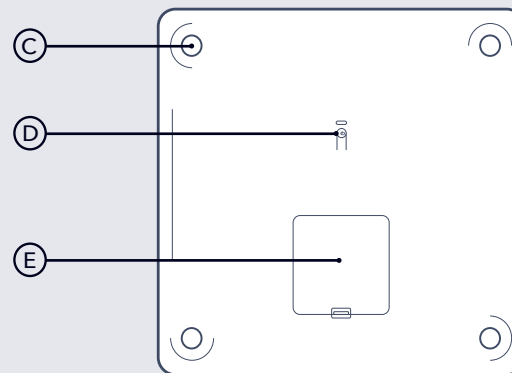
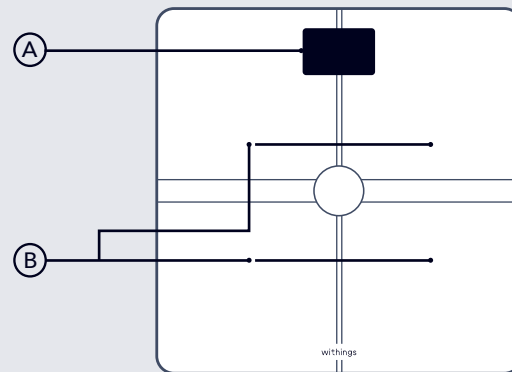
Use only a cloth dampened with water to clean your scale. The use of any detergents or cleaners may cause irreversible damage to the coating that protects the electrodes (B) used to measure your body fat.

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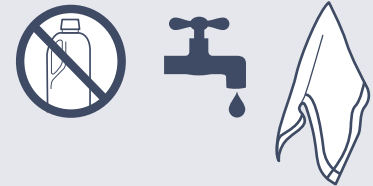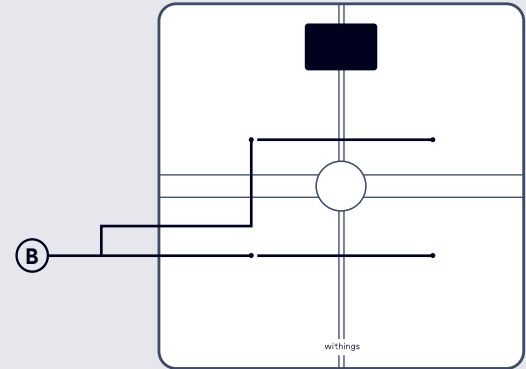12

This product is not suitable for people wearing a pacemaker or any other internal medical device.

- FR** N'utilisez pas ce produit si vous portez un stimulateur cardiaque ou autre appareil médical.
- DE** Benutzen Sie dieses Produkt nicht, wenn Sie einen Herzschrittmacher oder ein anderes medizinisches Gerät tragen.
- ES** Este producto no es adecuado para personas que llevan un marcapasos u otro aparato médico implantado.
- IT** Non utilizzare questo prodotto se porti uno stimolatore cardiaco o ogni altro dispositivo medico interno.

This product is not a medical device and should not be used to diagnose or treat any medical condition.

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13

IC: 11411A-WBSO5

This device complies with Industry Canada's licence exempt RSSs. Operation is subject to the following two conditions:

- (1) This device may not cause interference; and
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

Radiation Exposure Statement

The product complies with the Canada portable RF exposure limit set forth for an uncontrolled environment and are safe for intended operation as described in this manual. The further RF exposure reduction can be achieved if the product can be kept as far as possible from the user body or set the device to lower output power if such function is available.

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Garantie | Garantie | Garanzia | Garanzia

14

US Withings™ Two (2) Years Limited Warranty - Withings Body+

Withings SA, 2 rue Maurice Hartmann, 92130 Issy-les-Moulineaux, («Withings») warrants the Withings-branded hardware product («Withings Product») against defects in materials and workmanship when used normally in accordance with Withings' published guidelines for a period of TWO (2) YEARS from the date of original retail purchase by the end-user purchaser («Warranty Period»). Withings' published guidelines include but are not limited to information contained in technical specifications, safety instructions or quick start guide. Withings does not warrant that the operation of the Withings Product will be uninterrupted or error-free. Withings is not responsible for damage arising from failure to follow instructions relating to the Withings Product's use.

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This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) This device must accept any interference received, including interference that may cause undesired operation.

The crossed out wheeled bin symbol on your product, battery and/or accessories means that the item is classified as electrical and electronic equipment. Such items should not be mixed with general household waste, and must be taken to dedicated collection points at the end of their working life for proper treatment, recovery and recycling. This way you help save valuable resources and promote recycling of materials.

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QIG_WBS05_85

withings

Withings Body+ | Body Composition Wi-Fi Scale
© 2018 Withings SA. All rights reserved.