OXILOM

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Taille de la housse: 219 x 157 x 88 mm

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Production introduction

MassEye 2 massager is a product that will provide you relaxation and wellness. Designed from technical studies and massages, it rubs your temples and other acupressure points in the ocular region for a real sensation of tranquility. The method used combines pneumatic massage, soft vibration and hot compression in order to relieve eye fatigue and to fully distress.

MassEye 2 has an optional music function that can be activated during the massage. You can choose between soothing sounds by defaults or your favorite songs by creating your own playlist. You can download your music from your computer to the eye massager. This eye massager, your haven of peace, has a sleeping mask shape and an ergonomic and elegant design. It is convenient to carry with its collapsible form and is comfortable to wear.

Users

MassEye 2 massager can be used by anyone who works on a computer or watches their television or tab screen for a long time but also by anyone who stays up late or makes important ocular efforts on a daily basis. For the people who have trouble to fall asleep, MassEye 2 massager helps relax and sleep. It allows the users to relieve from stress, migraine and hormonal headaches and stimulates bloodstream around the eyes so they seem significantly well-rested and less ringed.

Features

- Smart pneumatic massage: massages temples and other acupressure points in the ocular region through air pressure.
- Vibrating massage: imitates soft hammering, shiatsu and other traditional massage techniques on the important acupressure points around the eyes.
- Hot pressure: improves bloodstream around the eyes thanks to its top quality infrared system.
- Collapsible shape easy to carry with its case.
- 2 buttons: one operating button and one music button for an ease of use.
- LCD display, figures and functions in white to be read better.
- Removable and washable cloth to avoid direct contact with the product and make it more hygienic to use.

Product content:



- 1.One MassEye 2 device
- 2.One USB cable
- 3.One case
- 4.One instruction manual
- 5.One protection cloth

Product Picture


Front of product

1. Adjustable buckle
2. Elastic band
3. On/off and switching mode button
4. Music button: on/off, adjusting music volume, switching song
5. Charging USB interface
6. LCD display : modes, duration and level of battery display

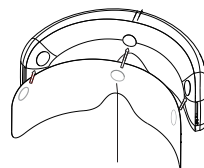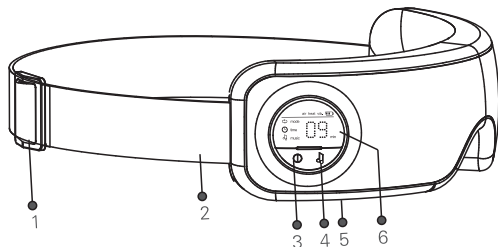Velcro tape

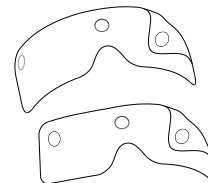Removable and washable cloth

Operating mode

Preparation before use:

1. Integrated rechargeable battery

Press briefly “” to turn the eye massager on and check the battery level icon on the LCD display. If the battery level indicator flickers, the battery is low. In that case, connect the product with an USB cable to an AC adaptor or to a computer. It is recommended to use MassEye 2 when it is fully charged.

2. Use of MassEye 2

Join the cloth to the Velcro tapes inside the eye massager according to the diagram. Put the MassEye 2 massager on and adjust the elastic band to make sure to be in a comfortable position.

Buttons and symbols meaning

1.Buttons

“” : On/off, switching mode, long press to turn it off

“” : Music function: on/off, adjust the volume and switch music

2.Icons

“” : Functional operation

“” : Massage time left

“” : Music on

“” : Battery level, indicates also that the battery is charging.

“air” : Indicates that air pressure function is on.

“heat” : Indicates that hot compression function is on.

“vib.” : Indicates that vibrating function is on.

Product function

1. Press “**⏻**” to turn the product on, LCD display lights up and Mode 1 starts by default as well as the music.

2. Press normally “**⏻**” to switch message mode. Whenever you press this button, message mode is switched. Default mode is Mode 1.

The 5 modes are as follow:

Mode 1: Air pressure + hot compression

Mode 2: Every functions

Mode 3: Air pressure + vibration

Mode 4: Air pressure

Mode 5: Hot compression

Press “**⏻**” for 2 seconds in order to switch off the eye massager.

3. Method of use for music function “**🎵**” :

- Press “**🎵**” once to turn on the music.
- Press “**🎵**” to switch up the music when the music is already on.
- Press “**🎵**” twice to adjust the volume when the music is on. The volume by default is the highest one. Volume switches up from the highest one, to the medium one then to the lower one.

Press “**🎵**” for a few seconds to switch off.

Compatible with music downloading from a computer

Connect MassEye 2 USB cable to a computer. Transfer your MP3 music to the eye massager. It only supports MP3 music. Thanks to this function, you can create your own soothing playlist with your favorite songs.

Note: MassEye 2 switches off when you connect it to a computer.

Cleaning and maintenance

1. Wipe the product with a slightly wet cloth or towel to clean MassEye 2. Do not use flammable or explosive product.
2. Keep the product in a dry and airy place.
3. Do not expose it directly to the sunlight or in high temperatures or in a damp place.

Removable and washable cloth maintenance

1. Wipe sweat stain and other stains off the removal and washable cloth with a slightly wet cloth after using the product.
2. The cloth can be washed. Please pay attention to the items below prior to cleaning:
 - a) The cloth can only be washed by hands. Warm water is better to use. Do not scrub it vigorously while cleaning.
 - b) Do not wring for drying, hang it up to dry air.
 - c) Do not use decolorant.
 - d) Do not iron it.

How to charge your battery

1. Low battery if you hear a buzzer ticks three time in a row or when the battery level icon is empty “**☐**” when you use the eye massager. In that case, it is necessary to use the USB cable to charge the battery. The 4 level battery icon bars flickers while charging. The battery level icon displays “**100**” when the battery is fully charged. It takes approximately 2 hours for fully charging the battery.
2. Full charged battery can be used until 4 to 5 days (15 minutes per day).
3. Rechargeable battery life is around 500 charges if not more, or to say about 3 years.

Security measures

1. Please read and understand entirely the instruction manual before getting the product to work.
2. Please clean the eye massager prior to using.
3. Children or any person without expression ability are forbidden to use the product.
4. Those who have sensitive skin or are insensitive to heat are forbidden to use the product.
5. Pregnant women, people with heart or artery or cerebrovascular disease, malignant tumor, or people who see doctors for major health problems are forbidden MassEye 2 without first checking their doctor and get their agreement.
6. Stop using if you feel uncomfortable or if your skin feels unusual.
7. Do not use the product or use it with your doctor's agreement if you have dermatological or ocular problems like those who did eye surgery or with glaucoma, cataract, retinal detachment or any other eye disease.
8. Those wearing eye contact lenses can only use the product when they take them out.
9. Do not use the product in high-temperature, inflammable or damp place, or a place with electromagnetic radiation.
10. Do not disassemble the product without permission. Make sure the product do not break nor get in contact with water for a proper functioning.

Product specification

Name: Eye Massager

Model: MassEye 2

Power supply: 5V DC 1A 5 (USB charger)5W

Lithium battery: 3.7 V, 1,450 mAh

Size of case: 219 x 157 x 88 mm

Weight: Around 660 g

Information on the guarantee

In case of warranty claim, please contact your reseller. If you send them back the product, please indicate the defect and enclose a duplicate of your receipt.

The guarantee conditions are as follow:

1. A two-year guarantee is provided on the product starting the purchase date that must be proved by the receipt or its invoice in case of guarantee.
2. The problems caused by material or manufacturing defects are freely repaired during the warranty period.
3. Excluded from this guarantee :
 - Any damage resulting inconsistent treatment like instruction manual non respect or water.
 - Any damage resulting repairs or intervention from the buyer or any unauthorized third party.
 - Any damage during the transportation from the manufacturer to the consumer or during the return from customer service.
 - Any available part as accessory and subject to normal wear.
 - Any problem resulting natural causes.

Liability for consequential damages, directly or indirectly, caused by the product, is excluded, even when the damage on the product is admitted as a case of warranty.

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Oxilom is a trademark and the product imported by: