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## USER'S GUIDE

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7

When using any electrical products, especially if children are present, certain basic safety precautions should always be followed, including the following:

#### DANGER To avoid the risk of electric shocks:

- Do not use this device in the bath or in the shower.
- Do not handle the USB-C charger with wet hands.
- Do not place the USB-C charger or drop into water or any other liquid.
- Check the condition of the USB-C charger power cord before first use, then regularly during the device life.

#### WARNING To prevent against the risk of fire, burn, or injury to persons:

- Do not connect this device to a voltage system other than the one specified on the device.
- Do not use this device if you notice that a cord or a plug is damaged.
- Do not direct the water under the tongue, inside the ear, nostrils or other delicate areas. This device is capable of producing pressure that could cause serious injury in these areas.
- Use the device as indicated in the instructions for use or follow the recommendations of your dentist, periodontist or pharmacist.
- Use only tips and accessories recommended by NEOPULSE®.
- Do not drop the device
- Do not insert objects into the device holes.
- Keep the USB-C charger away from heated surfaces.
- Do not use outdoors or operate the device, where aerosol products are being used or where oxygen is being administered.
- Do not use concentrated insoluble essential oil in this device. The use of these products will degrade product performance and longevity. The use of these products may cancel the warranty.
- Remove any oral jewelry items before use.
- Do not use the device if you have an open wound on the tongue, in the mouth or in the nose.
- If your doctor has advised you to take antibiotics before dental care, consult your dentist, periodontist or pharmacist before using this device or any other oral hygiene aid.
- Give specific instructions to children (from 6 years) and people with special needs for proper use of this device and supervise them.
- This device does not contain any userreplaceable components, and does not require maintenance.
- Fill the tank only with lukewarm water or another solution recommended by a healthcare professional. Be careful that it does not contain dyes that could damage the system.
- If the product is opened or taken apart for any reason, the warranty will be invalidated.

8

## 2 YEARS WARRANTY

NEOPULSE warrants that this product is free from material and manufacturing defects for a period of two years from the date of purchase.

Please contact NEOPULSE: [conso@neopulse.fr](mailto:conso@neopulse.fr)

## PRODUCT DESCRIPTION

1. Tip
2. Power on/off button
3. Pressure control selection button
4. LED indicators, 6 pressure levels
5. Battery charge indicator
6. USB-C charger plug
7. Removable tank
8. Tank filling lid

#### Accessories :

- USB-C charging cable
- Travel pouch
- Classic tip (x2): perfect for daily use (to replace every 6 months)
- Orthodontic tip (x2): perfect to clean orthodontic appliance (braces) and for general use (to replace every 3 months)

#### Pressure settings:


You can adjust the pressure on the dental water jet by successive presses on the pressure control button (3). The LED indicators (4) shows the pressure level. 6 pressure levels, or 2 levels per mode : 2 for Normal, 2 for Soft and 2 for Pulse.

**Normal pressure:** ideal for those who are used to dental water jet and who are looking for power for better efficiency:  
 - Normal 1 (7.6 bars): more intense blue light  
 - Normal 2 (6.1 bars): less intense blue light

**Soft pressure:** ideal for novices to get used to the pressure of the water, for children (+6 years) and people with sensitive teeth:  
 - Soft 1 (4.75 bars): more intense blue light  
 - Soft 2 (1.7 bars): less intense blue light

**Pulse pressure:** ideal for gum massage:  
 - Pulse 1 (6.9 bars): more intense blue light  
 - Pulse 2 (4.13 bars): less intense blue light

**Recommendation: Start with the lowest setting and increase the pressure according to your own preferences.**

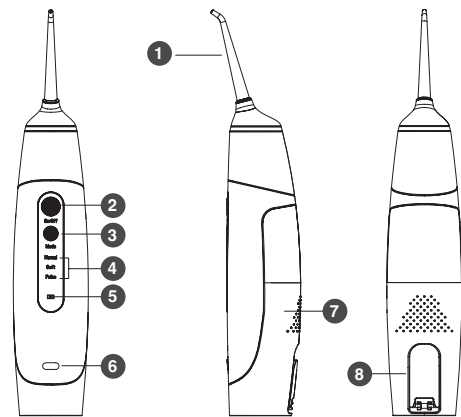9

10