# SOEHNLE

BEDIENUNGSANLEITUNG	8
USER MANUAL	11
MODE D'EMPLOI	14
ISTRUZIONI PER L'USO	17
GEBRUIKSAANWIJZING	20
INSTRUCCIONES DE MANEJO	23
MANUAL DE INSTRUÇÕES	26
BRUGSANVISNING	29
BRUKSANVISNING	32
INSTRUKCJA OBSŁUGI	35
NÁVOD NA OBSLUHU	38
NÁVOD K POUŽÍTÍ	41
ИНСТРУКЦИЯ ПО ЭКСПЛУАТАЦИИ	44
KULLANMA KILAVUZU	47
NAUDOJIMO INSTRUKCIJA	50
УПЪТВАНЕ ЗА УПОТРЕБА	53
دليل الاستخدام	56

Shape Sense Connect 200

Interpretation der Messwerte • Interpretation of measured values • Interprétation des valeurs de mesure • Interpretazione dei valori • Interpretatie van de meetwaarden • Interpretación de los valores medidos • Interpretação dos valores medidos • Tolkning af målingerne • Interpretation av mätvärden • Interpretace naměřených hodnot • Interpretácia meraných hodnôt • Interpretacja wyników pomiaru • Интерпретация измерений • Ölçüm değerlerinin yorumlanması • Matavimo vertių paaiškinimas • Интерпретиране на измерените стойности • شرح قيم القياسات

Display	1	2	3	4	5	6	7	8	9
	-			😊			+		
<b>BMI</b>	< 17,4	17,5 - 17,9	18,0 - 18,4	18,5 - 20,4	20,5 - 22,9	23,0 - 24,9	25,0 - 27,4	27,5 - 29,9	≥ 30
	< -12 %	- 12 %	- 7 %	-3 %	0 %	+ 3 %	+ 7 %	+ 12 %	> + 12 %
	< - 6 %	- 6 %	- 4 %	-2 %	0 %	+ 2 %	+ 4 %	+ 6 %	> + 6 %
	< - 6 %	- 6 %	- 4 %	-2 %	0 %	+ 2 %	+ 4 %	+ 6 %	> + 6 %

Alter/Age	♀						♂					
	-	😊	+	++	😊	😊	-	😊	+	++	😊	😊
17-18	< 15	≥ 15-25	≥ 25-33	> 33	> 52,5	> 28	< 8	≥ 8-18	≥ 18-24	> 24	> 58,5	> 38
18-30	< 20	≥ 20-29	≥ 29-36	> 36	> 49,5	> 28	< 8	≥ 8-18	≥ 18-24	> 24	> 57	> 38
30-40	< 22	≥ 22-31	≥ 31-38	> 38	> 47,5	> 28	< 11	≥ 11-20	≥ 20-26	> 26	> 56	> 38
40-50	< 24	≥ 24-33	≥ 33-40	> 40	> 46,5	> 28	< 13	≥ 13-22	≥ 22-28	> 28	> 55	> 38
50-60	< 26	≥ 26-35	≥ 35-42	> 42	> 46	> 28	< 15	≥ 15-24	≥ 24-30	> 30	> 54	> 38
60+	< 28	≥ 28-37	≥ 37-47	> 47	> 45	> 28	< 17	≥ 17-26	≥ 26-34	> 34	> 53	> 38

Alle Angaben in %, All figures in %, Toutes les données en %, Tutte le indicazioni sono in %, Alle waarden in %, Todos los datos en %, Todas as informações em %, Alle tal i %, Alla uppgifter i %, Všechny údaje jsou v %, Všetky údaje sú v %, Wszystkie dane w %, Все показаны в %, Tüm bilgiler % cinsindedir, Visi skačai pateikiami procentais (%), Всички данни в %, % كفة البيانات بالنسبة المئوية %

## BMI

Body-Mass-Index, Body Mass Index, Indice de masse corporelle, Indice di massa corporea (BMI), Body-Mass-Index, Índice de masa corporal, Índice de massa corporal, Body Mass Index, Body-Mass-Index, Index tělesné hmotnosti, Index telesnej hmotnosti, Wskaźnik masy ciała, Body-Mass-Index, Vücut Kitle İndeksi, Kūno masės indeksas (KMI), Body Mass Index مؤشر كتلة الجسم

Körperfettanteil, body fat content, taux de graisse, percentuale di grasso, lichaamsvetaandeel, porcentaje de la grasa corporal, percentagem de gordura corporal, fedtprocent, andelen kroppsfett, test zsirtartalmát, zawartość tłuszczu w ciele, podíl tělesného tuku, доля внутреннего жира, vücuttaki yağ oranını, Kūno riebalų procentinė dalis, Процент на телесната мазнина, نسبة الدهون في الجسم


Körperwasseranteil, body water content, taux d'eau dans le corps, percentuale acqua corporei, lichaamswateraandeel, porcentaje del agua corporal, percentagem de água corporal, vandprocent, andelen kroppsvatten, test vızirtartalmát, zawartość wody w ciele, podíl vody v těle, доля тканевой жидкости, ücuttaki su oranını, Kūno vandens procentinis, Процент на телесната вода, نسبة الماء في الجسم

Muskelanteil, muscle proportion, pourcentage de muscles, massa muscolare, spierenandeel, porción de los músculos, percentagem de músculo, muskler andelen, muskellandel, izom arány, masy mięśniowej, podílú svalů, доля мышечной, vücuttaki kas oranını, raumenų procentinė dalis, Процент на мускулите, نسبة العضلات

# A
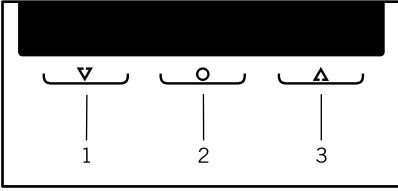
# B

- 
- 
- 
- 

# C

- 
- 
- 
- 
- 
- 
- 

# D

- 
- 
- 
- 
- 
- 
- 
- 
- 
- 

# E

- 
- 
-

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**EN** This device may be used by children over the age of 8 and persons with reduced physical, sensory, or mental abilities or with a lack of experience and knowledge only if such are supervised or have been instructed in the safe use of the device as well as the resulting risks. Cleaning and maintenance shall not be made by children without supervision, Children shall not play with the device.

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**International:** [soehnle.com](http://soehnle.com)

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Apple, iPhone, iPad, iOS, Apple Health and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Android, Google Fit, Google Play and the Google Play logo are trademarks of Google Inc.

The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

## A. Operating elements

1. Minus (▼)
2. Confirm (○)
3. Plus (▲)

Sensor touch:

Only touch the keys softly!

## B. Preparation

1. Insert batteries (3 x 1.5 V AAA).

**Note: Immediately** (as long as „0.0“ is still being displayed) place the body analysis scale on a level surface and wait until the scale powers down automatically.

Then begin with data entry. Otherwise, it might be possible that the weight of the scale would erroneously be added to the initial weighting process.

The same might occur during regular operation if the scale is held in the hand prior to the weighting process.

2. For all measurements, place scale level on a hard surface (not on carpet).
3. Cleaning and care: Clean only using a slightly damp cloth. Do not use solvents or abrasives. Never submerge scale in water.
4. A country-specific option to switch from kg/cm to st/in or lb/in is possible by pushing the key on the bottom of the scale to switch units

## C. Data input

1. Press ○ button.
2. Use ▼/▲ buttons to select memory slot and confirm with ○ . → P
3. Use the ▼/▲ buttons to set your height and confirm with ○ . → H
4. In the next step, set your age and confirm with ○ . → A
5. Set gender and confirm with ○ . → S
6. To finish, set the activity level and confirm with ○ . **(Hold down key for 3 sec.)** → AC

AC	Explanation	Hrs per week
1	Hardly any physical activity	
2	Slight, active, sitting and standing activity	
3	House and garden work, mainly standing, sitting only occasionally	
4	Athletes and people with plenty of movement	Intensive training 5-10 hrs/week
5	Heavy manual labourer, highly active athletes	Intensive training at least 10 hrs/week

### Athlete mode:

By setting the personal activity level, the respective calorie requirement can be precisely determined.

7. Now **immediately** place the body analysis scale on the floor, wait until the screen displays „0.0“ and then step on the scale with bare feet.

The initial measurement is indispensable and serves to store the data for later automatic recognition of persons. If no person steps on the scale, data input must be repeated. The first body analysis is performed after a person steps on the scale.

This scale is equipped with a Time-Out mode. It will automatically switch off after approx. 40 seconds if no keys are pressed.

**Alternative:** You can also perform this step conveniently via the Soehnle Connect app, in which case the data is transferred to the scale. The first measurement is necessary, however. Once you have entered the data on the scale for the recognition of persons as described under C, the Soehnle Connect app will automatically adopt this.

## D. Body analysis

1. Place scale on a solid, level floor and wait until scale switches off automatically.  
Step on the scale with bare feet. Stand still.
2. Body weight and the other calculated data are displayed.
3. If the determined values cannot be clearly allocated to one user, the possible memory slots will be displayed.

If the calculated weight can be allocated to more than 2 users, other possible memory slots will be displayed after 3 seconds.

Use the ▼/▲ keys (▼= left, ▲= right) to confirm the memory slot displayed there.

The scale measures the individual body values by means of bioimpedance analysis (body fat percentage, body water, muscle mass), calculates the basic metabolic rate per day (kcal) and the BMI (Body Mass Index). The 9-segment status bar shows the respective deviation from the normal value (see overleaf).

The scale switches off automatically.

**Note:** If needed, the user can step off the scale as soon as the first value appears on the display.

### Interpretation of measured values

In order to be better able to assess your measured values, at the front of this operating manual you will find an overview table, from which you can take the commonly accepted optimal values.

The optimal range is marked in the table with a 😊. Values outside this range are either too low (–) or too high (+/+/+).

Here is an example:

A female user is 35 years old and is shown a body fat content of 28.2 % from the scale. For the age group of 30-40 year old women, this is a value in the optimal range (22-31 % → 😊).

### Classification of weight according to BMI (source: World Health Organisation, 2016)

BMI (kg/m <sup>2</sup> )	Classification
< 18.5	Underweight
18.5 - 24.9	Normal weight
25.0 - 29.9	Overweight
30.0 - 34.9	Obesity level I
35.0 - 39.9	Obesity level II
> 40.0	Obesity level III

## E. Notifications

1. Replace battery (Lo = low battery)
2. Overloading of scale (max. 180 kg)
3. The measured value is above (H) or below (L) the threshold value specified below.

	L	H	L	H	L	H
♀	< 10 %	> 60 %	< 5 %	> 90 %	< 5 %	> 90 %
♂	< 4 %	> 60 %	< 5 %	> 90 %	< 5 %	> 90 %

## F. Other add-on functions via the Soehnle Connect app

You can connect this body analysis scale with your smartphone (iOS® or Android™) via **Bluetooth®**. Use the Soehnle Connect app for this to make use of more functions among other things:

1. Difference to last weighing
2. Difference to the desired weight

## G. Installation of the Soehnle Connect app

Install the Soehnle Connect app on your smartphone or tablet as follows:

### Android™ smartphones / tablets

If you use an Android device, download and install the Soehnle Connect app from the Google Play Store. To do so, scan the QR code or open the app with the name “Play Store” and look for “Soehnle Connect app”. Then load and install the free app.

### iPhone® / iPad® (iOS®)

If you use an Apple iPhone or iPad, download and install the Soehnle Connect app from the Apple App Store. To do so, scan the QR code or open the app with the name “App Store” on your iPhone and look for “Soehnle Connect app”. Then load and install the free app.

### System requirements for using the Soehnle Connect app and full functionality:

Android™: 4.4.4 or later

iOS®: iOS 9 or later

Bluetooth®: 4.0 or later

### Compatibility

The Soehnle Connect app is compatible with Google Fit® and Apple Health®.

You can find more information on the Internet at [www.soehnle.de](http://www.soehnle.de). Or contact our consumer service.

## H. Connect scale with smartphone / tablet

1. Activate **Bluetooth®** on your smartphone or tablet and start the Soehnle Connect app.
- 2.a. When you start the Soehnle Connect app for the first time, you are automatically guided through the process for connecting a Soehnle appliance. Make sure that the scale is switched on (step onto the scale briefly).
- 2.b. Once you have already used the Soehnle Connect app, go to the connection in the Soehnle Connect app on the “main menu” and press “my devices” and then also “Add device”.

## I. Technical data

**Bluetooth®:** 4.1

Frequency band: 2.4 GHz

Maximum transmission power: 10 mW

Max. 200 kg / 100 g

Max. 440 lb / 0.2 lb

Max. 31 st / 0.2 lb

Age: 17-99 years

Batteries required: 3 x 1.5 V AAA

### Guarantee

Leifheit AG offers a 3 year guarantee for this quality product as of the purchase date (or in case of orders, as of receipt of the goods). Guarantee claims must be asserted within the guarantee period immediately after occurrence of the defect. The guarantee includes product quality.

The following are excluded from the guarantee:

- (1) Wear caused by use or other natural wear;
- (2) Damage caused by improper use or handling (e.g. impact, shock, drop),
- (3) Damage caused by non-compliance with these operating instructions;
- (4) Battery, or rechargeable battery

In case of a guarantee claim, Leifheit, at its own discretion, either repairs the defective components or replaces the product. If a repair cannot be performed and an identical product is no longer available in the product range for the purpose of replacement, you will receive a replacement product of the closest equivalent value. Reimbursement of the purchase price is not possible in case of a guarantee claim. This guarantee furthermore does not grant any claims for damages.

To assert the guarantee, please contact the seller where you have acquired the product and present the defective product and the purchase receipt (copy). This guarantee is valid worldwide.

Your statutory rights, in particular warranty rights, continue to apply and are not limited by this guarantee.

### EC conformity

Soehnle herewith declares that this device is in compliance with the basic requirements and the other pertinent regulations

of the Directive 2014/53/EU. You can find the EU declaration of conformity at [www.soehnle.com](http://www.soehnle.com).

### Battery disposal, EC Directive 2008/12/EC

Batteries are not a part of your regular household waste. You must return batteries to your municipality's public collection or wherever batteries of the respective type are being sold.

Batteries that contain hazardous substances feature this sign:

Pb = Battery contains lead

Cd = Battery contains cadmium

Hg = Battery contains mercury

### Disposal of electric and electronic devices, EC Directive 2002/96/EC

This product is not to be treated as regular household waste but must be returned to a collection point for recycling electric and electronic devices. Further information is available from your municipality, your municipality's waste disposal services, or the retailer where you purchased your product.

### Consumer service

The following contact partners are happy to help you in case of questions and comments:

**Germany** (08 00) 5 34 34 34  
Monday to Friday 08:30 to 12:00

**International** [soehnle.com](http://soehnle.com)

### Legal information

Apple, iPhone, iPad, iOS, Apple Health and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries.

Android, Google Fit, Google Play and the Google Play logo are trademarks of Google Inc.

The **Bluetooth®** word mark and logos are registered trademarks owned by Bluetooth SIG, Inc.

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LEIFHEIT AG