# BRAUN


## Sensian™ 7

Non-contact thermometer

BNT400
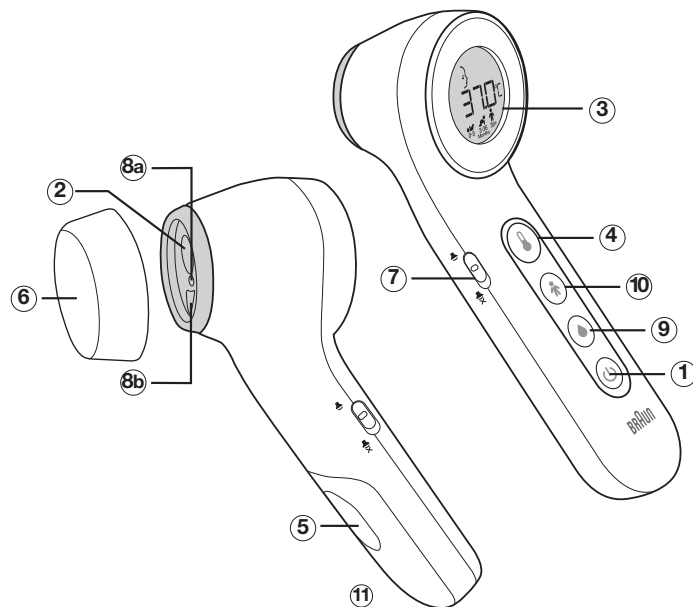
<b>GB</b>	English	2
<b>DE</b>	Deutsch	14
<b>ES</b>	Español	27
<b>FR</b>	Français	40
<b>GR</b>	Ελληνικά	53
<b>IT</b>	Italiano	66
<b>NL</b>	Nederlands	79
<b>PT</b>	Português	92
<b>TR</b>	Türkçe	105

Certain trademarks used under license from The Procter & Gamble Company or its affiliates.

Age Precision™, PositionCheck™ and Sension™ are trademarks of Helen of Troy Limited and/or its affiliates.

**Kaz Europe Sàrl**  
 A Helen of Troy Company  
 Q-Center, Chaux 4  
 CH-1030 Bussigny - Switzerland

© 2022, All rights reserved.

[www.BraunHealthcare.com/uk\\_en/](http://www.BraunHealthcare.com/uk_en/)

Patents: [www.BraunHealthcare.com/us\\_en/patents](http://www.BraunHealthcare.com/us_en/patents)

Made and printed  
 in China

**MD REF** BNT400WE

P/N: A007796R0  
 29NOV22

## Product description (See page 1)

- |                                    |  |
|------------------------------------|--|
| ① Power button                     | ② Infrared temperature scanner                                 |
| ③ LCD screen                       | ④ Temperature button   |
| ⑤ Battery door                     | ⑥ Protective scanner cap                                       |
| ⑦ Silent mode switch               | ⑧a PositionCheck™ aiming light                                 |
| ⑧b PositionCheck™ proximity sensor | ⑨ Surface temperature mode                                     |
| ⑩ Age Precision™ Technology button | ⑪ Fahrenheit to Celsius switch<br>(inside battery compartment) |

Thank you for purchasing the Braun Sensian™ 7 Non-contact forehead thermometer (BNT400). This thermometer is a high quality product incorporating the latest technology and tested in accordance with international standards. With its unique technology, the product can provide a stable, accurate reading with each measurement. The thermometer performs a self-test every time it is switched on for accurate measurement.

**Please read these instructions carefully before using this product and keep the instructions and the thermometer in a safe place.**

## Intended use and contraindications

This No touch + touch thermometer is a handheld, battery-powered, infrared thermometer intended for the intermittent measurement of human body temperature in a home-use environment on people of all ages, except pre-term and small-for-gestational-age babies, using the center of the forehead as the measurement site.

These thermometers are not intended for clinical use in a professional environment and are for home usage only.

This thermometer is not intended to interpret hypothermic temperatures. If the device displays a temperature of 36.4 °C (97.5 °F) or less with an LED backlight color of green or white, and the individual is exhibiting atypical symptoms or behaviors, contact your physician or health care professional.

This thermometer is not intended for pre-term babies or small-for-gestational-age babies. Pre-term is defined as babies with a calendar age of less than 37 weeks. Small-for-gestational-age is defined as a baby born at 37 weeks or later, with a weight below the 10th percentile for babies of the same gestational age.

Use of this thermometer is not intended as a substitute for consultation with your doctor.

## Warnings and precautions

Keep out of reach of children under 12 years.

**NEVER** use the thermometer for purposes other than those for which it has been intended. Please follow the general safety precautions when using on children.

**NEVER** immerse the thermometer into water or other liquids (not waterproof). For cleaning and disinfecting please follow the instructions in the "Care and cleaning" section.

**DO NOT** store this thermometer in temperature extremes below -25 °C or over 55 °C (below -13 °F or over 131 °F) or in excessive humidity (above 95% non-condensing relative humidity).

**If thermometer is stored in a location that is cooler or warmer than where it will be used, allow it to acclimate to the room temperature for 10 minutes before taking a measurement.**

**DO NOT** use the thermometer if there are signs of damage on the scanner or on the thermometer itself. If damaged, **DO NOT** attempt to repair the product.

**NEVER** scratch the sensor with a sharp object when cleaning.

**NEVER** insert a sharp object into the scanner area or any other open surface on the thermometer.

This thermometer consists of high-quality precision parts. **DO NOT** drop the instrument. Protect it from severe impact and shock. **DO NOT** twist the instrument or the measuring sensor.

Temperature elevation may signal a serious illness, especially in neonates and infants, or in adults who are old, frail, or have a weakened immune system. Please seek professional advice immediately when a temperature elevation occurs on persons that are:

- Neonates and infants under 3 months (Consult your physician immediately if the temperature exceeds 37.4 °C [99.4 °F])
- Individuals over 60 years of age (Fever may be blunted or absent in older individuals)
- Individuals having diabetes mellitus or a weakened immune system (e.g., HIV positive, cancer chemotherapy, chronic immunosuppressant treatment, splenectomy)
- Individuals who are bedridden (e.g., nursing home patient, stroke, chronic illness, paraplegia, quadriplegia, surgical recovery)
- A transplant recipient (e.g., liver, heart, lung, kidney).

**DO NOT** allow children to take their temperatures unattended.

Please consult your physician if you see symptoms such as unexplained irritability, vomiting, diarrhea, dehydration, changes in appetite or activity, seizure, muscle pain, shivering, stiff neck, pain when urinating, etc., even in the absence of fever.

Even in the absence of fever, those who exhibit a normal temperature may still need to receive medical attention. People who are on antibiotics, analgesics, or antipyretics should not be assessed solely on temperature readings to determine the severity of their illness.

**DO NOT** modify this equipment without the authorization of the manufacturer.

## Why Braun Sensian™ 7?

### Measurement in under 2 seconds

The innovative infrared technology allows measurement of forehead temperature in under 2 seconds whether you use the contact or non-contact option.

### Accurate and reliable

Due to the unique technology, the Braun Sensian™ 7 thermometer captures the heat naturally given off by the forehead to calculate body temperature value. Get the same professional accuracy whether touching your child's forehead or holding up to 2.5cm away.

### Easy to use

The Braun Sensian™ 7 Non-contact forehead thermometer is non-invasive. A measurement can be taken even while a child is sleeping.

The Braun Sensian™ 7 thermometer is safer to use on a child than a rectal thermometer and is easier to use than other methods.

### Safe and hygienic

Non-contact option helps minimize spreading of germs.

Completely safe for use on children and adults.

## How does Braun Sensian™ 7 work?

The Braun Sensian™ 7 Non-contact forehead thermometer measures infrared energy radiated from the skin at the center of the forehead just above the eyebrows. This captured energy is collected through the lens and converted to an oral equivalent value.

**The Braun Sensian™ 7 thermometer has been clinically tested and proven to be safe and accurate when used in accordance with its operating instruction manual.**

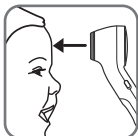
## Age Precision™ Technology

Clinical research shows the definition of fever changes as newborns grow into children and children grow into adults.<sup>1</sup>

Select the appropriate age setting, take a temperature, and the display will glow green, yellow, or red to indicate whether the temperature is in a normal, fever, or high fever range, respectively.

The device is only intended for the measurement of human body temperature. The device also has a surface mode. It is not intended to diagnose or monitor a serious medical condition or disease.

<sup>1</sup> Herzog L., Phillips S. Addressing Concerns About Fever. *Clinical Pediatrics*. 2011; 50 (#5): 383-390.

## Temperature taking hints

It is important to know each individual's normal temperature when they are well. This is the only way to accurately diagnose a fever. Take multiple readings when healthy to determine normal temperature.

A child's normal temperature can be as high as 37.7 °C (99.9 °F) or as low as 36.1 °C (97.0 °F). Be sure to note this unit reads 0.5 °C (0.9 °F) lower than a rectal digital measurement.

**An individual must be inside for 30 minutes before taking a measurement.**

**Note: The thermometer and individual should be in the same ambient temperature for at least 10 minutes.**

**ALWAYS** hold the thermometer and the forehead steady when taking a reading. **DO NOT** move the thermometer until you hear the final beep. **DO NOT** take a measurement while or immediately after nursing a baby.

Individuals should not drink, eat, or be physically active before/while taking the measurement. Remove hats and wait 10 minutes before taking a measurement.

Before taking a measurement, remove dirt or hair from the forehead area. Wait 10 minutes after cleaning before taking measurement.

**ALWAYS** take the temperature exactly as directed. Temperature results may vary if positioned in the wrong location.

For individuals measuring their own temperature, it is recommended to use the contact option instead of non-contact.

In the following situations it is recommended that three temperatures in the same location be taken and the highest one taken as the reading:

- Newborn infants in the first 100 days.
- Children under three years of age with a compromised immune system and for whom the presence or absence of fever is critical.
- When the user is learning how to use the thermometer for the first time until he/she has familiarized himself/herself with the instrument and obtains consistent readings.

### General precautions:

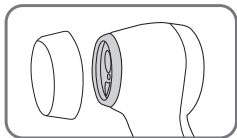
- Remove any sweat prior to measuring with a dry cloth.
- Avoid any cooling or warming cloths on the forehead for at least 30 minutes prior to measurement.
- **DO NOT** take temperature measurements over scar tissue, open sores, or abrasions.
- Keep the measurement sensor and lens clean. Avoid directly touching the sensor or lens.
- **ALWAYS** store the thermometer with the protective cap in place to prevent dirt and scratches from damaging the device.
- Make sure to remove the protective cap before taking a measurement and to put the cap back on after using the thermometer.

## How to use your Braun Sensian™ 7

**Note: The thermometer and Individual should be in the same ambient temperature for at least 10 minutes.**


### 1. Remove protective cap

### 2. Power on

Press and release the power button once. Backlight will come on and the start-up sequence begins.

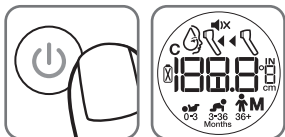NOTE: Make sure to remove the protective scanner cap before taking a measurement.

### 3. Select age range

Select age with Age Precision™ Technology button.

Device will start cycling through the age ranges, highlighting the 3 ranges you may select: 0-3 months (👶), 3-36 months (👦), or 36+ months (👤).

Press the Age Precision™ Technology button to toggle through the age ranges until you select the age appropriate for the individual.

#### 3.1 To set a preferred age range: (Locking an age setting)

To save a preferred age range, so that it will be automatically selected when the thermometer is turned ON, follow these steps:

A. Make sure the thermometer is OFF.

B. Press and hold the power button (⏻) and Age Precision™ Technology button (⊕) for 4 seconds. The unit will turn on, displaying a yellow backlight and the 3 age range icons will flash.

C. Press the Age Precision™ Technology (⊕) button. The age range selected will flash. Press the Age Precision™ button until your desired age range is flashing on the display.

D. Once the desired age range is selected, stop pressing the button and wait 3 seconds. The green backlight will illuminate, and the selected age range will turn solid. After 4 seconds, the thermometer will turn off.

E. The preferred age range will be automatically selected when the thermometer is turned ON.

Note: With a preferred age range, an alternate age range can be selected during each use, by pressing the Age Precision™ button when the age range icon appears, after the thermometer is turned ON.

F. To remove the preferred age range, open the battery door, remove the batteries, wait 30 seconds and re-insert the batteries.

**IMPORTANT:** To ensure accurate readings, you must always take a temperature reading in the age range that corresponds to the individual's age.

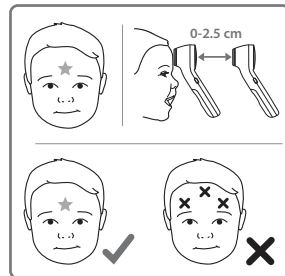
### 4. Position correctly with PositionCheck™

Position thermometer on or up to 2.5 cm away from the center of the forehead, just between the eyebrows.

When the thermometer is ready and on or 2.5 cm away from the forehead a horizontal line of dashes ("----") will appear on the screen.



For non-contact readings, the yellow aiming light will show you where you are aiming.

If the eyebrow area is covered with hair, sweat or dirt, please clean the area beforehand to improve the reading accuracy. It is important to hold the thermometer and the forehead steady during measurement. Movement will impact the temperature reading.


### 5. Take temperature

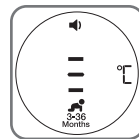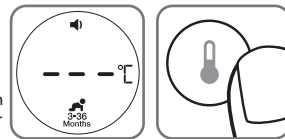Three dashed lines mean the thermometer is ready to take a temperature. Press the temperature button. The display will show a dashed line animation while it takes a reading. Once the reading is complete a confirmation beep will be heard and the appropriate fever light color will be displayed on screen.

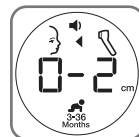If the device is positioned too far away from the forehead, it will prompt you to move closer by displaying the characters "0-2 cm". Slowly move the device toward the forehead until the dashed line animation starts and a reading is displayed.

## 6. Read temperature

If temperature is in the normal range, a long beep will sound when the measurement is complete. If the temperature is in the fever range, ten short beeps will sound to indicate the measurement is complete.

Age range	Green Normal temperature	Yellow Fever	Red High Fever
0*-3 Months	≥35.8 – ≤37.4 °C (≥96.4 – ≤99.4 °F)		>37.4 °C (>99.4 °F)
3-36 Months	≥35.4 – ≤37.6 °C (≥95.7 – ≤99.6 °F)	>37.6 – ≤38.5 °C (>99.6 – ≤101.3 °F)	>38.5 °C (>101.3 °F)
36 Months - adult	≥35.4 – ≤37.7 °C (≥95.7 – ≤99.9 °F)	>37.7 – ≤39.4 °C (>99.9 – ≤103.0 °F)	>39.4 °C (>103.0 °F)

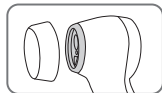
\* 0 months is defined as having a gestational age of greater than or equal to 37 calendar weeks at birth.

## 7. To turn off

Press the power button to turn off. Device will also shut off automatically after 30 seconds of no use.

## How to measure food and liquid temperature


### 1. Remove protective cap

### 2. Power on


Press and release the power button once. Backlight will come on and the start-up sequence begins.

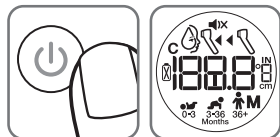NOTE: Make sure to remove the protective scanner cap before taking a measurement.

### 3. Select the surface temperature mode

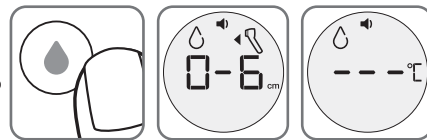Press and release the surface temperature mode button. A beep will sound, the surface mode icon (👉) will be displayed, the screen backlight will turn blue and three dashes will show on the screen indicating the thermometer is ready to take a temperature.

### 4. Position

Position device up to 6 cm away from the liquid or food.

**Note:** When measuring food or liquid temperature the thermometer must be aimed directly at the food or liquid. **DO NOT** aim at sides of container holding the liquids or food. **DO NOT** touch liquid or food with the thermometer.

### 5. Take temperature

Press the temperature button (you can press and hold the button or press and release it.) The display will show a dashed line animation while it takes a reading. Once the reading is complete, the display will show the temperature reading.

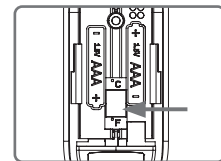
### 6. Read temperature

Remove the device and read the temperature.


## How to change temperature between °C and °F

1. Open battery compartment and remove the batteries.
2. °C / °F switch is accessible in the battery compartment.
3. Slide switch to °C or °F to set preferred temperature scale.
4. Replace the batteries and close the battery door.

## Changing the sound mode

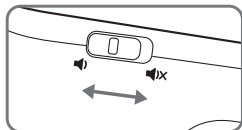The Braun Sensian™ 7 thermometer allows you to silence the beeps on the thermometer to avoid waking a sleeping child. To activate the silence feature, slide the switch to the silent mode position and take a temperature reading. An icon will appear on the screen to indicate that the thermometer is in silent mode.


## Calibration

This product has been calibrated at the time of manufacture. If used according to the use instructions, periodic calibration is not required. This device is not intended as a substitute for regular check-ups by your doctor, please consult your doctor if you have any doubt about the temperature reading.

## Errors and troubleshooting

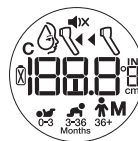Error message	Situation	Solution
	When 20% of battery life is left, the display flashes the low battery warning symbol; however, the device can still work until the battery life has 0% left.	Replace batteries.
	If the steady battery icon is the only symbol shown on the display, the device cannot work. The battery should be replaced immediately.	Replace batteries.
	This message displays when the measured temperature is lower than 34 °C (93.2 °F) or higher than 43 °C (109.4 °F) or when the room temperature is outside the operating range of 15 °C – 40 °C (59 °F – 104 °F).	Re-measure the temperature, carefully following the instructions in "How to use" section.

System error – self-check display flashes continuously and will not be followed by the ready beep and the ready symbol.

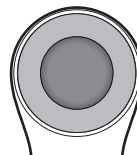
If error persists,

If error still persists,

Wait 1 minute until the thermometer turns off automatically, then turn on again.

... reset the thermometer by removing the batteries and putting them back in.

... please contact Customer Service.

Blank display. Thermometer does not have power.

Please check if the batteries have been loaded correctly. Also check polarity (<-> and <-> ) of batteries.

Contact customer service if thermometer still does not function.

## Care and cleaning

Use an alcohol swab or cotton swab moistened with alcohol (70% isopropyl) to clean the thermometer casing and the measuring probe. Ensure that no liquid enters the interior of the thermometer. Never use abrasive cleaning agents, thinners or benzene for cleaning and never immerse the instrument in water or other cleaning liquids. Wait 10 minutes after cleaning before taking a temperature measurement. Make sure to replace the protective cap after use to prevent scratches or damage from occurring to the thermometer.

Never insert a sharp object into the scanner area or any other open surface on the thermometer.

## Replacing the batteries

The Braun Sensian™ 7 thermometer comes with 2 AAA batteries. Replace with 2 new AAA batteries when the flashing battery symbol appears on the LCD screen. To change the batteries, slide open the battery cover and remove batteries. Replace the batteries being sure to align properly as indicated inside the battery compartment. Remove the battery from the product if it is not required for extended periods of time in order to avoid damage to the thermometer resulting from a leaking battery.

Remove the battery from the instrument if it is not required for extended periods of time in order to avoid damage to the thermometer resulting from a leaking battery.

Only discard empty batteries. They should not be disposed of in the household waste, but at appropriate collection sites or at your retailer.

## 2 Year Limited warranty

Please read all instructions before attempting to use this device. Please retain the receipt as proof of and date of purchase. The receipt must be presented when making any claim within the relevant guarantee period. Any claim under the guarantee will not be valid without a proof of purchase.

Your device is guaranteed for Two Years (2 Years) from date of purchase.

This guarantee covers defects in materials or workmanship that occur under normal use; defective devices meeting these criteria will be replaced free of charge.

The guarantee DOES NOT cover defects or damage resulting from abuse or failure to follow the user instructions. The guarantee becomes void if the device is opened, tampered with, or used with non-Braun branded parts or accessories, or if repairs are undertaken by unauthorized persons.

Accessories and consumables are excluded from any warranty.

For support requests, please visit [www.BraunHealthcare.com/uk\\_en/](http://www.BraunHealthcare.com/uk_en/) or find service contact information at the end of this owner's manual.

This Guarantee is applicable to Europe, Russia, Middle East, and Africa only.

UK Only: This does not affect your consumer statutory rights

## Inquiry

In Europe, any serious incident (e.g., death, life-threatening injury, surgical intervention, etc.) that has occurred in relation to this thermometer should be reported to Kaz Europe (see Helen of Troy consumer service numbers) and to the competent authority of your member state.

REF, LOT, and Serial numbers of your product must be provided as they are essential to record and follow up your inquiry or claim. The LOT and Serial numbers are located in the battery compartment. The production date is given by the LOT number and can be deciphered as explained below:

The first 3 numeric digits in the LOT number represent the day of the year of manufacture. The next 2 numeric digits represent the last two numbers of the calendar year of manufacture and the letter(s) at the end designate the manufacturer of the product. (e.g.: LOT No.: 12313tav this product was made on the day 123, year 2013 at the manufacturer code tav.)

## Product specifications

Type:	Non-contact forehead (BNT400)
Measuring range:	34 °C – 43 °C (93.2 °F – 109.4 °F)
Resolution:	0.1 °C (0.1 °F)
Laboratory Accuracy:	0.2 °C for 35 °C – 42 °C (0.4 °F for 95 °F – 107.6 °F) outside that range ±0.3 °C/0.5 °F (Ambient Temperature: 15 °C to 40 °C (59 °F to 104 °F) This thermometer displays a calculated oral equivalent estimate.
Display:	Liquid Crystal Display, 4 digits plus special icons
Acoustic:	Audio: Normal temp range = Green temp range: 1 long beep for 0.4 second duration Fever = Red or Yellow temp range: 10 short beeps for 0.2 second duration
Operating temperature:	15°C – 40 °C (59 °F – 104 °F)
Automatic Switch-off:	Approx. 30 seconds after last measurement has been taken
Weight:	100 gr (with batteries), 82.5 gr (w/o batteries)

Long term storage ranges

Storage/transport temperature: -25 °C – 55 °C (-13 °F – 131 °F)

Humidity: 15 – 95% non-condensing

Battery: (2) AAA Batteries - at least 500 measurements

Pressure: 700-1060 hPA (0.7-1.06 atm)

Reference Body Site: Oral Equivalent

This medical device carries the CE mark and is manufactured in conformity with RoHS Directive 2011/65/EU and other applicable directive and/or regulations as specified in the EU Declaration of Conformity. This infrared thermometer meets requirements established in ASTM Standard E1965-98 for the thermometer system. Full responsibility for the conformance of the product to the standard is assumed by Kaz Europe Sàrl, Switzerland.

ASTM laboratory accuracy requirements in the display range of 37 to 39°C (98 to 102°F) for IR thermometers is ±0.2°C (±0.4°F), whereas for mercury-in-glass and electronic thermometers, the requirement per ASTM Standards E667-86 and E1112-86 is ±0.1°C (±0.2°F).

Equipment with type BF applied parts

See Instruction for use

Operating temperature

Storage temperature

Operating and storage humidity

Caution

Batch code

Serial number

Catalogue number

Medical device

Manufacturer

Date of manufacture

Authorized representative in the European Community

Internally Powered Equipment

Continuous Operation

**IP22** Protected against solid foreign objects of 12.5 mm diameter and greater. Protected against vertically falling water drops when the device is tilted up to 15 degrees.

MEDICAL ELECTRICAL EQUIPMENT needs special precautions regarding EMC.

For detailed description of EMC requirements please refer to [www.hot-europe.com/support](http://www.hot-europe.com/support).

Portable and mobile RF communications equipment can affect MEDICAL ELECTRICAL EQUIPMENT.

Please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at your local retailer or at appropriate collection points provided in your country.

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## Helen of Troy Consumer Service Numbers

### Algeria

☎ +213 34 18 71 42

### Austria

☎ +43 (0)1 360 277 1225

### Bahrain

☎ +973 17228887

### Belgium/Luxemburg

☎ +32 (0) 2 620 01 01

### Bosnia

☎ +387 (33) 902 911

### Bulgaria

☎ +359 2 492 87 05

### Croatia

☎ +385 (01) 3444 856

### Cyprus

☎ +35722337593

### Cyprus Northern region

☎ +35722337593

### Czech Republic

☎ +420 (0) 2 25 43 97 69

### Denmark

☎ +45 35 15 80 40

### Finland

☎ +358 (0)9 81 71 00 14

### France

☎ +33 (0) 1 85 14 80 95

### Germany

☎ +49 (0) 21 173 749 003

### Greece

☎ +30 211 180 94 56

### Hungary

☎ +36 (06) 1 429 2216

### Iceland

☎ +354 555 3100

### Ireland (The Republic Of)

☎ +353 (0) 1 525 1809

### Israel

☎ +1 800250221

### Italy

☎ +39 02 3859 1183

### Jordan

☎ +962 64168003

### Kuwait

☎ +965 184 0011

### Lebanon

☎ +961 (01) 512002

### Netherlands

☎ +31 (0) 78 201 8001

### Norway

☎ +47 23 50 01 20

### Poland

☎ +48 22 512 39 02

### Portugal

☎ +351 21 060 8045

### Qatar

☎ +974 4 4075048

### Romania

☎ +40 213 520 167

### Saudi-Arabia

☎ +966 920001414

### Slovakia

☎ +421 2 5011 2115

### Slovenia

☎ +386 (01) 888 86 74

### South Africa

☎ ++27 (0) 11 844 6190

### Spain

☎ +34 913 754 176

### Sweden

☎ +46 (0)8 5199 3097

### Switzerland

☎ +41 (0) 22 567 5200

### Tunisia

☎ +216 98 780237

### Turkey

☎ +90 (0) 212 217 66 60

### UAE

☎ +971 4 349 4014

### United Kingdom/ Northern Ireland

☎ +44 (0) 207 949 0115

**Brand:** Braun  
**Category:** Thermometers  
**Model:** **BNT400WE**  
**Artwork Part #:** A007796R0  
**Die Line Part #:** NA

**Subject:** Owner's Manual  
**Region:** EMEA  
**Flat Size:** IN: **W8.5 x H5.5**  
MM: **W215.9 x H139.7**  
**Folded Size:** IN: **W4.25 x H5.5**  
MM: **W107.95 x H139.7**

**Scale:** 1/1  
**Material:** 120g wood free paper (uncoated)  
for Cover,  
80g wood free paper (uncoated)  
for interior

**Page count:** 120  
**Revision:** 10

**Date:** **29NOV22**

**Release Date:** 09NOV22

**Rerelease Date:**

**Colors:**

Dielines (Do not print)

Spot Colors

Cyan  
0%Magenta  
0%Yellow  
0%Black  
100%**Special Instructions:**

Quality Requirement of Artwork and Quality  
**Clarification Process of Artwork Printing:**  
**Meet Eng-QS-06&02**

Helen of Troy  
Creative Services  
Marlborough, MA 01752 USA  
+1 508 490 7000