
Panasonic

Operating Instructions (Household use) Rechargeable Oral Irrigator Model No. EW-DJ86

Please check here for Instructions
and Troubleshooting.

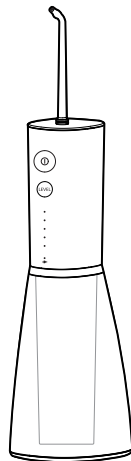Series 800

English	3	Português	47	Česky	89
Deutsch	10	Svenska	54	Slovensky	96
Français	18	Dansk	61	Magyar	103
Italiano	26	Norsk	68	Română	111
Español	33	Suomi	75	Türkçe	118
Nederlands	40	Polski	82		

Panasonic

Operating Instructions (Household use) Rechargeable Oral Irrigator Model No. EW-DJ86

Contents

Safety precautions.....	5	Removing the built-in rechargeable	
Important information	7	battery.....	8
How to use	7	Specifications	9
Battery life	8		

Thank you for purchasing this Panasonic product.

Before operating this unit, please read these instructions completely and save them for future use.

WARNING

- Do not immerse the USB power adaptor and USB cable in water or clean them with water.
Doing so may cause electric shock or fire due to a short circuit.
- Always unplug the adaptor from a household outlet when cleaning the USB power adaptor and USB cable.
Failure to do so may cause electric shock or injury.
- Oral irrigators can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.
- Oral irrigators shall be used only with cold or warm water, or specific solutions as defined in the manufacturer's instructions.
- Children shall not play with the appliance.
- Do not use anything other than the supplied USB power adaptor or USB cable for any purpose. Also, do not use any other product with the USB power adaptor or USB cable.
Doing so may cause burns or fire due to a short circuit.
- The following symbol indicates that a specific detachable power supply unit is required for connecting the electrical appliance to the supply mains. The type reference of power supply unit is marked near the symbol.

Safety precautions

To reduce the risk of injury, loss of life, electric shock, fire, and damage to property, always observe the following safety precautions.

Explanation of symbols

The following symbols are used to classify and describe the level of hazard, injury, and property damage caused when the denotation is disregarded and improper use is performed.

 DANGER	Denotes a potential hazard that will result in serious injury or death.
 WARNING	Denotes a potential hazard that could result in serious injury or death.
 CAUTION	Denotes a hazard that could result in minor injury or property damage.

WARNING

► Preventing accidents

- Do not store within the reach of children or infants. Do not let them use it.
 - Doing so may cause an accident or injury due to accidental ingestion of removable parts.

► About the USB power adaptor and USB cable

- Do not place the USB power adaptor or USB cable over or near water filled sink or bathtub.
- Do not use the appliance when the USB power adaptor or USB cable is damaged, or when the fitting into the household outlet is loose.
- Do not damage, modify, forcefully bend, pull, or twist the cord. Also, do not place anything heavy on or pinch the cord.
 - Doing so may cause electric shock or fire due to a short circuit.
- Do not connect or disconnect the USB power adaptor or USB cable to a household outlet with a wet hand.
 - Doing so may cause electric shock or injury.

- Do not use in a way exceeding the rating of the household outlet or the wiring.
 - Exceeding the rating by connecting too many power plugs to one household outlet may cause fire due to overheating.
- Use the specified USB power adaptor and USB cable.
 - Failure to do so may cause burns or fire due to a short circuit.
- Fully insert the power plug of the USB power adaptor.
 - Failure to do so may cause fire or electric shock.
- Be sure to wipe off any water droplets from the socket before charging.
 - Failure to do so may cause electric shock, injury or fire.

► In case of an abnormality or malfunction

- Immediately stop using if there is an abnormality or malfunction. Unplug the USB power adaptor from the household outlet when using.
 - Failure to do so may cause fire, electric shock or injury.
- <Abnormality or malfunction cases>
 - The main unit, USB power adaptor or USB cable is deformed or abnormally hot.
 - The main unit, USB power adaptor or USB cable smells of burning.
 - There is abnormal sound during use or charging of the main unit.
 - Immediately request inspection or repair to an authorized service centre.

► This product

- This product has a built-in rechargeable battery. Do not throw into fire, apply heat, or charge, use, or leave in a high temperature environment.
 - Doing so may cause overheating, ignition, or explosion of the built-in rechargeable battery.
- Do not modify or repair.
 - Doing so may cause fire, electric shock, or injury.
 - Contact an authorized service centre for repair (battery change, etc.).
- Never disassemble except when disposing of the product.
 - Doing so may cause fire, electric shock, or injury.

WARNING

► **Cleaning**

- Always disconnect the USB plug from the main unit when cleaning the appliance with water.
 - Failure to do so may cause electric shock or injury.
- Regularly clean the power plug and the USB plugs to prevent dust from accumulating.
 - Failure to do so may cause fire due to insulation failure caused by humidity.
- Unplug the USB power adaptor from the household outlet and wipe with a dry cloth.

CAUTION

► **Preventing damage to teeth and gums**

- People who are unable to operate should not use this appliance.
- People with weak oral senses should not use this appliance.
 - Doing so may cause injury or damage to teeth and gums.
- People who may have periodontal disease, have teeth being treated, or are concerned about symptoms within their mouth should consult a dentist before use.
 - Failure to do so may cause injury or damage to teeth and gums.

► **This product**

- Never use for a purpose other than oral cleaning.
 - Doing so may cause an accident or injury when the stream of water is pointed toward the eyes, nose, ears, tongue, or throat.
- Do not fill the tank with water warmer than 40 °C.
 - Doing so may cause burns.
- Do not share the nozzle with your family or other people.
 - Doing so may cause infection or inflammation.

► **Note the following precautions**

- Do not drop or subject to shock.
 - Doing so may cause injury.
- Do not allow needles or any other foreign objects to adhere to the power plug of the USB power adaptor or the USB plug.
 - Doing so may cause electric shock or fire due to a short circuit.

- Do not charge the appliance when the socket is wet.
 - Doing so may cause burns or a fire.
- Do not wrap the cord around the USB power adaptor when storing.
 - Doing so may cause the wire in the cord to break with the load, and may cause fire due to a short circuit.
- Disconnect the USB power adaptor or the USB plug (Type-C) by holding onto it instead of the cord.
 - Failure to do so may cause burns or electric shock due to breakage to the wire of the USB cable, or result in fire due to a short circuit.
- Disconnect the adaptor from the household outlet when not charging.
 - Failure to do so may cause electric shock or fire due to electrical leakage resulting from insulation deterioration.
- After use, clean the appliance and the nozzle and dry them.
 - Failure to do so may cause damage to your health due to propagation of mould and bacteria.

Handling of the removed battery when disposing

DANGER

- The rechargeable battery is exclusively for use with this appliance. Do not use the battery with other products.
- Do not charge the battery after it has been removed from the product.
- Do not throw into fire or apply heat.
- Do not hit, disassemble, modify, or puncture with a nail.
- Do not let the positive and negative terminals of the battery get in contact with each other through metallic objects.
- Do not carry or store the battery together with metallic jewellery such as necklaces and hairpins.
- Do not use or leave the battery where it will be exposed to high temperatures, such as under direct sunlight or near other sources of heat.
- Never peel off the tube.
 - Doing so may cause overheating, ignition, or explosion.

WARNING

- After removing the rechargeable battery, do not keep it within reach of children and infants.
 - The battery will harm the body if accidentally swallowed. If this happens, consult a physician immediately.
- If the battery fluid leaks out and comes in contact with the eyes, do not rub the eyes, and thoroughly rinse with clean water, such as tap water.
 - Failure to do so may cause eye injury. Consult a physician immediately.

CAUTION

- If the battery fluid leaks out and comes in contact with the skin or clothes, rinse off with clean water, such as tap water.
 - Failure to do so may cause inflammation.

Important information

- Please use the appliance together with toothbrush. It is hard to remove tartar or dental calculus attached to the tooth surface just by a water jet, so a toothbrush is required to brush the teeth. We recommend rinsing the inside of your mouth with the water jet after brushing your teeth.
- When you first begin using the oral irrigator, you may experience slight bleeding of the gums even though your gums may be healthy. This is because your gums are receiving stimulation by the oral irrigator for the first time. The bleeding should stop after one or two weeks of use. Bleeding that persists for more than two weeks may be indication of an abnormal condition of your gums, such as periodontal disease. In this case, stop using the oral irrigator and consult your dentist.
- The USB cable provides the overheat protection function. If charging is abnormal or the USB plug (Type-C) heats up, disconnect the power supply, check the socket, and remove any water or foreign objects from the socket.
- Do not put salt, mouthwash or other chemicals into the water tank. Doing so may cause malfunction.
- Do not dry the main unit with a dryer or a fan heater. Doing so may cause malfunction or deformation of the parts.
- Store the appliance in a place with low humidity after use. Leaving the appliance in a bathroom may cause malfunction.
- Charge the appliance correctly according to these Operating Instructions.
- This appliance contains a battery that is only replaceable by skilled persons. Please contact an authorized service centre for the details of repair.

How to use

Please refer to the Quick Start Guide for how to use the appliance.

Battery life

Approximately 3 years. The battery may have reached its operating life if the usage frequency decreases drastically, even after a full charge. (The operating life of the battery may differ depending on usage or storage conditions.)

Removing the built-in rechargeable battery


Remove the built-in rechargeable battery before disposing of the appliance.

The battery is to be disposed of safely.

Please make sure that the battery is disposed of at an officially designated location if there is one.

This figure must only be used when disposing of the appliance, and must not be used to repair it. If you dismantle the appliance yourself, it will no longer be waterproof, which may cause it to malfunction.

- Remove the appliance from the USB cable and USB power adaptor, press the power switch to turn on the power and then keep the power on until the battery is completely discharged.
- Remove the bottom cover.
- Perform steps ① to ⑥ to disassemble the main unit and remove the battery with a screwdriver.
- Insulate the ⊕ and ⊖ terminals by applying tape to them. (To prevent short circuit of terminals).

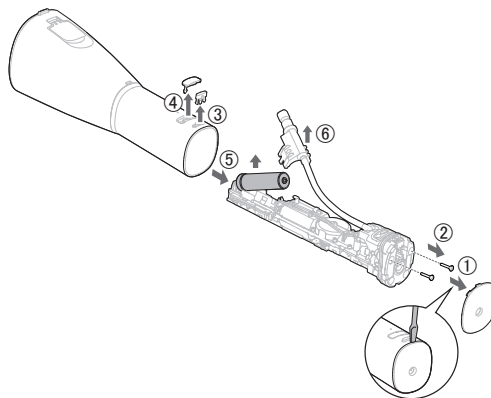For environmental protection and recycling of materials

This appliance contains a Lithium-ion battery.

Please make sure that the battery is disposed of at an officially assigned location, if there is one in your country.

Specifications

Power source	See the name plate on the product (Automatic voltage conversion)
Charging time	Approx. 1 hour
Operating time	LEVEL 5: Approx. 10 minutes* (*at 20 °C when fully charged)
Tank volume	Approx. 200 mL (For approx. 60 seconds of operation when the water pressure level is set to 5)
Dimensions	Main unit: 212 (H) x 68 (W) x 81 (D) mm (Not including nozzle)
Mass	Main unit: Approx. 275 g (Including nozzle)
Airborne Acoustical Noise	71 (dB (A) re 1 pW)

This product is intended for household use only.

For more information about the energy efficiency of the product, please visit our website, www.panasonic.com, and enter the model number in the search box.

Authorized Representative in EU:

Panasonic Marketing Europe GmbH
Panasonic Testing Centre
Winsbergring 15, 22525 Hamburg, Germany

Importer:

Panasonic Marketing Europe GmbH
Hagenauer Strasse 43, 65203 Wiesbaden, Germany

Manufactured by:

Panasonic Corporation
1006, Oaza Kadoma, Kadoma City, Osaka, Japan

Disposal of Old Equipment and Batteries Only for the European Union and countries with recycling systems.

These symbols on the products, packaging, and/or accompanying documents mean that used electrical and electronic products and batteries must not be mixed with general household waste.

For proper treatment, recovery and recycling of old products and used batteries, please take them to applicable collection points in accordance with your national legislation. By disposing of them correctly, you will help to save valuable resources and prevent any potential negative effects on human health and the environment.

For more information about collection and recycling, please contact your local authority. Penalties may be applicable for incorrect disposal of this waste, in accordance with national legislation.

Note for the battery symbol (bottom symbol)

This symbol might be used in combination with a chemical symbol. In this case, it complies with the requirement set by the Directive for the chemical involved.

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This image shows a blank sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

Panasonic Corporation
<https://www.panasonic.com>
© Panasonic Corporation 2024

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(Guangzhou) Co., Ltd.