LANAFORM

# Evolutionary Scale

2 IN 1 SCALE

HOLD

TARE

PRECISION

VARIATION

REMOVABLE  
TRAY

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

EVOLUTIVNA TEHTNICA

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
VAGA ZA PRAČENJE  
RAZVOJA


# 1
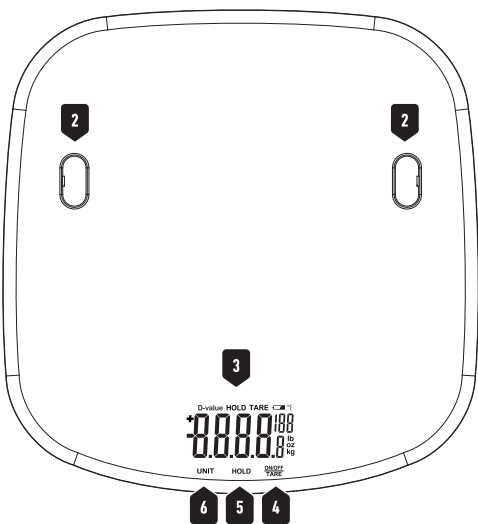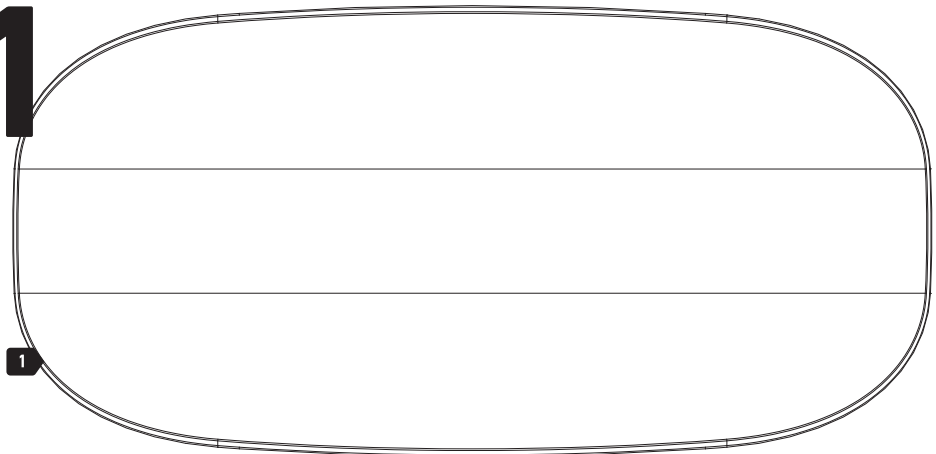
# 2
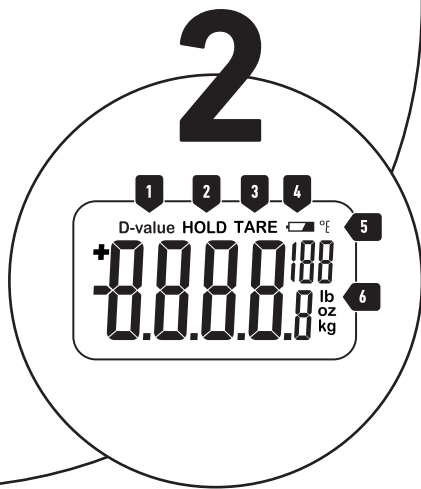
# 3
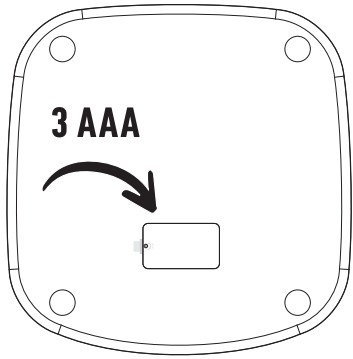
## ENGLISH

### CONTENTS OF THE PACKAGING

- 1 Evolutive Scale with detachable baby tray
- 2 silicone plugs
- 1 instruction handbook
- 3 AAA batteries

### PRECAUTIONARY STATEMENT

▲ **Please read all the instructions before using your product, particularly these basic safety instructions.**

■ The photographs and other images of the product in this manual and on the packaging are as accurate as possible, but do not claim to be a perfect reproduction of the product.

#### ▲ Danger

- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
- Never leave the baby unsupervised when using the Evolutive Scale.
- Do not immerse the device in water or any other liquid.
- To avoid electrocution, do not try to repair the device yourself.
- When not in use, the scale switches off automatically after a few minutes of being switched on. To switch off the scale more quickly, press the on/off button for about half a second.

#### ▲ Warning

- Keep the Evolutive Scale away from heat sources such as stoves, radiators, etc.
- Do not use the device near wet areas, such as the bath, sink etc.

#### ▲ Precautions

- Only use the Evolutive Scale for its intended purpose, as described in these instructions.
- The scale is not suitable for people weighing over 150kg. In the event of overloading, the "EEEE" Error message will be displayed on the screen. Overloading may cause damage to the scale.

- Do not use this device in a room in which aerosols (sprays) are used or in a room in which oxygen is administered.
- When using the scale, keep it away from sources of electromagnetic waves such as microwave ovens, mobile telephones and wireless networks, as this may affect the measurement.
- In the event of electrostatic failure, install the batteries again to resume normal performance.
- Never leave anything on the scale platform when not in use.

### INTRODUCTION

Thank you for choosing the LANAFORM Evolutive Scale. The accurate and easy to use digital Evolutive Scale allows you to weigh your baby to an accuracy of 5 grams. Once your baby becomes too big for the tray, the baby scale can be converted into scale for children and adults.

### DESCRIPTION OF THE DEVICE ■ 1

- 1 Baby tray
- 2 Silicone plugs
- 3 Screen
- 4 ON/OFF button – TARE weight function
- 5 HOLD function
- 6 Choice of weight-measuring unit (kg/lb/oz)

#### Screen ■ 2

- 1 D value
- 2 HOLD
- 3 TARE
- 4 Battery
- 5 Temperature
- 6 Unit of measurement

#### Features

- Scale capacity: up to 150 kg/330 lb/5921 oz
- Weight measurement division:  $\pm 0.005\text{kg}$  (0.05kg-20kg) |  $\pm 0.05\text{kg}$  (20kg-150kg)
- Weight measuring accuracy: 0-5kg:  $\pm 60\text{g}$ , 5-20kg:  $\pm 150\text{g}$ , 20-50kg:  $\pm 250\text{g}$ , 50-100kg:  $\pm 300\text{g}$ , 100-150kg:  $\pm 400\text{g}$
- LCD screen
- Low battery alert: the "Lo" icon appears when the battery is low
- Automatically switches off
- Batteries: 3 x AAA batteries (supplied)
- TARE weight function (deducts the initial weight)
- HOLD function (stabilises the weight measurement)
- Function D value (comparison with the last weight measured)
- Temperature function
- Screen size: 66 (L) x 34 (W) mm

- Device size: 285 (L) x 290 (W) x 25 (H) mm
- Tray size: 550 (L) x 269 (W) x 60 (H) mm
- Net weight of device: 1.5 kg
- Storage environment: -20°C to 60°C (humidity: 25%-85%)
- Operating area: 0°C to 40°C (humidity: 20%-95%)

### PREPARATION PRIOR TO USE

Open the cardboard packaging.

Leave the table and the detachable baby tray, leaving the protective film behind.

Now remove all the protective film.

Check that the device is not damaged.

### AREA OF USE

Always place the Evolutive Scale on a hard, flat and level surface.

### INSTRUCTIONS FOR USE

#### Inserting and replacing batteries ■ 3

- 1 Remove the cover of the battery compartment located underneath the scale.
- 2 Remove the plastic protection inserted above the 3 AAA batteries or insert 3 new AAA batteries if necessary.
- 3 Replace the cover.

#### UNIT SELECTOR FUNCTION

- Weight  
The scale is set to the "kg" unit on delivery. Short press the UNIT button to select the kg/lb/oz weight whenever you require it.
- Temperature  
Long press the UNIT button to switch the temperature unit between C°/F°.

#### NORMAL WEIGHING FUNCTION

Before putting the baby on the scale, make sure that the surface of the scale is completely dry and not slippery. You can also place a towel on the scale.

- 1 Press the "ON/OFF" button to switch on the scale.
- 2 0.0kg (lb or oz depending on the weight unit selected) will be displayed on the screen.
- 3 Put the baby gently on the scale, taking care of its head. If the baby moves too much, select the "Hold" function which will set the baby's weight when they are more still.
- 4 The weight measured is displayed for 5 seconds followed by the D value icon and the difference between the actual weight and the previous weight is displayed twice in rotation.

ⓘ **Note: The difference may be positive in the event of weight gain or negative in the event of weight loss.**

- When the baby is too heavy for the baby scale, remove the detachable baby tray and use the table as shown above.

## WEIGHING FUNCTION WITH A TOWEL

- Press the "ON/OFF" button to switch on the scale.
- 0.0kg (lb or oz depending on the weight unit selected) will be displayed on the screen.
- Place a towel on the scale. The weight of the towel will be displayed.
- When the "Hold" icon appears and stops flashing, press the "Tare" weight button.
- The "Tare" icon and "0.0 kg" (lb or oz depending on the weight unit selected) will be displayed on the screen again.

ⓘ **Note: Leave the towel on the scale. If you take it off the scale, the screen will show a negative weight.**

- Put the naked baby gently on the scale. Repeat the operation with the "Hold" button as explained above in the "normal weighing" section, if necessary.

## CLEANING

Do not wash the device and do not submerge it in water or any other liquid.

To clean it by hand, used a damp cloth or sponge and a small quantity of liquid detergent, if necessary.

Do not use abrasive detergents.

## STORAGE

Store the device in a dry, safe place. Do not put any object on the scale when they are not in use. Remove the batteries if the device is not to be used for a long time.

## ⚠ BREAKDOWNS AND SOLUTIONS

ⓘ **Nothing is displayed on the screen.**

- ⓘ Contact failure at battery level or batteries are not charged.
- ⓘ Make sure that the batteries have been installed correctly or change them.

## ⚠ ADVICE ON WASTE DISPOSAL

The packaging consists entirely of environmentally-friendly materials that can be taken to your local waste sorting centre to be used as secondary materials. Cardboard can be placed in a paper collection container. Packing film must be taken to your local sorting and recycling centre.

When you are no longer using the device, dispose of it in an environmentally responsible manner and in compliance with the statutory requirements.

## LIMITED WARRANTY

LANAFORM guarantees that this product will be free from any material and workmanship defects for a period of two years from the date it was purchased, other than as specified below.

The LANAFORM warranty does not cover damage caused by normal wear and tear on this product. In addition, the warranty on this LANAFORM product does not cover damage caused as a result of any abusive, inappropriate or incorrect use, accident, attachment of any unauthorised accessory, modification of the product or any other situation of whatever kind, beyond the control of LANAFORM.

LANAFORM shall not be held liable for any type of incidental, consequential or special damage.

Any implied warranties of the product's fitness for purpose are limited to a period of two years from the original purchase date on presentation of a copy of the proof of purchase.

Upon receipt, LANAFORM will repair or replace your device, as appropriate, and return it to you. The warranty may only be used at the LANAFORM Service Centre. Any maintenance and repair work carried out on this product by any person other than the LANAFORM Service Centre will invalidate this warranty.

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