
NOOANCE
PARIS

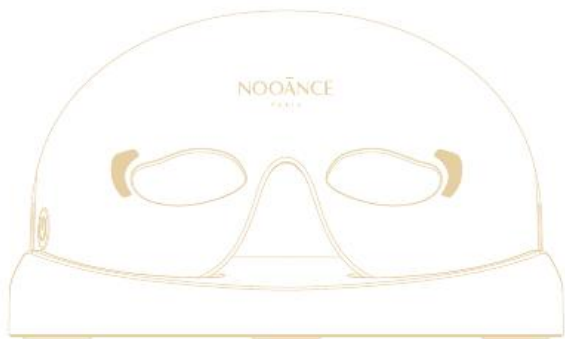MASQUE LED YEUX
LED EYE MASK

X200

fr- Français	02-21
en- English	22-41
de- Deutsch	42-61
es- Español	62-81
it- Italiano	82-101
ar- العربية	102-121
pl- Polski	122-141
nl- Nederlands	142-161
cs- Čeština	162-181
sv- Svenska	182-201
da- Dansk	202-221
pt- Português	222-241
ro- Română	242-261
bg- български	262-281
el- Ελληνικά	282-301
hr- Hrvatski	302-321
tr- Türkçe	322-341

SOMMAIRE

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Table of Contents

Recommendations for Use	23
Safety Instructions	
Warnings	
Precautions	
Cleaning	
Contraindications	
Potential side effects	
Disclaimer	29
Product Representation	30
Customer Service	30
The Science of Photobiomodulation	30
Product Information	31
Device Specifications	
Contents of the box	
Product Overview	33
Usage	34
Before First Use	
First Use and Preparation	
User Guide	35
Frequency of Use	
Storing	
Warranty	
Air travel	
Disposal	
Troubleshooting	38
Symbols	39

Recommendations for Use

Safety instructions

- Read all instructions carefully before using this device.
- All instructions and this manual are available on the www.nooance-paris.com website or upon request to customer service at: bonjour@nooance.com
- DO NOT charge or use the device near water or if it is at risk of coming into contact with water. Do not immerse the device.
- DO NOT use the device if it is not working properly, has been dropped, or damaged. DO NOT use the device if the timer is faulty. For warranty and repairs, contact customer service.
- DO NOT attempt to repair the device. Repairing or opening the device will void your warranty. The device contains no user-serviceable parts. This device contains non-replaceable batteries.
- DO NOT perform maintenance or disassemble the product. In case of malfunction or damage, contact customer service at bonjour@nooance.com
- DO NOT use the device if it is damaged or if the cable is damaged. If any part of the cable is damaged, stop using it. Contact customer service at bonjour@nooance.com to replace the cable.
- The device can be charged using either the included USB-C cable or the charging base. In both cases, a Safety Extra-Low Voltage of 5V DC must be used to charge (not provided).

- DO NOT charge with a damaged or defective power adapter.
- DO NOT excessively bend or mistreat the USB charging cable. You could damage the device.
- Please refer to the 'Box Contents' table for additional instructions regarding each component of the product.

⚠ Warnings

- DO NOT let children manipulate or play with the device. Cleaning and charging the device should not be done by children.
- This device can be used by individuals with reduced physical, sensory, or mental abilities, or those lacking experience or knowledge, provided they are properly supervised or have been instructed on the safe use of the device and understand the risks involved.
- DO NOT use the device if it is wet. If the device has come into contact with water, do not touch it without unplugging it first.
- STOP using the device and consult your doctor if you experience any adverse reactions.
- STOP using the device if you experience symptoms of low blood pressure (dizziness or light-headedness) during use.
- STOP using the device if it overheats. Turn off the device immediately.

- STOP session if the device feels hot or if your skin emits an uncomfortable warmth.
- DO NOT point the LED light directly into your eyes or another person's eyes. You may experience temporary visual impairment.
- DO NOT use the device with unauthorized accessories.
- DO NOT use the device for purposes other than those specified in this manual.
- DO NOT use the device excessively – Respect the recommended conditions of use and frequency of use.
- DO NOT use the device if your skin is sunburned. Wait until your skin heals before using it.
- DO NOT use the device on painful areas. Consult your doctor if you experience pain during use.
- DO NOT use photosensitizing active ingredients or abrasive products (such as scrubs) before using the mask.
- No modifications to this equipment are allowed.
- Contact your doctor if you have any medical questions regarding the use of this device.

Precautions

- DO NOT scratch or damage the light-emitting surface.
- Your warranty will be void if you do not use the device in accordance with the instructions in the manual.
- DO NOT use the device outdoors to avoid damage.
- DO NOT place the device in direct sunlight.
- DO NOT sit, lie down, or step on the device to avoid breaking it.
- DO NOT use or store the device near hot surfaces.
- Store the device according to the protocol described in this manual.
- The device is designed for cosmetic use only by healthy individuals. Please contact your doctor if you are unsure of your health condition.
- The device's battery is not removable by the user. If the device needs to be disposed of, do not burn it. Dispose of the device in accordance with local regulations.
- Inspect the device and do not use it if you see signs of damage, such as battery leaks. Contact our customer service.
- Consult your doctor before using the device if you have a fever, acute inflammation, severe circulatory problems, light sensitivity, anaesthesia, or desensitized areas.
- We strongly recommend exercising caution when taking medications. If in doubt, consult your doctor.

Note: For hygiene reasons, we advise not sharing your device with others without disinfecting the part in contact with the skin using 70% alcohol.

Cleaning

- Clean the device with a soft, dry cloth.
- DO NOT use solvents or abrasive cleaners other than 70% alcohol on the device.
- The device should not be cleaned or charged by children

The mask is reusable, and it needs to be cleaned after each use. Follow the instructions below to clean your device.

- 1) Before cleaning, please turn off the device.
- 2) Use a clean cotton cloth moistened with water to clean the casing and ensure that no liquid enters the interior of the device. Stubborn stains can be removed by wiping gently with 70% alcohol.
- 3) Visually check the cleanliness of the device and goggles, make sure there are no dust, impurities, or any other pollutant, otherwise, please repeat the above step again.
- 4) Drying: please dry the device in the air naturally.

Note: Never immerse the device in water or other cleaning liquids to avoid damage to the device.

Contraindications

- DO NOT use the device if you are receiving medical treatment for conditions involving photosensitivity or if you are taking medication that causes photosensitivity.
- DO NOT use the device if you have precancerous lesions, active cancer or a cancerous lesion.
- DO NOT use the device on areas with phlebitis, thrombophlebitis, or varicose veins.
- DO NOT use the device if you have a known allergy or sensibility to silicone. If in doubt about allergies, contact your doctor before using the device.
- DO NOT use the device if you have a medical history of seizures or epilepsy triggered by light.
- DO NOT use the device on individuals under 18 years of age.
- DO NOT use the device if you are pregnant, planning to become pregnant, or breastfeeding.
- DO NOT use the device if you experience headaches caused by light.
- DO NOT use the device to target rosacea, moles, warts, open wounds, or any other skin condition without medical advice.
- DO NOT use the device on depigmented spots or areas of skin affected by vitiligo.
- DO NOT use the device if you have open wounds.
- DO NOT use the device on bacterial infections as this could worsen the symptoms.
- DO NOT use the device if you have had recent surgery or other surgical procedures.

- DO NOT use the device if you are hypersensitive to light. Perform a light sensitivity test before using the device (refer to section “Before first use”)

Potential side effects

- Immediately after your session, the skin may appear red. The skin should return to normal within 1 to 2 hours.
- Temporary headaches may occur after light session. Consult your doctor if headaches persist for an extended period and contact the customer service at: bonjour@nooance.com
- If unexpected side effects, such as itching, occur within 48 hours following the first session of using the device, consult your doctor before continuing the exposure and contact the customer service at: bonjour@nooance.com

Disclaimer

- This product is not intended to diagnose, treat, or prevent any disease. NOOANCE makes no claims or warranties regarding this product's ability to treat any mental, physical, or dermatological condition. Consult a qualified healthcare professional for any condition requiring medical attention.
- This product is not intended for use in the diagnosis of diseases and other conditions, nor for the cure, mitigation, treatment, or prevention of diseases, and is not intended to alter the structure or function of the body.
- Aesthetic, non-medical device intended for home use.

Product Representation

NOOANCE continually improves its products using the latest technologies. As a result, the actual product may differ slightly from the product shown in the images on the packaging and in this user guide.

Customer Service

For more information on using or maintaining the device, or to obtain additional information, please contact Customer Service via email at bonjour@nooance.com.

The Science Of Photobiomodulation

According to the literature, the skin has the remarkable ability to absorb light energy. This energy is associated, in the scientific literature, with increased cellular energy production (ATP), supporting cells in carrying out their natural functions. This process is also associated with supporting fibroblasts in their natural production of collagen and elastin, and with natural cell renewal processes.

The NOOANCE EYE Mask translates this science through 200 LEDs with optimal precision (± 2 nanometres), with two wavelengths selected for their complementary properties as documented in the scientific literature:


- 100 red LEDs at 633 nm
- 100 near-infrared LEDs at 830 nm

Product Information


Device Specifications

Product name	LED Eye Mask X200
Model	EY-80,EY-80A, EY-80B, EY-80C, EY-80D, EY-80E
LED wavelengths	633 nm $\pm$ 2 nm (Red) 830 nm $\pm$ 2 nm (Infrared)
Max output optical power density	50mW/cm ²
Energy delivered	30J/cm ²
Lamp classification	Exempt group (Risk group according IEC 62471: 2006)
Operating temperature	0°C to 40°C 32°F to 104°F
Storage and transport temperature	- 5°C to 55°C 23°F to 131°F
Operating humidity	30% to 95% without condensation
Storage and transport humidity	20% to 95% without condensation
Charging supply / base	5V~2A
Charging port / base	USB-C
Charging cable / base	USB-C to USB-C
Charging cable length	1 meter / 3.28 feet
Battery	3.7V, 700mAh
Battery type	Lithium-ion
Charging time	120 minutes/ 2 hours
Recommended use	3 * 10 minutes per week
Timer	10 minutes
Dimension of the device	189.8*102*107 mm
Weight of the device	0.365 kg / 12.89 oz

Contents Of The Box

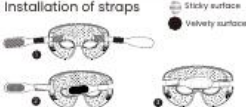Name	Package contents	Additional information
EYE MASK		● DO NOT clean the appliance when it is switched on.
Charging Base		● Only for charging your device EYE MASK.
Adjustable Straps		● Installation of straps 
Removable eye protections		● Use the ocular protection to avoid exposing your eyes to the light. ● The red light of the LEDs may appear orange after a few seconds of use as the eye adapts to this brightness. This is completely normal.
USB Cable		● Use the USB-C cable to charge the eye mask or to connect it to the charging base.
Manual		● Read this manual carefully before using your mask.

Product Overview

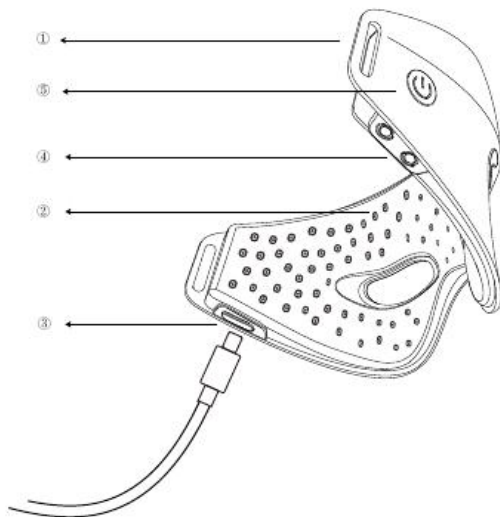- 1) Adjustable Straps Connection Port
- 2) Led -Emitting Surface
- 3) USB Cable Connection Port
- 4) Charging port for the charging base
- 5) Power ON/OFF Button

Usage

Before First Use :

We strongly recommend that users carry out a light sensitivity test before first use.

1. Place the device on the inside of the forearm.
2. Switch on the device to start your session.
3. Turn the mask off after 5 minutes (the timer is set for 10 minutes by default).


Visually inspect the treated area. If your skin feels excessively warm or itchy during the session, stop immediately. If the area turns red and this redness persists for more than 2 hours, your skin may be hypersensitive to light. Discontinue use and contact customer service at: bonjour@nooance.com.

First Use and Preparation :

1) Unpack the package and check off all the parts described in section « Contents of the box ».

Your product may be partially or fully charged upon arrival. To verify the battery level, simply plug it : a steady indicator light confirms a sufficient charge, whereas a blinking light indicates that the device requires further charging.

2) Use the USB-C cable provided in the package and connect it to a 5V Safety extra-low voltage adaptor to plug the charging base. Put the mask on the base to charge it. When fully charged, the mask indicator lights will stop flashing and change to a constant white light.

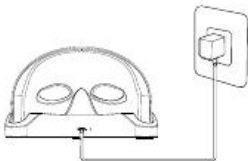
3) When the mask battery level is less than 10%, the indicator lights flash white. To recharge the device, follow the instruction above.

User Guide

Preparation: Before starting the session, please thoroughly cleanse your skin.

Insert the Removable eye protection: Align the narrow edge of the eye protection with the eye socket of the mask's inside, and push it into the outside of the mask. The broader edge of the eye protection should be facing the user's skin.

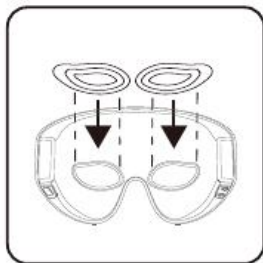Connect The Adjustable Straps

Place the mask over your face: Apply the mask to your face, and adjust the straps to ensure that the mask sits well on your face / desired position.

Tip: To optimise the results of your Nooance mask, we recommend that you use one of our aqueous serums recommended on our website www.nooance-paris.com. At the end of the session, you can also apply a moisturising cream or serum.

Frequency Of Use

We recommend applying the mask for 10 minutes, 3 times a week. For best results, use consistently as directed. Individual results may vary.

Storing

- 1) Avoid direct sunlight on the device while storing.
- 2) It is recommended to store the device in its original packaging.
- 3) Store in an environment with a temperature of $-5-55^{\circ}\text{C}$, a humidity of 20-95%RH and an atmospheric pressure of 700hPa-1060hPa.

Warranty

- The device comes with a 24-month guarantee for the original purchaser, subject to proof of purchase from an authorised dealer. This guarantee covers normal use with approved accessories. In the event of a fault, contact our customer service department.
- Please note that the guarantee does not cover damage caused by accidents, incorrect use, use of non-recommended accessories, or poor electrical connection

Air Travel

The device is allowed in the cabin. The main airlines refuse to accept lithium batteries in the hold. If in doubt consult your airline.

Disposal

DO NOT dispose of the device with household waste. The device should be recycled. The device contains a Lithium-ion battery. For protection of the environment, please safely dispose of the device following your local guidelines.

Troubleshooting

If the actions suggested fail to fix the problem, or the problem reoccurs, and your device is within the warranty period please contact our customer service by writing to bonjour@nooance.com.

⚠ DO NOT attempt to investigate any further.

⚠ DO NOT attempt to modify the device in any way. The device has NO parts you can repair yourself.

Problem	Cause	Solution
The device/base doesn't switch on	The device is not fully charged, because the charging cable was not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	The device is not fully charged, because the charging cable was not properly inserted into the device.	Insert the charging cable properly.
	The power outlet didn't work during charging.	Try a different 5V adapter or power outlet.
The device doesn't switch on	The on/ off button has not been pushed 3 full seconds.	Push the on/off button for 3 seconds straight.
The device/base switches off unexpectedly.	The battery level is too low.	Charge the device.

Symbols

Symbol	Explanation
	Indicates that the device is CE marked and authorized for sale in Europe in accordance with applicable safety and conformity assessment standards (refer to the EU Declaration of Conformity, available upon request, for the list of relevant standards).
	The Waste of Electrical and Electronic Equipment (WEEE), indicates the device may not be discarded as household waste (2012/19/EU).
	Indicates that the electromagnetic radiation from the device is below the limits specified by the Federal Communications Commission and the manufacturer has followed the requirements of the Supplier's Declaration of Conformity authorization procedures.
	This warning symbol on the device label is a warning to read the enclosed documents before use and take all safety precautions required.
	Read operator's manual.
	Date of Manufacture.
	This symbol means keep the device dry.
	This symbol indicates the serial number of the device.
	This symbol refers to the batch code of the device.
	This symbol indicates the USB connector.

	This symbol means that the device uses a Lithium-ion Polymer battery.
	This symbol means keep the device away from direct sunlight.
	This symbol means the upper and lower operating temperatures that the device can be safely used at.
	This symbol indicates the upper and lower operating humidity that the device can be safely used at.
	This symbol means that the cardboard packaging that the device came in can be recycled.
	Protected against solid foreign objects of 12.5 mm and greater, and Protection against vertically falling water drops when ENCLOSURE tilted up to 15°.
	Unique Device Identifier.
	Indicates that the device is UKCA marked and authorized for sale in the United Kingdom in accordance with applicable safety and conformity assessment standards (refer to the UKCA Declaration of Conformity, available upon request, for the list of relevant standards).
	Indicates that the device is tested according to RoHS requirements.
	To identify equipment meeting the safety requirements specified for class III according to IEC 61140.
	Manufacturer

	Country of manufacture
	Importer
	EU Representative

Notice version A0 – 2026/04/11

Legal manufacturer
NOOÂNCE
179 boulevard Haussmann, 75008 Paris, France
bonjour@nooance.com

Country of manufacture
Made in China

Importer
NOOÂNCE
179 boulevard Haussmann, 75008 Paris, France
bonjour@nooance.com

EU Representative
Company: NOOÂNCE
Address: 179 boulevard Haussmann, 75008 Paris, France
Email: bonjour@nooance.com

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AB Temsilcisi