absolute

by  EPI LADY®

Dispositif d'épilation au laser
 Laser hair removal device
 Laser-Haarentfernungsgerät
 Dispositivo de depilación láser
 Dispositivo laser per la depilazione
 Laser ontharing apparaat

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WARNINGS

-  Please observe the following instructions when using the appliance:
 - This sign anywhere in the brochure indicates safety and risks information, which should be read carefully and understood fully.
 - This appliance is exclusively for home use.
 - This appliance is intended for removal of hairs from the human body. Any other use, in particular removal of hair from the head, is absolutely forbidden.
- 
 - Use only the USB cable and / or adapter supplied with the device.
 - Do not drop the appliance.
 - Do not switch on the appliance if the appliance shows visible signs of damage.
 - Do not maintain the appliance while it is switch on.
 - This appliance can be used by adult persons from 18 years old. Children shall not play with the appliance.
 - This appliance can be used by adult persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance.
- 
 - Although the device is equipped with a safety capacitance sensor, do not point the device towards the eyes, yours or anyone else's. Do not look directly into the device.
 - Never open the appliance. The manufacturer's guarantee does not cover any damage caused by mishandling.
 - Keep the device out of reach of children.
 - Never submerge or wet the appliance.
-  **ATTENTION-** Do not use this appliance near bathtubs, showers, basins or other vessels containing water.

SAFETY CHART

 Check if you can use the appliance. Women 18 and older; use below the nose line

1. Check your Skin Tone:

Match your skin tone of the area(s) you wish to treat to the skin tones shown below to determine if you can safely use this device. The appliance is intended for use only on Type V (dark brown) to Type VI (black) skin tones.

2. Check your Hair Color:

The device works best on naturally dark hair (black, black brown, medium brown, or light brown). The device will work but will not be as effective on white, gray, blond or red hairs.

3. Where you can and cannot use this appliance:


<p>WHERE TO USE THIS APPLIANCE ✓ : The appliance can be used on the following areas of the body:</p> <ul style="list-style-type: none"> ✓ Women's Facial Hair (cheeks and below) ✓ Underarms ✓ Hands and Feet ✓ Bikini Line ✓ Arms and Legs ✓ Back 	<p>WHERE NOT TO USE THIS APPLIANCE ✗ : The appliance is not safe for use on the following areas of the body:</p> <ul style="list-style-type: none"> ✗ Eye Areas (eyes, eyebrows, and eyelashes) ✗ Ears ✗ Nipples ✗ Genitals ✗ Anus
--	---

4. Who should not use the appliance? Contraindications.


DO NOT use the appliance:

-  On or around the eyes, eyebrows, or eyelashes, anywhere on or around the ears, nipples, genitals, or around the anus. These areas may have more sensitive skin, a darker skin tone, and using the appliance there may cause skin injury or skin pigmentation changes.
-  On tattoos, dark brown or black spots (such as dark freckles, birthmarks, moles or warts). Tattoos or dark spots may absorb too much light energy which could cause skin swelling, temporary or permanent blisters or discoloration of the skin or whitening of the skin.
- If you have epilepsy or light-induced seizures or migraines.
- If you have a history of keloidal scar formation, a known sensitivity to light (photosensitivity) or are taking medication that makes the skin more sensitive to light, including non-steroidal anti-inflammatory agents, (e.g., aspirins, ibuprofens, cetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulfonyluraes, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha-Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A, Accutane and/or topical retinoids.

- If you have any skin condition in the treatment area, including acne, sores, psoriasis, eczema, burns, herpes simplex, wounds or a rash. These conditions can make the skin more sensitive to the light treatment and could make the treatment uncomfortable or cause skin swelling, blisters or infection. Wait for the affected area to heal before using. If you are not sure if your skin condition has healed enough to use the device, then please consult with your doctor.
- If you have been treated with Accutane® (isotretinoin) within the past 6 months. This can make the skin more sensitive to the light treatment and could make the treatment uncomfortable or cause skin swelling, blisters or infection.
- Do not use the device if you are allergic to hydroquinone or other skin-whitening substances.
- DO NOT use if you are under the age of 18. The effects of using the device in this age group are unknown.
- DO NOT use on any area where you may want your hair to grow back.
- DO NOT use on a man's face, jaw or neck.
- DO NOT use if you are pregnant, trying to become pregnant (conceiving) or breast-feeding. The device has not been tested on pregnant or breast-feeding women. Hormonal changes could increase skin sensitivity and the risk of skin damage.

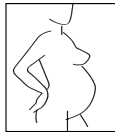
Not safe skin tones


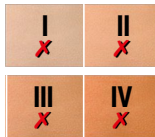
DO NOT use the appliance if your skin tone matches one listed as "Not Safe". If you use the appliance on medium brown, dark brown or white skin tones - all indicated as "Not Safe" for treatment - the skin may absorb too much light, resulting in burning, discoloring, or scarring.

HOW TO USE YOUR APPLIANCE

How it works: the laser beam impacts the hair and the follicle ("root") under the skin surface. On impact, the beam heats the follicle, weakening it and impeding its future growth. To be more effective, treat each hair follicle during its period of active growth. The hairs grow in cycles, it will be necessary to carry out several treatments. The hair falls off after a certain time, otherwise you can depilate it once the treatments are done.


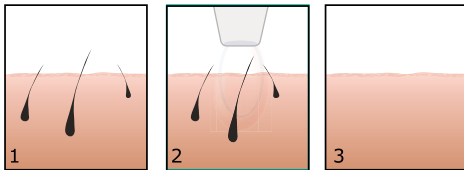
1. Ensure your skin is clean with no oily products and completely dry before to proceed.
2. For the first time; you should test the compatibility of the device with your skin. The purpose of the test is to find out whether your skin reacts naturally well to laser treatment.
3. Choose a concealed area of your body, such as a thigh.
4. **Do not** remove hair using an epilator/tweezers/wax **before** the treatment. It is better if the hair is visible.
5. If the hair is long, it should be trimmed to 1-2 mm, using scissors or a hair trimmer.
6. Use the power button and turn on the device. The **green** indicator light will flash. Aim the Laser-beam opening to the treated hair.
7. Attach the Laser-beam opening to your skin. The laser automatically will operate and will treat the hair for 1 second.
8. When you place the laser beam window against your skin, the device beeps and the **green / red** LED flashes, and you will feel a slight sensation of heat on the affected area for the second time needed for treatment. If you remove the laser beam window from the area to be treated too early, the device will alert you by emitting vibrations.
9. After 5 minutes check the treated area.
10. If the treated spot appears red or irritated, please wait 24 hours and repeat the test.
11. Otherwise consult your doctor.
12. After you made sure that your skin reacts well to laser, you can use it as follows:
 - Turn the device on in the same manner as for the test.
 - Place the device against the skin on the hair to be treated and treat the hair the same manner as for the test.
 - Always turn off the device after use.
13. After the 8 weeks treatment, the follicle will be weakened and fall down. It can take some times, depending on each body and skin (from 3 to 14 days). Between 2 treatments, if you feel uncomfortable with the remaining hair, do not hesitate to shave or epilate them. It will not affect the laser efficiency and result of the treatment.
14. Use the device once or twice a week for 8 to 12 weeks. In order to optimize your treatment, we recommend that you perform three laser shots on the hair to be treated continuously.
15. Using the device may cause slight tingling sensation. This means that the beam "hit" the follicle. However, if the use of the device is accompanied by pain, a burning sensation, a burn, swelling of the skin, or any other unusual reaction, stop using the device. Once again, do not hesitate to shave or epilate the hair. It will not affect the laser efficiency and result of the treatment.
16. Then make spot repairs or additional treatments if necessary. When the device is used according to the instructions in this manual, it is a device for personal use that provides permanent hair removal. The results vary among people and depend on which areas of the body are treated and how well the treatment instructions are followed.

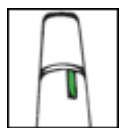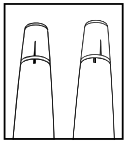BATTERY OPERATION

Attention:

1. Turns blinking red when the device is connected to the power adaptor and is charging.
2. Turns orange when the device is connected to the power adaptor and ready to use.
3. When the device is connected to the power adaptor and is charging, the device will not operate and you cannot treat yourself.
4. Turns green-red alternately when the device is operating and targets the laser to the hair.
5. The port is located at the bottom of the operating handle. One end of the cable connects to this port, while the other connects to the USB port on the power adaptor

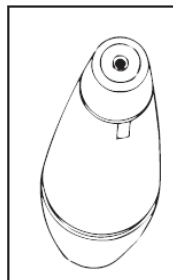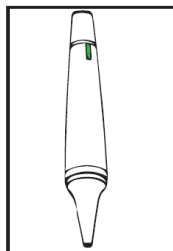CLEANING

- Clean the appliance after each use.
- Clean the treatment window and the outside surface of the system using ONLY a dry soft cloth.
- Hold the device with the laser-beam opening facing downwards, and using circular motions, clean the opening. Do not insert the swab into the device, but only clean the outer rim which touched the skin.
- Do not wet the appliance.
- Store the appliance in a cool and dry place AWAY from the reach of children.

Operation and storage	
Temperature (C°)	Operation 10° - 40° Storage -40° - 70°
Relative humidity	Operation 15%-90% Storage 10%-90%
Atmospheric Pressure (kPa)	1000-700

TROUBLESHOOTING


- No indicator lights are on: Use the power button and turn on the device and check the device is charged.
- No Impulse: The skin contact surface is not in contact with the skin in the treatment area. Place correctly the device for the skin contact surface to be in contact with your skin.

WARNING

- (1) Never let children use the appliance.
- (2) Never plug the appliance with wet hands.
- (3) Do not attempt to repair your device.
- (4) Always switch the appliance off before putting it down, even if it is only for a moment.
- (5) DO NOT use perfumed products for several hours before or after a treatment.
- (6) Always unplug the appliance after use.
- (7) Never submerge or wet the appliance/adaptor.
- (8) Do not insert metal objects into the air grilles to avoid electric shock.

GUARANTEE

The warranty shall apply in the country in which the product was sold via an authorized dealer. If the product should become defective within the warranty period, please contact the dealer where you bought the product. The warranty does not include damage to the product by accident or misuse, abuse, alteration to the product or use inconsistent with the technical and/or safety instructions required. This warranty shall not apply if the product has been dismantled or repaired by a non-authorized person.

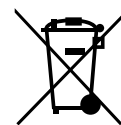This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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