YBRUSH

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FasTeesh SAS, 49, Avenue Rosa Parks, 69009 Lyon

E-mail : contact@y-brush.com

Tel : +33 7 83 47 71 58

Https: www.Y-Brush.com

USER GUIDE (EN)

IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

READ ALL INSTRUCTIONS BEFORE USING

DANGER – To reduce the risk of electrocution:

1. Do not use while bathing.
2. Do not place or store product where it can fall or be pulled into a tub or sink.
3. Do not place in or drop into water or other liquid.
4. Do not reach for a product that has fallen into water. Unplug immediately.

WARNING

To reduce the risk of burns, electrocution, fire, or injury to persons:

1. Close supervision is necessary when this product is used by, on, or near children or invalids.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.
4. Keep the cord away from heated surfaces.
5. Never use while sleeping or drowsy.
6. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

1. SAVE THESE INSTRUCTIONS

Operating instructions

Please read the operating instructions carefully before using the device for the first time. Toothbrushes may be used by children and persons with reduced physical, sensory or mental abilities or insufficient experience or knowledge, provided that they are supervised or instructed in the safe use of the appliance and understand the potential hazards. Cleaning and maintenance by the user must not be performed by unsupervised children. The device should be used only with the power supply or USB charger provided with the product and the brush sold by Y-Brush. Connect the charger only to an AC outlet (100–240V). The product should only be operated with the safety extra-low voltage corresponding to the marking on the product. Children should not play with the product.

» Introduction

Y-Brush is a device designed and developed by oral health specialists to facilitate dental brushing for both children and adults. It allows you to brush your teeth in a simple, effective way, and in a few seconds thanks to a patented system.

This device is composed of a flexible brush that fits perfectly in the jaw. Thanks to the nylon fibers lining the brush and the vibrations of the handle, the teeth are brushed simultaneously and efficiently.

Remember that the recommendations in terms of oral health are to brush your teeth at least twice a day, in the morning and in the evening. An effective brushing helps to prevent the consequences related to a bad oral health. (periodontal, cardiovascular diseases, etc.) It helps to remove plaque, but is not intended to remove tartar deposits. It is recommended to see a dental surgeon every six months before using this device.

» Directions for use

+ Charging Device

- Before each use of the Y-Brush, make sure it is fully charged. In normal use, the indicator lights on the handle are lit in white. They will light in red when the battery level is too low and in green when the Y-Brush is plugged into the power supply.

- If the battery level of the handle is too low:

- Connect the charger to a 110-220V power outlet or use the USB charger
- Plug the charging cable into the part provided for this purpose on the handle.

An illuminated halo around the on/off button lights up cyclically and indicates that the device is being charged. If the halo stays on green continuously, it indicates that the handle is fully charged. A period of 8 hours is sufficient to fully charge the handle. It will then have an autonomy of up to 3 months at a rate of 2 brushes per day.

- The handle is ready to be used. It is imperative to close the cover that protects the connection of the recharging cable on the handle side to guarantee its watertightness.

- The handle cannot be used during charging.

» Warning

+ The Brush

- The appliance and its brush have the sole function of brushing teeth, they should not be used for anything else. Do not attempt to swallow all or part of the brush.

- The handle must not be inserted in the mouth without a brush.

- Make sure that the nylon fibers of the brush reach all areas of all teeth once in the mouth. If any area are not reached, brush them separately.

USER GUIDE (EN)

- Be aware of the risk of pinching when the brush is attached to the handle.
- The brushes are individual and have a maximum service life of 6 months under normal use. If there is any doubt as to whether the brush belongs to a third party and/or can be used by a third party, change the brush to avoid the risk of cross-contamination.
- Replace the brush after severe impact (even if there is no visible sign of damage).

+Power Supply & Connection

- Use the device at room temperature (0–40°C). Avoid prolonged exposure of the product to sunlight or any other heat source. Do not use the charger near water (bathtub, washbasin, etc.) or outdoors. Never immerse or rinse the charger.
 - Never touch the power cable with sharp or pointed objects. Never touch the plug with wet hands, and never unplug by pulling on the cable. Regularly check the cable for cracks or damage.
 - Be sure to close the cover that protects the handle charging plug connection so that it remains watertight. The handle is only splash-proof.
 - Do not place in a microwave oven or any other household appliance.
 - After use, rinse the brush head thoroughly, shake off excess water and store it on the brush holder with the connector facing up.
 - Do not dispose of the product in the household garbage. The handle has batteries and should therefore be taken to a collection point for electronic products. The charger must also be taken to such a collection point. Do not attempt to remove the batteries.

+Other important instructions

- CAUTION: Not suitable for children under 4 years old. The elements, in both the adult and child version, contains small parts that may cause a choking hazard. Therefore, they must not be used by a child without adult supervision.
 - Make sure that the use of the device is not unpleasant and does not impair dental health. If in doubt, contact a dentist. If you have dental problems, wear braces, have dental implants or crowns, consult your dentist or physician before using the product.
 - Do not open the product. Any warranty is void if the product is opened. Improper repair is dangerous for the user.
 - Never use a product that has been opened or damaged.
 - Do not wash the product in the dishwasher.
 - People with pacemakers or other electronic implants should consult their physician beforehand.
 - The product should not be used for commercial purposes.
 - Keep the packaging out of the reach of children.
 - The manufacturer is not liable for any failure to observe these operating instructions, for repairs that are not carried out in accordance with the state of the art, or for improper use of the device.

USER GUIDE (EN)

- In the event of a serious incident, please contact the manufacturer immediately (contact below).
- Do not use the device while showering or bathing.
- Dry the handle with a towel if it's wet.

FasTeesH SAS, 49, Avenue Rosa Parks, 69009 Lyon

E-mail : contact@y-brush.com

Tel : +33 7 83 47 71 58

Https: www.Y-Brush.com

. ATTENTION .

Please note that the warranty does not apply if water is found in the handle or if the handle has not been used according to these instructions. It is therefore important to close the connector cover and not to leave the handle within reach of a sink, shower or bath.

- This device complies with part of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesirable operation.

- Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

- Note: This equipment has been tested and found to comply with the limits for a class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.
- This digital apparatus does not exceed the class B limits for radio noise emissions from digital apparatus set out in the interference-causing equipment standard entitled: "Digital Apparatus" ICES-003 OF THE Canadian Department of Communications.
- The device complies with Industry Canada License-exempt RSS standard(s).
- Operation is subject to the following two conditions:
 - This device may not cause interference.
 - This device must accept any interference, including interference that may cause undesired operation of the device.

USER GUIDE (EN)

2. CONTENTS


1 BRUSH | 1 HANDLE | 1 USB CHARGER

WITH ACCESSORIES (OPTIONAL):

1 STORAGE STAND | 1 TOOTHPASTE APPLICATOR (KID OR ADULT)
1 TRAVEL CASE | 1 BOX OF TOOTHPASTE TABLETS | 1 TOILETRY BAG

3. FEATURES

Take a tour around your new sonic toothbrush

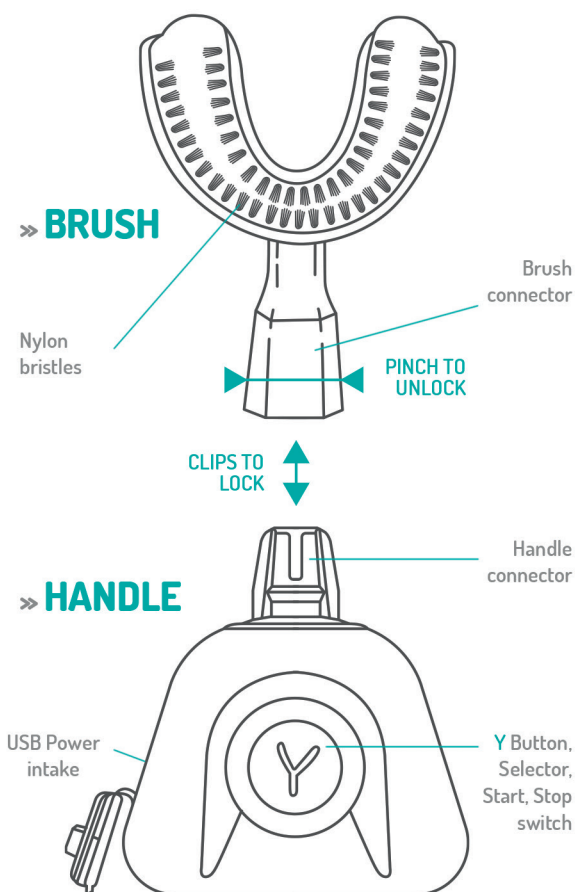4. FIRST BRUSHINGS

Y-Brush is a new automated toothbrushing device, different from an electric toothbrush. For this reason, you will need **1 to 2 weeks of learning and adaptation**, and the feeling of freshness may be different at start. Follow the guide !

» **WEEK 1** . Start with **Y-Brush**

Start with our **Soft Mode**. It is a double cycle of **15 seconds** that will allow you to familiarize with the new feeling of your brushing ritual. We recommend you to alternate the use of this product with your regular brush for the first days.

With a new design comes new habits:

See next pages to get more insights on how to use your Y-Brush for the first time.

» **WEEK 2** . Switch to **10 Seconds**

Try our **Normal Mode**. It is a double cycle of **10 seconds** a little more intense than our Soft Mode. The Normal mode is the recommended cycle for your daily oral care routine.

» **WEEK 3** . Try the **Intense mode**

Brush all your teeth thanks to our **fastest mode**. With its double cycle of **5 seconds**, you will complete your oral care ritual in only 10 seconds when you are in a hurry.

Do you need help ? Write us at customers@y-brush.com

Our services will make their best to ease your brushing experience. We will answer within 24H (Monday to Friday).


5. INSTRUCTIONS

4 steps to fast & effective toothbrushing

» 5.1. Toothpaste

If you don't have a **Toothpaste Applicator** ^(optional), put a dab of toothpaste on your tongue and run it over all your teeth.

Y-Brush Applicator: Spread a thin line of toothpaste along the entire length of the brush. You can also use our **Toothpaste Tablets** ^(optional). In all cases, **wet the brush thoroughly after you put your toothpaste.**

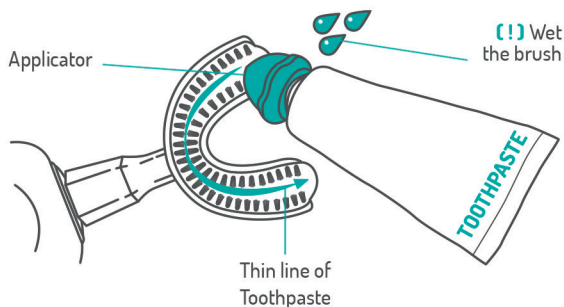» 5.2. Vibration modes


Press the "Y" button to **power ON** the brush. Each new pressure will change the duration of your cycle of brushing. The LED lights around the switch indicate your cycle selection:

1 Light: Soft 15s | **2 Lights:** Normal 10s | **3 Lights:** Intense 5s


» 5.3. Do the twist

Insert the brush on the upper teeth and hold the button. Until the vibrations stop, chew gently & quickly, and twist the handle from left to right. Repeat this operation for the bottom teeth.

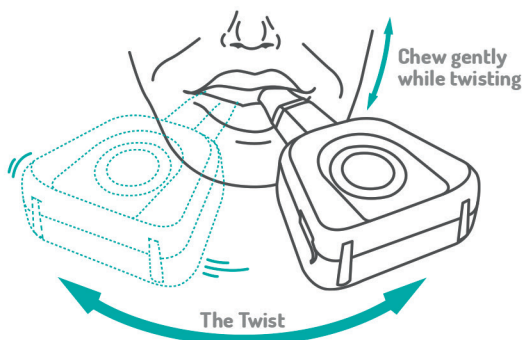» 5.4. Cleaning & Storage

When you have finished to brush all your teeth, pinch the base of the brush to unlock it and **rinse** your brush with **water**. Feel free to use a finger to remove the excess of toothpaste between the filaments. Let it dry on its **stand** ^(optional) or head down on the edge of a glass in a dry and ventilated room.

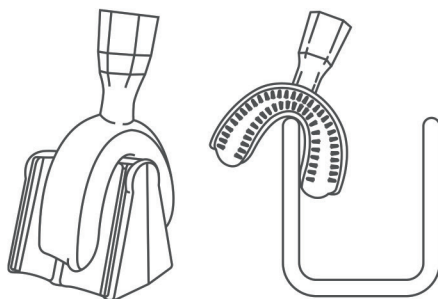6 . Frequently Asked Questions

Why is it important to chew and rotate the brush ?

The chewing gesture allows to brush teeth on their whole height as well as the gums, and the rotating motion the last molars. It is therefore very important to follow this procedure.

I don't have the same feeling of freshness as with a classic toothbrush.

It is very important to follow the procedure (gestures, etc.) and to give 1 to 2 weeks to Y-Brush to get used to this new feeling.

The brush filaments get crushed quickly.

When you chew too hard, the filaments, usually at the bottom and at the molar level, can be damaged quickly. Try to chew less hard but just as quickly as before for a regular use of the device.

I bleed slightly after using the Y-Brush, is this normal ?

Don't worry, this can happen with any brush. Y-Brush is removing plaque from bottom to top, but also rubs your gingivas. Make sure all your teeth are well inserted in the brush every time you use Y-Brush. If you keep bleeding after a week of use, please consult your dentist.

Still need some help ?

Our customer service is here for you

customers@y-brush.com

Find more information on our support page

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Discover our video tutorials
by scanning this code, or visit
y-brush.com/support

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[@ybrush_france](https://www.instagram.com/ybrush_france)