# Y BRUSH

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WATER FLOSSER

FR

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EN

USER GUIDE

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## IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

### READ ALL INSTRUCTIONS BEFORE USING.

**DANGERS**—To reduce the risk of electrocution:

1. Do not use while bathing.
2. Do not place or store product where it can fall or be pulled into a tub, shower stall or sink.
3. Do not place in or drop into water or other liquid.
4. Do not reach for a product that has fallen into water. Unplug immediately.
5. Do not handle charger with wet hands when charging.

**WARNING**—To reduce the risk of burns, electrocution, fire, or injury to persons.

1. A product should never be left unattended when plugged in.  
Exception: Contact lens disinfectors, toilet seat assemblies, and similar products.
2. Close supervision is necessary when this product is used by, on, or near children or invalids.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.
5. Keep the cord away from heated surfaces.
6. Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
7. Never use while sleeping or drowsy.
8. Never drop or insert any object into any opening or hose.
9. Unplug this product before filling. Fill (reservoir) with water only unless otherwise specified by manufacturer. Do not overfill (or specify filling instructions)
10. This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance (5V).

11. This appliance contains battery that are non-replaceable.
12. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge under supervision.  
Cleaning and user maintenance shall not be made by children without supervision. Children shall not play with the appliance.
13. This appliance shall be used only with cold or warm water(max 40°C)
14. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

## SAVE THESE INSTRUCTIONS

### FCC STATEMENT

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications to this unit not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

## Safety Measures

Basic safety precautions should be observed during the use of electrical products, especially when children are present.

Please read all safety guidelines before using this device.

**Warning: To reduce the risk of electric shock and other hazards, please observe the safety guidelines below:**


1. Check if the charger or power cord are damaged before initial and during continued use of this device. If either are damaged, do not attempt to use the device. Instead, contact after-sales services for a replacement.
2. Do not throw the device into a fire or other heated environments.
3. Do not immerse the device into water, as this may cause damage.
4. Do not attempt to plug or unplug the charger while your hands are wet. When plugging or unplugging the device, do so using the charger, not the power cord.
5. Do not modify, damage, or destroy the charging equipment. Do not forcefully stretch, bend, or twist the power cord. Do not place heavy objects on the power cord or stick the power cord between other objects.
6. Do not drop or insert foreign objects into any of the device openings.
7. Do not place the device, charger, or power cord on heated surfaces.
8. Do not place the device in a hot or humid location, nor any locations exposed to direct sunlight.
9. Before cleaning and maintenance, ensure that the device is turned off and not charging.

**Warning:** To reduce the risk of burns, electric shock, fires, and other injury hazards, please observe the safety guidelines below.

1. Do not use this device in any manner other than those outlined in the manual.
2. Do not use this device while a Tip is not installed.
3. Persons (including children) without the full mental capacity or relevant experience and knowledge must clean and maintain this device under the supervision of a guardian.
4. Place this device outside of the reach of children and infants. Do not allow this device to be used as a toy.
5. Users suffering from periodontitis or those who have undergone oral surgery within the past 2 months should consult a doctor before using this device.
6. Only the Tips and accessories included with, sold with, or recommended by this device should be used.
7. Do not disassemble or swap out parts from this device.
8. To prevent burns, do not fill the tank or attempt to clean your mouth with water over 40°C.
9. Do not use this device for any purpose other than oral hygiene.
10. To prevent injury, do not press the Tip forcefully into your gums or teeth.
11. Do not use this device in outdoor locations where aerosols (sprays) or oxygen supplies are in use.
12. Remove all cosmetic teeth accessories before using this device.
13. Do not use iodine, salt water, or concentrated essential oils that are insoluble in water with this device. Doing so will reduce its performance and shorten its service life.
14. Do not fill the device with pure mouthwash liquid. If you wish to use a mouthwash, dilute it in accordance with the proper proportions.

Store this manual for future reference!

## I. Pack Content

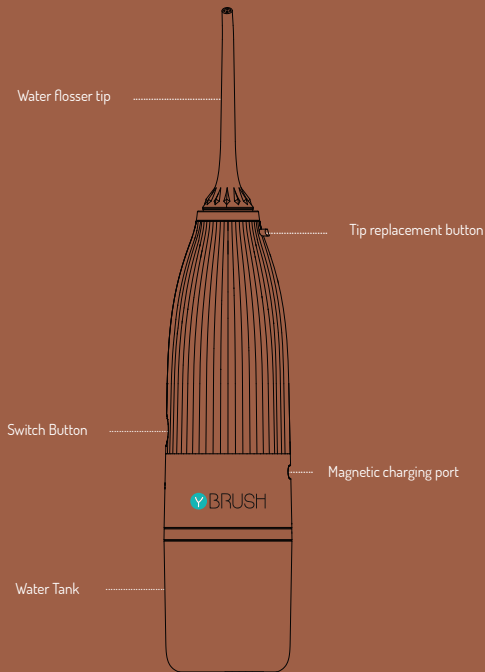
## II. Accessories and Specifications

### Product Accessories:

Water flosser tips x2

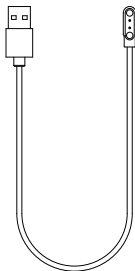Power Cable x1

### Product Specifications:

Product Name		Y-BRUSH Oral Irrigator
Oral Irrigator	Product No.	FC3840
	Battery Rated Voltage	3.7V ===
	Power Rating	5W
	Waterproof Rating	IPX7
	Charge Time	2-3H
	Input	5V === 1A
	Adaptor Rating Value	Input: 100-240V-50/60Hz 0.3A Output: 5V === 1A 5W
Place of origin	China	


## III. Product Use

**Before Starting: Very Important!** Ensure the device is fully charged before using it for the first time by connecting the USB cable to the charging port and turning on the power. This device cannot be used while charging.

### Step 1: Fill the Water Tank

Remove the water tank by rotating it clockwise and pulling downwards, then reattach the tank by aligning it with the snap joint and rotating counterclockwise.

Note: Ensure that the device is turned off before filling the water tank.
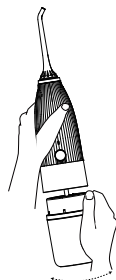
### Step 2: Install the Tip

Insert the tip into the hole on the top of the controller. You should hear a "click" sound on successful installation. Rotate the tip to adjust the direction of water outflow. If you need to replace the irrigator, press and hold the uninstall button, then remove the irrigator.


Installation

Replacement

### Step 3: Select Mode

This device can only operate in a single mode.
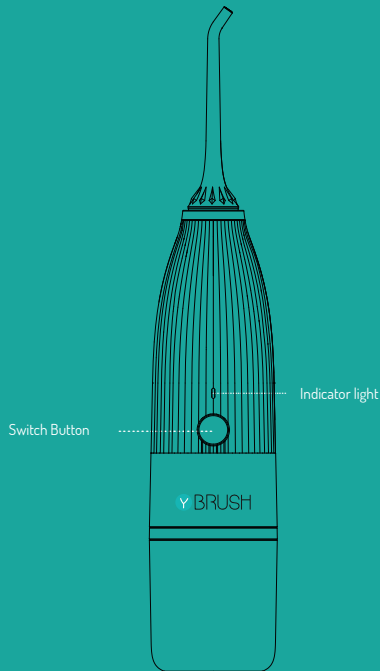
### Step 4: Before Use

Before using the device, lean slightly over the sink, extend the Tip into your mouth and aim it at the teeth, gums, or between the teeth. Rotate the tip to adjust the direction of the water outflow so that it is perpendicular to the part of your mouth you wish to clean. Ensure that your mouth is open enough so that any water flows right out.

#### Note:


You may experience itching or minor gum bleeding during use. This is normal and should typically stop after 1-2 weeks.

### Step 5: Use

Press and hold the switch button for 1 second to start the device. To pause the device during use, tap the switch button. Press and hold the button again to resume operation. This device is equipped with a 120-second shutdown timer. Automatic shutdown will only trigger if the water tank is empty. If you wish to keep using the device, tap the switch button to turn it off, then restart the device after refilling the water tank. Turn off the device after you are finished using it.


#### Usage tips :

For best results, start with the molars (posterior teeth) and start cleaning towards the front teeth. While cleaning, move the tip 90° along the gum line. Stay between each tooth for a short period. Close your lips slightly to avoid the water from spraying out, but leave enough space so that water can flow freely from your mouth to the sink.
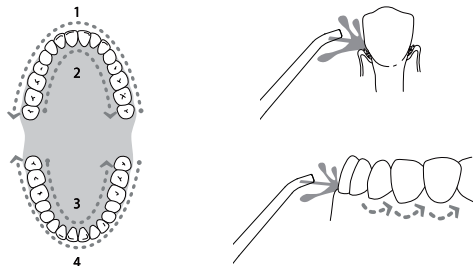
### Step 6: After Use

After use, open the water tank and pour out any remaining liquid. Turn on the device and drain any liquid left until the irrigator no longer sprays out any liquid, then turn off the device.
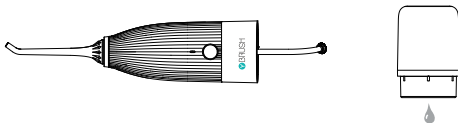
### Step 7: Charging

The indicator light is white, and will remain on while the device is in use. When the device is on low charge, the indicator light will flash for 5 seconds, after which the device will turn off. When the device is turned on while on low charge, the indicator light will flash for 3 seconds and then the device will turn off, which indicates that charging is required. While the device is charging, the indicator light will gradually darken over time. While the device is fully charged, the indicator light will remain on.

## IV. Cleaning and Maintenance

### Cleaning the Body:

1. Clean this device only with clean water or a neutral detergent. Do not use corrosive or abrasive detergents. Do not use water above 50 C to clean any part of this device, as this may lead to damage or reduced performance.
2. Do not soak the oral irrigator in water. If necessary, clean the surface of the device by wiping with a soft cloth dipped in mild, non-corrosive detergent.

### Cleaning the Tank:

1. Remove the water tank.
2. Rinse the inside of the tank with clean room temperature water. Do not clean the tank with any type of detergent.
3. Reinstall the water tank after it has been cleaned.

### Removing Hard Water:

1. Depending on the mineral content of the water used, hard water deposits may build up on the device. Failure to remove hard water buildups will reduce product performance.
2. Point the Tip at the sink, turn on the device, then let it run until all water has been used. Fill the water tank with clean warm water, then repeat the above steps to clean the tank. This cleaning process should be repeated every 1-3 months to maintain product performance.

### Maintenance:

1. Ensure that the battery is fully charged before initial use.
2. Ensure that the battery is fully charged before storing the device for an extended period.
3. In order to prolong battery life, charge the device immediately when the battery is low.

## VI. Troubleshooting

Problem	Cause	Solution
There is a leak between the tip and the body	The seal of the Tip holder is damaged, or the tip was not installed in the proper position	Contact the vendor for repairs, or reinstall the tip in the proper position
There is a leak in the water inlet	The seal is damaged, or there are impurities on the inlet	Contact the vendor for repairs
There is a leak in the water tank	The seal is damaged	Replace the seal, or contact the vendor for repairs
There is not enough water pressure	Not enough charge	Try again after charging the device
The device cannot start	Not enough charge	Try again after charging the device or contact the vendor for repairs
The device is not releasing any water	The device is being used for the first time, or has not been used in an extended period	Add a suitable amount of water to the tip installation hole and try again after reinstalling the tip.

If you were not able to solve the problem with the troubleshooting instructions above, please contact the vendor you purchased the device from or visit a designated after-sales repair location.

## VII. After-Sales Services

Thank you for purchasing a Y-BRUSH Water Flosser. Please keep the below in mind.

1. When you purchase this device, ask the vendor to provide a valid sales certificate to serve as your warranty.
2. In accordance with national regulations, this device includes a 2-year free warranty period for failures not attributable to user error. This warranty does not include Tips. If the warranty period is exceeded, the device can be sent to designated after-sales locations for paid maintenance.
3. Charges and service fees may still be applied during the warranty period in the below circumstances.
  - The device was damaged due to improper use, maintenance, or storage by the consumer.
  - The device was damaged by a service agent that has not undertaken the "three guarantees".
  - The sales voucher is not valid or already two years over the expiration date.
  - The device was damaged during commercial use.
  - The device was damaged by a force majeure event (such as lightning, earthquake, or fire).

Email : [customers@y-brush.com](mailto:customers@y-brush.com)

N° : + 33 7 83 47 71 58

Site : [www.y-brush.com](http://www.y-brush.com)

We reserve the right to final interpretation in case of any typographical errors or misunderstandings found in the contents of this manual. Device upgrades and updates will be incorporated into a new version of the manual without notice. The appearance and color of the product is subject to the packaged model, which shall take precedence over any illustrations of descriptions found in this manual.

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---

Find more about our products  
on: [www.y-brush.com](http://www.y-brush.com)

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### ANY QUESTIONS?

Write us at:  
[customers@y-brush.com](mailto:customers@y-brush.com)

Our customer service is committed to respond within 24 hours to give you the best possible experience.

[y-brush.com](http://y-brush.com)