Oral-B®

**sonic
complete™**

OxyJet® Center

powered
by **BRAUN**

Deutsch	4, 9, 51
English	10, 15, 51
Français	16, 22, 51
Español	23, 29, 51
Português	30, 36, 51
Italiano	37, 43, 51
Nederlands	44, 50, 51

Internet:
www.oralb.com
www.braun.com

Type 3721

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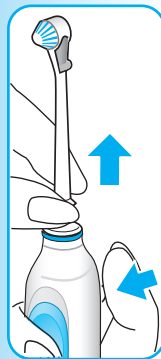
Thank you for purchasing this fine product from **Oral-B®**, the brand more dentists use themselves.

The **Oral-B® Sonic Complete™ OxyJet® Center** is a dentist-inspired oral care system caring for the health of your mouth by removing plaque, stimulating gums and cleaning harmful, odour-causing bacteria from your tongue.

The **Oral-B® Sonic Complete™** toothbrush offers customised brushing modes («Clean», «Soft» and «Massage») targeting the needs of your mouth. It also features multi-angled, multi-textured bristles that help clean tooth surfaces, between teeth and along the gumline.

The **OxyJet®** irrigator mixes air into water in the form of micro-bubbles. These micro-bubbles in the water stream attack plaque bacteria. At the same time, the irrigator removes remaining food particles from between your teeth and massages your gums. Always use the irrigator after brushing your teeth.

Please read the use instructions carefully before first use.

Important

The appliance must be mounted in a way that it cannot fall into water. The plug should always be accessible. Periodically check the cord for damage. A damaged or non-functioning appliance should no longer be used. Take it to an Oral-B Braun Service Centre.

Description

Toothbrush

- a Charger base
- b Charge indicator
- c Brushing modes
- d Handle (rechargeable)
- e On/off button
- f Brushhead
- g Storage tower for additional brushheads

Irrigator

- h Jet compartment
- i Irrigator jet
- j Jet switch
- k Jet release button
- l On/off switch
- m Water pressure control dial
- n Handle switch (water stop and go)
- o Irrigator handle
- p Water container (600 ml)
- q Cord storage

Specifications

Power supply: see bottom of housing
Power intake: 24 W

Using your Oral-B® Sonic Complete™ toothbrush

Oral-B® Sonic Complete™ has the following features designed to offer you individual mouth care:

Premium bristling technology



- 1) CrissCross™ bristles
Advanced multi-textured bristles to clean tooth surfaces, between teeth and along the gumline.
- 2) Interdental bristles
Raised to reach between teeth for interdental cleaning.
- 3) Power Tip® bristles
Extra-long bristles at the tip to help clean hard-to-reach areas.
- 4) Indicator® bristles
Blue colour fades halfway with use, signaling the need to replace the brushhead.


Different brushing modes for individual mouth care

- «Clean»
for a thorough cleaning of your teeth and gums.
- «Soft»
for the gentle care of your tongue and other sensitive areas.
- «Massage»
for a gentle stimulation of your gums.

Timers to encourage thorough brushing

- Professional Timer: A short stuttering sound at 30-second intervals reminds you to brush the four quadrants of your mouth equally for a more thorough clean.
- 2-Minute Timer: A longer stuttering sound indicates when the dentist recommended 2-minute brushing time has been reached.

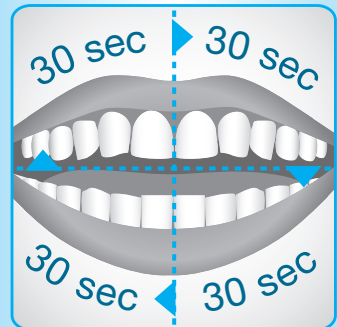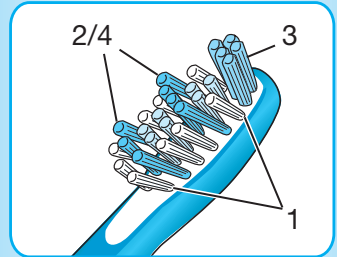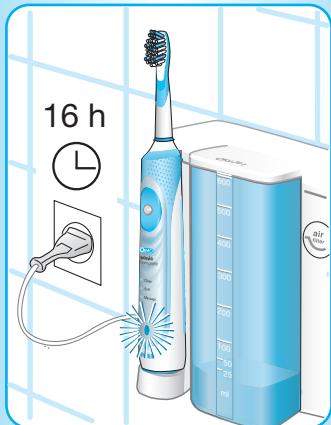
Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and can be used in the bathroom without hesitation.

- When you place your handle (d) on the charger base (a) and plug the base into an electrical outlet, the blue charge indicator (b) will light up. While charging, the light blinks. Once the handle is fully charged, the light will stop blinking and shine continuously. If you then use your brush and place it back on the charger, the charging light will again blink to let you know it is charging.
- A full charge lasts for approximately two weeks of regular brushing, so you can easily take your **Oral-B® Sonic Complete™** with you on holiday. Please note: Assumes the handle is used by one individual, twice a day for the dentist recommended 2-minute brushing time.
- Your **Oral-B® Sonic Complete™** is designed to consistently provide you with an effective clean. If the battery has run too low to provide this clean, your brush will stop.
- We recommend keeping your brush on the charger base when not in use to maintain a full battery charge. If your rechargeable battery does run down completely, it will take about 16 hours for your **Oral-B® Sonic Complete™** to fully charge.
- However, to maintain the maximum capacity of the rechargeable battery, refrain from placing the handle on the charger base and fully discharge it through regular use at least every 6 months.

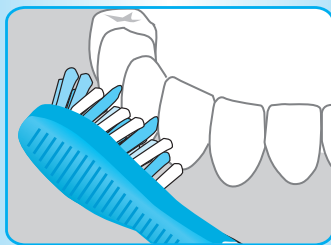Brushing technique

- Slide the brushhead on the handle. There will be a tiny space (2 mm) between the brushhead and the handle. This is to allow for a full range of motion for the many brush-head movements.
- Place the brushhead in your mouth before turning on your **Oral-B® Sonic Complete™**.
- Angle the bristles toward the gumline (see diagram).
- Hold the handle with a light grip and turn your **Oral-B® Sonic Complete™** on.
- Gently move the brushhead in a slightly circular motion. Allow the longer interdental bristles to reach between your teeth, while the CrissCross bristles focus on cleaning the tooth surface and gumline regions.
- After a few seconds, glide the bristles to the next section and pause while you gently circulate the bristles again. Try using only light pressure and let your **Oral-B® Sonic Complete™** do the brushing for you.
- Adjust your brushing, so that you keep as many bristles in contact with your teeth as possible. Don't forget to reach all areas, including the inner and chewing surfaces, and behind the back teeth.

Like a manual toothbrush, brushheads wear out with use and need to be replaced periodically.

Remember to floss daily.

Switching between brushing modes

- To start your **Oral-B® Sonic Complete™**, push the on/off button once, and your toothbrush automatically begins vibrating in the «Clean» mode.
- To shift from the «Clean» mode to the «Soft» mode, push the on/off button once more.
- Pushing the on/off button again will bring you from «Soft» to «Massage» mode.
- To turn your brush off from the «Massage» mode, push the on/off button again.
- If you wish to turn your **Oral-B® Sonic Complete™** off in the «Clean» or «Soft» mode, push and hold down the on/off button.

You may experience a slight tickling or tingling sensation when you first use **Oral-B® Sonic Complete™**. This will subside as your teeth and gums become accustomed to the **Oral-B® Sonic Complete™** cleaning experience. If you experience this initial tingling, you may want to brush in the «Soft» mode at first.

During the first few days of using the **Oral-B® Sonic Complete™**, it is possible that your gums may bleed slightly. This should stop after a few days. Should it persist, please consult your dentist or hygienist.


After use

After brushing, switch off the **Oral-B® Sonic Complete™** and remove the brushhead from the handle. Rinse the brushhead under water. Wipe the handle dry after each use. Periodically wipe the charger base clean with a damp cloth. Do not use alcohol to clean any parts of the **Oral-B® Sonic Complete™**.

Clean

Soft

Massage

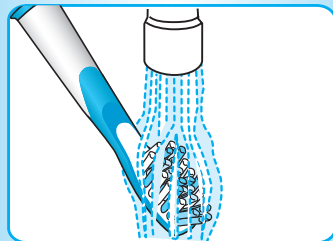Using your OxyJet® irrigator

Before using for the first time, fill the water container (p) and allow the water to run through completely by holding the irrigator handle (o) over a washbasin.

1. Attach one of the irrigator jets (i) to the handle until it locks into place. Fill the water container (p) with water and add mouthwash, if desired.

2. Turn on the irrigator by pressing the on/off switch (l). (When doing so, the handle switch (n) should be in the lower position = water stop).

3. Select the water pressure with the pressure control dial (m) (1 = gentle, 5 = strong). We recommend that you start with setting 1, especially if your gums tend to bleed or if children are using the irrigator.

4. Bend over a washbasin and guide the jet to your teeth and gums. Slide the handle switch (n) up to start the water stream. Keep your mouth slightly open to allow the water to run into the washbasin.








You can choose between two different jet functions by sliding the jet switch (j) down or up (before doing so, turn off the water stream by sliding the handle switch (n) down):

Rotating stream: for general cleaning and gum massage (jet switch in lower position). The rotation of the water stream and its micro-bubbles is designed for easy penetration below the gumline.

Mono stream: straight water stream for site-specific cleaning (jet switch in upper position).

Max. operating time: 15 minutes
Cooling off time: 2 hours

To ensure the jet's functions, a small sieve is incorporated in the lower end of the shaft. Should you notice a reduced water pressure, check the sieve for any deposits and remove them, e.g. by carefully tapping the shaft on a flat surface.

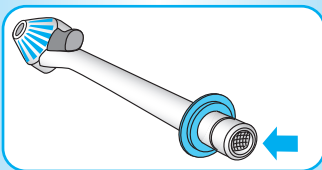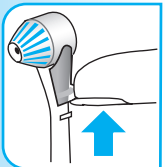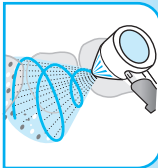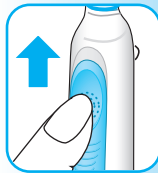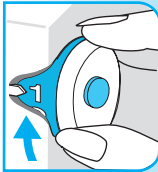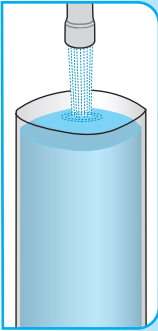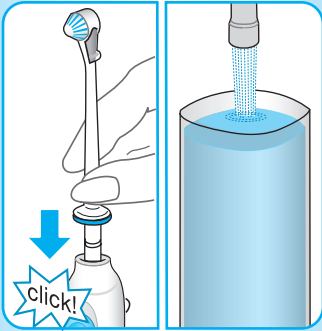After use

Always empty the water container (p) entirely to prevent bacteria from collecting in the standing water. If your dentist has recommended the use of a chlorhexidine solution with your irrigator, allow some water to run through the irrigator after use to prevent any deposit build-up.

Press the on/off switch (l) to turn off the irrigator. It will automatically turn off after 10 minutes if you forget to turn it off manually. In this case, to turn it on again, you have to press the on/off switch twice.

To take off the irrigator jet, press the jet release button (k). The jets can be stored in the jet compartment (h). Wipe the handle dry and place it back onto its base.

Cleaning

Regularly clean the housing with a damp cloth. The water container and the jet compartment are removable and dishwasher-safe.

Subject to change without notice.

This product conforms to the European Directives EMC 89/336/EEC and Low Voltage 73/23/EEC.

This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Braun Service Centre or at appropriate collection points provided in your country.

Guarantee

We grant 2 years guarantee on the product commencing on the date of purchase. Within the guarantee period we will eliminate, free of charge, any defects in the appliance resulting from faults in materials or workmanship, either by repairing or replacing the complete appliance as we may choose.


This guarantee extends to every country where this appliance is supplied by Braun or its appointed distributor.

This guarantee does not cover: damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The guarantee becomes void if repairs are undertaken by unauthorised persons and if original Braun parts are not used.

To obtain service within the guarantee period, hand in or send the complete appliance with your sales receipt to an authorised Oral-B Braun Customer Service Centre.

For UK only:

This guarantee in no way affects your rights under statutory law.

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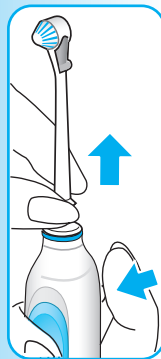Hong Kong

Audio Supplies Company,
Room 506,
St. George Building,
2 ICE House Street,
Hong Kong,
☎ (852) 2524 9377

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Tehran Bouran Company
145 Ghaem Magham Farahani,
Tehran
☎ +982 18 31 27 66

Ireland (Republic of)

Gavin's Electronics,
83-84, Lower Camden Street,
Dublin 2,
☎ 1800 509 448

Israel

S. Schestowitz Ltd.,
8 Shacham Str.,
Tel-Aviv, 495 17,
☎ 1 800 335 959

Italia

Servizio Consumatori Braun
Gillette Group Italy S.p.A.,
Via G. B. Pirelli, 18,
20124 Milano,
☎ 02 / 6678623

Jordan

Interbrands
Wadi Al-Sir, Adwiat alhekma
St. Mekadabi bld.
Amman
☎ +692 582 75 67

Kenya

Radbone – Clark Kenya Ltd.,
P. O. Box 40833,
Mombasa Road,
Nairobi,
☎ 2 82 12 76

Korea

Gillette Korea Ltd.
144-27 Samsung-dong,
Kangnam-ku,
Seoul, Korea,
☎ 080-920-6000

Kuwait

Union Trading Company,
Braun Service Center,
P. O. Box 28 Safat,
Safat Code 13001, Kuwait,
☎ 04 83 32 74

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Philippines

Gillette Philippines Inc.,
Corporate Corner Commerce
Avenues
20/F Tower 1,
IL Corporate Centre
1770 Muntinlupa city
☎ 02-771071 02-06/-16

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South Africa (Republic of)

Fixnet After Sales Service,
159 Queen Street, Kensington South,
P.O. Box 751770,
Johannesburg 2094,
Kensington South,
☎ 11 615 6765

St. Maarten

Rupchand Sons n.v. (ram's),
P. O. Box 79, Philipsburg,
Netherlands Antilles,
☎ 52 29 31

St. Thomas

Boochoand's Ltd.,
31 Main Street,
P. O. Box 5667,
00803 St. Thomas,
US Virgin Islands,
☎ 340 776 0302

Suomi

Gillette Group Finland Oy,
P.O. Box 9,
Niittykatu 8, PL 9,
02200 Espoo,
☎ 09-45 28 71

Sverige

Gillette Group Sverige AB,
Dept. Sweden, Stockholm Gillette
Räsundavägen 12,
Box 702, 16927 Solna,
☎ 020-21 33 21

Syria

Ahmed Hadaya Company
Hadaya building
Ain Keresh
Unisyria, P. O. Box 35002,
Damascus,
☎ 963 011-231433

Taiwan

Audio & Electr. Supplies Ltd.,
Brothers Bldg., 10th Floor,
85 Chung Shan N Rd., Sec. 1,
Taipei (104),
☎ (886) 02 2523 3283

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Yemen (Republic of)

Saba Stores for Trading,
26th September Street,
P. O. Box 5278,
Taiz,
☎ 4-25 23 88

For latest updates, please refer to
www.Braun.com.

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Gillette de Venezuela S.A.,