
IDItalian  
*design*

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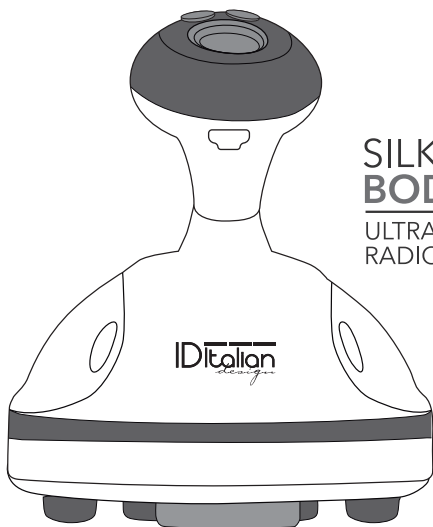ULTRASONIC &  
RADIOFREQUENCY

## EN PRECAUTIONS

Please read the instructions carefully before using the appliance.

To reduce any risk of accident, take the following precautions:

**1.** Always unplug the appliance immediately after use.

**2.** Never use the appliance while you are in the bath.

**3.** IEC 60417 - 5582: Do not

place or store the appliance anywhere where it could fall or be pushed towards the bath or sink.

**4.** Do not immerse or drop into the water or any other liquid. You could get an electric shock.

**5.** If the appliance falls into the water, unplug it immediately. Do not place

your hand in the water.

**6.** Keep the appliance away from water.

**7.** An appliance should never be left unattended while it is plugged into the mains.

**8.** This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

**9.** Use this appliance only for its intended purpose, as indicated in this manual. Do not use any accessory that is not recommended by the manufacturer.

- 10.** Never use this appliance if the power cord or plug is damaged, if it is malfunctioning, if it has been dropped or has fallen into water. Take the appliance to your nearest service centre for inspection or repair.
- 11.** Keep the power cord away from sources of heat or hot surfaces. Do not coil the power cord around the appliance.
- 12.** Do not use the appliance if you are drowsy.
- 13.** Do not drop the appliance or insert any object in the openings or tube.
- 14.** Disconnect the instrument with power supply before cleaning it.
- 15.** Before cleaning the instrument, make sure the power is off.
- 16.** Clean the residue of gels with a soft damp cloth or paper towl immediately after each use.
- 17.** Keep the instrument dry & clean.
- 18.** Store the device in the protective case after being cleaned.
- 19.** This care instruments suitable for private use only, do not share with your family or friends.
- 20.** Do not use an extension cable with this appliance.
- 21.** Preserve the device in a cool, dry place.
- 22.** If your skin is under nursing, or under any unhealthy condition, please consult your doctor before using the device.
- 22.** Do not use the appliance in case of failure and/or

malfunction. Contact directly the technical service.

**23. DO NOT USE** if you are sensitive to light. Test the device on your arm or leg prior to using on your face if you think you may be sensitive. Low-level light may cause pain and redness. If pain and redness last more than 24 hours, see your doctor.

**24. DO NOT** look directly into the light. **DO NOT USE** over the eye lids. The light is bright and you may feel brief pain or headache.

**25. DO NOT** fall asleep while using the device. You could burn yourself.

**26.** This product must not be used by children or by other people without assistance or supervision

if their physical, sensory or mental condition limits safe use of the appliance. Children must be watched to make sure that they do not play with the appliance.

## EN PRODUCT STRUCTURE

MODE, GEAR SELECTION KEY

SMART  
DISPLAY

TURN ON/OFF,  
FEATURES BUTTON

MAIN BODY

USB CHARGE PORT

ENERGY MASSAGE HEAD

ENERGY RED LIGHT

## EN FEATURES

1. 300.000 times/s ultrasonic high frequency vibration, instantly wake up the cells skin inside, improve cell metabolism, promote fat block decomposition, dredge lymph nodes.
2. Powerful Radiofrequency shaping

function, heat effect on the dermal layer of the skin, and enhance fat-burning.

**3.** Wavelength 625nm energy red with high purity, powerful photosynthetic energy and uniform energy density characteristics, activate skin cells, promote cell metabolism.

**4.** Mechanical vibration to make collagen fiber contraction, accelerate the rate of fat decomposition, effective reduction of stubborn fat.

**5.** Unique zinc alloy massage head, feel extraordinary, equipped with LCD HD Smart display, easy to operate mode.

**6.** With user-friendly design, sophisticated internal structure and perfect light extravagant appearance, take the ultimate fitness experience for you.

**7.** USB charge, use lasting, more energy efficient and environmentally friendly.

## EN SPECIFICATION

- Vibration Frequency: 300 KHZ.
- Operating Voltage: 8.4 V.
- Operating Current: 0.4 - 1.4 A.
- Charging voltage: 5.0 V.
- Power: 12W.
- Charging Time: 2.5 h.
- Continuous Use Time: 1 h.

## EN HOW TO USE

(waist, thighs, arms, etc) and spread evenly.

touch the skin to work.

parts from down to up with circle way, repeat 2 - 3 times.


- 1.** After bathing, apply slimming cream / oil on the required thin parts (waist, thighs, arms, etc) and spread evenly.
- 2.** Press On/Off button, select the appropriate mode, energy massage head touch the skin to work.
- 3.** Hold the handle tightly, uniform press massaging the required thin parts from down to up with circle way, repeat 2 - 3 times.
- 4.** After using, please clean the massage head with soft cloth or wet wipes in time.

**5.** Placed the machine in dry ventilated place

**6.** Please use the attached USB cable to charge.

## EN FUNCTION INTRODUCTION

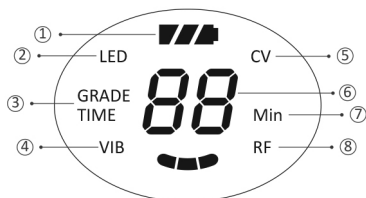1. Battery level
2. Red LED Therapy
3. Grade and time
4. Vibration Mode
5. Ultrasonic
6. Time
7. Minutes
8. Radio-frequency

### USING TIME OF EVERY MODE DOESN'T EXCEED 15MINS.

#### Brighten Skin Mode

Energy LED red light + Radiofrequency (red light flash slowly, flash quickly, steady on, three modes).

#### Burn fat Mode

Energy LED red light + ultrasonic CV.

#### Burst fat Mode

Energy LED red light + RF + Mechanical vibration VIB (red light flash quickly + intermittent mechanical short vibration, red light flash slowly +

intermittent mechanical long vibration, red light steady on + mechanical long vibration, three modes).

#### Shaping Mode

Mechanical vibration VIB + Ultrasonic CV (intermittent mechanical short vibration, intermittent mechanical long vibration, mechanical long vibration, three modes).

## EN TIPS

1. Boot to brighten skin mode, can change working mode by pressing function selecting key
2. Please use the matched adapter to charge.
3. One key to turn on, press 2s to turn off with "di di" sound; When power is not enough, will with "di di" sound, please charge in time.
4. If many hair on using part, please shave firstly, then use the machine.
5. Using effects is better after bathing.
6. Using the machine with slimming cream / oil together, enhance the effect.
7. No exceed 3 times one week.
8. This product is suitable for hand-arm, abdomen, leg, thigh fat and other fat easy accumulation part. (Doesn't use for: head, neck, collar bone, elbow, wrist, Knee, ankle, foot, spine and other parts).

## EN CAUTIONS

1. Pregnant women and patients with skin diseases don't use this product.
2. Avoid using within an hour after meal.
3. Don't overuse this product, avoid overheating, causing machine loss.
4. Don't use when bath or steam bath, avoid water to enter, causing machine loss.
5. This product is not waterproof, please use soft cloth or wet wipe to clean; Please keep machine clean if no use long time.
6. If feel unwell when use the machine, or skin appear abnormal, please stop to use and ask doctor.
7. Please placed the machine in dry ventilated place.
8. Please keep out of reach of children.

## EN CHARGING TIPS

1. For a better experience of this product, when used at the first time, please charge.
2. This product use a USB cable to charge and need 2,5 hours for fully charged.
3. Press the button for each time, the indicator light display red show flash, red flash, long red light, light off, 4 level modes.
4. Low power, if battery icon on the led display flash and with a beep

sound, please charge.

5. When you press the button without reaction, please charge it, but still don't work after charging, please check the notes or troubleshooting to solve the problem.
6. During charging, the screen displays the battery charging prompted. When power to the saturation, full battery tips.
7. After charging, place it in a dry and clean place for safe keeping.
8. Fully charge can use about 1 hour. About 20 times each 3 - minute operation.

## EN TROUBLE SHOOTING

### How to know when electricity is finished?

- Display battery icon flashing, with "di di" sound, when the final power consumption is finished, will automatically shut down.

### How to know when charge full?

- Charge on the case of turning on, battery icon will show incremental pattern, charge on the case of turning off, display will show "CH".

### No vibration?

- Charge.  
- Please contact Customer Service to repair.

**Can work when charge?**

- No, will lead to leakage accident, reduce product life, must be charged full, then can use.

**EN ENVIRONMENT**

Disposal of old Electrical & Electronic Equipment (Applicable throughout the European Union and other

European countries with separate collection programs. This symbol, found on your product or on its packaging, indicates that this product should not be treated as household waste when you wish to dispose of it. Instead, it should be handed over to an applicable collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent negative consequences to the environment and human health, which could otherwise be caused by inappropriate disposal of this product. The recycling of materials will help to conserve natural resources.

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