
ID **Italian**  
*design*

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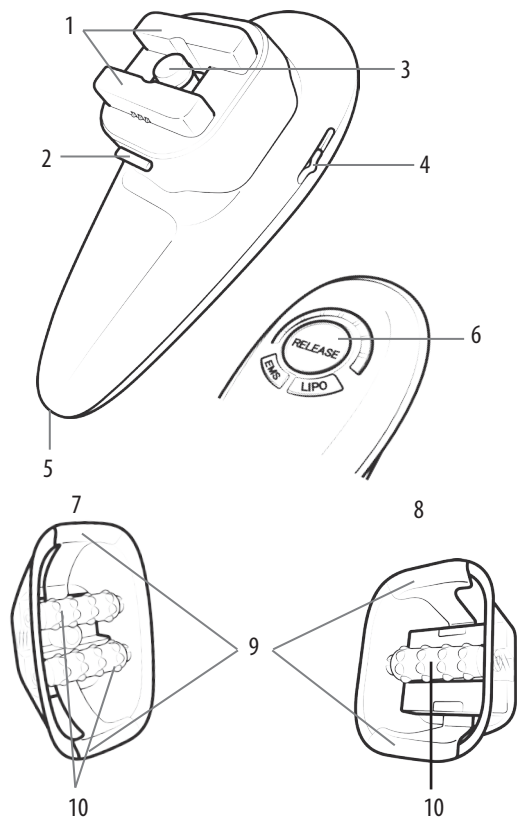7. Small suction cup (For use in the arm area, etc.).
8. Large suction cup (For use on the belly, legs, etc.).
9. EMS polarity plates (electrostimulation).
10. Roller.

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## EN PARTS SPECIFICATIONS

1. Hooking area of the suction cup.
2. Suction cup release button.
3. Filter.
4. Suction mode button.
5. USB charging Port.
6. Emergency stop button (temporarily stops the suction function when press).

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## EN PRECAUTIONS

Please read the instructions carefully before using the appliance.

To reduce any risk of accident, take the following precautions:

1. Always unplug the appliance immediately after use.
2. Never use the appliance while you are in the bath.
3. IEC 60417 - 5582: Do not place or store the appliance anywhere where it could fall or be pushed towards the bath or sink.
4. Do not immerse or drop into the water or any other liquid.
5. If the appliance falls into the water, unplug it immediately. Do not place your hand in the water.
6. Keep the appliance away from water.
7. An appliance should never be left unattended while it is plugged into the mains.
8. This appliance can be used by children aged from 8 years and above

and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.

**9.** Use this appliance only for its intended purpose, as indicated in this manual. Do not use any accessory that is not recommended by the manufacturer.

**10.** Do not use the appliance if you are drowsy.

**11.** Do not drop the appliance or insert any object in the openings or tube.

**12.** Disconnect the instrument with power supply before cleaning it.

**13.** If your skin is under nursing, or under any unhealthy condition, please consult your doctor before using the device.

**14.** Clean the residue of gels with a soft damp cloth or paper towel immediately after each use.

**15.** This care instruments suitable for private use only, do not share with your family or friends.

**16.** Preserve the device in a cool, dry place.

**17.** Before cleaning the instrument, make sure the power is off.

**18.** Do not use the appliance in case of failure and/or malfunction. Contact directly the technical service.

**19.** This product must not be used by children or by other people without assistance or supervision if their physical, sensory or mental condition limits safe use of the appliance. Children must be watched to make sure that they do not play with the appliance.

**20.** 2-3 minutes / area of use, maximum 30 minutes to cover the whole body.

**21.** Do not overuse, maximum 5 minutes in one area of use.

**22.** IMPORTANT: do not use on the face.

**23.** For people with loose muscles, this device can leave a bruised mark after use.

**24.** In case of bruising occurs, stop to use it immediately. Normally, it will gradually disappear after 1 week.

**25.** Do not push the unit with force against the body during use.


**26.** ATTENTION: After turning on the

unit, avoid contact with other beauty devices to avoid damage.

## EN BEFORE USE

Button functions and indications:

1. EMS level indicator (red, 5 levels).
2. LIPO level indicator (blue, 3 levels).
3. Suction emergency stop.
4. EMS button.
5. LIPO button.

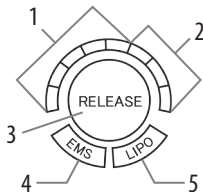**EMS:** Electrical pulse current that stimulates muscle contractions.

**LIPO:** 2 suction modes, squeezing the muscle directly:

- Suction mode 1: Suction and release, Intermittent suction in cycle.
- Suction mode 2: continuous suction.

### EMS button operations

Press and hold the "EMS Button" for 3 seconds to activate the EMS function with the related indication illuminated in RED. The EMS levels will change in the following sequence as the EMS button is pressed further. In initial use, start with level 1. Gradually increase EMS levels as the body becomes accustomed to using EMS. Level 1 (by initially pressing the EMS button) -> Level 2 -> Level 3 -> Level 4 -> Level 5 -> EMS OFF.

Wearing EMS feels like a light electric shock with muscle contractions, which is a normal phenomenon. The body will gradually get used to EMS.

### LIPO button operations

Hold down the "LIPO Button" for 3 seconds to activate the suction function with the related indication illuminated in BLUE. The suction levels will change in the following sequence by further pressing the LIPO button. In initial use, start with level 1. Gradually increase levels as the body becomes accustomed to using suction. Level 1 (by initially pressing the LIPO button) -> Level 2 -> Level 3 -> SUCTION OFF.


**NOTE:** In the initial use after purchase, or if it is not used for a long time, you must fully charge the unit before use.

Use at an ambient temperature of 25°C or 77°F.

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
1. Choose the appropriate suction cup size to fit the area of use.
2. Mount the suction cup in the hooking area of the main unit. Make sure to fix the cup firmly to the main unit to prevent it from falling out of or any air leakage.

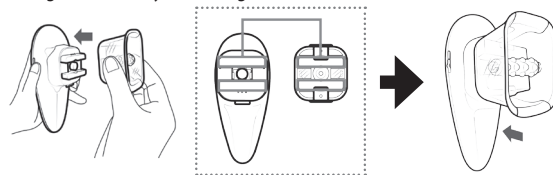**Removing the suction cup from the main unit.**

While holding down the cup release button, remove the cup from the main unit.
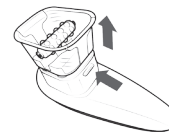
## EN PRECAUTIONS FOR USE

**IMPORTANT:** do not use on the face.

To avoid bruising of the skin, be sure to follow the following precautions:

- In case of bruising occurs, stop to use it immediately. Normally, it will gradually disappear after 1 week.
- Do not push the unit with force against the body during use.
- Press the button "RELEASE"; to release the cup from the body when you feel the suction power is too strong.

**NOTE:** Begin with level 1 on initial use.

**NOTE:** If the suction power is reduced check that the cup is correctly installed as otherwise air leaks may occur.

**It is advised to apply cream, lotion, body oil in the area of use to have a smooth movement on the body.** In addition, the EMS together with the rollers will facilitate its absorption.


**NOTE:** After use, the skin may become red, which is normal. After 1 day of non-use, the red color will disappear.

**Time of use**

- 2-3 minutes / area of use, maximum 30 minutes to cover the whole body.
- Do not overuse, maximum 5 minutes in one area of use.
- Keep a slow and continuous movement during use. You should not stay in one position for a long time.

**EN HOW TO USE**

1. Select and mount the suction cup on the main unit.
2. Press and hold the EMS button or the LIPO button for 3 seconds to activate the unit. After powering on, the corresponding LEDs will light up.
3. Use the EMS button or the LIPO button to select the appropriate function or levels.
4. At the first use, level 1 is suggested to start EMS or LIPO function.
5. When the suction mode is selected, be sure to check the selected suction mode:  
**A. Mode 1 intermittent suction:**  
 In initial use, Mode 1 is recommended to start.  
**B. Modo 2 continuous suction :**  
 Press the "release button" to release the cup from the body when you feel the suction power is too strong.
6. During use, avoid contact with water.
7. Firmly place the cup against your body and keep a slow movement. (2-3 minutes / for each area of use). During use, keep moving your body. Simply hold the unit firmly against your body to prevent it from falling, without the need to push hard. If you feel that the suction power is too strong, press the button "RELEASE" to release the cup from the body.
8. After use, press the EMS or LIPO button for 3 seconds to turn off the unit.

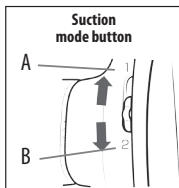**NOTA:** 10 minutes after power on, the unit will automatically shut down.  
**ATTENTION:** Using EMS feels like light electric shock with muscle contractions, which is normal phenomenon. The body will get used to EMS gradually.

**BELLY - WAIST.**

Slow movement from waist to thigh.  
 Great movement between left and right.

**ATTENTION: No use on the navel.**

Up and down movement between the lower chest and the lower waist.

**LEGS**

Up and down movement on the back of the calf.  
 Up and down movement in the thigh.


**HIP**

Move from the inside of the knee to the upper hip area.

**INNER SIDE OF ARM**

Move between the elbow and the armpit.

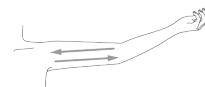**NOTE:** For sensitive areas or weakly toned muscle areas, use your hand to slightly stretch the skin and muscle to prevent the suction from hurting them

**EN CLEANING**

1. Make sure the unit is turned off.
2. Remove the used suction cup from the main unit.
3. Use a pliers to replace the old filter with a new one.  
 Replacement frequency: once a month.  
 Filter replacement when you feel that the suction power is lower.
4. Wash the clean suction cup with tap water.
5. Lightly dampen a clean cloth to keep the main unit clean.  
**Do not use detergent to clean the main unit or suction cup, to avoid damage.**  
**IMPORTANT:** This unit is NOT waterproof.

**EN PRODUCT SPECIFICATIONS**

Rated voltage: DC 5.0V.

Rated current: 0.5A.

Power: about 3W, under load.

The Charging time is about 2-3 hours. Charging time may vary due to different voltage remaining in the battery and under different circumstances of use.

Use time: About 30 minutes.

## EN STORAGE

After use, be sure to turn off the unit and remove the adapter from the socket.

Store the unit away from areas with high humidity or high temperature. If it is not going to be used for a long time, thoroughly clean the main unit and the suction cups for storage.

## EN ENVIRONMENT

[EN] Disposal of old Electrical & Electronic Equipment (Applicable throughout the European Union and other European countries with separate collection programs.

This symbol, found on your product or on its packaging, indicates that this product should not

be treated as household waste when you wish to dispose of it.

Instead, it should be handed over to an applicable collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent negative consequences to the environment and human health, which could otherwise be caused by inappropriate disposal of this product. The recycling of materials will help to conserve natural resources.

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**EN PROBLEM SOLVING**

Phenomenon	Possible reason	Countermeasure
Can't turn on?	Did you press the EMS or LIPO Button long enough?	Press and hold the EMS or LIPO button for 3 seconds
	Not enough battery left?	Recharge the unit with the correct method
Not enough battery after charging?	Not enough charging time?	Recharge the unit with the correct method
	Is the battery life exhausted?	Return to the specified service center for battery replacement
	Is the adapter not properly connected to the sockets	Plug the adapter securely back into the sockets
Can't load?	Is the USB connector not properly connected to the charging port?	Plug the USB adapter back into the charging port securely
Tingling feeling when wearing EMS?	Does the cup not make good contact with the skin?	Make sure the suction cup touches the skin well, especially when moving it by the body.
	Is the EMS level used too high?	Start from level 1. Only when the body is used to EMS increase gradually the EMS level.
	Is the skin too dry?	Apply a layer of lotion / cream to the skin before use.
Nothing feels when using EMS?	Is the skin too dry?	Apply a layer of lotion / cream to the skin before use.
	Is the EMS level used too low?	Gradually increase the EMS level.
	The used area is not sensitive to EMS sensation?	Please try another position
	There is a Difference in body sensation?	N / A
Suction power too low?	Does the cup not make good contact with the skin?	Make sure the cup touches the skin well, especially when moving the cup over the body, to avoid air leakage.
	Not enough battery left?	Recharge the unit properly.
	Cup is not properly mounted on main unit?	Make sure the cup is securely mounted on the main unit.
	Is the cup dirty?	Clean the glass only with tap water.
	Is the filter dirty?	Replace the filter
Can't move the suction comma during use?		Use a lower suction power.
	Is the suction feeling too strong?	Press the RELEASE button, to remove the cup from the body Stretch loose skin a bit during use.
	Is the skin too dry?	Apply a layer of lotion to the skin before use.

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