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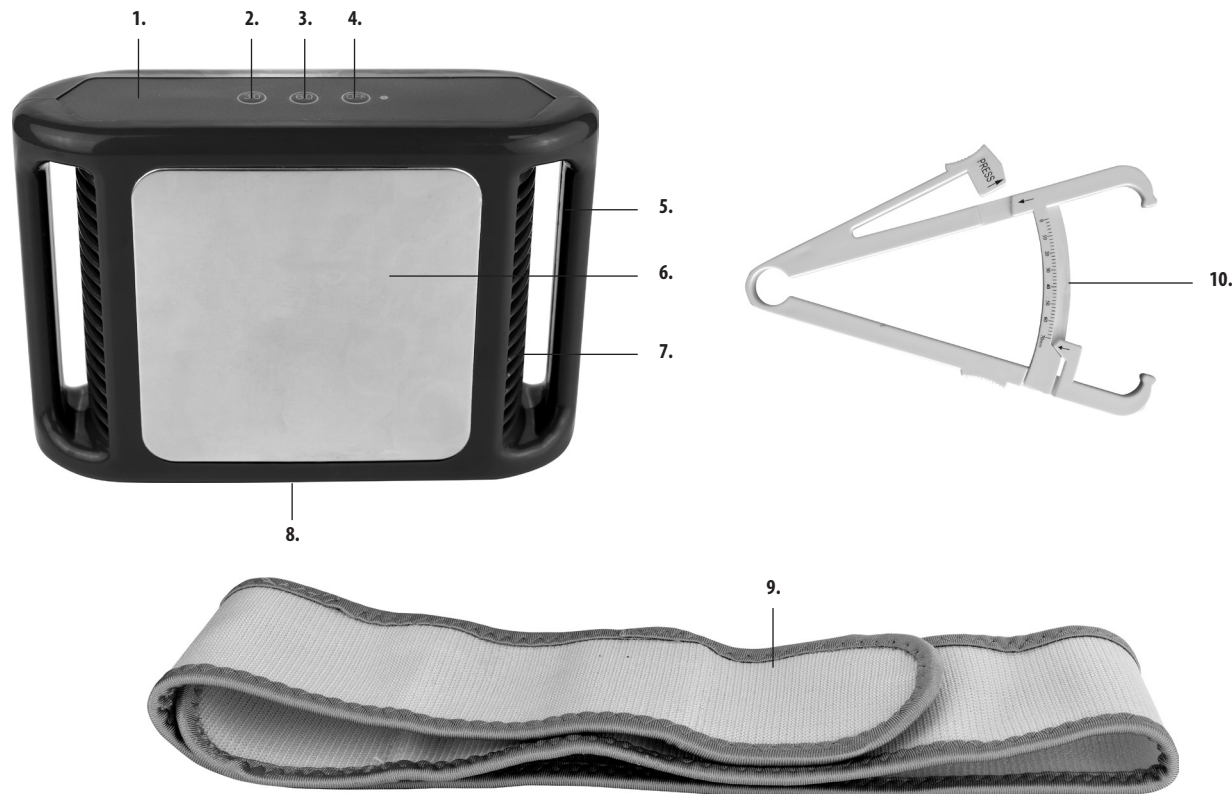
EN PRECAUTIONS

Please read the instructions carefully before using the appliance. To reduce any risk of accident, take the following precautions:

1. Unplug the appliance from outlet when it not in use or before cleaning.
2. IEC 60417 - 5582: Do not place or store the appliance anywhere where it could fall or be pushed towards the bath or sink.
3. Do not immerse or drop into the water or any other liquid.
4. If the appliance falls into the water, unplug it immediately. Do not place your hand in the water.
5. Keep the appliance away from water.
6. An appliance should never be left unattended while it is plugged into the mains.
7. Use this appliance only for its intended purpose, as indicated in this manual. Do not use any accessory that is not recommended by the manufacturer.

8. Never use this appliance if the power cord or plug is damaged, if it is malfunctioning, if it has been dropped or has fallen into water. Take the appliance to your nearest service centre for inspection or repair.
9. Keep the power cord away from sources of heat or hot surfaces. Do not coil the power cord around the appliance.
10. Do not use an extension cable with this appliance.
11. Do not drop the appliance or insert any object in the openings or tube.
12. Never use a sponge or abrasive cloth on the unit.
13. This appliance is for household use only, not professional or industrial use.
14. This appliance should be stored out of reach of children and/or the disabled.
15. Do not try to open the control housing under any circumstance and never attempt to repair the equipment yourself.
16. This device is provided with a universal AC-DC adapter. This means connecting the power cord to any 100V-200V or 220V-240V outlet.
17. DO NOT connect the AC power adaptor to the AC outlet before connecting the adaptor to the appliance, as this may damage the AC-DC adaptor.
18. IMPORTANT Damage may occur if you use an adaptor other than the one supplied to power this appliance or if you use the adaptor of this appliance to power other electrical devices. If smoke, a burning smell, or extreme heat comes out of the AC-DC adapter, you may damage this appliance with a defective AC-DC adapter.
19. WARNING The AC-DC adapter may become hot or warm when in use. Be sure not to cover the adapter and keep it away from your body.
20. IMPORTANT! Unplug the power adapter or switch off the AC outlet to minimize power consumption when the appliance is unused.
21. Do not use this product if you are pregnant, under 18, immobile, or anyone who cannot clearly understand the instructions and operate the controls.
22. Do not use this product if you have a heart-related condition, epilepsy, fever, inflammation, dermatitis (acne, eczema, herpes), warts, open or recently stitched wounds, severe or progressive diseases (cancer, AIDS, hepatitis, blood diseases, etc.), blood clotting problems or mental illness.
23. Do not use this device on the head, eyes, neck (throat), breasts, genital area, and inside the mouth.
24. IMPORTANT! Do not use this appliance without a protective pad or cloth material between the freezing plate and the skin. DO NOT APPLY DIRECTLY TO THE SKIN SURFACE.

25. It is unsafe to use the fat freezer, its accessories, detachable parts, and materials other than as described in the instructions for use.
26. Do not clean or perform any maintenance on the fat freezer while the appliance is in use.
27. Do not position the appliance to disconnect the power cable from the adapter. Keep the device in the correct position with the CONTROL PANEL buttons facing upwards. This will allow quick access to the OFF button and quick access to removing the DC plug from the appliance.
28. IMPORTANT!: Never use the fat freezer twice (in succession) on the same body area in a single session. You must allow at least one week between sessions on the same site of the body.
29. This device is not intended to diagnose, treat, cure, or prevent any disease, illness, or medical condition. It is not a weight loss device.

EN ENVIRONMENTAL CONDITIONS

Transport and storage

Ambient temperature: -20°C to 60°C (-4°F to 140°F).

Relative humidity: 0% to 85%.

Atmospheric pressure: 700hPa to 1060hPa (700hPa to 1060hPa).

Normal Use

Ambient temperature: 20°C to 25°C (68°F to 77°F).

Relative humidity: 30% to 85%.

Atmospheric pressure: 700hPa to 1060hPa.

EN POSSIBLE SIDE EFFECTS

Temporary but rare side effects include bruising, numbness, redness, and skin irritation where the fat freezer has been placed. These side effects are rare and only temporary. If skin irritation occurs after treatment, discontinue use immediately. If symptoms persist, contact a physician.

EN BEFORE USE

- Ensure the appliance is working correctly by switching it on briefly before use.
- Make sure that the appliance is correctly connected.
- Avoid placing the appliance near a strong magnetic field or static electricity.
- Avoid placing the appliance near high-frequency equipment.
- Avoid placing the appliance in direct contact with the skin.

EN USE

1. Wrap one end of the Velcro belt around the handle at one end of the device to secure the belt properly.
 2. Apply a protective pad to the area of your body you wish to treat. You can use this device on your arms, legs, buttocks, back, and abdomen. Avoid the breast and genital area - AVOID DIRECT CONTACT WITH THE SKIN!
 3. Carefully place the device (metal plate) on the protective pad to avoid direct contact with the skin.
 4. Wrap the Velcro belt tightly around the body part and attach the Velcro ends together to secure the appliance firmly to the body part to prevent it from moving during treatment.
- WARNING!:** To avoid strangulation caused by cutting or obstruction of blood circulation, DO NOT tighten the strap around the body part excessively. The Velcro belt holds the fat freezer securely to prevent it from slipping and is NOT intended to obstruct blood circulation.
- WARNING!** If you feel numbness in the part of your body where the device is working, loosen the strap.
5. Plug the AC/DC adapter/power cord into the device and then into an electrical outlet.
 6. Press the 30- or 60-minute setting button to turn on the fat freezer.
- Note:** 30 minutes are used to tone the area to be treated.
Note: 60 minutes are used to remove fat cells from that area.
Note: A green light comes on after pressing the "30" or "60" button to confirm that the appliance works properly.
7. After the time has elapsed, switch off the fat freezer and disconnect it from the power cord.
- NEVER** use the appliance twice (in succession) on the same body area in a single session. It would help if you allowed at least one week between sessions on the same site of the body. Once your session has been completed (either a 30-minute- or 60-minute), please make sure you drink at least 250ml of water to aid in the lipolysis process. This will allow your body to quickly and adequately eliminate lipids.
8. OFF button: Pressing this button will turn the device off. Press this button to interrupt the timed operation (30 or 60 minutes).
 9. Gently massage the body area with circular movements for at least 2 minutes after the treatment.
- WARNING!** Stop using the appliance if you experience unusual pain. Immediately press the "OFF" button to stop the machine and disconnect it from the power supply.

WARNING! Stop using the appliance if you notice any abnormality (e.g., the liquid inside the device, inability to switch on as intended, or unusual sound with the internal fan).

EN AFTER USE

- Let the corporate stand for 10 minutes after switching it off and unplugging it. This allows the metal plate to return to average (room) temperature.
- Clean the separate with a soft cloth.
- You can use an organic solvent (e.g., antiseptic solution) to clean the device.
- **WARNING!** After using the device, you MUST wait one week before using it again in the same area.
- Keep the original box and packaging for storage and transport.

EN CLEANING AND MAINTENANCE

- Dirt and dust can accumulate and adversely affect the fat freezer's performance and decrease the appliance's longevity.
- Use a soft cloth to clean the appliance after each use. The entire appliance can be cleaned.
- Periodically, you may use a dust removal blower to clean the fan area of dust build-up.
- DO NOT attempt to open the appliance, as this may damage the delicate electrical components of the device.
- The protective pads are the only disposable accessories used with the grease freezer.
- The protective pads should be discarded after use. They can be disposed of in your household waste receptacle/container in the same manner as you would dispose of standard household waste.
- There are no replaceable parts that require installation by the consumer to repair this appliance.
- Please do not leave the appliance where it may be affected by vibrations or shocks.
- Do not leave the appliance in a dusty, saline, or sulphuric environment. These elements may adversely affect the performance of the fat freezer and shorten its service life.

EN TROUBLESHOOTING

It does not turn on:

1. Ensure the appliance is properly connected to an electrical outlet, the power cord is connected correctly to the AC-DC adapter,

the AC-DC adapter is connected to the appliance, and the power cord is connected to an electrical outlet.

2. Make sure that the electrical outlet has power.

It does not turn off:

1. If you press the "OFF" button to stop the fat freezer and it does not turn off, immediately disconnect the AC-DC adapter from the machine.
2. Then disconnect the power cord from the power outlet.

EN ENVIRONMENT

Disposal of old Electrical & Electronic Equipment (Applicable throughout the European Union and other European countries with separate collection programs.)

This symbol, found on your product or on its packaging, indicates that this product should not be treated as household waste when you wish to dispose of it. Instead, it should be handed over to an applicable collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent negative consequences to the environment and human health, which could otherwise be caused by inappropriate disposal of this product. The recycling of materials will help to conserve natural resources. For more detailed information about the recycling of this product, please contact your local city office, household waste disposal service or the retail store where you purchased this product.

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