SCH400

www.philips.com/AVENT

Philips Consumer Lifestyle BV
Tussendiepen 4, 9206 AD Drachten, The Netherlands


Trademarks owned by the Philips Group.
©2014 Koninklijke Philips N.V. All rights reserved

www.philips.com

100% recycled paper
100% papier recyclé

4213.354.3583.4

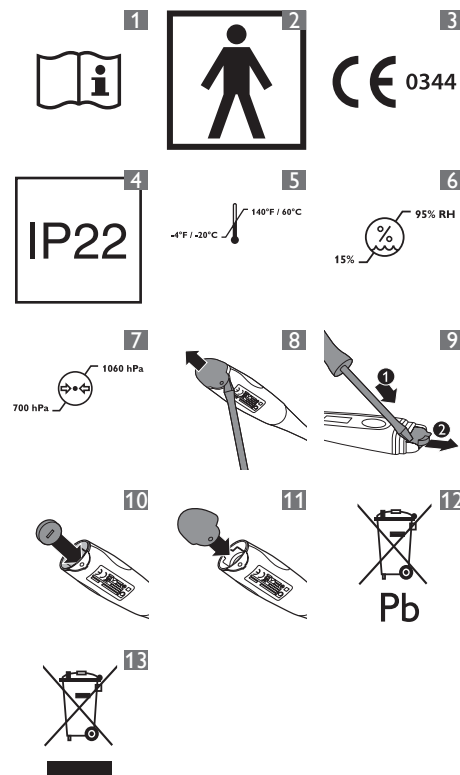ENGLISH

Important

Congratulations on your purchase and welcome to Philips! To fully benefit from the support that Philips offers, register your product at www.philips.com/welcome. Read this user manual carefully and save it for future reference.

Warnings

- The digital thermometer is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the digital thermometer by a person responsible for their safety.
- Never sterilise the digital thermometer in boiling water.
- **Risk of swallowing batteries!** The digital thermometer contains a button battery. If the button battery is swallowed, it can cause severe internal burns in just 2 hours and can lead to death.
- Battery compartment cover of the digital thermometer must be replaced properly if removed. If loose, it may choke your child.
- These products contain items with sharp edges and small parts. These products are not toys. Do not allow children to play with them. When not in use, keep these products out of the reach of children.
- Never leave a child unattended with these products.
- Clean all products before and after each use. See section 'Cleaning'.
- Do not clean any product with aggressive solvents or chemicals.
- Do not store or leave any product in direct sunlight or near any source of heat.

Cautions for digital thermometer

- The digital thermometer is suitable for children and adults of all ages.
- The digital thermometer is suitable for use at room temperatures between 10° and 35°C.
- Do not drop the digital thermometer or expose it to shocks.
- Do not bend or disassemble the tip of the digital thermometer.
- Do not use the digital thermometer if the display shows the battery symbol. This means that the battery is low and therefore measurements will be inaccurate.
- Do not bite the metal tip while you measure temperature orally. Simply place the tip under the tongue.

Symbols on digital thermometer

- 1 This symbol means: Read this manual (Fig. 1).
- 2 This symbol is to identify a type BF applied part complying with IEC60601-1 (Fig. 2).
- 3 This symbol means: Conforms to EC Directives. CE stands for 'Conformité Européenne'. 0344 is the number of the notified body (Fig. 3).
- 4 This symbol means: protected against access to hazardous parts with a finger and against vertically falling water drops (Fig. 4).
- 5 This symbol means: storage temperature (Fig. 5).
- 6 This symbol means: relative humidity (Fig. 6).
- 7 This symbol means: atmospheric pressure (Fig. 7).

Electromagnetic fields (EMF)

The Philips AVENT digital thermometer complies with all applicable standards and regulations regarding exposure to electromagnetic fields.

Digital thermometer (SCH410)

The digital thermometer measures your child's body temperature rectally (in the rectum), axially (in the armpit) or orally (in the mouth).

Note: Test in bath water for 10 seconds.

Note: Results may vary based on the individual and the method used.

Note: To ensure the best accuracy, we advise you to measure the temperature rectally.

Note: Oral and axial measurements are usually 0.5 to 1°C lower than actual body temperature.

Children's temperature

The temperature of a child is normally between 36.5 and 37.5°C. Keep in mind that a child's temperature varies depending on age, activity, time of day, etc. Remember that some illnesses do not have high temperature as a symptom. Temperature measurement should not be a substitute for consulting your physician. When reporting your child's temperature to a professional, state the measurement method.

Rectal measurement

For rectal measurement of a baby's temperature, proceed as follows:

- 1 Remove the protection cap from the digital thermometer. Make sure that the metal probe of the flexible tip is still firmly attached to the digital thermometer.
- 2 Lubricate the tip of the digital thermometer with a water-soluble jelly.

Note: Do not use petroleum jelly (e.g. vaseline).

- 3 Press to switch on the digital thermometer.
- 4 Lay the baby on its stomach or its side.
- 5 With one hand, separate the baby's buttocks. With the other hand, gently slide the tip of the digital thermometer no more than 12mm (1/2 inch) into the rectum. If you detect any resistance, stop immediately.

Hold the baby firmly so that it cannot turn while the digital thermometer is inserted.

Never push the flexible tip too far into the rectum.

- 6 Keep the digital thermometer in place until the measurement is finished (when the °C indication stops flashing and the thermometer beeps).

See section 'Audible signals' at the end of this chapter for the meaning of the beeps.

Note: If for some reason the temperature cannot be measured, the display shows 'Lo' or 'Hi'. This means that the temperature is outside the temperature range. Switch off the digital thermometer, check if it is correctly positioned and switch it on again.

- 7 Remove the digital thermometer and read the display.
- 8 Press to switch off the digital thermometer.

Note: Switch off the digital thermometer after use. The digital

thermometer switches off automatically after 1 min. if the temperature measured is lower than 32°C, or after 10 min. if the temperature measured is higher than 32°C.

- Measuring the temperature of older children**
- For older children and adults, the digital thermometer can also be used to measure temperature orally (in the mouth) or axially (in the armpit).
- **How to measure the temperature orally:** Place the tip of the thermometer under your child's tongue.

Note: Make sure your child seals its lips.

- **How to measure the temperature axially:** Place the tip of the thermometer in your child's armpit. Press its arm against its body and hold its elbow to keep the thermometer in place.

Note: Make sure your child's armpit is dry.

Do not use the digital thermometer for oral or axial temperature measurements on babies.

After using the digital thermometer rectally or in the armpit, do not use it for oral measurement for hygienic reasons, not even after cleaning.

Audible signals

1 beep per second: Temperature is normal (<37.8°C).
3 beeps per second: Possible fever (>37.8°C).

Note: The thermometer only indicates fever. It does not beep when your child's temperature is below normal.

Features

Last temperature recall

The digital thermometer has a memory function. When the digital thermometer is switched on, the display automatically shows the last measured temperature for approx. 2 seconds. The display then shows the normal body temperature of 37°C.

Battery-low indication

Batteries can last for about 2 to 3 years. If the display shows a battery symbol , the battery is low. Replace the battery (see section 'Replacing the battery').

Error function

In case of malfunction, the display shows 'Err'. Proceed as follows if this happens:

- 1 Press  to reset the digital thermometer.
- Replacing the battery**

Batteries may explode if exposed to high temperatures or fire.

Never attempt to recharge the battery.

If the display shows a battery symbol, the battery is low. Replace the battery of the digital thermometer as described below.

- 1 Replace the battery of the digital thermometer using a small flat-blade screwdriver. See illustrations fig. 8, 9, 10 and 11.
 - 5 Dispose the battery properly (see chapter 'Environment').
- Cleaning**

Be careful not to damage the parts when you clean the digital thermometer.

- 1 Use warm water and some washing-up liquid.

Do not use abrasive cleaning agents, anti-bacterial cleaning agents or chemical solvents. Excessive concentration of cleaning agents may eventually cause plastic components to crack.

Do not use cracked parts.

Storing

- 1 Put the protection cap on the digital thermometer after use and store it in a dry place.

Technical specifications

Digital thermometer	
Battery	1.55VDC (LR41 or SR41)
Battery life	> 100 operating hours
Measurement range	32°C to 42.9°C
Accuracy	± 0.1°C for 35°C to 42°C range
Compliance	Conforms to MDD93/42/EEC
Storage temperature	-20°C to 60°C / -4°F to 140°F
Humidity level	15% to 95% RH
Atmospheric pressure	700hpa - 1060hpa

Nasal aspirator

- 1 Hold the baby upright when you use it.
- 2 Squeeze the aspirator bulb and place the tip gently into the baby's nostril and tilt the tip to the side.
- 3 Release the bulb slowly to let the suction draw mucus out of the baby's nose.
- 4 Repeat the process with the other nostril.

Never use the nasal aspirator in the baby's mouth or ears.

Do not insert the aspirator into the baby's nose too far.

Never point the tip directly in the centre of the nostril.

Cleaning

- 1 Disassemble all parts of the nasal aspirator and wash them in warm water with washing-up liquid.
- 2 Sterilise the nasal aspirator in boiling water for 10 min. and dry it afterwards.

Finger toothbrush

- 1 Place the finger toothbrush on your fingertip to massage the baby's gums gently.

Cleaning

- 1 Sterilise with a steriliser or boiling water. Let it cool down and dry.

Brush and comb

- 1 Brush/comb the baby's hair from root to tip.

Tip: Carefully keep a hand on the roots while brushing or combing to avoid pulling the baby's scalp.

Cleaning

- 1 Clean with warm water and some washing up liquid, rinse them thoroughly and let them dry.

Nail cutter and scissors

Be careful not to cut the nails too close to the baby's skin.

Tip: We advise you to cut the baby's nails after a bath when the nails are soft, or when the baby sleeps.

Cleaning

- 1 Wipe the nail cutter and the scissors with a damp cloth and dry them.

Note: You can sterilise the nail cutter and the scissors by wiping their cutting edges with alcohol.

Emery boards

- 1 Gently glide the emery board across the baby's nails.

Pouch

Save all items in the pouch, and store them in a dry place.

Environment

- Non-rechargeable batteries contain substances that may pollute the environment. Do not throw away non-rechargeable batteries with the normal household waste, but dispose of them at an official collection point for batteries. Always remove the non-rechargeable battery before you discard and hand in the digital thermometer at an official collection point (Fig. 12).
- Do not throw away the digital thermometer with the normal household waste at the end of its life, but hand it in at an official collection point for recycling. By doing this, you help to preserve the environment (Fig. 13).

Guarantee and support

If you need information or support, please visit www.philips.com/support, or for the digital thermometer, read the separate worldwide guarantee leaflet.

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