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Dear Customer,

Thank you for purchasing this product. To be completely satisfied with your product, we recommend reading this user manual carefully.

If you have any difficulties or need more information, you can find our help centre by scanning the QR code. If you cannot find the answer to your question, you can also contact our After Sales Service: support.uk@terraillon.com

PRECAUTIONS

The following precautions should always be taken before using your TB251 blood pressure monitor, hereinafter referred to as “the device”.

GENERAL INFORMATION

1. Read the operating instructions carefully and keep them for future reference. If you pass this device on to others, please also pass on this manual.
2. Use this device only in accordance with the instructions in this user manual.

PRODUCT INFORMATION

- The blood pressure monitor is a digital monitor for measuring arterial pressure and heart rate with an arm circumference of 22 cm to 42 cm. It is intended for adult use only.
- This device uses the oscillometric method. This non-invasive method is used to measure systolic pressure, diastolic pressure and heart rate. This technique detects variations in the brachial artery and interprets them digitally to facilitate reading.
- This device meets the requirements of ISO 81060-2, ensuring reliability and accuracy.
- The blood pressure monitor and cuff are suitable for use in the patient environment. That is, for indoor use at home.
- This device must be installed and started up in accordance with the information provided in this user manual.
- You can use this device to take your own measurements — no third party operator is required.
- Do not use this device for the diagnosis or treatment of a health condition or disease. Contact a healthcare professional in case of repeated abnormal measurements. Do not change your medication without the advice of your doctor or healthcare professional.
- Only use accessories approved by the manufacturer. The use of unapproved accessories may damage the device, injure users or provide incorrect data.
- Please use the device in the environment provided in the user manual. Otherwise, the performance and lifetime of the device will be affected and reduced.
- Similarly, performance can be affected by extreme temperature, humidity and altitude.

- Do not drop your device or expose it to strong shocks or vibrations.
- No modification of this equipment is permitted.
- Report any malfunctions or unexpected events to the manufacturer.
- Keep the device, cuff and hose, charging cable and batteries out of reach of infants, children and elderly persons as they may present a choking or strangulation risk if used incorrectly.
- Keep small parts out of reach of children as they may swallow them. Contact a doctor immediately if ingested accidentally.

This device may only be used for the intended purpose described in this manual. The manufacturer cannot be held liable for incidental, consequential or special damage caused by misuse or abuse.

CONTRAINDICATIONS

1. The device uses the oscillometric method to measure arterial pressure. The arm used for the measurement must have normal blood circulation. It is not suitable for people with restricted or impaired blood circulation.
2. Avoid taking measurements on people who are unable to communicate clearly, such as newborns, infants, children, or unconscious patients.
3. The device is not suitable for use on women who are or may be pregnant; people with diabetes; patients with implanted electronic devices such as cardiac stimulators, defibrillators, pacemakers, etc.; patients with pre-eclampsia, ventricular extrasystole, atrial fibrillation, peripheral arterial disease or kidney disease. Patients undergoing intravascular treatment or arteriovenous anastomosis; individuals who have undergone mastectomy.
4. The device is not suitable for measurements on people with illnesses that cause uncontrollable movements (shaking, chills, etc.).

If you are in any of these four situations or suffer from an illness, please consult your doctor before using the device.

- Similarly, if you are taking medication, consult your doctor to determine the appropriate time to measure your arterial pressure.
- Do not use this device if you are allergic to polyester, nylon or plastic.

USEFUL TO KNOW:


- The cuff meets the requirements of ISO 10993-5, ISO 10993-10 and ISO 10993-23.

However, some sensitive people may have allergic reactions.

PACKAGE CONTENTS

Blood pressure
monitor
x1

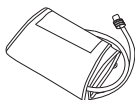Arm cuff
22-42 cm
x1

Carry pouch
x1


Instruction
manual
x1

AAA batteries
x4

PRODUCT OVERVIEW

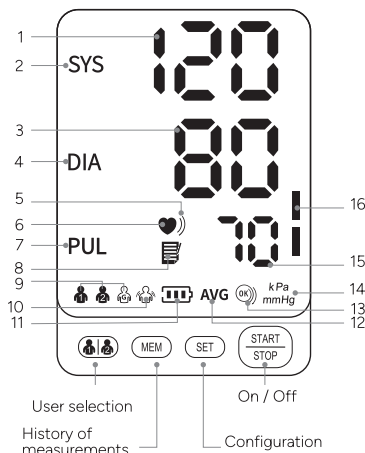READING THE DATA

1. Systolic pressure value: arterial pressure in the heart contraction phase
2. Systolic pressure symbol
3. Diastolic pressure value: arterial pressure in the heart relaxation phase
4. Diastolic pressure symbol
5. Irregular heartbeat symbol: irregular heartbeat detected
6. Pulse symbol: heartbeat detection during measurement
7. Heart rate symbol (pulse)
8. Memory: stored measurements
9. User symbols
10. "Stay still" symbol
11. Battery symbol
12. Average for the last three measurements
13. Warning for incorrect cuff positioning
14. Blood pressure unit
15. Heart rate (pulse): number of heart beats per minute
16. Coloured voltage level indicator (WHO — World Health Organisation classification)

PREPARATION BEFORE USE


Before use, make sure that the device is operating safely and in good working order. Do not use if it is damaged in any way. The continuous use of a damaged device may cause injury, improper results or serious danger.

I. POWERING UP THE PRODUCT

INSERTING BATTERIES

1. Open the battery cover according to the method shown in the figure.
2. Insert four AAA batteries into the battery compartment, observing the + and - polarity.

CAUTION: Remove the batteries if you do not intend to use the device for a long period of time (more than 3 months)

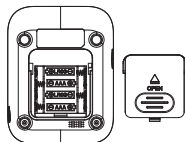USB-C CONNECTION FOR POWER

In addition to the batteries, power can also be provided by connecting to an external 5V DC, 1A power supply via a USB-C port. (Power cord is not included.)

NOTE: The USB-C power connection is suitable for temporary use when you do not have suitable AAA batteries.

The adapter must comply with the IEC 60601-1 standard, and the specifications must comply with the following requirements: Input: AC 100~240V 50/60 Hz, Output: DC 5V 1.0A. Other AC adapters may vary in terms of output voltage and polarities and may pose a risk to your life and damage the device.


External power supply port

II. USER SELECTION

This monitor offers three user profiles: two users with a memory capacity of 199 records each, and a guest mode.

When the device is switched off, press the  button corresponding to the user to choose: user 1, user 2, guest.

Note: The guest's data will not be saved.

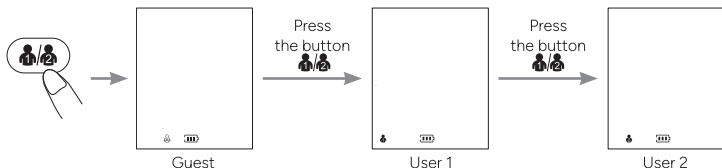
III. UNIT DISPLAY SETTINGS

There are two units for arterial pressure: mmHg and kPa.

The default unit is mmHg.

When the device is switched off, press the button ^{START}/_{STOP} device is switched off, press the button for approx. 5 seconds to access the unit selection. Press the button MEM to toggle between mmHg and kPa, then press the button ^{START}/_{STOP} to confirm the selection.

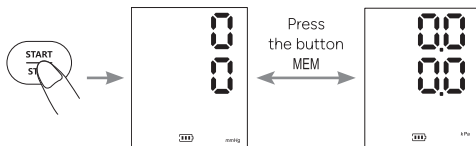TAKING MEASUREMENTS

I. PLACING THE CUFF


The cuff is placed on the left arm.

However, if your doctor has diagnosed you with poor blood circulation on this arm, or if you are in one of the situations explained in the "Warnings" section on page 24, you can use the right arm.


1. Remove all jewellery, such as watches and bracelets.

Roll or push up your sleeve to expose the skin. Make sure that your sleeve does not squeeze your arm too tightly.

2. Connect the cuff to the monitor by inserting the air Intake into the Air Jack.

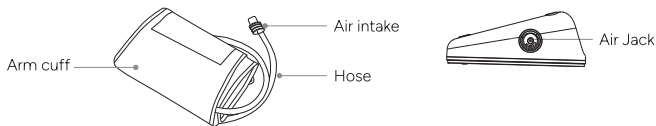
3. Place your hand in the cuff loop. Pull on the cuff until it reaches the upper part of your arm.

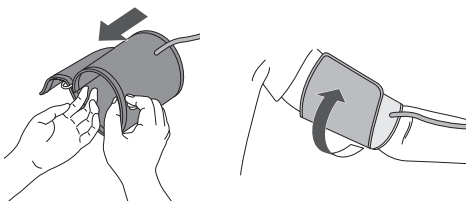4. Make sure that:

- the lower edge of the cuff is 2 to 3 cm above the elbow,
- the hose is outside, aligned with your middle finger.

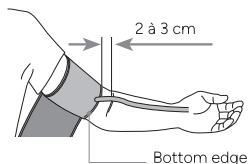5. Adjust your cuff.

Tip: your cuff is tight if you can insert one finger between your arm and the cuff.

Warnings:

- Do not place the cuff on:
 1. Extremities other than the arm.
 2. A wound as this can lead to additional injuries.
 3. The arm of a person whose arteries or veins are undergoing medical treatment, i.e. intravascular access, intravascular therapy or an arteriovenous (A-V) shunt, which could disrupt blood circulation and cause injury.
 4. The arm on the same side as a mastectomy (especially when lymph nodes have been removed). It is recommended to measure on the unaffected side.
 5. An arm to which an intravenous drip or blood transfusion is attached.
 6. An arm where another monitoring device is applied. One or both devices may temporarily stop working if you try to use them at the same time.
 7. Patient areas where skin is delicate or damaged.
- Frequently check the cuff location for irritation or to ensure there is no irritation.
- Do not bend, stretch, squeeze or deform the hose during the measurement, as the cuff pressure may increase continuously, which may prevent blood flow and cause injury.

Safety:

In the rare event that the cuff remains fully inflated during the measurement, loosen and remove the cuff immediately. Prolonged high pressure applied to the arm (cuff pressure >300 mm Hg or constant pressure >15 mm Hg for more than 3 minutes) may cause ecchymosis and discolouration of the skin.


II. SITTING IN THE RIGHT POSITION

To take a measurement, you should be relaxed and comfortable in a room with a pleasant temperature.

- Sit upright in a comfortable chair with your back resting on the backrest.
- Do not cross your legs and keep your feet flat on the floor.
- Your arm should rest comfortably on a table or flat surface, palm up.

- Your cuff should be placed on your arm at the same level as your heart.
- Relax as much as possible.

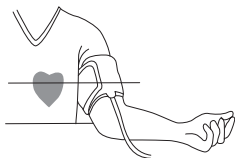III. MEASURING ARTERIAL PRESSURE

This product uses the oscillometric measurement method to detect arterial pressure. Before each measurement, the unit establishes a “zero pressure” equivalent to the air pressure. Then the product starts to inflate the cuff, while the unit detects the pressure oscillations generated by the pulses, beat by beat, which are used to determine the systolic and diastolic pressure as well as the pulse rate.

1. Press the button START STOP to switch on the device and measure the arterial pressure.
2. The blood pressure monitor will start to inflate for the measurement and display “00”.
3. Inflating and measuring

The value shown on the display will increase, before slowly decreasing, and the heart rate symbol will flash.

4. Display and storage of results

Once the measurement is complete, the systolic pressure, diastolic pressure and pulse will be shown on the display.

Consult your doctor if you obtain unexpected results.

Safety:

If you feel uncomfortable during the measurement, press the button START STOP measurement, press the button immediately to stop it.

If you experience any discomfort during a measurement, such as arm pain or other problems, press the button START STOP immediately to release the air from the cuff.

Recommendations:

- Do not measure your arterial pressure immediately if you feel like urinating, after exercising, bathing or eating or after drinking tea, coffee or smoking a cigarette. Rest for 20 to 30 minutes before taking the measurement.
- Rest for 5 minutes before the first measurement.
- Stay still and do not speak during the measurement.
- Wait at least 3 minutes between each measurement.
- Measure your arterial pressure at approximately the same time every day.

IV. TURNING OFF THE DEVICE

Press the button START STOP to switch off the blood pressure monitor. Otherwise, the monitor will switch off automatically after 1 minute.

Warnings:

- Avoid using your blood pressure monitor simultaneously on the same limb. Similarly, too frequent arterial pressure measurement can disrupt blood circulation, cause injuries and display inaccurate measurements.
- Please ensure that the device does not cause any prolonged impairment of blood circulation during operation.

- The device contains sensitive electronic components. To avoid measurement errors, avoid taking arterial pressure measurements in the vicinity of an interference signal emitted by a strong electromagnetic field or a fast transient/burst electric signal.
Examples: high-frequency (HF) surgical equipment; magnetic resonance imaging (MRI); scanner, computer-aided tomography (CT), etc.
- Wireless communications equipment such as wireless home network devices, mobile phones, cordless telephones and their base stations and walkie-talkies can cause interference that can affect the accuracy of measurements. A minimum distance of 1 foot (30 cm) from these devices must be maintained during a measurement.
- This device is portable, but not intended to be used during patient transport.
- This device is not suitable for continuous monitoring during emergencies or medical operations.
- This device is not AP/APG equipment. Do not use the device in the presence of flammable anaesthetics or in environments mixed with air or oxygen or nitrous oxide.

DATA MANAGEMENT

MEMORY FUNCTION

1. Each value measured is automatically stored under the previously selected user. This device can store up to 199 measurements for User 1 and User 2. When the memory is full, the old values are replaced by the new ones.
2. In off mode, press the button MEM once and the device will display the average value of the arterial pressure measurements of the last two or three readings. Press the button MEM again and the last measured value will be displayed. Press the button MEM and the last measured value will be displayed. Press the button again and the other values measured will be displayed one by one. Press the button SET to display the most recent measurements.
3. Clearing the memory

In off mode, press the button  to select the user whose measured values should be deleted.

Press the button  to switch off the device and press the button MEM once to activate the display. Next, press the button MEM for approx. 3 seconds to erase the memories of the selected user and the icon  will appear on the display.

DETECTION OF "CUFF" POSITION

The icon  will appear on the display when the cuff is wrapped correctly. If the cuff is too loose, the icon  loose, the icon will flash continuously to indicate a problem.

In this case, press the button  to stop the measurement and reposition the cuff correctly before restarting the measurement.

IMMOBILITY INDICATOR

The icon  flashes if you move or shake your arm during the measurement, which can lead to incorrect results. Please adjust your posture and start the measurement again.

DETECTION OF AN IRREGULAR HEARTBEAT

The heartbeat icon  flashes as soon as a pulse is detected during measurement.

The icon  will be displayed if an irregular heartbeat is detected.

WHO ARTERIAL PRESSURE INDICATOR READING

For information, please refer to the pressure classification table (mmHg unit) established by the WHO and the International Society of Hypertension (ISH) in 1999.

Colour display		BP classification:	Systolic blood pressure (mmHg)	Diastolic blood pressure (mmHg)
	Red	Severe hypertension	≥ 180	≥ 110
	Red	Moderate hypertension	160-179	100-109
	Red	Mild hypertension	140-159	90-99
	Orange	High normal value	130-139	85-89
	Green	Normal arterial pressure	120-129	80-84
	Green	Optimal arterial pressure	< 120	< 80

Reminder: This device enables you to monitor your blood pressure.

It is not intended to be a diagnostic device. Always contact your doctor for advice, diagnosis and treatment.

SPECIAL INDICATIONS

ERROR SYMBOLS

If the measurement is abnormal, one of the following symbols may appear. Please use the recommended measurement method.

Error	Cause/Solution
Er U	Pressure cannot reach 30 mmHg (4 kPa) in 12 seconds.
Er H	Inflation reaches 295 mmHg and deflates automatically after 20 ms.
Er 1	Pulse is not detected correctly.
Er 2	Too many disturbances (movement, conversation or magnetic disturbance during a measurement).
Er 3	The measurement result is abnormal.
Er 23	The SYS value is below 57 mmHg.
Er 24	The SYS value is greater than 255 mmHg.
Er 25	The DIA value is below 25 mmHg.
Er 26	The DIA value is greater than 195 mmHg.

ANOMALIES & TROUBLESHOOTING

Anomaly	Possible fault	Solution
No power supply	If the battery level is low	Replace batteries or insert power cable for power supply
	Batteries are upside down	Install batteries correctly
No pressurisation	The air intake is inserted incorrectly	Firmly insert the air intake
	The cuff is damaged	Please contact the retailer to replace it with a new cuff
Unable to measure due to display error	If the arm is moved when pressurising	Keep your arm and body still
	If you talk during the measurement	Do not make any noise during measurement
Slow pressurisation	The cuff is too loose	Please tighten the cuff
 If the arterial pressure can still not be measured after trying the solutions mentioned above, please contact the dealer. Do NOT attempt to disassemble the device yourself.		

SAFETY INFORMATION

	Caution
	TYPE BF APPLIED PART
	Symbol for marking of electrical and electronic equipment in accordance with the directive 2012/09 EU
	Refer to instruction manual
	Keep dry
	Low battery
	Keep out of direct sunlight
	Vertical upwards
IP21	2 protected against solid foreign bodies 12.5mm long $\varnothing$ and above; 1 protection against

CE 0123	CE
	Manufacturer
	Date of manufacture
SN	Serial number
LOT	Batch number
EU REP	Authorised representative in the European community
RoHs	RoHs Mark
MD	Medical appliance
UDI	Unique device identifier
	Importer in Switzerland
CH REP	Authorised representative in Switzerland

MAINTENANCE & CARE

CLEANING

- Clean the device and cuff with a soft, dry cloth. If necessary, use a damp cloth and a natural detergent. Do not use alcohol, benzene or other harsh chemicals.
- When cleaning, make sure that you never immerse any part of the blood pressure monitor in order to prevent liquid from entering the device.
- Do not clean the cuff in a washing machine or dishwasher!
- •Clean within 3 minutes each time. The number of consecutive cleans must not exceed three.
- The service life of the cuff may vary depending on washing frequency, skin condition and storage condition. Typical service life is 10,000 uses.

DISINFECTION

- The body of the device can be cleaned with a soft, clean cloth slightly dampened with 75% medical alcohol for disinfection.
- Avoid using disinfection methods such as high-temperature steam or ultraviolet radiation, as these can damage the device and reduce its service life.
- It is recommended that the display unit be disinfected before and after each use. The number of consecutive cleans must not exceed two.
- The blood pressure monitor must not be used again if liquid has penetrated and damaged the device or cuff.

STORAGE

- Make sure your device is clean and dry before storing it.
- Store the device and cuff in a clean, dry place and protect them from heat, lint, dust, direct sunlight and corrosive gases. Never place heavy objects on top of it. Avoid places where there is a high risk of vibration, shock, moisture or exposure to splashes, as this may damage the device.
- Remove the cuff from your device and fold it gently.
- Make sure that the rubber hose of the cuff is not crushed, stretched or kinked during storage.
- Store your device in its original cover or packaging. Store indoors at temperatures between -20°C and 55°C, 10% to 93% relative humidity, and 70 kPa to 106 kPa atmospheric air pressure. Make sure that it is stored in a well-ventilated area, free of corrosive gases.
- It may take up to 30 minutes for the device to warm up/cool down from the minimum/maximum storage temperature before it is ready for use.

MAINTENANCE

- No calibration is required.
- Do not attempt to repair the device yourself if it is not working properly. Repairs must only be carried out by authorised service centres.
- Do not use the device during maintenance or service.

WARRANTY AND ENVIRONMENTAL PROTECTION

This appliance is guaranteed against material and manufacturing defects for two years. During this period, the product will be replaced free of charge (proof of purchase must be presented if the product is under warranty). This warranty does not cover damage resulting from accidents, improper use or negligence. If you have a complaint, please contact the store where you purchased your product first.

Waste electrical products should not be disposed of with household waste.

Please recycle where facilities exist. Check with your local authority or retailer for recycling advice.

TECHNICAL CHARACTERISTICS

Model	ARM-30Y	
Display	LED screen	
Measurement method	Oscillometric measurement	
Measurement part	Upper arm	
Pneumatic pressure measurement range	0~295 mmHg (0~39.3 kPa)	
Maximum pressure protection	295 mmHg (39,3 kPa)	
Measurement range	Arterial pressure value	SYS : 57~255 mmHg (7,6~34,4 kPa) DIA : 25~195 mmHg (3,3~26,0 kPa)
	Pulse frequency	40~199 bpm
Precision	Arterial pressure value	±3 mmHg (±0,4 kPa)
	Pulse frequency	±5%
Low battery	When the voltage is below 4V ±0.1V, the device switches off automatically.	
Power source	4 x AAA batteries or 5V, 1A power supply	
Memory	2 users x 199 memories + guest mode	
Dimensions	140,8 mm (L) x 110 mm (l) x 61 mm (H)	
Screen size	73 mm (L) x 88 mm (l) (4,5 pouces)	
Cuff size	22~42 cm	
Weight	Approx. 265g (without batteries)	
Automatic Shut-Off	1 minute without use	

Degree of protection against electric shock	Type BF
Operating mode	Continuous operation
Protection against ingress of water or particles	IP21
Blood pressure monitor service life	5 years
Cuff service life	10,000 uses
Protection against electric shock	Internally powered ME equipment (battery-powered only) Class II ME equipment (when mains power is used)
Operating environment	Temperature: 5°C ~ 40°C Humidity: 15% ~ 90% RH Air pressure: 70 kPa ~ 106 kPa
Transport and storage conditions	Store in a clean and safe place. Storage temperature: -20°C ~ 55°C Humidity: 10% ~ 93% RH Air pressure: 70 kPa ~ 106 kPa Store in a well-ventilated area, away from corrosive gases.
Standard	The product has undergone clinical testing in compliance with ISO 81060-2
Essential performance	1. Measurement range (arterial pressure): - SYS: 57-255 mmHg - DIA: 25-195 mmHg - Pulse rate: 40-199 bpm 2. Accuracy: Arterial pressure: ± 3 mmHg (± 0.4 Kpa); Pulse rate: $\pm 5\%$
Power supply	The specified power supply must meet the following requirements: - Output voltage: DC 5V, - Output current: 1000 mA, - Class II - Complies with IEC 60601-1, - Must provide at least two layers of MOPP insulation between the AC input and DC output, in accordance with US and Canadian deviation requirements.

Declaration: "The ARM-30Y arm blood pressure monitor has been tested in accordance with the recommendations of technical report IEC TR 60601-4-2: Medical electrical equipment – Part 4-2: Guidance and interpretation – Electromagnetic immunity; performance of medical electrical equipment and medical electrical systems."

FAQS

What are systolic pressure and diastolic pressure?

Arterial pressure is produced by the heart beating as it contracts and relaxes. The maximum (systolic) and minimum (diastolic) arterial pressure can be measured.

Why is the arterial pressure value obtained at home different than that obtained in the hospital?

- The difference in arterial pressure between home and hospital measurements is approximately 20 mmHg - 30 mmHg (2.7 kPa - 4.0 kPa). This is because people tend to be more relaxed at home than in the hospital.
- Also, when the device is placed above the heart, the arterial pressure value tends to be much lower than it actually is. Make sure that the device is positioned at the same level as the heart.
- Pressure varies throughout the day, depending on time, physical effort and emotions. Furthermore, "white coat syndrome" often causes the pressure to rise when taken by a healthcare professional.
- The antihypertensive drug may have lost its efficacy. Please follow your doctor's instructions.
- The cuff may not be in the correct position. If the cuff is not placed correctly, no arterial pressure value will be obtained, and the arterial pressure value may be much higher than it is. Therefore, position the cuff correctly.
- The cuff is not tight enough. If the cuff is too loose, the compressive force may not be transferred to the artery and the arterial pressure value may be much higher than it is. Therefore, readjust and tighten the cuff further.
- The patient is not sitting correctly during the measurement. It is not advisable to squat, tilt, lean over or sit with legs crossed when taking arterial pressure due to increased abdominal pressure or the position of the arm below the heart.
- Please take the measurements in the correct position.

Why does my arterial pressure change during the day?

1. Arterial pressure is not stable, it varies during the day. It is therefore important that the measurements are taken under identical conditions: regularly and at a fixed time, in quiet conditions and in the same environment.
2. If you are taking medication, the measurements may vary further.

When can I get better measurements?

- It is best to take the measurements in the morning, just after urinating, or when your mind and body are stable. We recommend taking the measurements at the same time of day each time.

Why is the arterial pressure value measured different each time?

- During systole, the arterial pressure varies to a certain extent. For example, a person with a pulse rate of 70 beats per minute will see their arterial pressure vary from 100 to 80 mmHg per day. Since the arterial pressure changes constantly, it is difficult to obtain the correct arterial pressure value by measuring it only once. Therefore, the measurement should be taken two to three times. The first measurement will usually be higher due to nervousness or inadequate preparation; in the second measurement the nervous emotion will be slightly attenuated, so that the second measurement will typically be 5 mmHg to 10 mmHg (0.7 kPa to 1.3 kPa) lower the first time. This phenomenon is more obvious for people with high arterial pressure.

- When measuring continuously, please note that: There may be extravasation of blood because the arm is squeezed, causing blood from the fingertip to not flow smoothly. If you continue the measurement in the event of blood extravasation, you will not be able to obtain the correct value measured. Loosen the cuff, lift your hand over your head, grasp and stretch the left and right palms 15 times in a row. The extravasated blood can then be removed and you can continue measuring your arterial pressure.
- Cuff position and wrapping method. The value measured varies depending on the cuff size. In particular, if the cuff is wrapped around your elbow, you will not be able to obtain the correct measured value.
- Please use the correct cuff wrapping method to take the measurement. The arm circumference of the cuff is between 22 cm and 42 cm (centre of arm). If the model does not match, please purchase separately.

Are the results the same with the right arm?

- Arterial pressure may be measured on both arms, but the results may differ from person to person. We therefore recommend that you always measure your arterial pressure on the same arm.

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EMC GUIDANCE

The ME EQUIPMENT or ME SYSTEM is suitable for home healthcare environments.

Warning: Don't be near the active HF surgical equipment and the RF shielded room of an ME system for magnetic resonance imaging, where the intensity of EM disturbances is high.

Warning: Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.

Warning: Use of accessories, transducers and cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

Warning: Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the equipment TMB-1490-BHJ including cables specified by the

Technical description:

1. All necessary instructions for maintaining BASIC SAFETY and ESSENTIAL PERFORMANCE with regard to electromagnetic disturbances for the expected service life.
2. Guidance and manufacturer's declaration-electromagnetic emissions and Immunity.

Table 1

Guidance and manufacturer's declaration - electromagnetic emissions	
Emissions test	Compliance
RF emissions CISPR 11	Group 1
RF emissions CISPR 11	Class B
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations / flicker emissions IEC 61000-3-3	Comply

Table 2

Guidance and manufacturer's declaration – electromagnetic Immunity		
Immunity Test	IEC 60601-1-2 Test level	Compliance level
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Electrical fast transient/burst IEC 61000-4-4	±2 kV for power supply lines ±1 kV signal input/output 100 kHz repetition frequency	±2 kV for power supply lines Not Applicable 100 kHz repetition frequency
Surge IEC61000-4-5	±0.5 kV, ±1 kV differential mode ±0.5 kV, ±1 kV, ±2 kV common mode	±0.5 kV, ±1 kV differential mode Not Applicable
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% UT; 0,5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°. 0% UT; 1 cycle and 70% UT; 25/30 cycles; Single phase: at 0°. 0% UT; 250 / 300 cycle	0% UT; 0,5 cycle. At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°. 0% UT; 1 cycle and 70% UT; 25/30 cycles; Single phase: at 0°. 0% UT; 250 / 300 cycle
Power frequency magnetic field IEC 61000-4-8	30 A/m 50 Hz / 60 Hz	30 A/m 50 Hz / 60 Hz
Conducted RF IEC61000-4-6	3 V 0,15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80% AM at 1 kHz	3 V 0,15 MHz – 80 MHz 6 V in ISM and amateur radio bands between 0,15 MHz and 80 MHz 80% AM at 1 kHz
Radiated RF IEC61000-4-3	10 V/m 80 MHz – 2,7 GHz 80% AM at 1 kHz	10 V/m 80 MHz – 2,7 GHz 80% AM at 1 kHz
NOTE U _T is the a.c. mains voltage prior to application of the test level.		

Table 3

Guidance and manufacturer's declaration - electromagnetic Immunity							
Radiated RF IEC61000-4-3 (Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment)	Test Frequency (MHz)	Band (MHz)	Service	Modulation	Modulation (W)	Distance (m)	IMMUNITY TEST LEVEL (V/m)
	385	380-390	TETRA 400	Pulse modulation 18Hz	1.8	0.3	27
	450	430-470	GMRS 460, FRS 460	FM $\pm$ 5kHz deviation 1kHz sine	2	0.3	28
	710	704-787	LTE Band 13, 17	Pulse modulation 217Hz	0.2	0.3	9
	745						
	780						
	810	800-960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation 18Hz	2	0.3	28
	870						
	930						
	1720	1700-1990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4,25; UMTS	Pulse modulation 217Hz	2	0.3	28
	1845						
	1970						
	2450	2400-2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation 217 Hz	2	0.3	28
	5240	5100-5800	WLAN 802.11 a/n	Pulse modulation 217 Hz	0.2	0.3	9
	5500						
	5785						

Terraillon SAS France & Headquarters

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Terraillon Asia Pacific Ltd

4/F, Eastern Centre
1065 King's Road
Quarry Bay – Hong Kong
Tel: +852 (0)2960 7200
customerservice@terraillon-asia.com

Terraillon Corp USA

contact@terraillon.fr

Terraillon®
www.terraillon.com