Wrist Ambulatory Blood Pressure Monitor

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

Models: RIV-B10/RIV-B11/RIV-B12

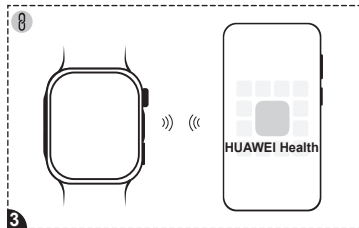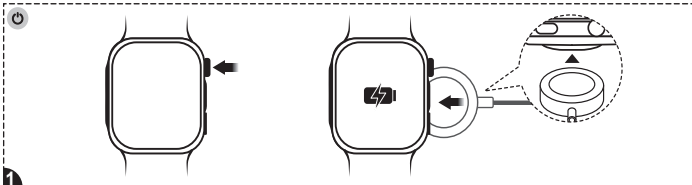Quick Start Guide

This device is intended for use in home healthcare settings. It is designed for adults to measure blood pressure at the wrist. If a defibrillator is used on a person wearing this device, its blood pressure measurement function may be affected. The Blood Pressure (BP) and Ambulatory Blood Pressure Monitoring (ABPM) features on this device have obtained CE certification in compliance with the Regulation (EU) 2017/745. The ECG app and Pulse Wave A-fib Prompt (PWAP) app have also obtained CE certification in compliance with the Regulation (EU) 2017/745. The device and two apps are three separate medical products, but work together as a single unit, which complies with Article 22 of Regulation (EU) 2017/745.

If you have any questions regarding this device, please seek assistance through official channels, such as the Huawei official website <https://consumer.huawei.com/en/>.

Warnings

- This device has not been validated for use on newborn babies or pregnant individuals (including those with pre-eclampsia).
- Do not use in emergency care settings.
- If your wrist circumference is less than 13.0 cm or more than 21.0 cm, blood pressure measurement accuracy may be affected.
- NEVER diagnose, treat yourself, or adjust medications based on the readings. ALWAYS consult your physician and take medicines as prescribed.
- Do not use if allergic to nylon, fluoroelastomer, or hydrogenated nitrile butadiene rubber (HNBR).
- Do not wear during intravenous infusions, blood transfusions, or when a catheter is inserted through the wrist. Doing so may restrict blood flow and cause injury.
- Do not use on a wounded arm or an arm undergoing medical treatment as this may cause further injury.
- If you have had a mastectomy or lymph node dissection, consult your physician before using on the affected arm.
- If you are undergoing intravascular intervention or treatment, or have an arteriovenous (A-V) shunt, consult your physician before using on the same arm. Using the device in these circumstances may temporarily restrict blood flow and cause injury.
- Ensure the device is not restricting blood flow by observing the arm during blood pressure measurements.
- Do not take blood pressure measurements too often as bruising may occur when your blood flow is restricted.
- Prolonged wear may cause friction or pressure that irritates your skin. Remove the device from time to time to let your wrist rest.
- The device may temporarily affect other electrical medical monitoring equipment used simultaneously on the same arm.


- Do not place or use near portable or mobile radio frequency (RF) communication equipment, high frequency surgical equipment, magnetic resonance imaging (MRI) equipment, or computed tomography (CT) scanners—this may affect accuracy.
- Do not use in an oxygen-rich environment or near flammable gases.
- Keep the charging cradle and power cable away from newborns, infants, and children to minimize strangulation risk.
- Keep this device away from newborns, infants, and children as it contains small parts that may be a choking hazard.
- Do not disassemble or tamper with this device, as this may damage it and/or cause injury.
- This device has a non-removable battery. Do not try to replace it, as doing so may damage the device and/or cause injury.

Precautions

- Keep this device clean. Stop using immediately and consult your physician if it causes skin irritation (e.g., from sweat or contamination).
- If you have an allergic reaction, stop using and consult your physician immediately.
- If the strap does not deflate after a blood pressure measurement, remove the strap.
- Do not use this device with other medical electrical equipment, as this may affect its functions and accuracy.
- Consult your physician before using if you have arrhythmias (e.g., atrial or ventricular premature beats, or atrial fibrillation), arteriosclerosis, malperfusion, diabetes, or kidney disease. Any of these conditions may affect blood pressure measurement accuracy.
- If you have atherosclerosis, arterial blockage, or if measurement conditions are not met, blood pressure readings on your wrist and upper arm may vary significantly. Adjust your position according to measurement requirements or consult your physician.
- Before using, consult your physician if you have excess wrist soft tissue, a weak pulse, or other factors that could affect measurement accuracy.
- Moving or trembling during measurement may affect accuracy.
- This device has been carefully tested. It meets required accuracy standards and is approved for use throughout its expected service life. Performance may deteriorate beyond that period.
- Users of this device should be intended operators (patients or assistants). No special expertise is required. Features should be used in accordance with this guide and with the device worn properly. Multiple users may use it at different times in guest mode; no special cleaning or disinfection is required before use.
- By default, this device runs and monitors data in the background and automatically uploads related logs.
- To prevent misoperation or damage, keep the device out of reach of children and pets.


1. Structural Diagram

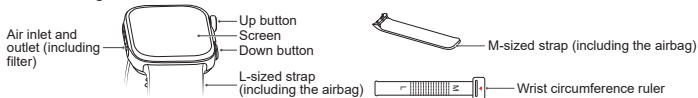
Notes: Press the Up button to power this device on, off, or restart it. Press the Down button to start a blood pressure measurement.

2. Selecting a Strap

Select the strap size based on the scale on the ruler.

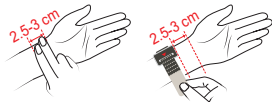
If the scale points to a number between 12 and 20, select the L-sized strap. (Ensure that the arrow on the inner side of the watch clasp points to the same number indicated by the arrow on the ruler.)

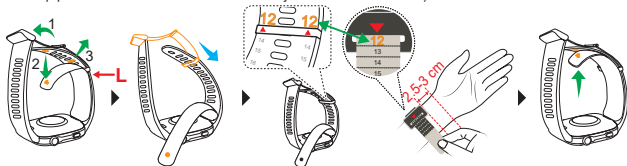Ensure that the arrow on the inner side of the watch clasp points to the same number indicated by the arrow on the ruler.


If the scale points to a number between 5 and 11, select the M-sized strap. (Replace the strap with the M-sized strap, and ensure that the arrow on the inner side of the watch clasp points to the same number indicated by the arrow on the ruler.)

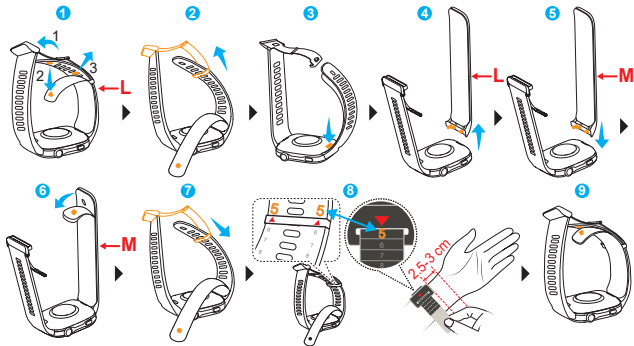Ensure that the arrow on the inner side of the watch clasp points to the same number indicated by the arrow on the ruler.

3. Pairing and Connection

1. Press and hold the Up button to power on the device or charge the device to power it on. (Line up the charger with the charging area on the back of the device, plug the cable into a power adapter, and connect the power adapter to a power source. Wait until the battery level appears on the screen.)

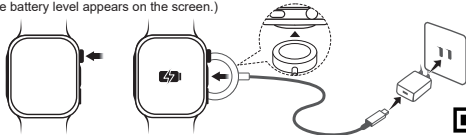2. Use your phone to scan the QR code, or search for "HUAWEI Health" in the AppGallery or on consumer.huawei.com to download the HUAWEI Health app.
3. Open the HUAWEI Health app, go to the devices screen, choose the devices you wish to add, and follow the onscreen instructions to complete pairing.

4. Blood Pressure Measurement


4.1 Measurement Requirements

Measurement results may be affected by position, posture, movements, and physiological conditions. To ensure the accuracy:


- Avoid eating, alcohol, smoking, and exercise for at least 30 minutes before a single measurement. Rest for 5 minutes before measuring and keep your body relaxed throughout. Avoid measuring when you feel emotional or stressed.
- Do not speak or move your fingers, arms, or body during measurement.
- If taking multiple blood pressure measurements in a row, rest your wrist for 1-2 minutes between each to allow arteries to return to their regular state.
- If your result falls outside the normal range, rest for 5 minutes and measure again or consult your physician.
- If your device gets wet, wipe it immediately and dry it before measuring. Failure to do so may cause water infiltration, leading to measurement failure or slow airbag deflation.
- When using at night, wear the device on your wrist as instructed, with wrist naturally positioned on one side of the body and no foreign objects pressing against the strap.
- For patients unable to lie flat at night, position the device at the same level as the heart. If this is not possible and the device must be used under special circumstances, consult your physician first.

During measurement, sit comfortably with feet flat on the floor, legs uncrossed, back and arm supported, and the device at heart level. Avoid pressing your chest. Extend your wrist and palm naturally; do not clench your fist. Maintain this posture throughout the measurement. Any deviation may affect the result.

Standard posture for a single measurement

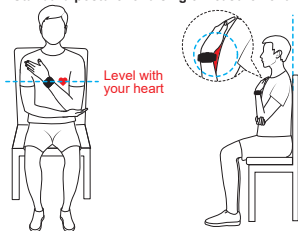Standard posture for ABPM/pre- and post-exercise measurements

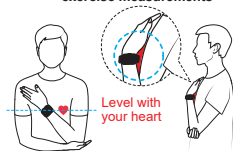Note: This device will automatically inflate and perform measurements while you're asleep, not requiring you to assume the standard posture.

Common incorrect postures


Suspended

Making a fist


Shrugging

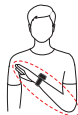Unsupported or everted

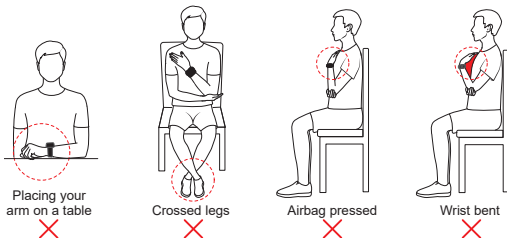4.2 Measuring Your Blood Pressure

If you're using the device's blood pressure monitoring app for the first time, follow the onscreen instructions to complete setup, including setting your wrist circumference.

Notes:

- During blood pressure measurement, press the Up button or tap the cancel button on the screen to exit the measurement. If the device does not automatically deflate, please unfasten the strap.
- After the measurement is done, the device deflates automatically and shows the result.
- Blood pressure cannot be measured when the battery is below 5% (approx. 3.4 V).
- In your results, SYS is systolic, DIA is diastolic, and pulse is your pulse rate.
- All measurement modes will display error alerts if an issue occurs. Please follow the onscreen instructions.


4.2.1 Performing a single blood pressure measurement

Go to the app list on your device, select the blood pressure app, and start a measurement.

Note: Measurements include high blood pressure (systolic), low blood pressure (diastolic), and pulse rate.

4.2.2 Performing ABPM

1. Go to the app list on your device, select ABPM, and set up a monitoring plan (start time, measurement frequency, and day/night measurement schedules).
2. Once the plan starts, the device will notify you when a measurement is due.

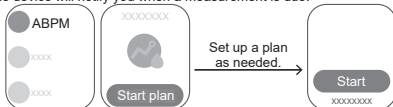Note: Measurements include 24-hour, daytime, and nighttime averages, nocturnal dip, valid measurements, maximum, median, minimum, overall average, standard deviation, morning surge, blood pressure load, and coefficient of variation.

4.2.3 Performing blood pressure measurements pre- and post-exercise

1. Go to the app list on your device, select the Pre- & Post-exercise BPM app, review the description, and follow the onscreen instructions.
2. On the measurement screen, start the measurement.
3. Once complete, begin your workout: select the exercise type, agree to the statements, and start.
4. When prompted, measure your blood pressure again within 120 seconds of finishing your workout. Once complete, your results and recommendations will appear on the screen.

Notes:

- This feature measures blood pressure at rest pre- and post-exercise. Do not use this feature while working out. It helps self-assess exercise-induced blood pressure changes, abnormalities, and potential health risks. Consult a doctor for medical advice if needed.
- Measurements include systolic and diastolic pressure pre- and post-exercise.

5. Blood Pressure Information

5.1 Principles

This device measures blood pressure via the oscillometric method. As the strap (airbag) pressure increases, the pressure sensor in the device strap monitors and analyzes pressure oscillations, automatically calculating your systolic and diastolic blood pressure and pulse rate.

5.2 Operating, Storage, and Transportation Requirements

Category	Operating Requirement	Storage and Transportation Requirement
Temperature	5°C to 40°C	-20°C to +55°C
Humidity	10% RH to 95% RH (non-condensing)	
Atmospheric pressure	70 kPa to 106 kPa	

Notes:

- This device transmits update packages and blood pressure data via Bluetooth. If the Bluetooth connection is lost, data sync will be interrupted. This device only supports a single device connection at a time. To switch to another device, re-store this device to factory settings first.
- This device is not intended for use in public settings.
- Blood pressure cannot be measured while charging.
- Avoid using on moving vehicles (e.g., cars, trains, or airplanes).
- Use in the recommended environment to ensure optimal performance and prevent damage. Readings may be inaccurate in extreme temperatures, humidity, or atmospheric pressure.
- If stored at maximum or minimum temperature limits, allow the device to acclimate to the specified operating temperature range for about 1 hour before use to prevent malfunctions.
- Store in a clean, secure environment. Avoid areas with chemical exposure, direct sunlight, high humidity, high temperatures, dust, corrosive gases, or salt. Also avoid strong impacts, collisions, rain, tilting, vibrations, or physical shocks.
- When using this device at a maximum operating temperature of 40°C, the highest temperature of the applied part may reach 41.5°C.

5.3 Measurement Irregularities and Solutions

Fault Type	Solution
Measurement timeout	If the fault persists after several attempts, contact an authorized Huawei service center.
Inflation outside the rated range	If the fault persists after several attempts, contact an authorized Huawei service center.
Inflation exception	Check whether the installed airbag is connected improperly and whether it is damaged.
Unsuitable airbag	Follow the instructions to measure your wrist circumference, install a strap (including the airbag), and set your wrist circumference correctly.

Slow deflation	Remove the device immediately and follow the onscreen instructions or contact an authorized Huawei service center.
----------------	--

5.4 GUI Abbreviations

Abbreviation	SYS	DIA	cm	in	ABPM	BP	bpm	avg
Full Name	Systolic blood pressure	Diastolic blood pressure	Centimeter	Inch	Ambulatory Blood Pressure Monitoring	Blood Pressure	beats per minute	averag

5.5 Static Pressure Test Mode

This involves dedicated personnel testing the static performance of the device through a stringent process, after entering the device's static pressure test mode.

5.6 Others

- **Storage function:** This device supports the storage and review of blood pressure data, with specific storage limits based on the type of measurement. For single blood pressure measurements, up to 100 records can be stored, while ambulatory blood pressure monitoring can store up to 300 records. Additionally, pre- and post-exercise blood pressure measurements can store up to 30 records. Once the storage capacity is reached, the earliest records are automatically replaced to make room for new data. Upon starting a new ambulatory blood pressure monitoring plan, all existing records from the previous plan will be cleared.
- **Deleting data:** Data deletion can only be performed by restoring the device to its factory settings. No other options are available for deleting data.

6. Maintenance and Care

6.1 Maintenance

- **Airbags and filters** are consumables. When used daily, they are designed to function reliably for six months from the date of first use. If an airbag or filter is damaged due to non-human causes during the warranty period, visit an authorized Huawei Customer Service Center for warranty support.
* **Tips:** Keep the airbag and filter clean and dry. Proper usage and maintenance ensure optimal performance and comfort.
- Do not attempt maintenance during blood pressure measurement. **Warning:** Replacing an original component with a non-Huawei component may cause measurement errors.
- This device is equipped with a non-removable battery. Do not try to replace or remove it. Doing so may damage this device and/or cause injury.
- The blood pressure-related software is pre-installed on the device. The device supports over-the-air updates (when paired with a phone). The software updates whenever the device updates.
- The pressure measurement system has been verified; no need to calibrate it yourself. If required, Huawei can provide circuit diagrams, component lists, illustrations, calibration details, and other repair-related information to Huawei-authorized maintenance personnel. For assistance, please contact an authorized Huawei service center.

6.2 Cleaning and Maintenance

- Keep the device clean and dry. If the device body or strap (including airbag) is dirty, use a soft, dry cloth to clean it. For stubborn dirt, use a cloth dampened with clean water or 75% rubbing alcohol (fully wrung out). Avoid leaving alcohol on the airbag or soaking it for extended periods.
- Remove the device before prolonged water immersion (e.g., swimming, diving) or fast-flowing water (e.g., car washing). Also, remove before entering high-temperature and/or high-humidity environments (e.g., hot showers, hot springs, saunas/steam rooms). Keep away from shower gel, shampoo, sunscreen, hand sanitizer, and other household chemicals. If contact occurs, clean and dry promptly with a soft cloth.

Precautions

- When cleaning the device and strap (including airbag), ensure no liquid enters the device.
- Do not immerse the device body or strap (including airbag) in water, and do not use heat to dry.
- Do not use organic solvents (e.g., hand sanitizer, soapy water, other acidic solvents), strong detergents, or strong disinfectants to clean the device or strap (including airbag).
- Do not use sharp or hard objects (e.g., steel wool) to wipe the device or its accessories.

7. More Information

7.1 Structure and Composition

This device consists of a device body, a strap (with an airbag inside), and a wireless charger.

7.2 Intended Purpose

This device is a digital self-monitor intended for use in measuring blood pressure, 24-hour ambulatory blood pressure, and pulse rate in adult population with wrist circumference ranging from 13.0 cm to 21.0 cm.

7.3 Contraindications

No contraindications have been found.

7.4 Clinical Benefit

Provides 24-hour ambulatory blood pressure monitoring throughout your daily activities, as well as pre- and post-exercise blood pressure measurement, to help manage hypertension by:

- Assisting healthcare professionals in identifying patterns of hypertension.
- Aiding healthcare professionals in assessing the effectiveness of hypertension treatment.
- Evaluating exercise-induced blood pressure change to facilitate early identification of cardiovascular risks.

7.5 Data Transmission

This device transmits 2.4-GHz RF signals for device updates and data transmission. Do not use it in places where the use of RF is restricted. After Bluetooth is disconnected, blood pressure measurements continue to work as normal but the software version will fail to update and data cannot be synced.

7.6 Technical Specifications

Name	Wrist Ambulatory Blood Pressure Monitor
Models	RIV-B10/RIV-B11/RIV-B12
Software Component Release Version	V1
Bluetooth	Bluetooth 6.0 and compatible versions
Wireless Communication	Frequency range: 2.4 GHz (2,402–2,480 MHz) Modulation: GFSK, pi/4-DQPSK, and pi/8-DPSK Effective radiated power: < 20 dBm
Hardware Configuration	Frequency ≤ 266 MHz, SRAM 11.5 MB, NORFLASH 16 MB × 2
Type of Protection Against Electric Shock	Internal power supply
Applied Parts	Airbag and bottom cover (Classification type: Type BF)
Level of Protection Against Penetration by Harmful Liquids or Particles	IP68 (The device provides the level of protection according to the IP68 rating under the IEC 60529 standard. It is dust-resistant, and can be safely immersed in water at a depth of 1.5 meters for 30 minutes.)
Level of Safety When Exposed to Flammable Anesthetic Gasses Mixed with Air, Oxygen, or Nitrous Oxide	Non-AP/APG type
Running Mode	Continuous running
Power Supply	Rechargeable lithium-ion polymer battery (rated voltage: 3.83 V) Operating time: When fully charged and with ambulatory blood pressure monitoring enabled, the device can perform a full cycle (24-hour) of blood pressure monitoring. Typical battery lifespan: A battery on a full charge will last for at least 500 cycles when being used at room temperature (25°C). The battery level is shown as a percentage. Simply swipe down on the watch home screen to check the battery level.
Expected Service Life	Device body: 3 years
Measurement Method	Oscillometric method
Pressure Measurement Range	0–300 mmHg (0–40 kPa)
Pulse Rate Measurement Range	40–180 times/min
Accuracy	Pressure: ±3 mmHg (±0.4 kPa); Pulse rate: ±5%

Blood Pressure Measurement Range	Diastolic blood pressure: 40–130 mmHg (5.3–17.3 kPa) Systolic blood pressure: 60–230 mmHg (8.0–30.7 kPa)
Inflation Method	Automatic pressure pump
Deflation Method	Automatic fast deflation
Performance Efficiency	Single blood pressure measurement time < 180 seconds
Resolution	Blood pressure display: 1 mmHg or 0.1 kPa; Pulse rate display: 1 time/min
Range of Wrist Circumference	M-sized strap: 13.0–16.0 cm; L-sized strap: 16.1–21.0 cm Please select the appropriate strap (including the airbag) for your wrist circumference. If your wrist circumference exceeds 21.0 cm, the accuracy of blood pressure measurements may be affected.
This device complies with IEC 80601-2-30, and has been clinically tested in accordance with ISO 81060-2:2018 & AMD1:2020 & AMD2:2024. The clinical trial included: A resting blood pressure test on 85+ subjects; an ambulatory blood pressure test on 35+ subjects; an exercise stress test on 35+ subjects. For all three tests, the average deviation of systolic and diastolic pressures did not exceed ± 5 mmHg, and the standard deviation did not exceed 8 mmHg.	

7.7 Packaging List

Device (with an L-sized strap (including the airbag)) × 1	M-sized strap (including the airbag) × 1	Wrist circumference ruler × 1
Quick start guide/Warranty card × 1	Water-resistance notification card × 1	Charger (with a 0.77 m shielded cable) × 1
Note: Actual packaging may vary.		

8. Explanations of Symbols

	Type BF Applied Part
	Caution
IP68	Indicates that the product is dust-proof and water-resistant.
	Indicates that the product is capable of being recycled.
	Short-distance radio technology, Bluetooth certification mark.

CE 0197	CE mark. Indicates that the device complies with EU 2017/745.
	Indicates that this device is endorsed by the Forest Stewardship Council, a non-profit NGO committed to promoting sustainable forest management.
	Indicates that you should recycle the object when you dispose of it.

	Non-ionizing electromagnetic radiation.
	Medical device
	Follow this guide carefully when using this device.
	Authorized representative in the European Union.
	Manufacturer
	The temperature range that the device can safely tolerate.
	The humidity range that the device can safely tolerate.
	The atmospheric pressure range that the device can safely tolerate.
	Serial number
	Model number

	Unique device identifier
	Date of manufacture
	How to wear the device properly
 13.0 cm ~ 21.0 cm (5.12 in ~ 8.27 in)	Supported wrist circumferences
	Remove the packaging and dispose of it in an appropriate recycling bin according to local laws and regulations.
	Recyclable non-corrugated fiberboard (board)
	Importer

Note: The two digits after the 16-digit characters on the back of the product indicate the manufacturing year and month. The first digit is the last digit of the year, the second digit indicates the month, and the months from 1 to 12 are represented by 1 to 9 and A, B, and C, respectively. For example, 68 corresponds to August 2026.

In accordance with Huawei's device power supply standards, the charging cradle's serial number (SN) consists of the BOM code and production line details. For further assistance, visit Huawei's official website.

9. Environmental Protection Statement

- Do not dispose of this device and its accessories in a garbage can.
- Please comply with all relevant laws and regulations regarding the disposal of this device and its accessories. Support recycling efforts.
- This device and its accessories shall be disposed of according to local environmental protection regulations.

10. Declaration of Electromagnetic Compatibility

This device complies with the IEC 60601-1-2:2014/AMD1:2020 Electromagnetic Compatibility (EMC) standard. For details, please visit <https://consumer.huawei.com/en/>. While other devices may meet national emission requirements, ensure sufficient spacing between this device and other devices to prevent interference and maintain normal operation.

Essential performance

- Within the temperature range of 5°C to 40°C and the relative humidity range of 10% to 95% (non-condensing), the maximum error for the measurement of the cuff pressure at any point of the nominal measurement range shall be less than or equal to ± 3 mmHg (± 0.4 kPa) or $\pm 2\%$ of the reading, whichever is greater.
- Blood pressure measurements done under laboratory conditions (for device testing) must be reproducible to within 3.0 mmHg (0.4 kPa).

Warning

- Install and use this device in accordance with the EMC guidelines provided.
- Portable and mobile radio frequency (RF) communications equipment may affect the performance of this device. If improper operation occurs, which may indicate a loss or degradation of essential performance, please contact an authorized Huawei service center.
- Do not use accessories or cables other than those supplied with this device. Doing so may increase electromagnetic emissions or decrease electromagnetic immunity, potentially compromising normal operation.
- Do not stack this device or place it adjacent to other equipment. If such use is unavoidable, monitor the device closely to ensure normal functionality.
- This device may experience interference from other devices, even if those devices comply with applicable national emission standards.

Copyright © Huawei Device Co., Ltd. 2026. All rights reserved.

THIS DOCUMENT IS FOR INFORMATION PURPOSE ONLY, AND DOES NOT CONSTITUTE ANY KIND OF WARRANTIES.

Trademarks and Permissions

The *Bluetooth*® word mark and logos are registered trademarks owned by *Bluetooth SIG, Inc.* and any use of such marks by Huawei Technologies Co., Ltd. is under license. Huawei Device Co., Ltd. is an affiliate of Huawei Technologies Co., Ltd. Other trademarks, product, service and company names mentioned may be the property of their respective owners.

Privacy Policy

To better understand how we protect your personal information, please see the privacy policy at <https://consumer.huawei.com/privacy-policy>.

Blasting Caps and Areas

Turn off your mobile phone or wireless device when in a blasting area or in areas posted turn off “two-way radios” or “electronic devices” to avoid interfering with blasting operations.

Petrol stations and explosive atmospheres

In locations with potentially explosive atmospheres, obey all posted signs to turn off wireless devices such as your phone or other radio equipment. Areas with potentially explosive atmospheres include fuelling areas, below decks on boats, fuel or chemical transfer or storage facilities, areas where the air contains chemicals or particles, such as grain, dust, or metal powders.

Operation and safety

[Warning] Please read all of the safety information carefully before using your device to ensure its safe and proper operation and to learn how to dispose of your device properly.

- Using an unapproved or incompatible power adapter, charger or battery may cause fire, explosion or other hazards.
-  To prevent possible hearing damage, do not listen at high volume levels for long periods.
- Ideal temperatures: 5 °C to 40 °C for operating, -20 °C to +55 °C for storage.
- Avoid placing metal objects on the surface of the wireless charging cradle, as they could cause the equipment to heat up and malfunction. Only use the wireless charging cradle in environments with an ambient temperature ranging from 0°C to 45°C.
- If the ambient temperature is too high, the wireless charger may stop working. If this occurs, you are advised to stop wireless charging.
- To keep the wireless charger in good condition, avoid using it inside a vehicle.
- Adult supervision is advised if children are to use or come into contact with this product.
- Please consult your doctor and the device manufacturer to determine if operation of your device may interfere with the operation of your medical device.
- Switch off your wireless device when requested to do so in hospitals, clinics or health care facilities. These requests are designed to prevent possible interference with sensitive medical equipment.
- Switch off your wireless device whenever you are instructed to do so by airport or airline staff. Consult the airline staff about the use of wireless devices on board the aircraft, if your device offers a ‘flight mode’ this must be enabled prior to boarding an aircraft.

- Some wireless devices may affect the performance of hearing aids or pacemakers. Consult your service provider for more information.
- Pacemaker manufacturers recommend that a minimum distance of 15 cm be maintained between a device and a pacemaker to prevent potential interference with the pacemaker. If you are using a pacemaker, hold the device on the side opposite the pacemaker and do not carry the device in your front pocket.
- Keep the battery away from fire, excessive heat, extremely low air pressure and direct sunlight. Do not place it on or in heating devices. Do not disassemble, modify, throw, or squeeze it. Do not insert foreign objects into it, submerge it in liquids, or expose it to external force or pressure, as this may cause it to leak, overheat, catch fire, or even explode.
- Observe local laws and regulations while using the device. To reduce the risk of accidents, do not use your wireless device while driving.
- To prevent damage to your device's parts or internal circuits, do not use it in dusty, smoky, damp or dirty environments or near magnetic fields.
- Do not use, store or transport the device where flammables or explosives are stored (in a gas station, oil depot, or chemical plant, for example). Using your device in these environments increases the risk of explosion or fire.
- Dispose of this device, the battery, and accessories according to local regulations. They should not be disposed of in normal household waste. Improper battery use may lead to fire, explosion, or other hazards.
- Do not attempt to replace the battery yourself — you may damage the battery, which could cause overheating, fire, and injury. The built-in battery in your device should be serviced by Huawei or an authorised service provider.
- If you experience skin discomfort when wearing the device, take it off and consult a doctor.
- If the battery leaks, ensure that the electrolyte does not make direct contact with your skin or eyes. If the electrolyte touches your skin or splashes into your eyes, immediately flush with clean water and consult a doctor.
- Do not place sharp metal objects, such as pins, near the microphone. The microphone may attract these objects and result in injury.
- Wireless devices may interfere with the airplane's flight system. Do not use your device where wireless devices are not allowed according to the airplane company's regulations.
- To ensure the best audio performance of your device's microphone, use a cloth to wipe the area of the microphone which has been exposed to water. Do not use the microphone until the wet areas are totally dry.
- Protect the wearable and its accessories from strong impact, strong vibration, scratches, and sharp objects, as

these may damage the device.

- Make sure that the wearable and its accessories stay clean and dry.
- Ensure that the power adapter meets the safety requirements of IEC/EN 62368-1 and has been tested and approved according to national or local standards.
- Your device meets requirements of IP68 (IEC 60529).
- Please report to the manufacturer and the competent authority of the Member State in which you are established about any serious incident that has occurred in relation to this device.
- This device contains a coin, button, or small-sized cell battery. Do not swallow the battery, as it constitutes a chemical burn hazard.
- If the coin, button, or small size cell battery is swallowed, it can cause severe internal burns in just 2 hours and can lead to death.
- Keep new and used batteries away from children. If the battery compartment does not close securely, stop using the product and keep it away from children. If you think batteries might have been swallowed or placed inside any part of the body, seek immediate medical attention.

Disposal and recycling information

The symbol indicates that products and batteries should be disposed of at designated local waste collection points when no longer usable. This will ensure that EEE waste is recycled and treated in a manner that conserves valuable materials and protects human health and the environment. For details, contact local authorities, retailers, waste disposal services, or visit <https://consumer.huawei.com/en/>.

Reduction of hazardous substances

Complies with EU REACH and RoHS requirements.

Full conformity details: <https://consumer.huawei.com/en/legal/certification/>.

UK PSTI Compliance

This product meets the relevant security requirements specified in Schedule 1 of the Product Security and Telecommunications Infrastructure (Security Requirements for Relevant Connectable Products) Regulations 2023. The full text of the PSTI Statement of Compliance is available at the following website: <https://consumer.huawei.com/en/legal/certification/>.

RF exposure requirements

SAR limits: 2.0 W/kg (head, 10 g), 4.0 W/kg (limbs, 10 g).

Operating distance: When wearing on the wrist, keep at least 1 cm away from the face.

Highest reported SAR value: Head SAR: 0.04 W/kg; Limbs SAR: 0.28 W/kg.

EU regulatory conformance

Complies with the Radio Equipment Directive (RED) 2014/53/EU and the Medical Device Regulation (MDR) 2017/745/EU. Full Declaration of Conformity: <https://consumer.huawei.com/en/legal/certification/>.

Frequency bands and power

Frequency bands in which the radio equipment operates: Some bands may not be available in all countries or all areas. Please contact the local carrier for more details.

Bluetooth 2.4GHz: 17 dBm, NFC 13.56 MHz: 42 dBμA/m at 10 m, Watch Wireless Charger 110.5 kHz-119 kHz, 135 kHz-140 kHz: 42dBμA/m at 10 m, 119 kHz-135 kHz: 66 dBμA/m descending 10 dB/dec above 119 kHz at 10m, 140 kHz-148 kHz: 37.7dBμA/m at 10 m.

Please go to **Settings > About > Regulatory information** on the device to view the E-label screen.

Restrictions in Norway

Use of this radio equipment is not allowed in the geographical area within a radius of 20 km from the center of Ny-Alesund, Svalbard.

CE0197

Lotus NL B.V.

Koningin Julianaplein 10, 1e Verd, 2595AA, The Hague, Netherlands.

Huawei Device Co., Ltd.

No.2 of Xincheng Road, Songshan Lake Zone, Dongguan, 523808 Guangdong,

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Copyright © Huawei Device Co., Ltd. 2026. Alle Rechte vorbehalten.
Copyright © Huawei Device Co., Ltd. 2026. Tutti i diritti riservati.
Copyright © Huawei Device Co., Ltd. 2026. Todos los derechos reservados.
Copyright © Huawei Device Co., Ltd. 2026. Wszystkie prawa zastrzeżone.
Copyright © Huawei Device Co., Ltd. 2026. Toate drepturile rezervate.
Copyright © Huawei Device Co., Ltd. 2026. Kaikki oikeudet pidätetään.