

# Smart Heated Belt


## A Welcome Guide




# Product Description

**The Vulpés "BellyBelt" is a smart heated belt that has been developed to help with belly, abdominal, and menstrual pain. Its warmth supports the relaxation and dissolution of painful cramps in the upper and lower abdominal area. The belt can be used for both, the abdominal as well as the back-pelvic area by placing the heating pad to the desired area.**


Vulpés BellyBelt is made of Teddy plush fabric with velvet texture to give comfort. Teddy Superflausch is a soft and warm lining with skin-friendly characteristics. The fabric is insulating, durable, and cuddly at the same time.

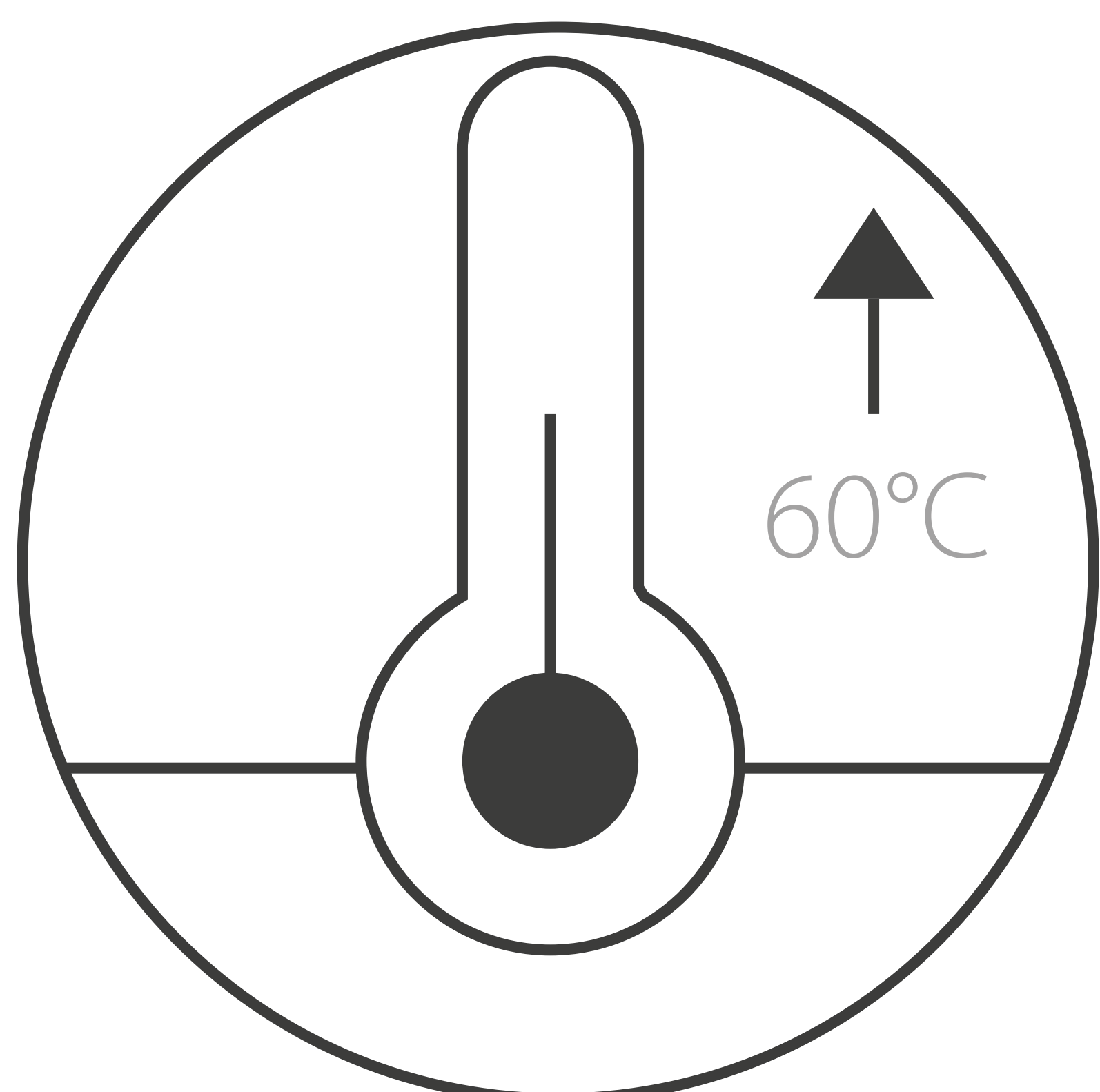


Electrically warming

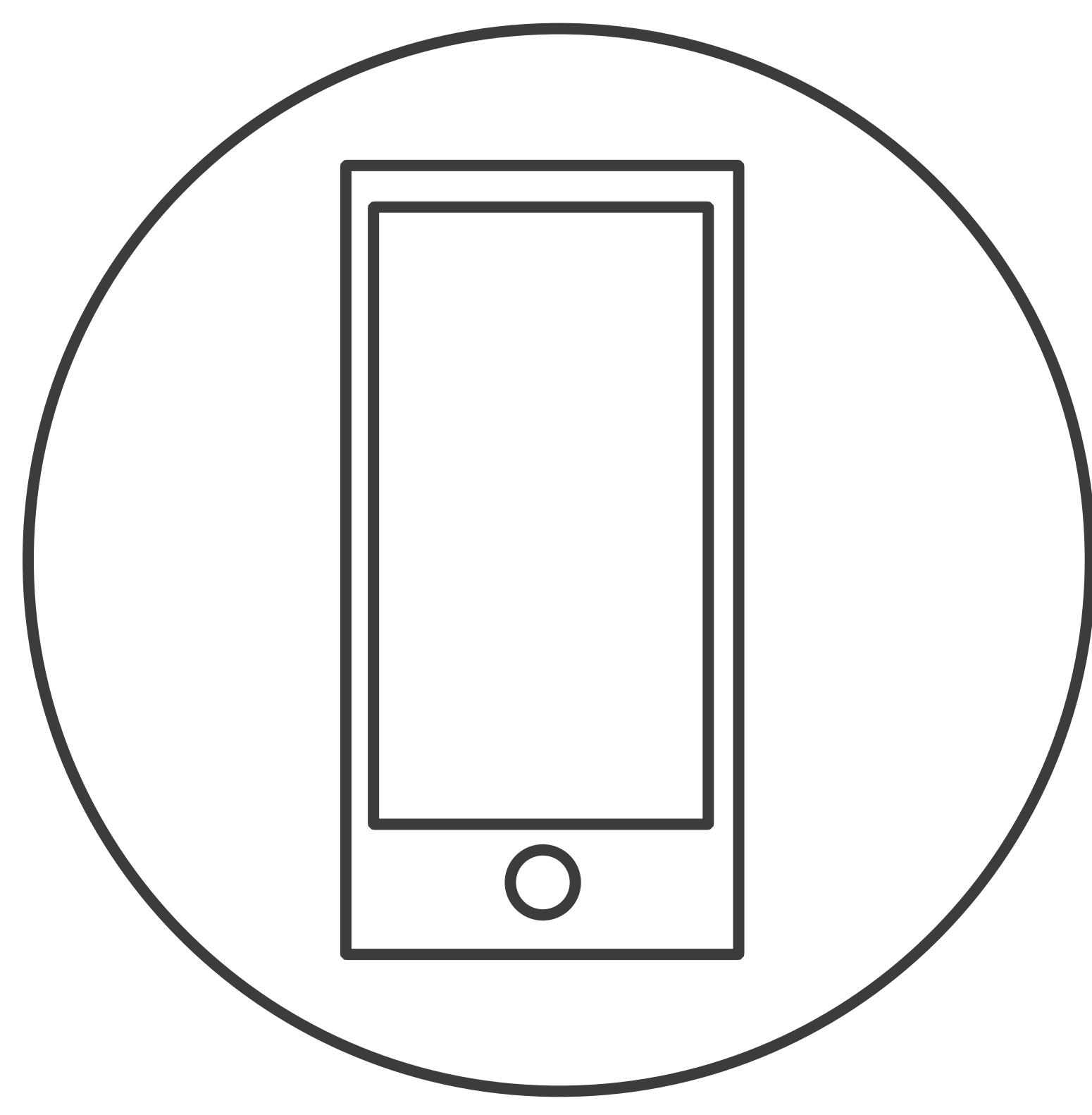


Smart Phone Controlled

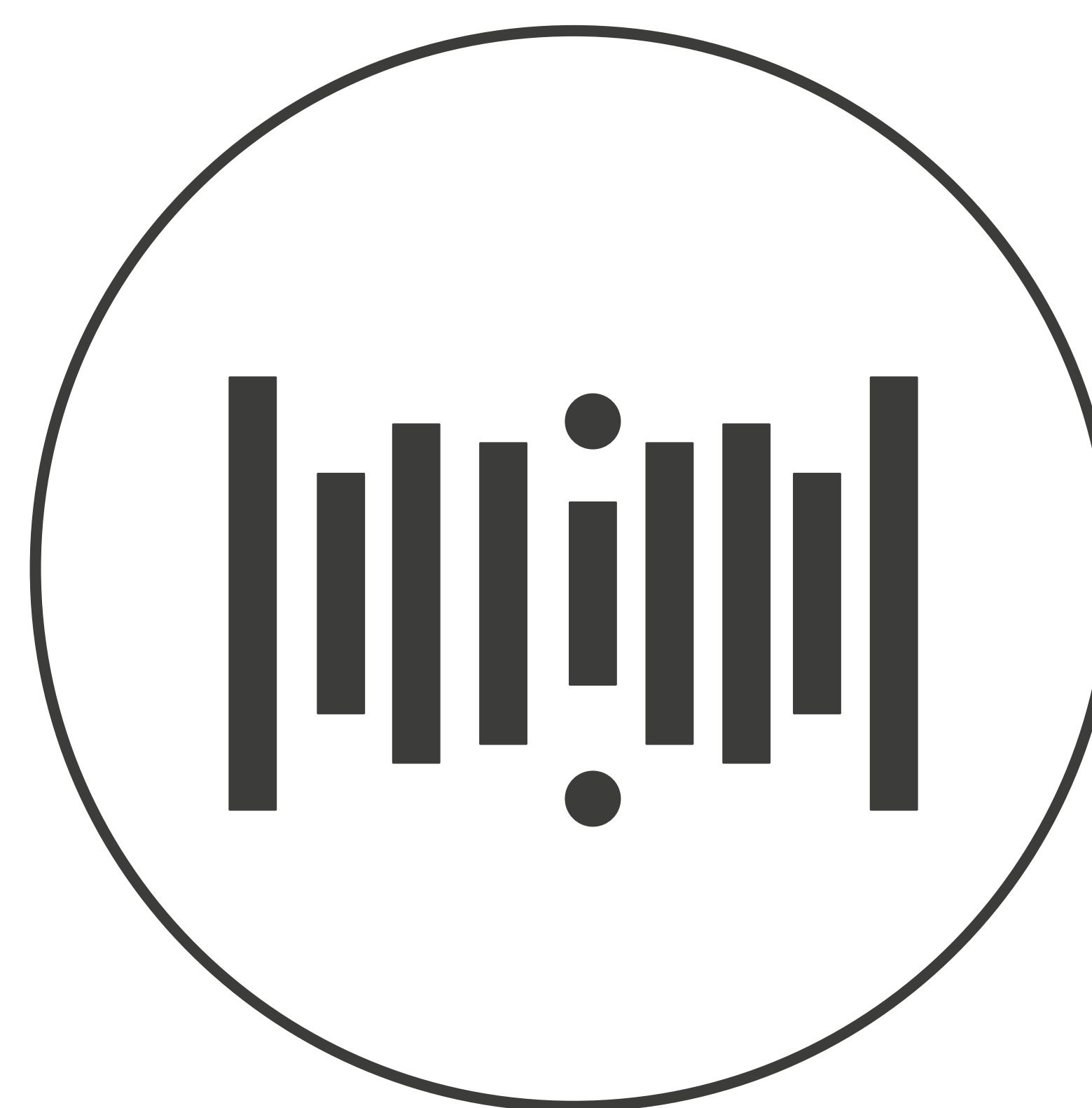


Voice Command iOS

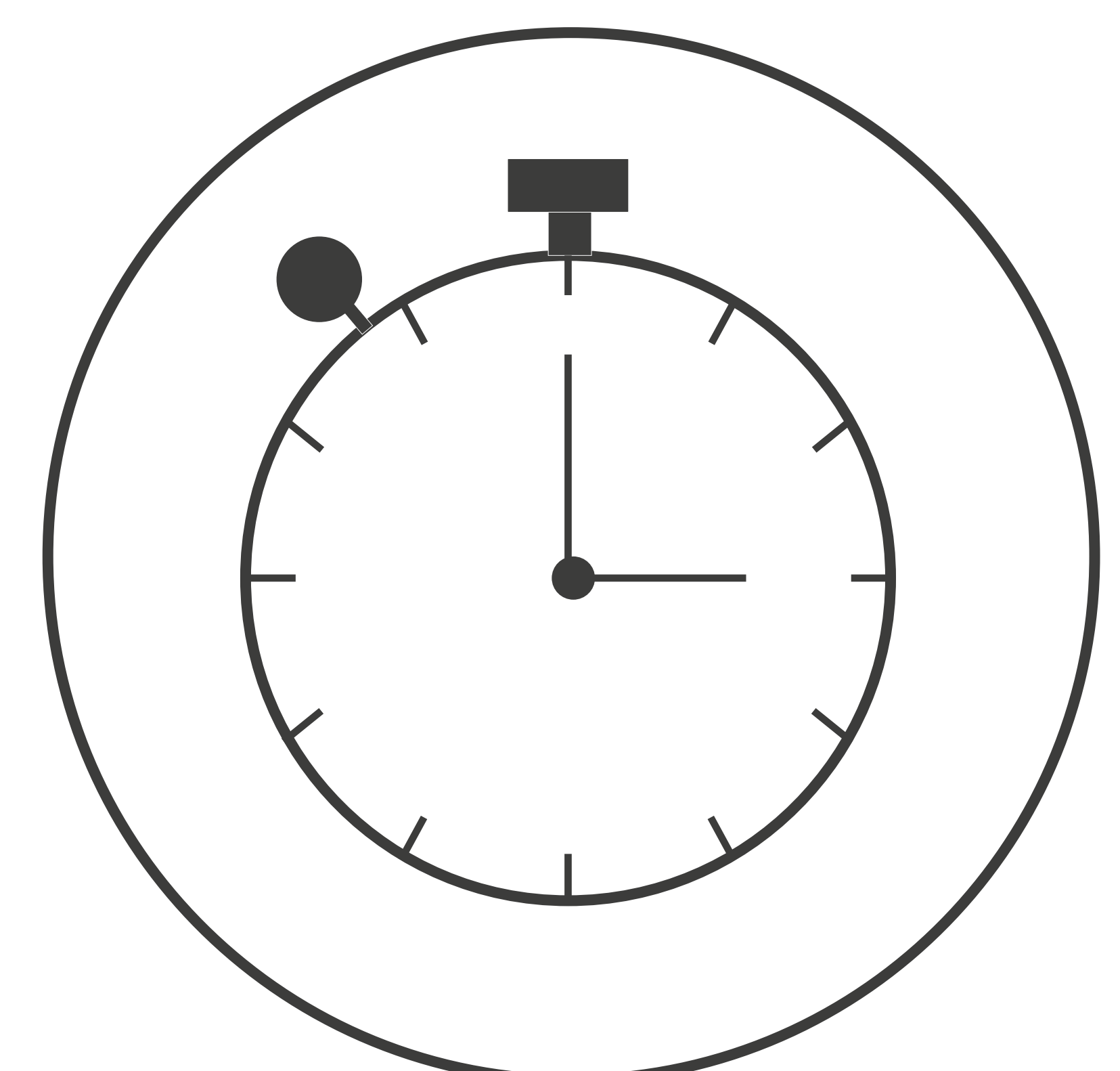


Auto-Timer

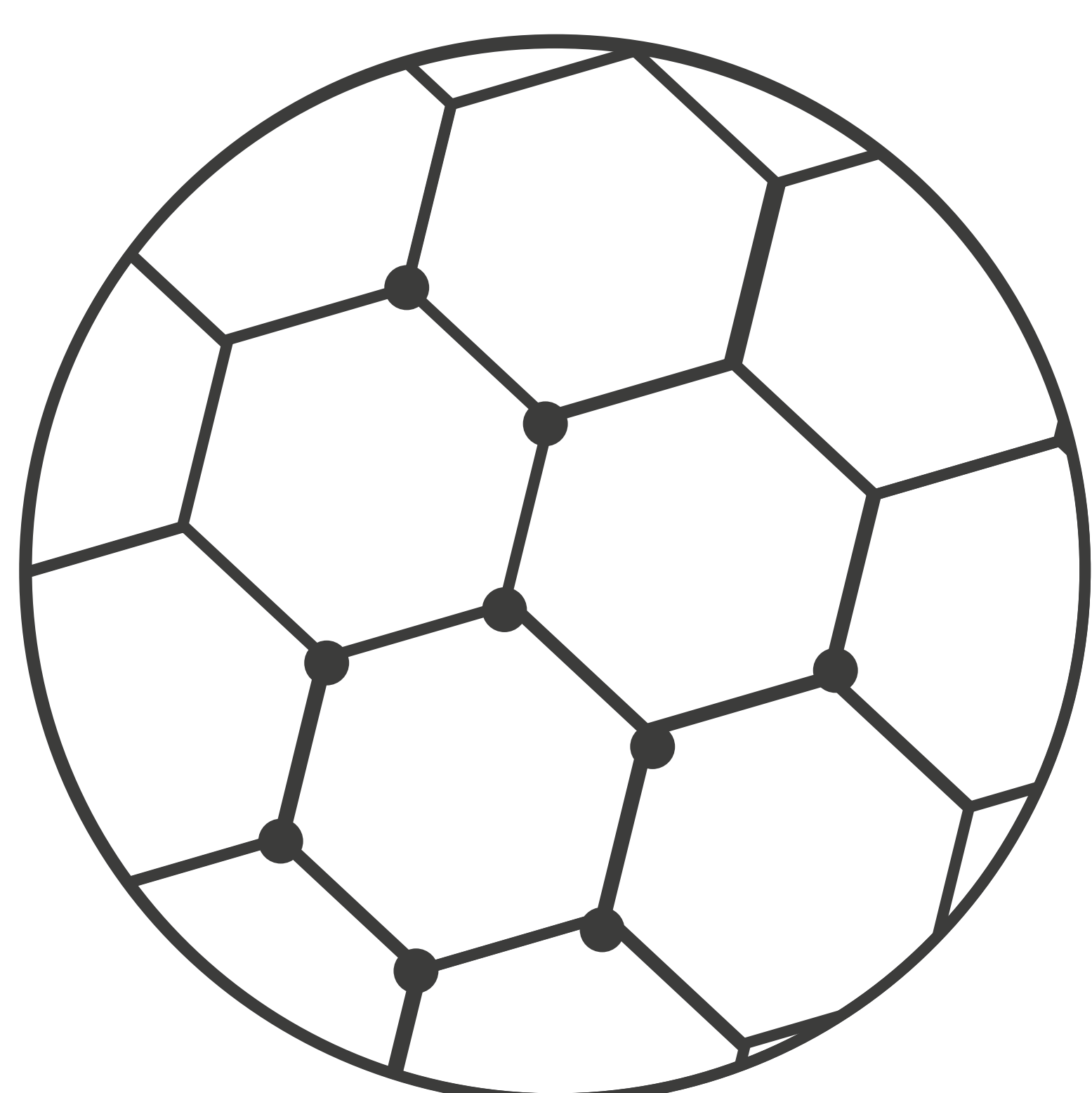


Graphene-Tech

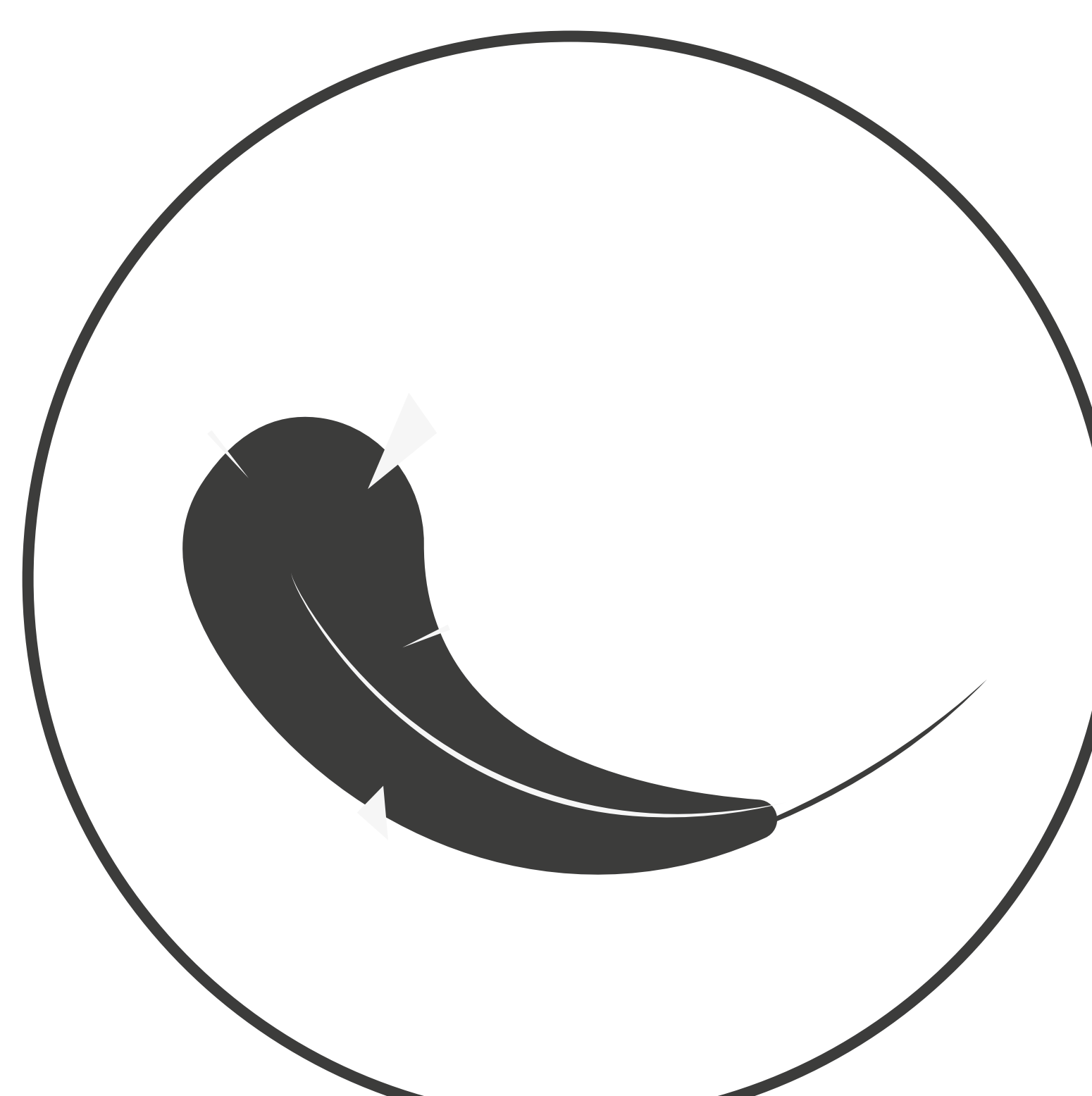

Ultra Soft

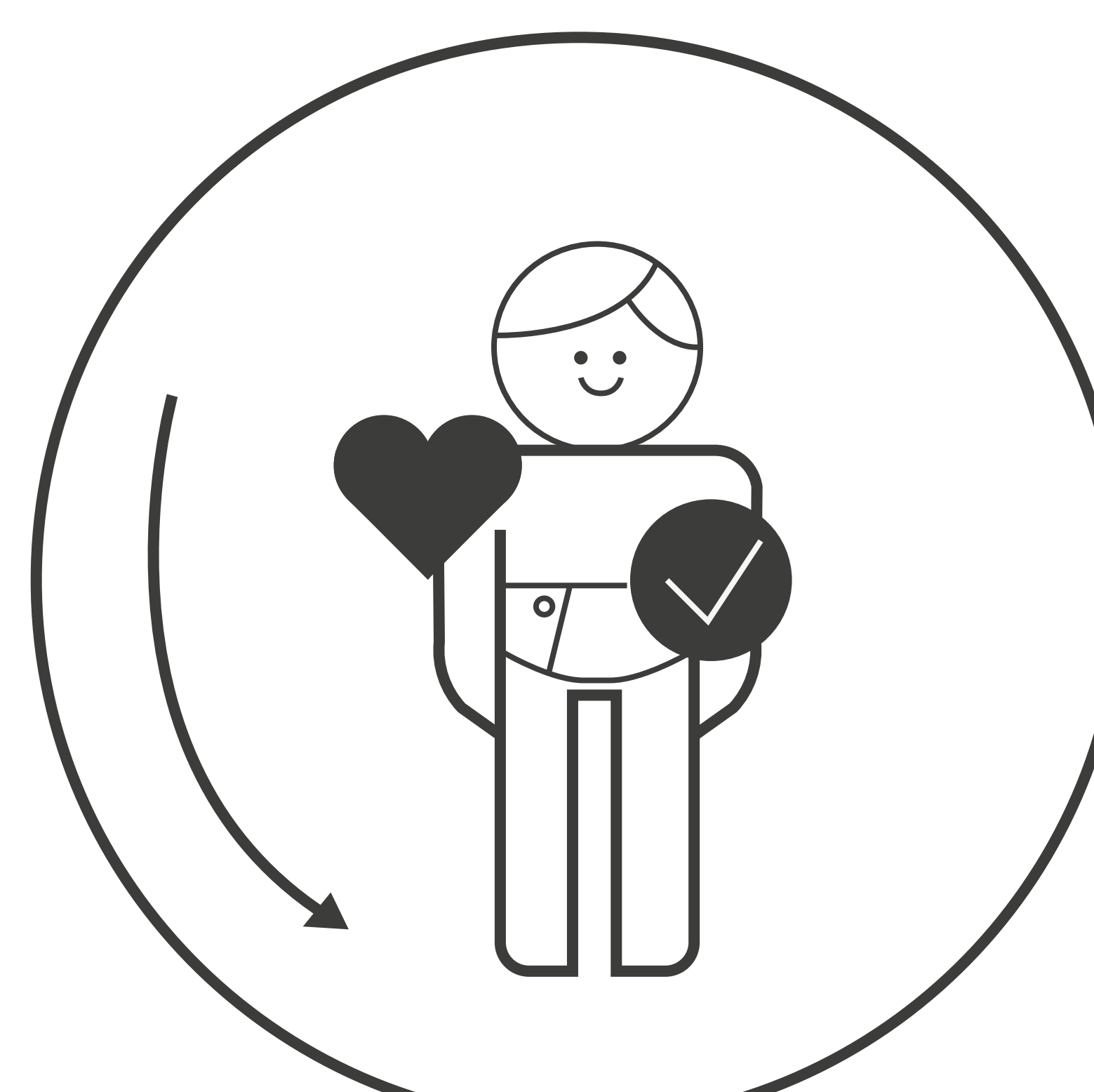


Ergonomic

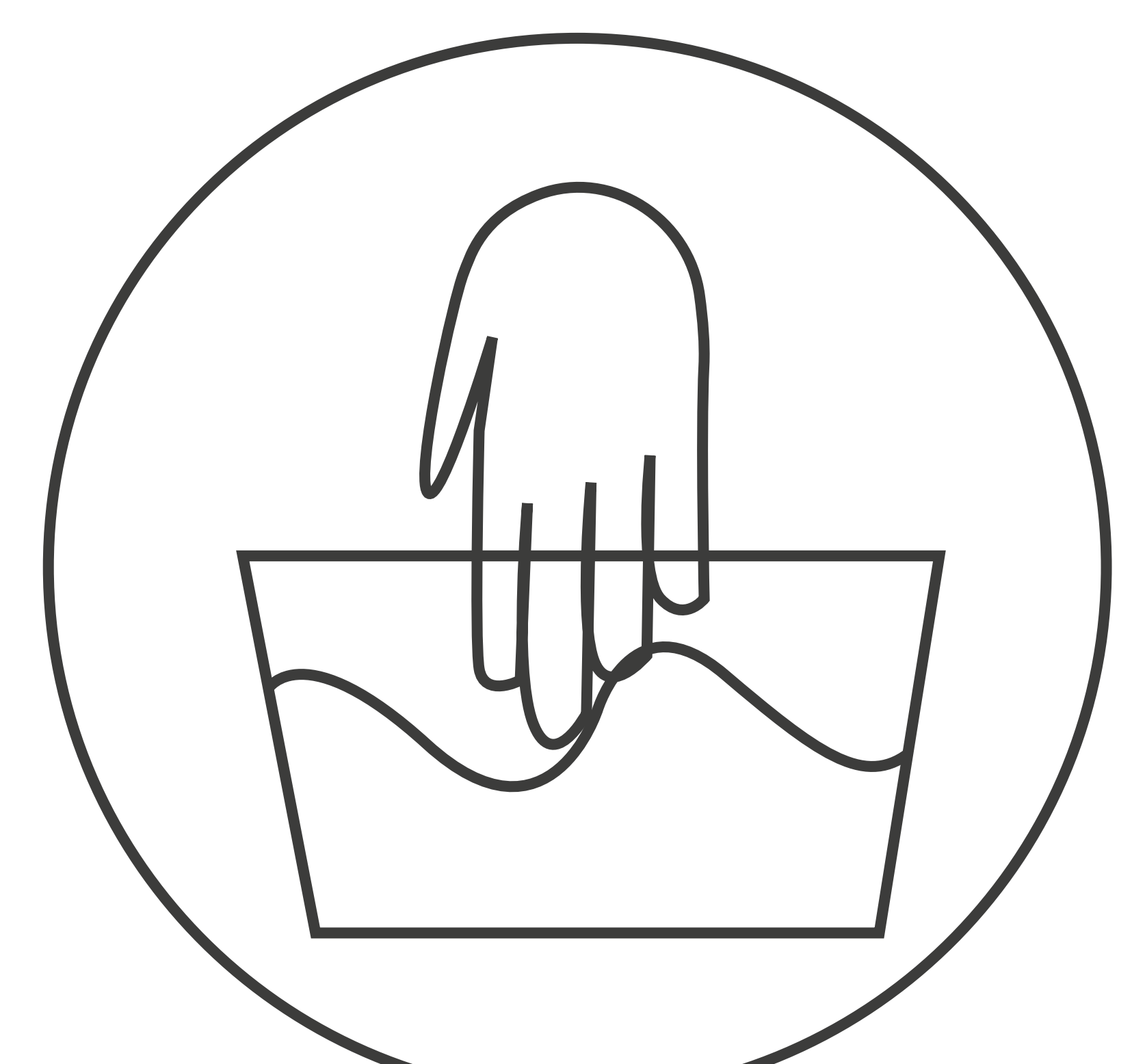


Hand wash




The Vulpés Smart Heated Belt can be used to support different conditions during the daily activities. The belt is suitable to use in daily activities, being at home, during work or on the go

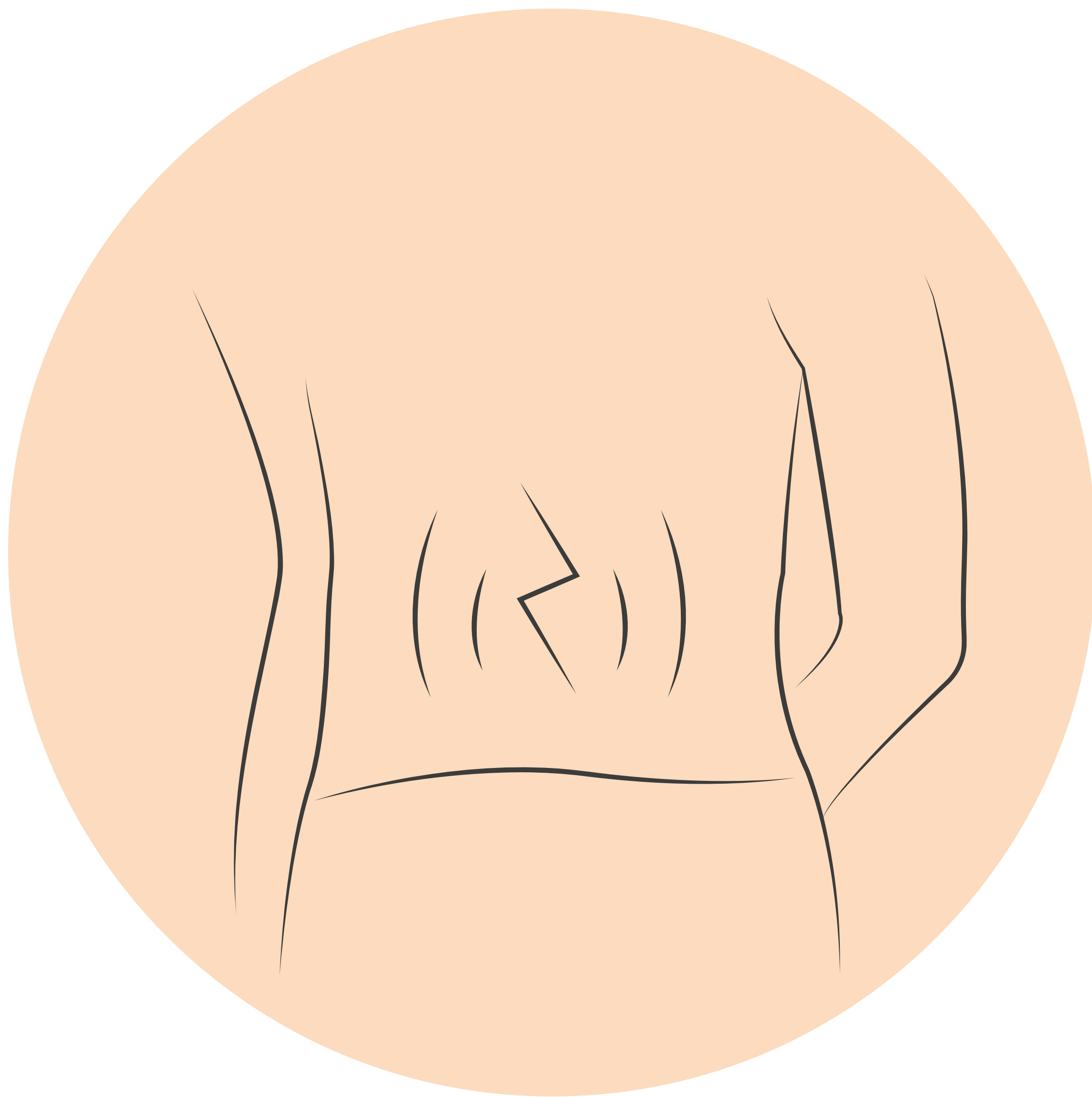


Endometriosis

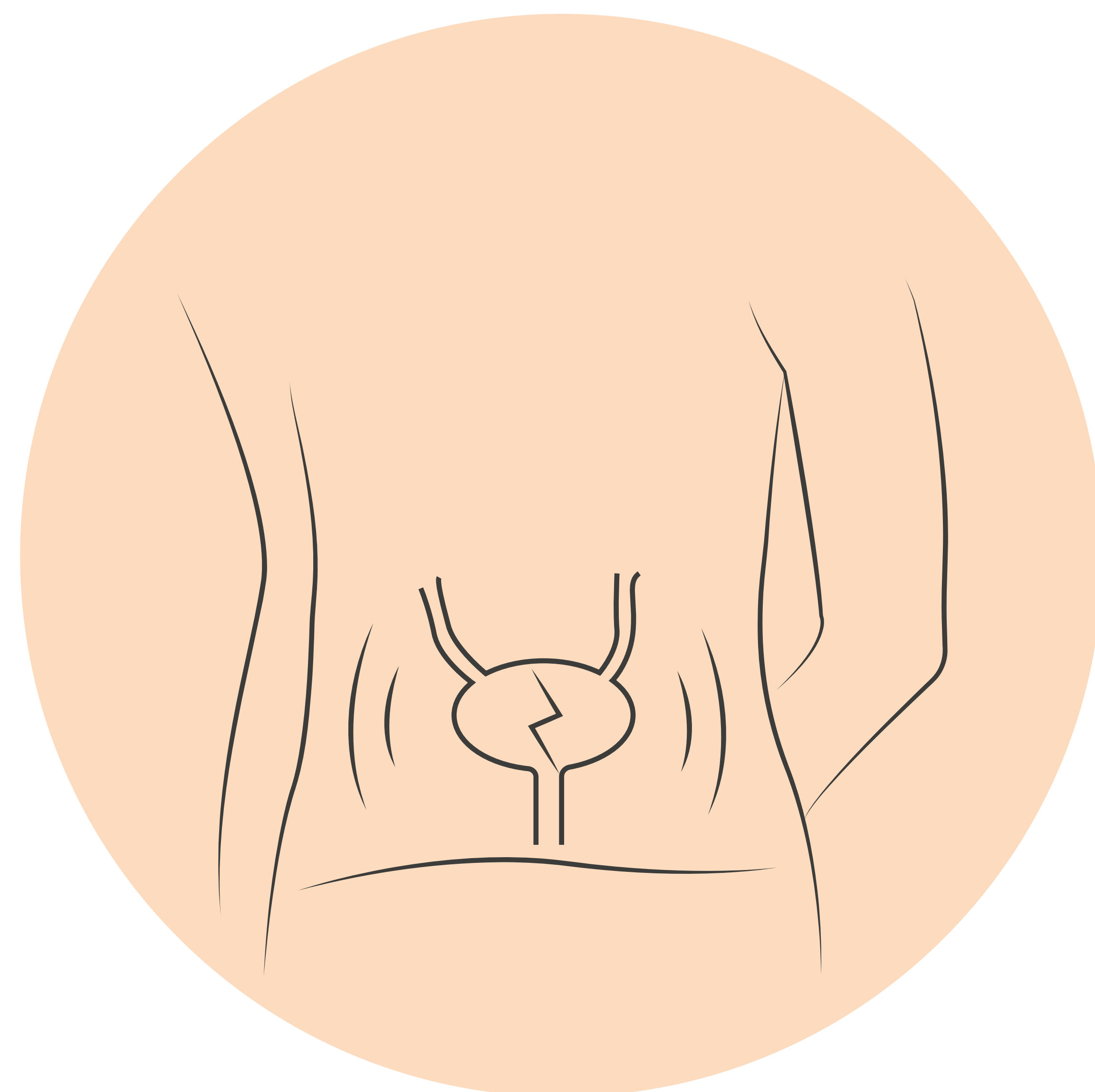


Cystitis

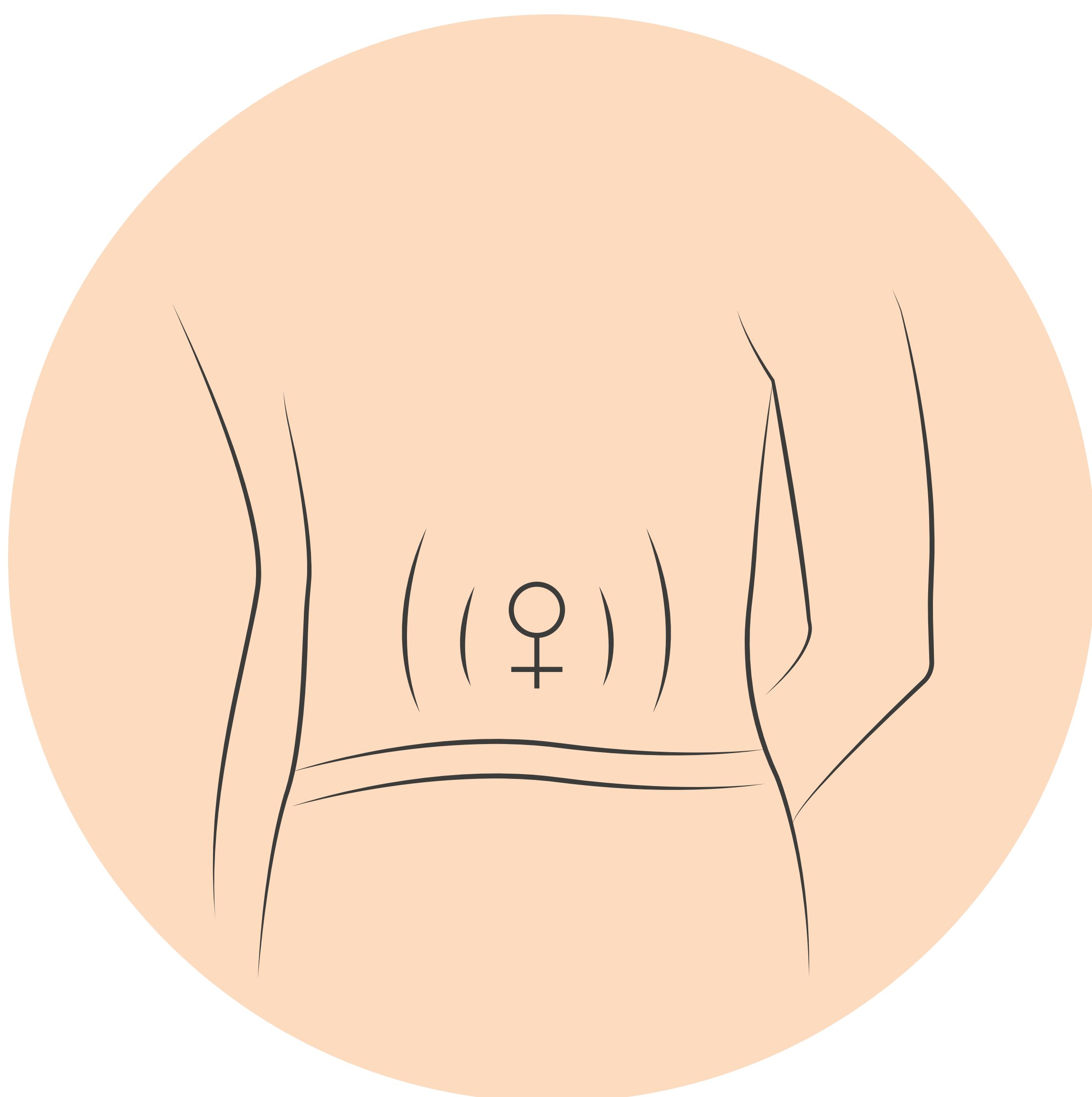


Menstruation

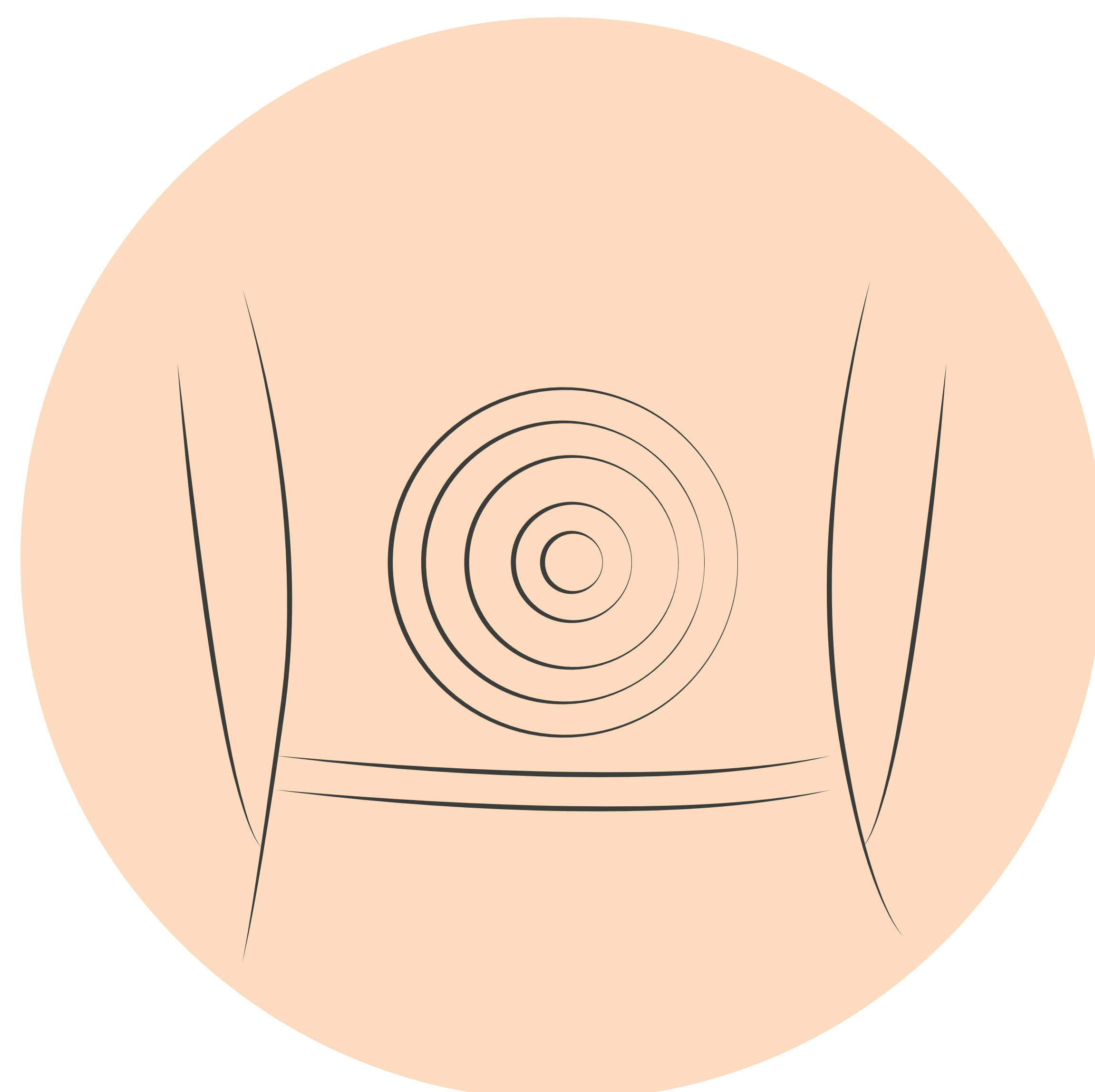

Muscle Spasm

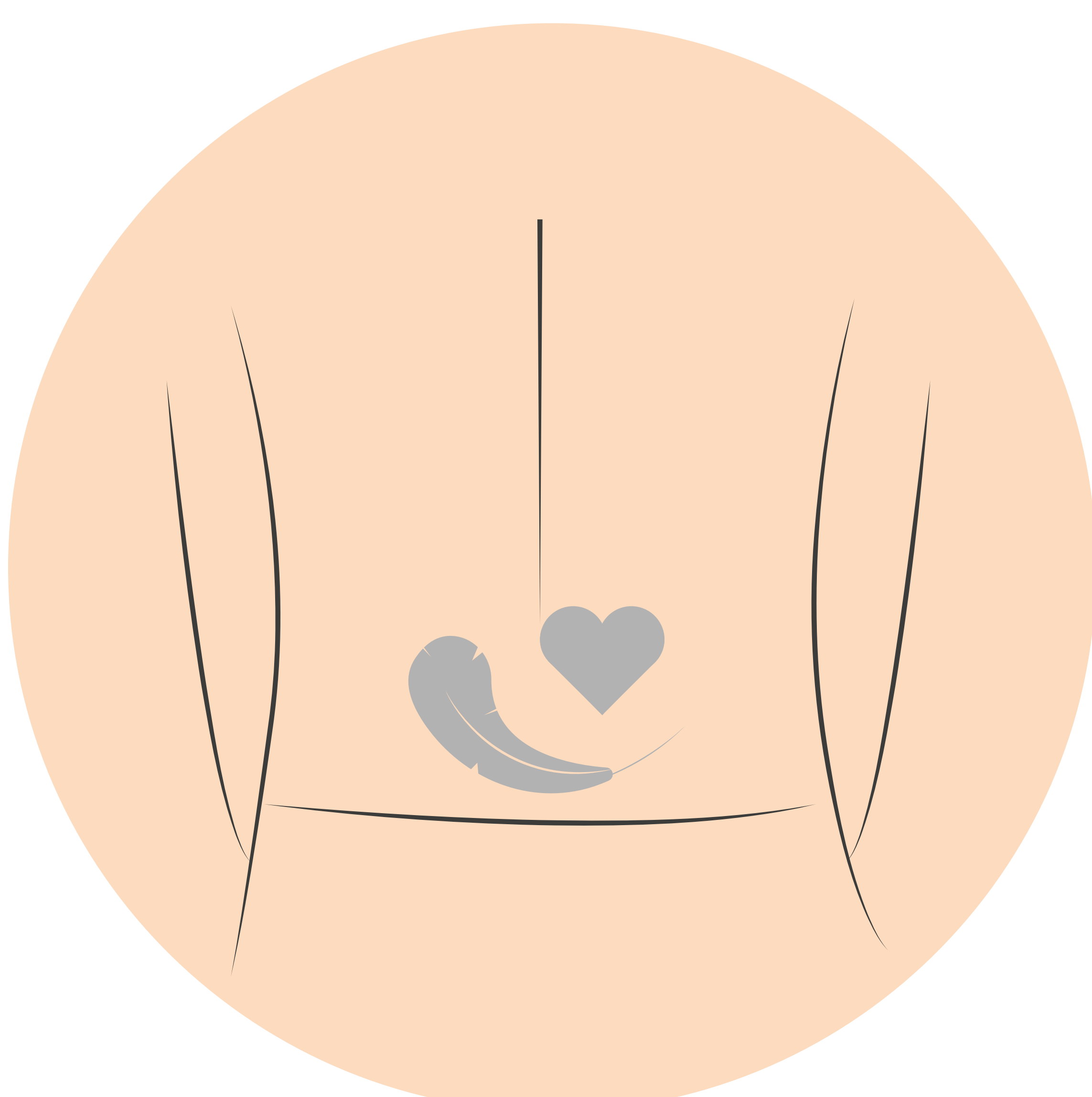

Relaxation



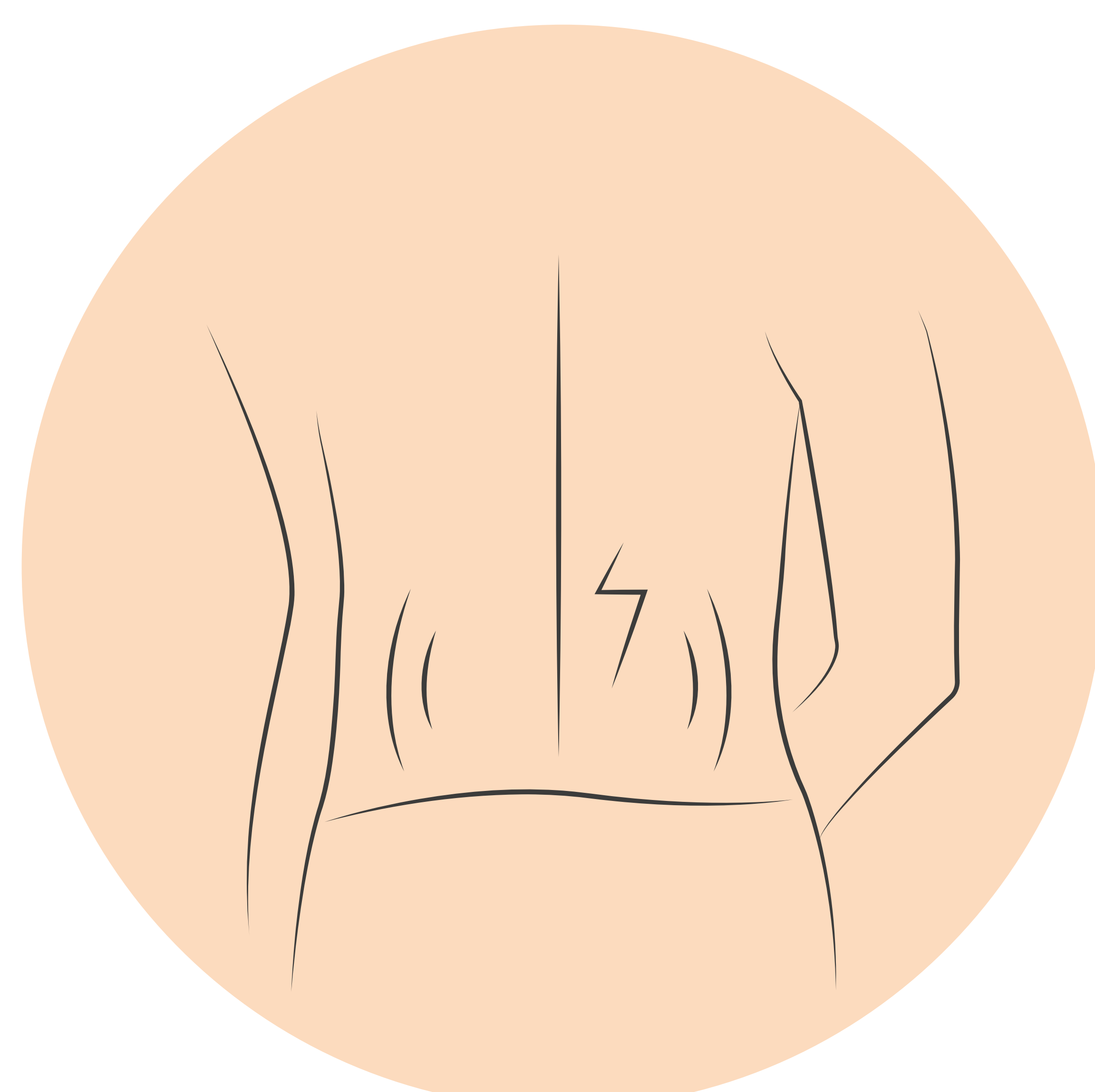


Lumbago

**The Vulpés Bellybelt is not a “medical product”. It is not promising any cure for illnesses or diseases. Please contact the healthcare personnel for medical advice or treatment options.**

# APPLICATIONS



# WHAT'S IN THE BOX ?

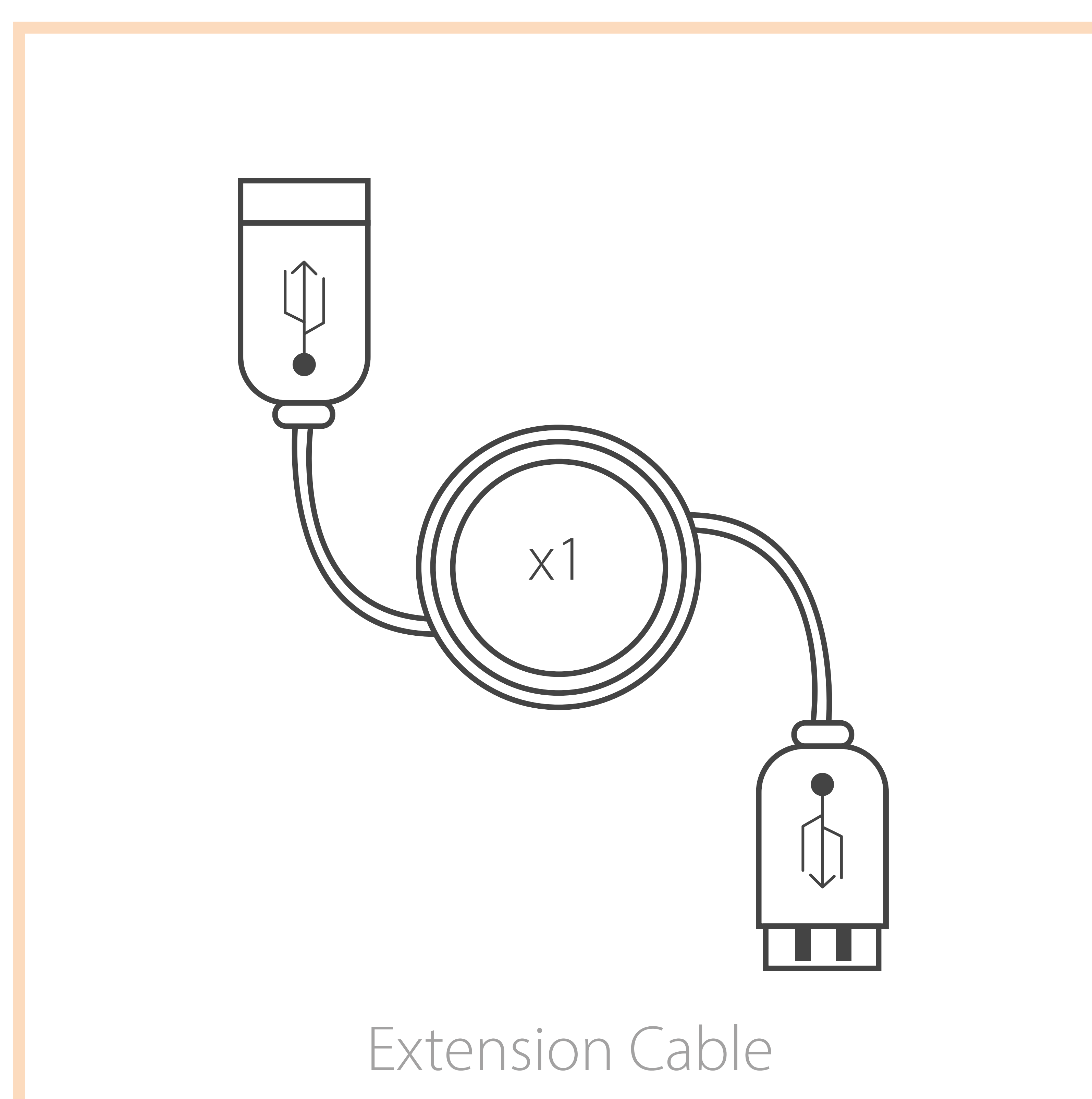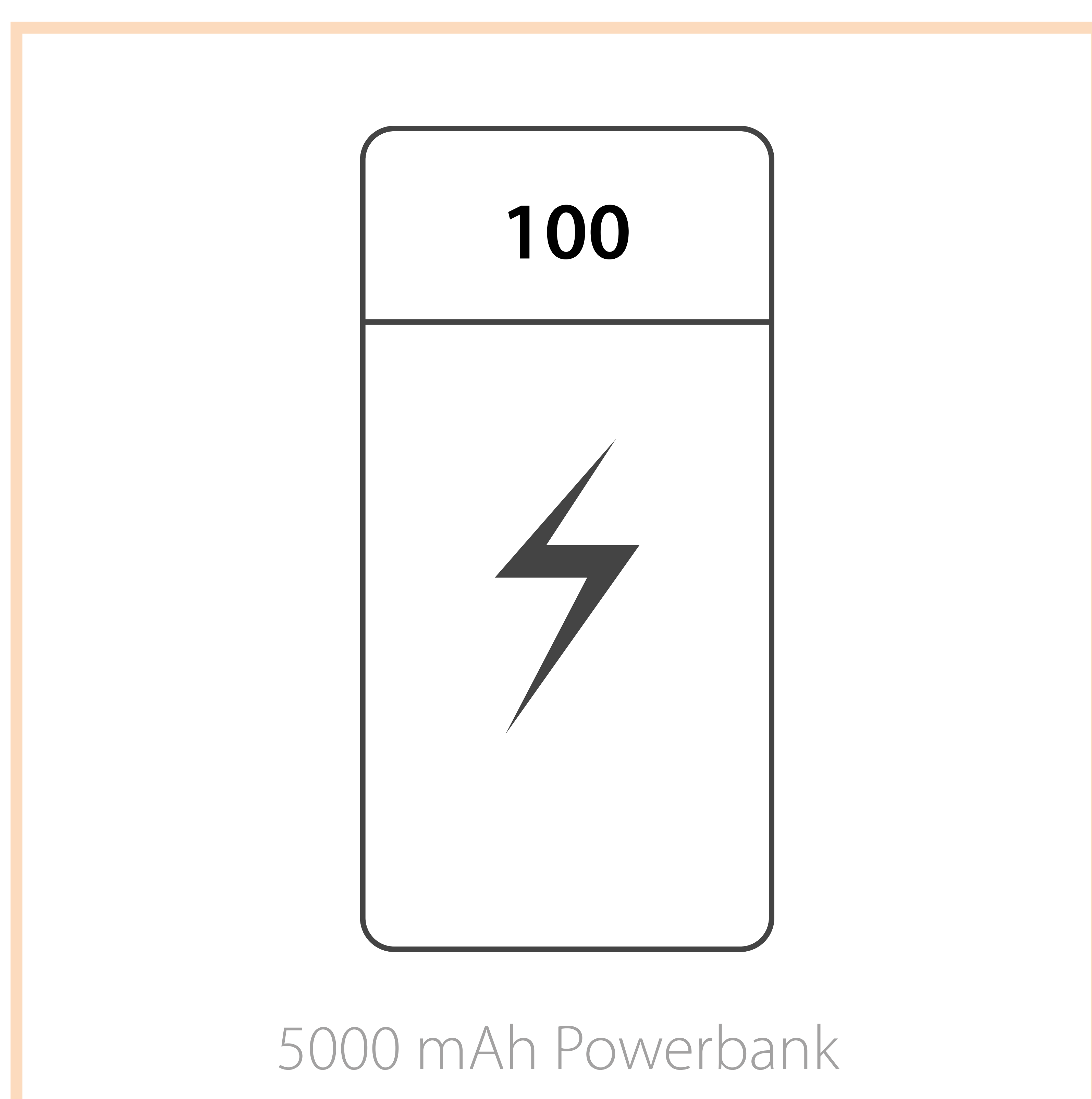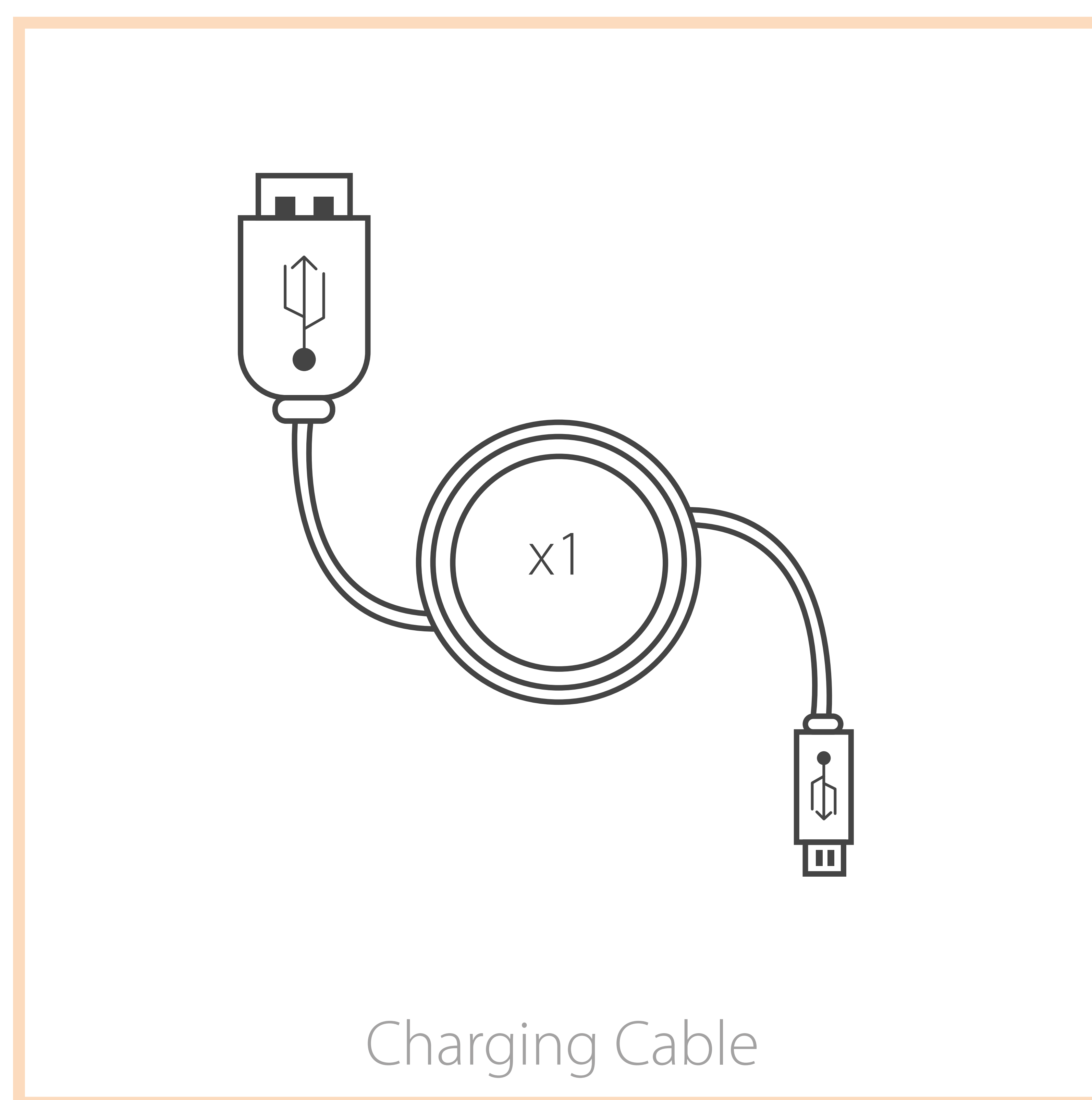

\* Vulpés BellyBelt Size S : 60 - 90 cm



\* Vulpés BellyBelt Size M: 85 - 120 cm

# PRODUCT

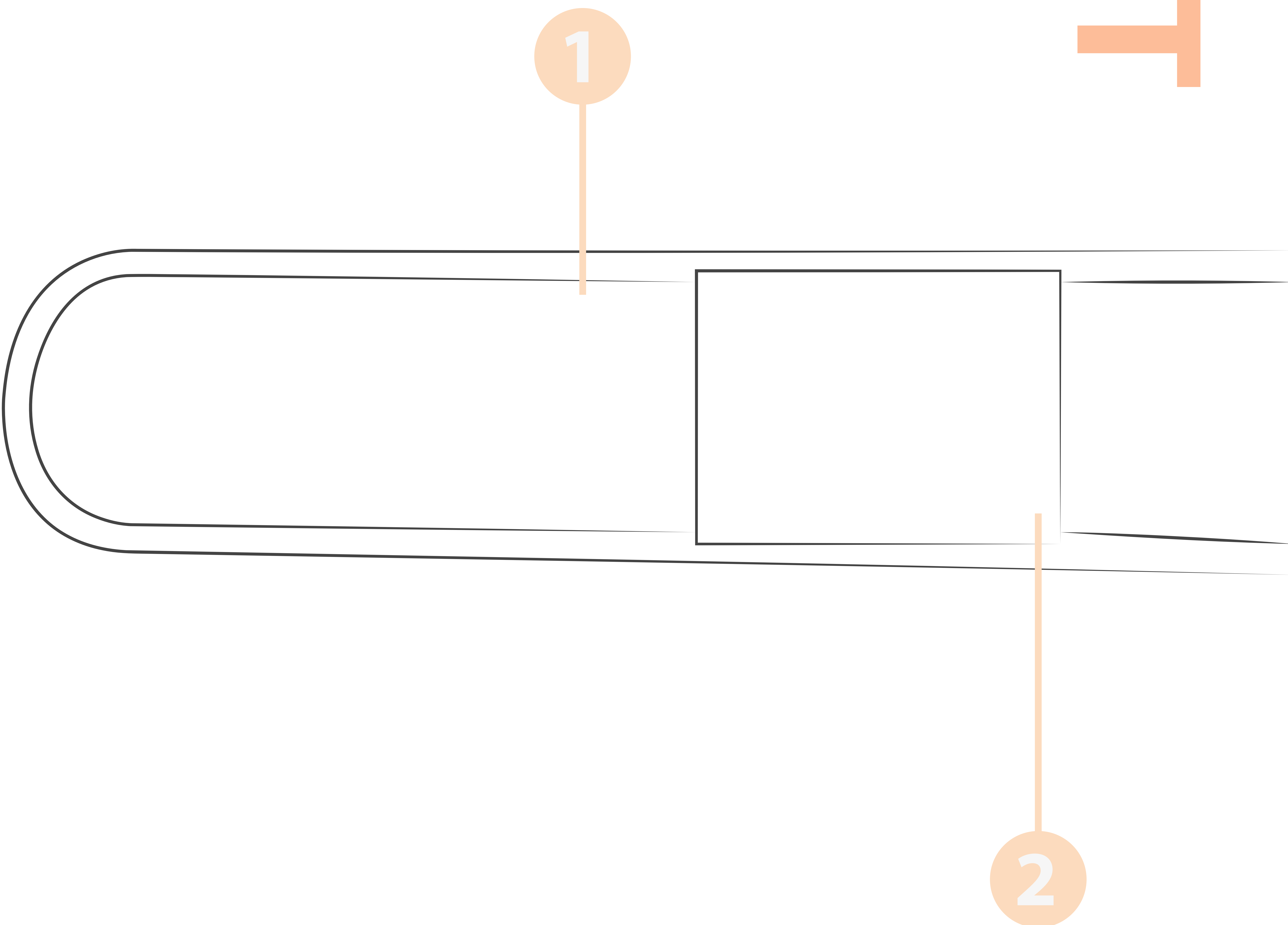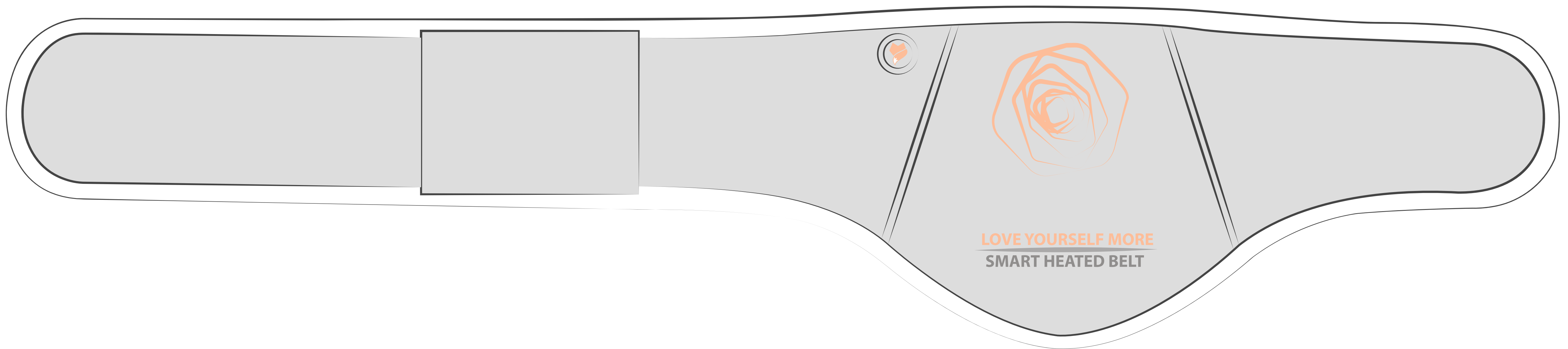

**1** Velcro fastener

**2** Battery pocket with integrated USB cable

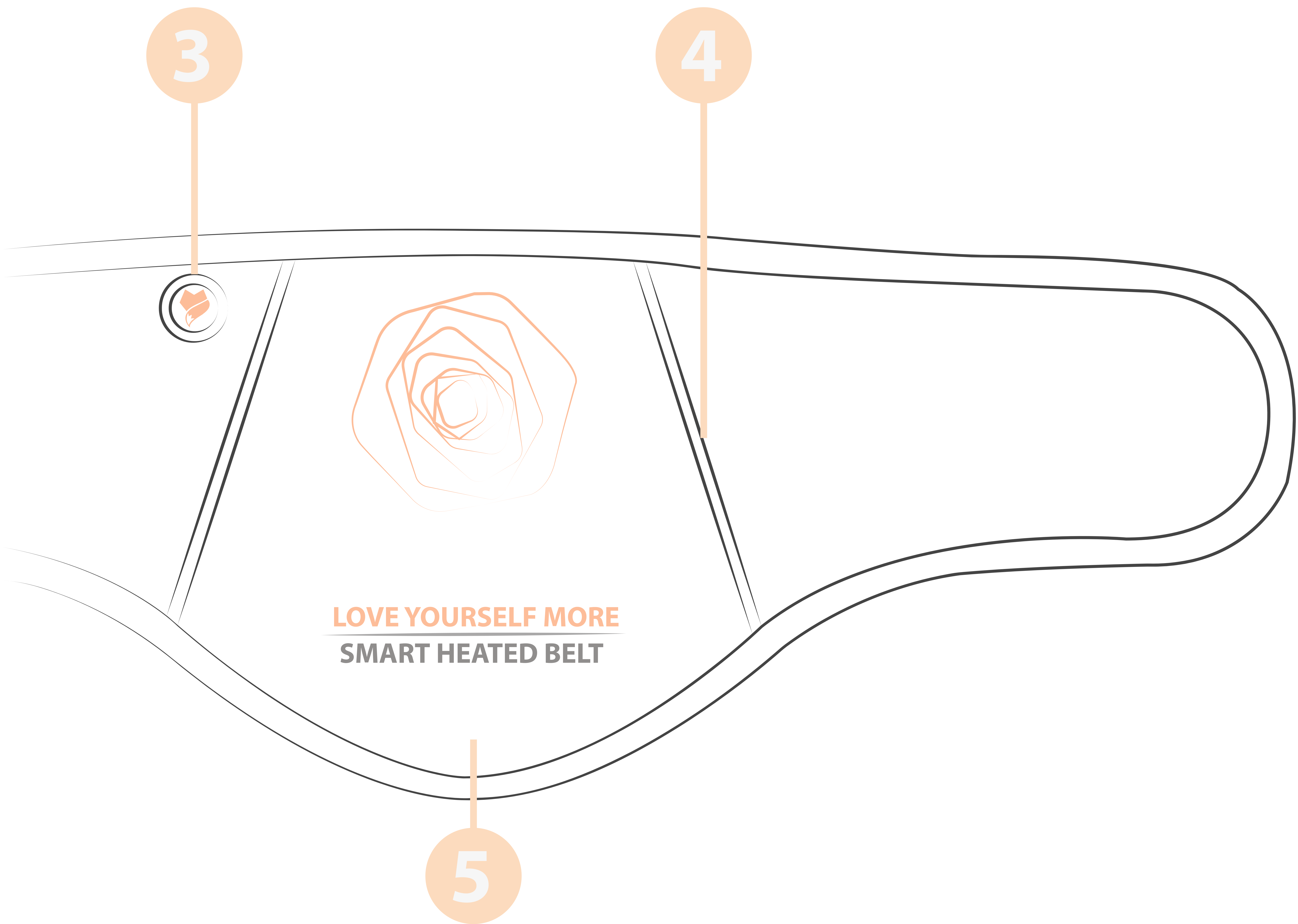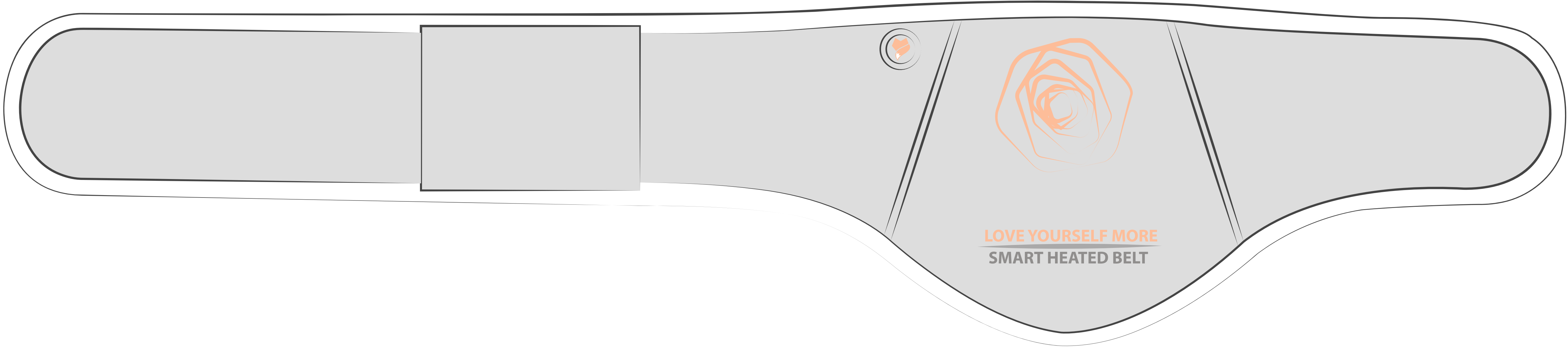


**3** Wireless, 3 mode controller

**4** Pocket for handwarming

**5** Heating area


# SET-UP

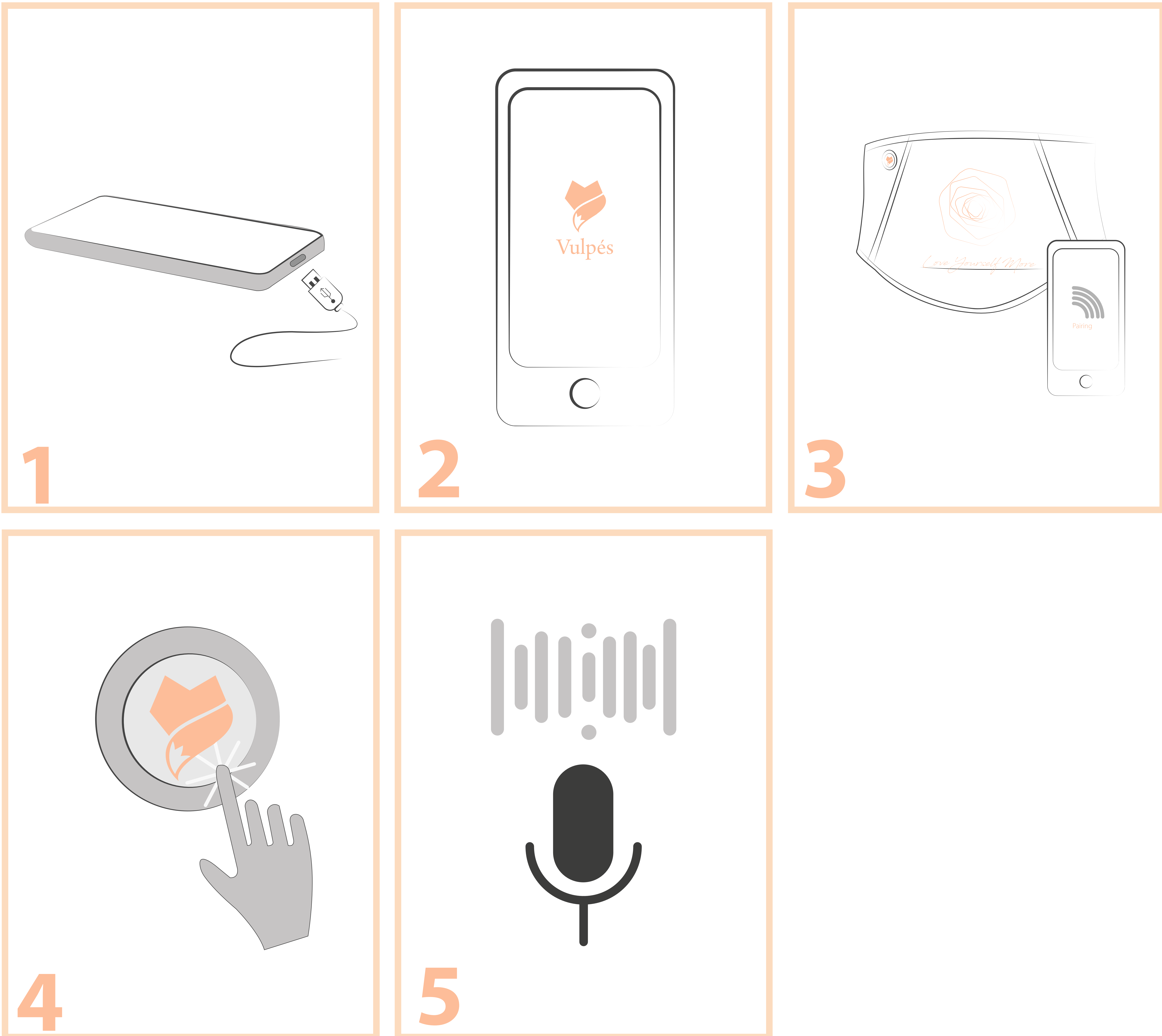

**1** Plug in the 5V / 2A power bank to the USB port which is located in the pocket and switch on the power bank.  
**Note: Please make sure before the activation of the abdominal belt that the powerbank is turned on. The powerbank turns off automatically after 30 seconds in case the belt has not been activated.**

**2** Download and install the Vulpés mobile application from the Apple App Store or from the Google Play Store.


**3** Pair your belt with the mobile application by following the steps inside the app.  
**Note: You have the option to control the heating settings via 3. Smartphone, 4. Button and 5. Voice**

**4** Press the controller to activate the heating settings. Push the controller button to adjust the desired heating state for **low (yellow)**, **medium (orange)** and **high (red)**. 

**5** Please follow the instruction in the App for Voice Command setup.  
**Note: The voice command feature is available for Apple iOS only.**



# USE

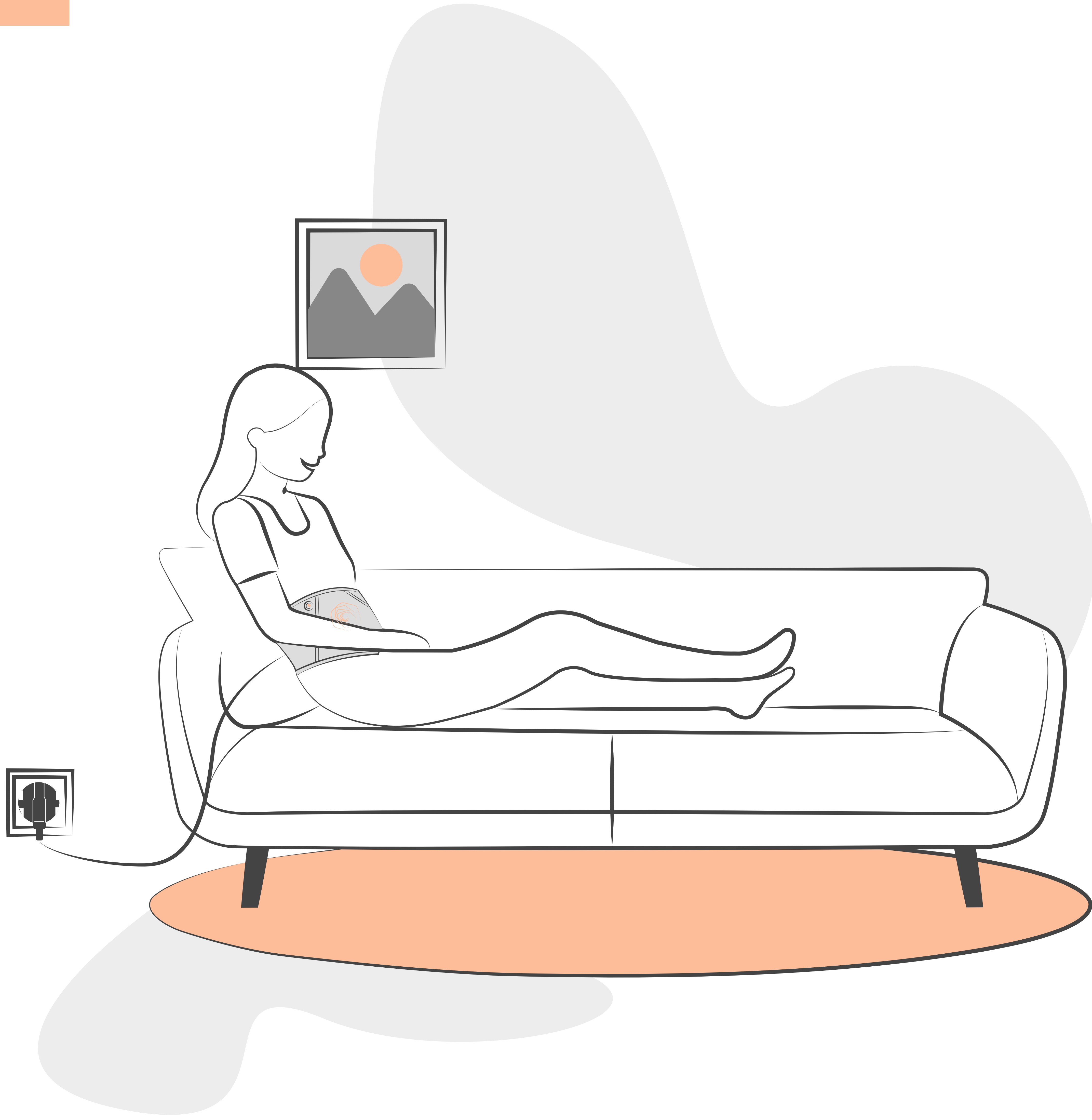

The Vulpés Belly Belt can be used stationary in the home or office environment by connecting the belt to the USB-A plug



The Vulpés Belly Belt can be used at the mobile state with the USB Powerbank



# PRODUCT





# USE


# POWER BANK

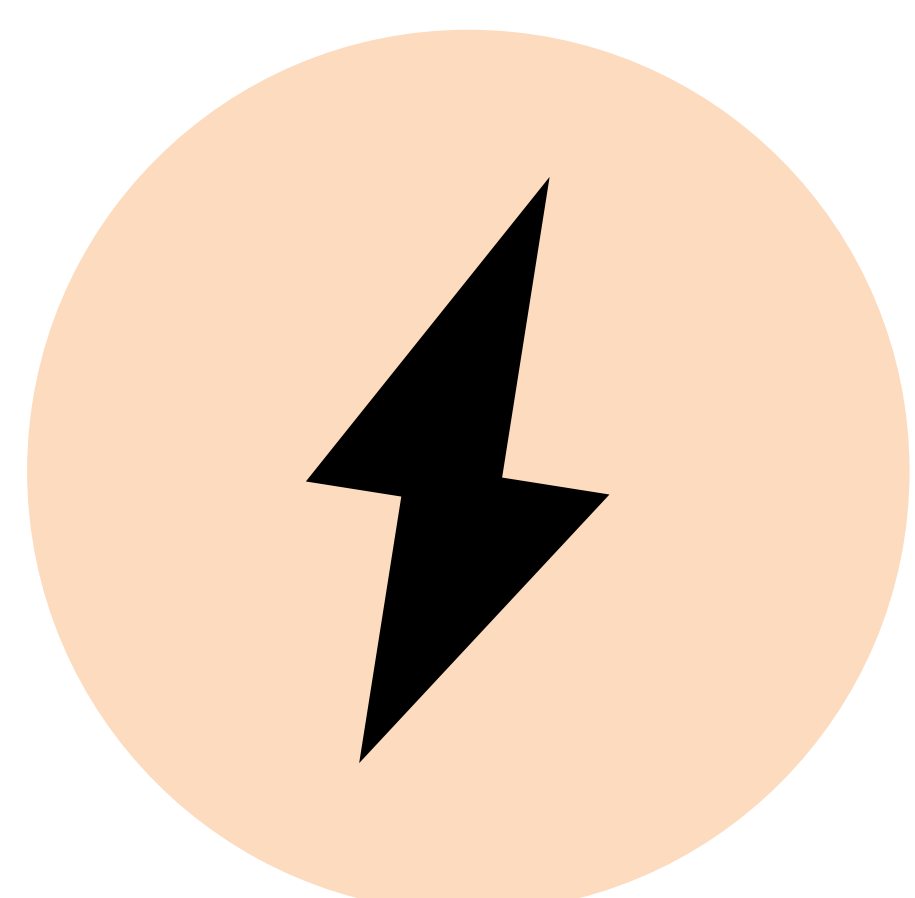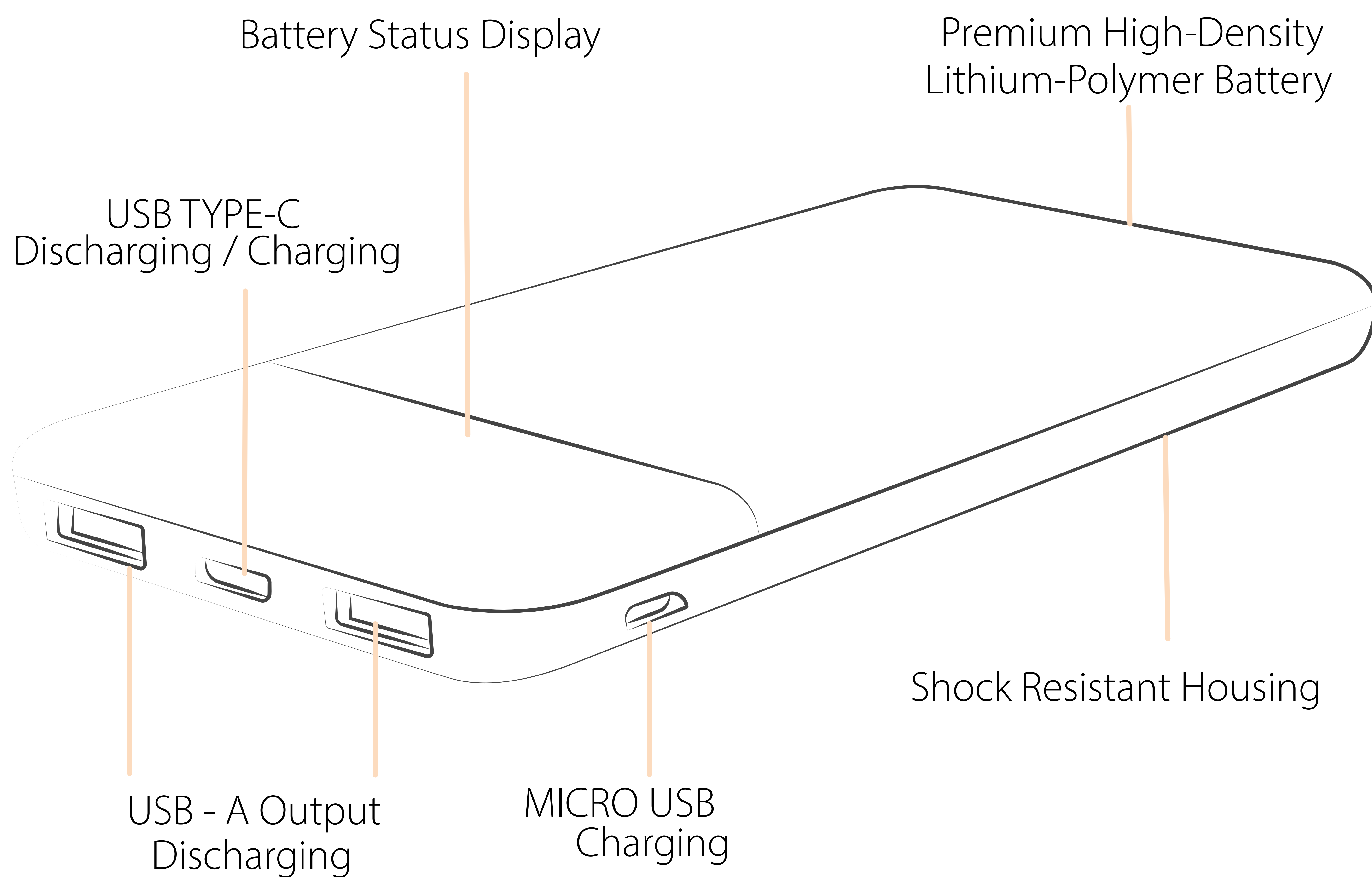


Battery Status Display



Light-weight

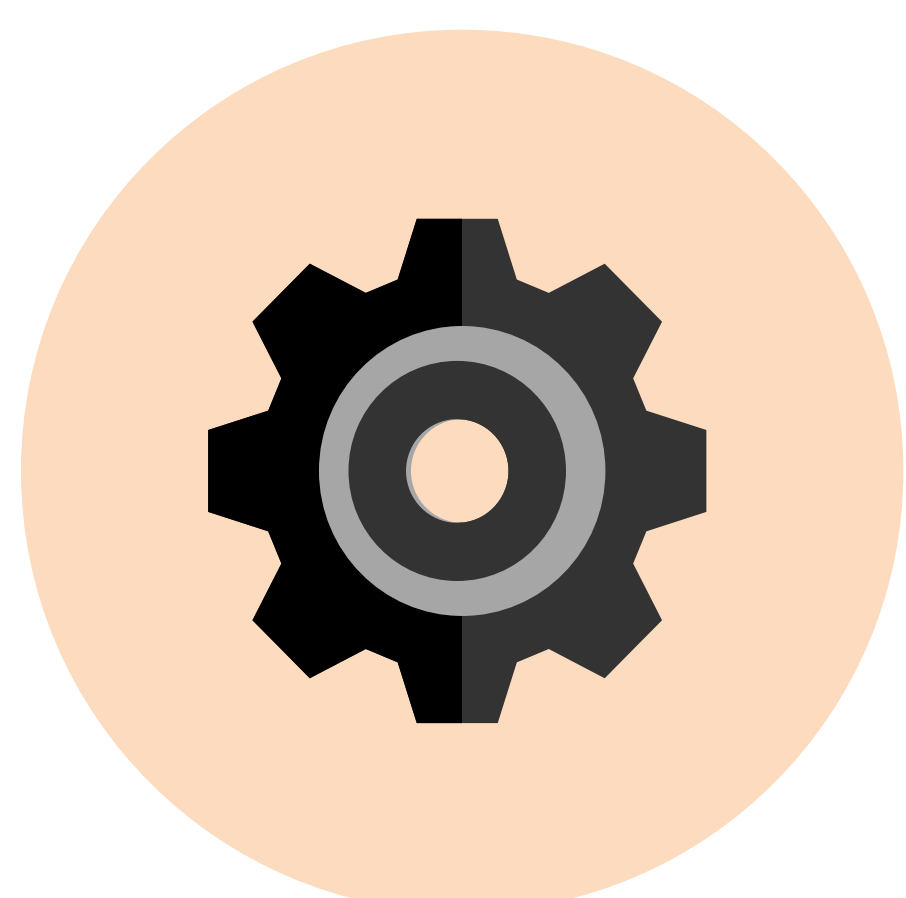


High-efficient

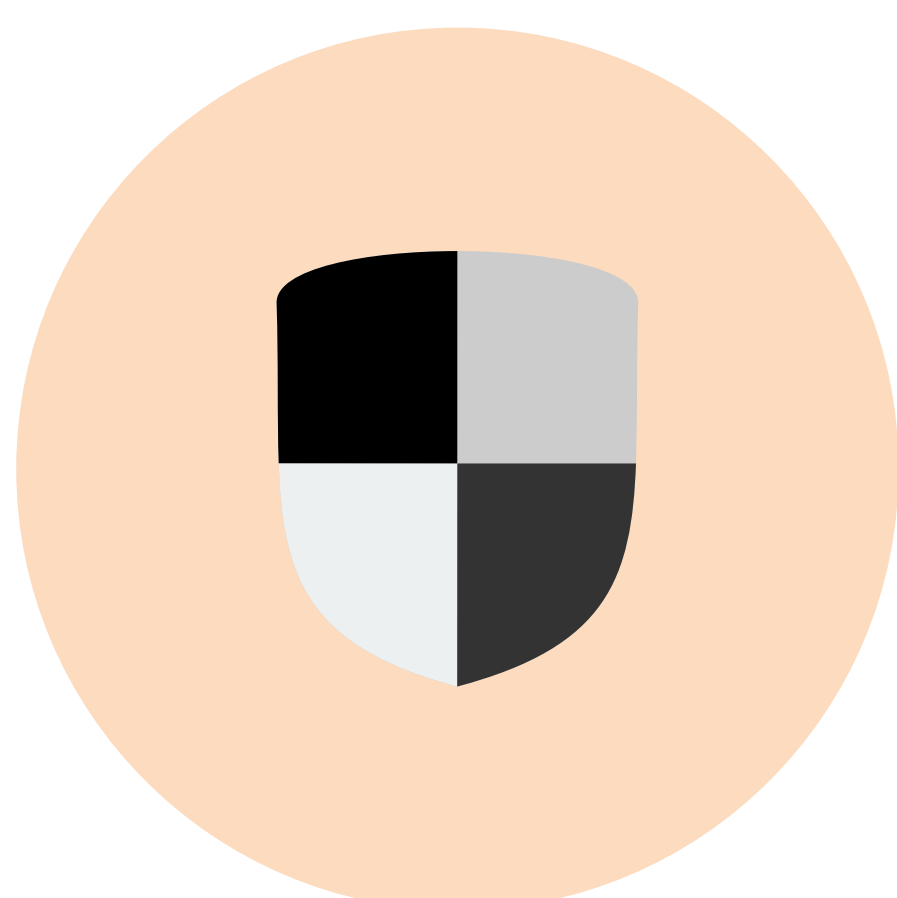

High safety




Lithium Polymer Battery

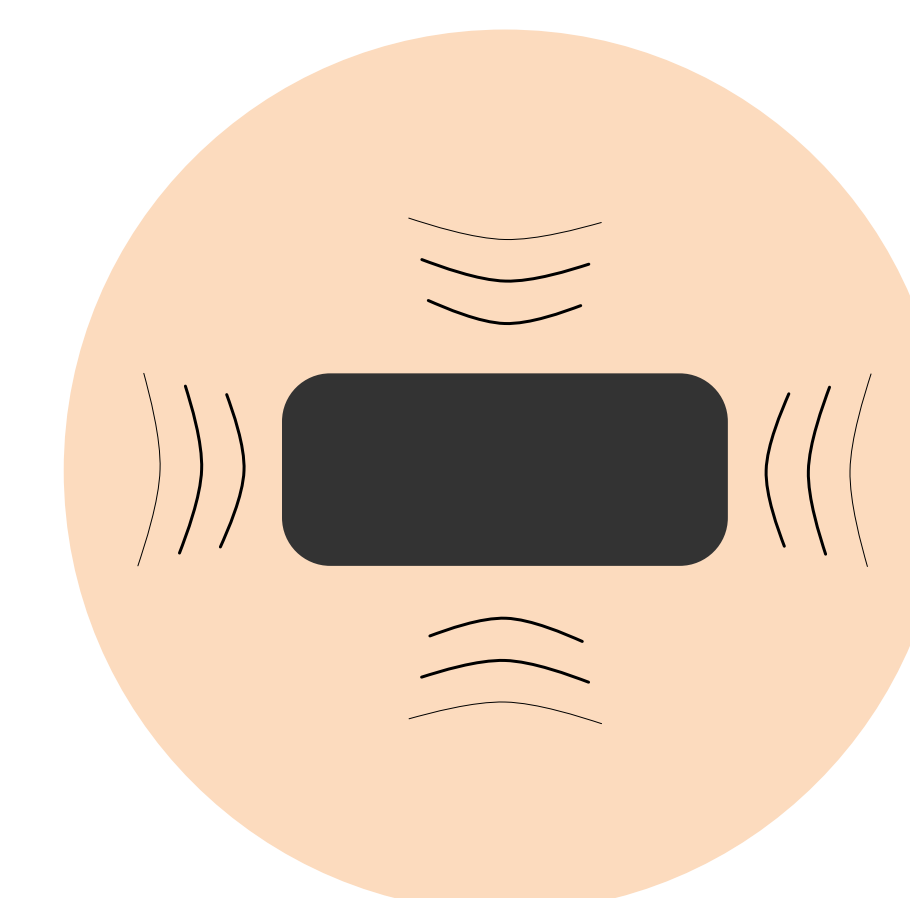


Shock resistant


## SAFETY FEATURES

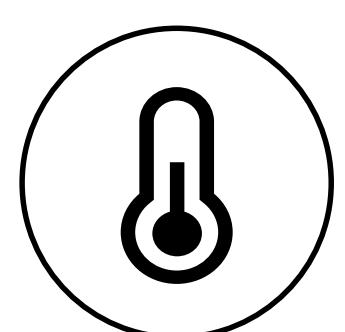


Temperature protection

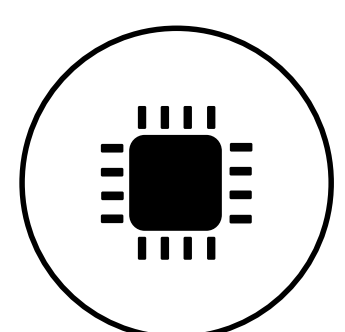


Overload and short-circuit protection

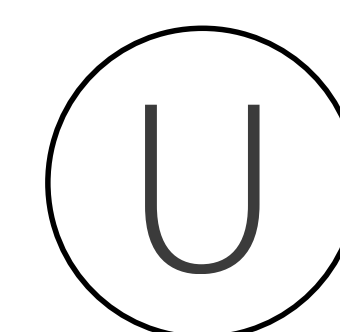


Input over-voltage protection

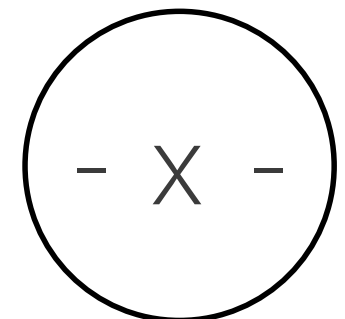


Short-circuit protection

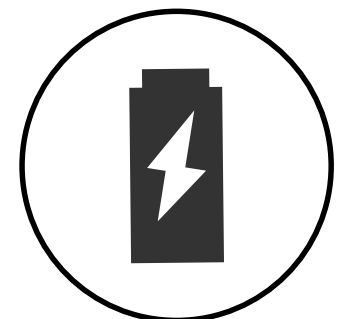


Overcharge and over-discharge protection

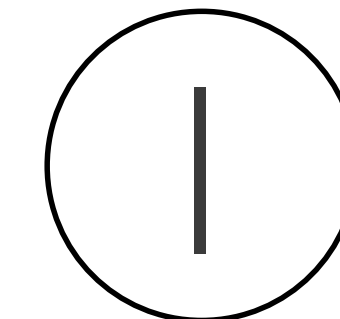


Output over-current protection

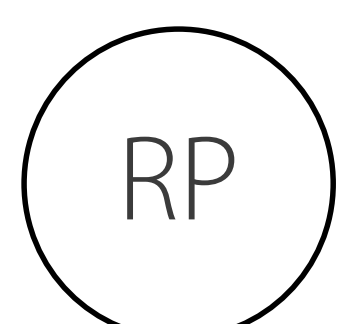

Reset protection

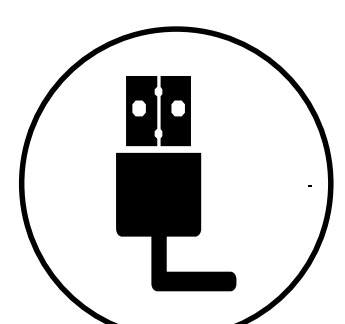

Input anti-reverse protection

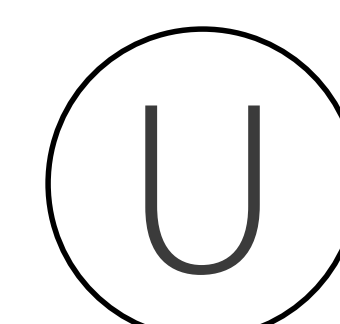

Output over-voltage protection

## SPECIFICATIONS

Battery type	Lithium Polymer	Input	Micro USB	DC 5.0V / 2.0	A
			USB Type C	DC 5.0V / 2.4	A
Capacity	5000 mAh	Output	USB Type A	DC 5.0V / 3	A
			USB Type C	DC 5.0V / 2.4	A
Net weight	110 g	Dimensions	100 x 63.3 x 13.5 mm		




**The product may not be used by people with reduced sensitivity to heat by individuals who might be vulnerable to injuries.**

(e.g. diabetics, persons with skin changes related to illness or persons with scarred skin in the area where the pad is to be applied, after taking any pain-relieving medication or alcohol)

**This product is not a “medical product”. It is not promising any cure for illnesses, diseases. Please contact the healthcare personnel for medical advice or treatment options.**

**People who experience any discomfort while wearing the Vulpés Belt, must remove the belt immediately and seek medical advice before resuming the use of the belt.**

**People who have skin allergies to synthetic, rubber or plastic materials that may be found in the belt, should monitor carefully their skin each time the belt is worn. In the event any irritation of the skin occurs, immediately stop wearing the belt, and seek medical advice.**

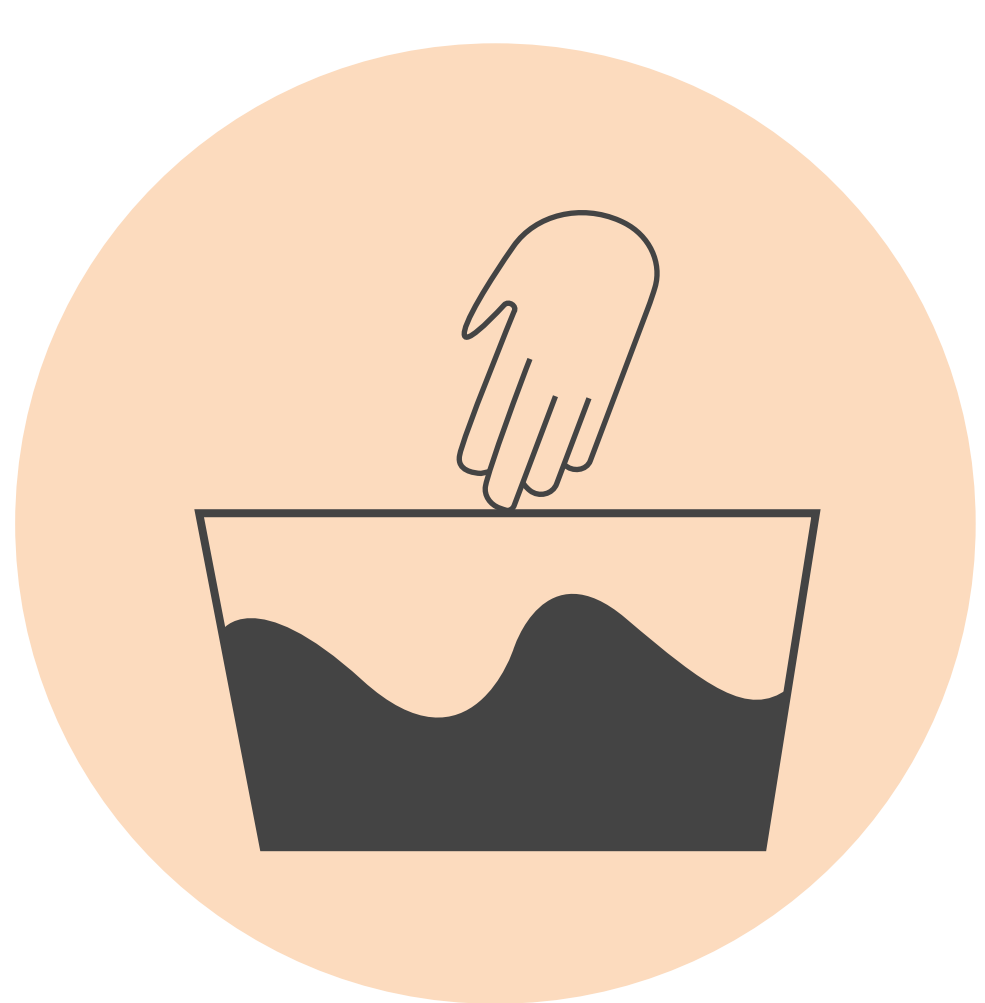

**Cleaning:** Please switch off the device before cleaning and remove the powerbank from the pocket. Please clean the device only by hand. Before you again attempt to use the thermal system, it must be completely dried out. During the drying process please prevent any contact with heaters.

Please prevent general contact with water during the operation.



**Children: Warning!** Keep out of the reach of children. Do not allow children to play with devices or accessories.



**Power Bank:** Please use always the **5V/2A power bank** output port. By using output ports with higher voltage you damage the control system permanently.

**Vulpés Power bank Description, Function and Intended Use:**

The product is made for charging or operating mobile devices and Vulpés specific products running with 5 V DC. The Power Bank comes with integrated Lithium Polymer battery, operating as energy storage.

We recommend to use Vulpés Powerbanks or powerbanks from brand suppliers. **Please carefully check the instructions of 3rd party powerbanks before using in connection with the Vulpés belt.**

# IMPORTANT



**Note: In case of any issues with the product please contact the Vulpés team or your local distributor.**



**Maintenance and Care:** In case of disfunction please let a professional do the maintenance and repair work. Please only use a dry and soft cloth for cleaning. Do not use detergents and chemicals.



**Safety Instructions - Please do not do following aspects:**

**Do not modify product, product parts or accessories.**

**Do not open the housing.**

**Avoid extreme conditions, such as extreme heat, direct sun light, as well as micro waves, vibrations and mechanical pressure.**

**Do not throw into fire.**

**NEVER let the product unattended when in use!**

**Do not short-circuit ports and circuits.**


**Dispose leaked, deformed or corroded by appropriate protectives.**

**Do not use in explosive ambiance.**

**When not using for longer time, charge every 3 months to keep the capacity.**



**Material damage by improper use:** Compare the specifications of product, power supply and peripherals. These must be identical. Use product, product parts and accessories only in perfect condition.

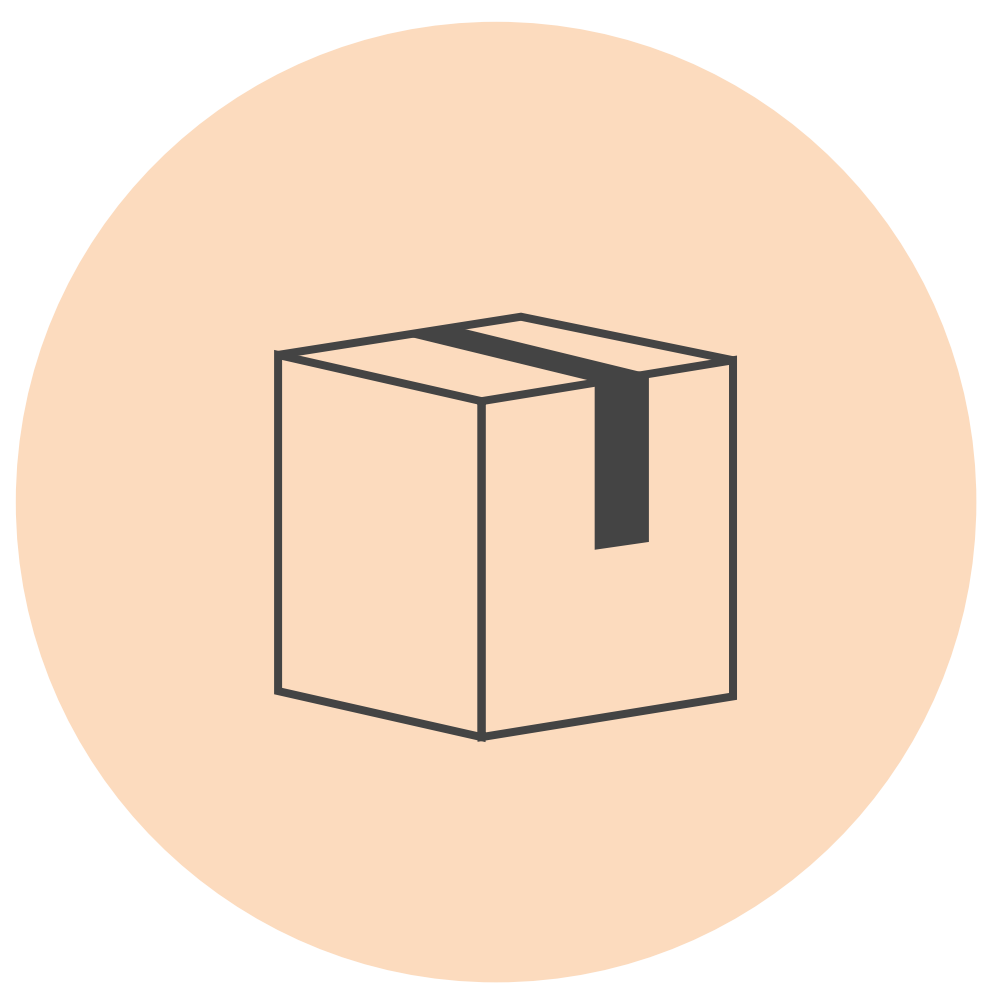

**Storage:** Please handle the device with care and keep it in a clean and dust-free place. **Please do not expose it to flames, drop, throw or try to bend or cut** the product.

**Please do not attempt to repair any electrical belt part or component. This will void the warranty and can result in temperature control malfunctions, which may cause minor to severe burning.**

**Please do not cut or amend the product. This will void the warranty and can affect the function of the product.**


Please do not cut the belt and keep it away from sharp objects. A cut of the belt can damage the product. In case the belt has been cut by accident please remove the product from the body and do not use it again.



**Power-off:** Please make sure that the device is in the off state when worn overnight and switch it off before falling asleep.

Users should be aware that different environmental conditions including high or low temperatures can change the characteristics of the Belt and may reduce the performance as well as the battery runtime.

Please do not activate the powerbank and the Vulpés belt in frozen state.

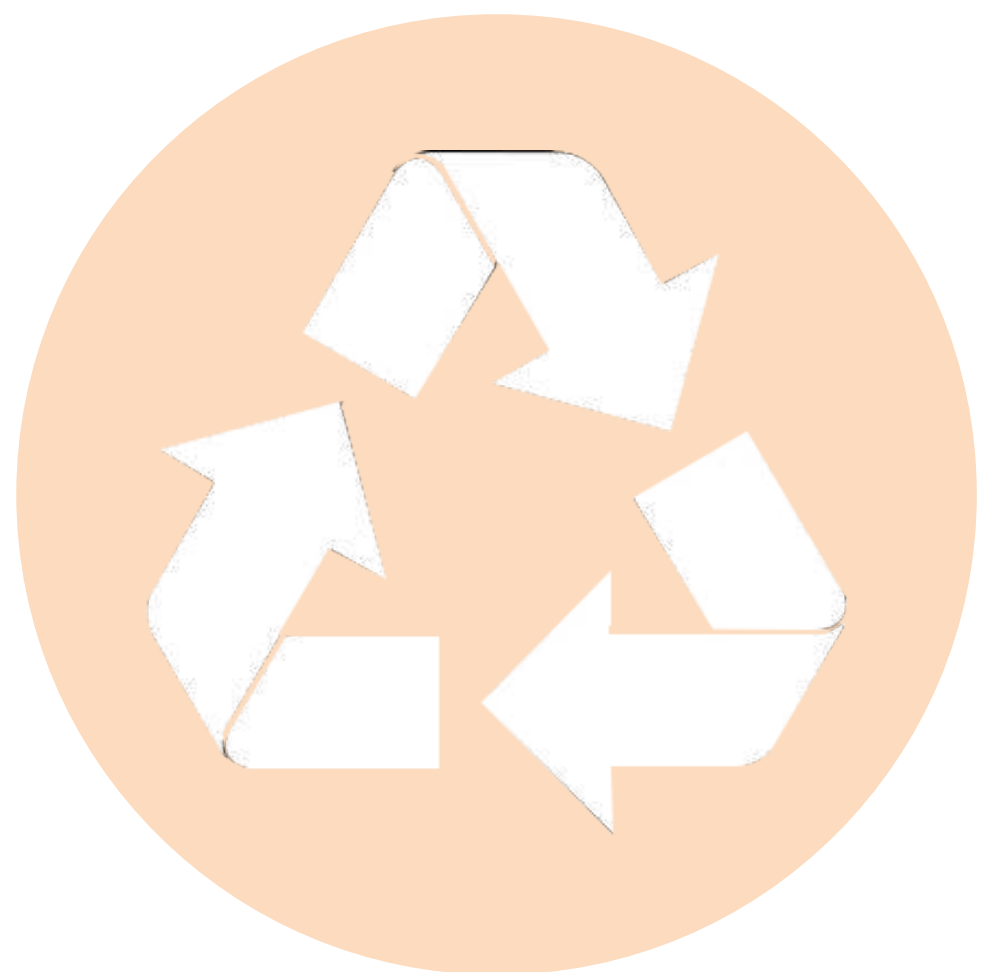

**Recycle:** According to the European WEEE directive, electrical and electronic equipment must not be disposed with consumers waste. Its components must be recycled or disposed apart from each other. Otherwise, contaminant and hazardous substances can pollute our environment. As a consumer, you are committed by law to dispose electrical and electronic devices to the producer, the dealer, or public collecting points at the end of the devices lifetime for free. Particulars are regulated in national right. The symbol on the product, in the user manual, or at the packaging alludes to these terms. With this kind of waste separation, application, and waste disposal of used devices you achieve an important share to environmental protection. **WEEE No. 28527803**



**Our Warranty:** Subject to the conditions of this Limited Warranty, Vulpés warrants this Product to be free from defects in design, material and workmanship at the time of its original purchase by a consumer. This Limited Warranty will last for a period of one (1) year as from the original date of purchase of the Product.

Please be aware that given the physio-chemical state of the Lithium-Polymer technology, the battery runtime might decrease after a specific time of use. The higher usage intensity accelerate the battery degradation process. **For the powerbank Vulpés does not provide any company warranty.**

For further details please visit the Vulpés website  
<https://www.vulpes-electronics.com/warranty>

# IMPORTANT

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Any Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

**Note:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

RF warning for Portable device: The device has been evaluated to meet general RF exposure requirement. The device can be used in portable exposure condition without restriction.